

# **Temperature Measurement and Control Catalog**











For over 45 years, Lake Shore Cryotronics has provided reliable and robust temperature measurement and control products for customers like you. Supporting your low-temperature research has been our primary pursuit since the beginning.

This Temperature Measurement and Control Catalog is designed to assist you, our valued customer, in finding the most appropriate solution for your specific cryogenic application.

In these pages, you'll see a number of solutions for measurement, control, and monitoring. Featured products include industry-leading Cernox® and silicon diode sensors, ultra-low-temperature controllers, an AC resistance bridge, and a superconducting magnet power supply. But you'll also see a number of new products. These additions are reflective of our commitment to continuous innovation in cryogenic product development.

In addition to product selection guides, the catalog includes a wealth of technical information, including technical data, performance characteristics, and a comprehensive reference section. The Appendices section is where you can learn more about thermometry, sensor characteristics, and tips for proper sensor installation, as well as common units and conversions charts, cryogenic reference tables, and other useful information.

Customers tell us that our catalogs, like our website and product documentation, are very useful in helping them identify their needs and improve their research. We hope this catalog is as valuable of a resource as Lake Shore's other materials.

If you cannot find the information that you need, please contact Lake Shore. Our global sales and service representatives are ready to answer your questions and to help you find the right solution for your application.

Thank you for your support and feedback as we continue to innovate and create products that advance your pursuit of science.

ikal A. Swarts

Michael S. Swartz President and CEO

The technical and pricing information contained herein is subject to change at any time.

CalCurve, SoftCal, Duo-Twist, Quad-Twist, Quad-Lead, and Rox are all trademarks of Lake Shore Cryotronics, Inc.

Cernox and Lake Shore Cryotronics are both registered trademarks of Lake Shore Cryotronics, Inc.

American Express and the American Express Box Logo are registered trademarks of American Express Company; Apiezon is a registered trademark of M & I Materials, Ltd.; CryoCable is a trademark of Omega Engineering, Inc.; Dacron, Kapton, Mylar, Teflon, and Vespel are registered trademarks of E. I. Du Pont de Nemours & Co.; Evanohm is a registered trademark of Carpenter Technology Corp.; Fair-Rite is a registered trademark of Fair-Rite Products Corp.; Gel-Pak is a registered trademark of GEL-PAK, LLC; Kester is a registered trademark of Litton Systems, Inc.; Kynar is a registered trademark of Atochem North America, Inc.; LabVIEW is a trademark of National Instruments Corporation; Lemo is a registered trademark of Lemo USA, Inc.; MasterCard and the Distinctive Interlocking Circles Design are registered trademarks of MasterCard International Incorporated; Ostalloy is a registered trademark of Umicore; Pyrex is a registered trademark of Corning Incorporated; Scotch is a registered trademark of 3M Company; Stay-Silv is a registered trademark of J.W. Harris Co., Inc.; Stycast is a registered trademark of Emerson & Cuming; Swagelok is a registered trademark of Visa.

All other trademarks or service marks noted herein are either property of Lake Shore Cryotronics, Inc., or their respective companies.



# Contents

	Introduction	5
	Sensors	9
F	Instruments	63
I	Accessories	143
	Appendices	163
	Customer Service	227



# Lake Shore Cryotronics, Inc. Company Overview

#### Over 45 years of cryogenic excellence

Leading researchers around the world trust Lake Shore for advanced solutions that drive the discovery and development of new materials for tomorrow's technologies. In electronics, clean energy, nanotechnology, and many other applications, Lake Shore provides the products and systems needed for precise measurements over a broad range of temperature and magnetic field conditions.

Serving the needs of the research community since 1968, Lake Shore has grown its product solutions to keep pace with evolving interests in scientific exploration, from the physics lab to deep space.

Lake Shore provides solutions to an international base of research customers at leading university, government, aerospace, and commercial research institutions and is supported by a global network of sales and service facilities.

Through our global technical service and sales teams, we foster a culture of collaboration and innovation, and a commitment to the pursuit of science. When you work with Lake Shore, you're dealing with a company that is run by scientists and engineers for the purpose of advancing the work of scientists in many research fields.

This ongoing interaction with the global research community spurs continual innovation of our product offering, advancements that will, in turn, enable researchers to explore new phenomena for new insight into cryogenic measurement.



### Our Mission

# Lake Shore is committed to our customers' pursuit of the science that benefits mankind.

We want to advance science by providing easy to use, high value and high performance products.

We support an international base of scientists and researchers as they explore and develop tomorrow's technologies. To us, "Advancing Science" is more than just a tagline. Helping customers around the world in their pursuit of science is the driving force behind our ongoing product research and development.

Lake Shore stands behind its customers, with exceptional, technically knowledgeable customer service and a three-year warranty on every product we sell.

Lake Shore is ISO 9001:2008 certified, a sign of our commitment to continuous improvement.



### Online Resources

#### The Lake Shore website and blog

#### Our website and blog offer:

- Software and firmware downloads
- Application notes
- Installation videos
- Calibration service information
- RFQ and pricing information
- Quality control and ISO documentation
- RoHS and WEEE certificates
- Product registration forms
- Terms and conditions
- Warranty information
- Shipping information
- Sales contact information
- Catalog request forms
- New product announcements
- Product roadmap news
- Lake Shore press releases
- Directions to Lake Shore

#### Plus see our website for our full product offering, including:

- Hall effect measurements systems
- Cryogenic probe stations
- VSM and AGM systems
- Terahertz materials characterization system
- Gaussmeters, Hall probes, Hall magnetic sensors, fluxmeters, and Helmholtz and search coils
- And a number of other solutions

#### Also, follow us at::









# Sensors

- 10 Sensor Selection Guide
- 14 Sensor Characteristics
- 20 Sensor Packages and Mounting Adapters
- **27 Temperature Probes**
- 32 Cernox® RTDs
- 36 DT-670 Silicon Diodes
- 41 GaAlAs Diodes
- 44 Germanium RTDs
- 47 Ruthenium Oxide (Rox™) RTDs
- 51 PT-100 Series Platinum RTDs
- 54 Rhodium-Iron RTDs
- 56 Capacitance Temperature Sensors
- 58 Thermocouple Wire
- 59 Cryogenic Hall Generators and Probes

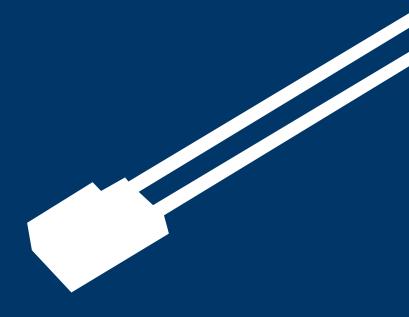

## Sensor Selection Guide

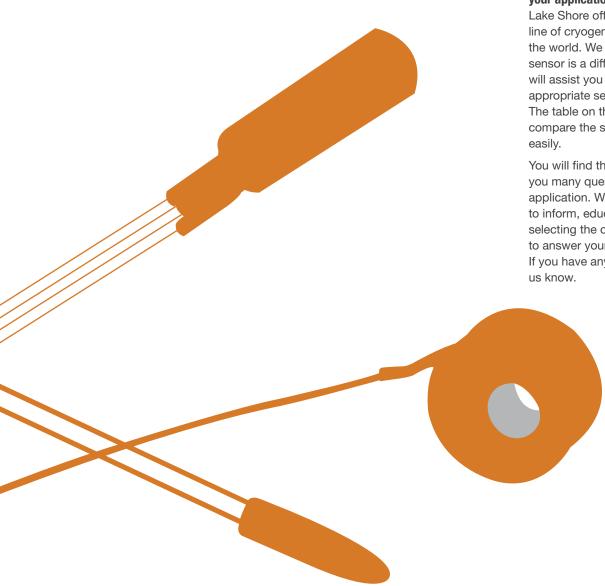

# How to select a temperature sensor for your application

Lake Shore offers the most comprehensive line of cryogenic temperature sensors in the world. We understand that selecting a sensor is a difficult procedure. This catalog will assist you in selecting the most appropriate sensor for your application. The table on the next page is designed to compare the sensor characteristics more easily.

You will find that our sales staff will ask you many questions regarding your application. We ask a lot of questions to inform, educate, and to assist you in selecting the correct sensor. We are here to answer your questions and concerns. If you have any specific needs, please let us know

Any one or several of the following environmental factors may be important to you in selecting a sensor:

- Temperature range
- Package size
- Fast thermal response time
- Fast electrical response time
- Heat sinking
- Small thermal mass
- Robustness
- Compatibility with harsh environments
  - Magnetic fields
  - Ionizing radiation
  - Ultra high vacuum (UHV)
  - Vibration/mechanical shock
  - Thermal shock
  - Temperatures above 323 K
- Easily measured signal
- Compatibility with sources of error
  - Thermal EMFs
  - Self-heating
  - Noise pickup
- High sensitivity
- High accuracy\*
- High repeatability—long and short term
- Low power dissipation
- Interchangeability
- Ease of use
- Low cost
- Available accessories
- Available instrumentation

#### Sensor overview

	Temperature range	Standard curve	Below 1 K	Can be used in radiation	Performance in magnetic field
Diodes					
Silicon	1.4 K to 500 K	×			Fair above 60 K
GaAlAs	1.4 K to 500 K				Fair
Positive temperature co	pefficient RTDs				
Platinum	14 K to 873 K	×		×	Fair above 30 K
Rhodium-iron	0.65 K to 500 K		×	×	Fair above 77 K
Negative temperature	coefficient RTDs				
Cernox®	0.10 K to 420 K		×	×	Excellent above 1 K
Germanium	0.05 K to 100 K		×	×	Not recommended
Ruthenium oxide*	0.01 K to 40 K	×	×	×	Good below 1 K
Other					
Thermocouples	1.2 K to 1543 K	×			Fair
Capacitance	1.4 K to 290 K				Excellent

<sup>\*</sup>RX-102B not recommended for use in magnetic fields

Unfortunately, you can't have it all in one sensor. The most stable and accurate temperature sensors are very large, have slow response times and are extremely fragile. The sensors with the highest sensitivity and resolution have the smallest range. Choosing the appropriate sensor for a particular application necessitates prioritizing the requirements for that application.

The sensors described in this catalog are manufactured for the rigors of cryogenic environments, and are designed with specific applications in mind. For much of its 45 year history, Lake Shore has focused on cryogenic sensors used for the precise measurement of temperatures from near absolute zero to well above room temperature.

As you continue through the Sensor section of the catalog, you will notice that information is presented in both graphical format as well as in more detailed specifications, pertaining to topics such as the sensor's highlights, typical magnetic field-dependent data, resistance, and sensitivity values.

Characteristics such as packaging are incorporated into each sensor's design with the customer in mind. To learn more about what package would be best for your application, please refer to the Sensor Packages and Mounting Adapters section. For more detailed information, see Appendix C.

<sup>\*</sup> The use of the terms accuracy and uncertainty throughout this catalog are used in the more general and conventional sense as opposed to following the strict metrological definitions. For more information, see Appendix B: Accuracy versus Uncertainty.



# Sensor Types

#### **Cernox®**

Cernox® sensors can be used from 100 mK to 420 K with good sensitivity over the whole range. They have a low magnetoresistance, and are the best choice for applications with magnetic fields up to 30 T (for temperatures greater than 2 K). Cernox® are resistant to ionizing radiation, and are available in robust mounting packages and probes. Because of their versatility, they are used in a wide variety of cryogenic applications, such as particle accelerators, space satellites, MRI systems, cryogenic systems, and research science.

#### Silicon diodes

Silicon diodes are the best choice for general-purpose cryogenic use. The sensors are interchangeable (they follow a standard curve) and are available in robust mounting packages and probes. Silicon diodes are easy and inexpensive to instrument, and are used in a wide variety of cryogenic applications, such as cryo-coolers, laboratory cryogenics, cryo-gas production, and space satellites.

#### **GaAlAs diodes**

GaAlAs diodes offer high sensitivity over a wide range of use (1.4 K to 500 K). They are useful in moderate magnetic fields, and offer many of the advantages of silicon diodes—easy to instrument, wide range, and robust packaging. They do not follow a standard curve. GaAlAs diodes are used in moderate magnetic field applications when instrumentation constraints (e.g., legacy installations, cost) prevent the use of Cernox®.

#### Germanium

Germanium RTDs have the highest accuracy, reproducibility, and sensitivity from 0.05 K to 100 K. They are resistant to ionizing radiation, but are not recommended for use in magnetic fields. Germanium RTDs are used mostly in research settings when the best accuracy and sensitivity are required. Germanium and Ruthenium oxide are the only two sensors that can be used below 100 mK.

#### Ruthenium oxide (Rox™)

Ruthenium oxide RTDs can be used to below 10 mK. Their unique advantage is that they have a low magnetoresistance and follow a standard curve (with the exception of the RX-102B). Their upper temperature range is limited to 40 K, and Cernox® are better in magnetic fields above 2 K. Rox™ sensors are often used for applications that require a standard curve in magnetic fields, such as MRI systems. Along with germanium, they are the only sensors that can be used below 100 mK.

#### **Platinum**

Platinum RTDs are an industry standard. They follow an industry standard curve from 73 K to 873 K with good sensitivity over the whole range. Platinum RTDs can also be used down to 14 K. Because of their high reproducibility, they are used in many precision metrology applications. Platinum RTDs have limited packaging options, but they are inexpensive and require simple instrumentation. They are widely used in cryogenic applications at liquid nitrogen temperatures or greater.

#### **Rhodium-iron**

Rhodium-iron temperature sensors can be used over a wide temperature range, and are resistant to ionizing radiation. Lake Shore RF-800s have excellent stability and are widely used as secondary temperature standards by many national standards laboratories.

#### **Capacitance**

Capacitance sensors are ideally suited for use as temperature control sensors in strong magnetic fields because they exhibit virtually no magnetic field dependence.

Small variations in the capacitance/ temperature curves occur upon thermal cycling. It is recommended that temperature in zero field be measured with another temperature sensor, and that the capacitance sensor be employed as a control element only.

#### **Thermocouples**

Thermocouples can be used over an extremely wide range and in harsh environmental conditions, and follow a standard response curve. Less accurate than other sensors, special techniques must be employed when using thermocouples to approach temperature accuracies of 1% of temperature. Thermocouples are used for their small size, extremely wide temperature range (exceeding high temperature limits of platinum RTDs), and simple temperature measurement methodology.



#### Lake Shore calibrations

Lake Shore offers complete calibration services from 50 mK to 800 K. Above 0.65 K, Lake Shore calibrations are based on the International Temperature Scale of 1990 (ITS-90). For temperature below 0.65 K, calibrations are based on the Provisional Low Temperature Scale of 2000 (PLT-2000).

Each scale is maintained on a set of germanium, rhodium-iron, and/or platinum resistance secondary thermometers standards. These secondary standards are calibrated at various national labs: NIST, PTB, and NPL. Working thermometers are calibrated against, and routinely intercompared with these secondary standards. For PLTS-2000 calibrations, working sensors are also compared to a superconducting fixed-point set and nuclear orientation thermometer.

Lake Shore offers sensor calibrations down to 20 mK. Our enhanced ultra-low temperature calibration facility includes dilution refrigerators, a nuclear orientation thermometer, and a superconducting fixed point set.

All calibration reports include:

- Certificate of calibration
- Calibration test data and data plot
- Polynomial fit equations and fit comparisons
- Interpolation tables
- Instrument breakpoint tables and data files

Lake S	hore offers thre	e classifications of calibration:
Good	Uncalibrated	Silicon diodes follow standard curve
		Platinum resistors follow standard curve
		Ruthenium oxide (Rox™) resistors follow standard curve (except RX-102B)
		GaAlAs diode, Cernox®, germanium, Rox™ RX-102B, and rhodium-iron sensors can be purchased uncalibrated but must be calibrated by the customer
Better	SoftCal™	An abbreviated calibration (2-point: 77 K and 305 K; or 3-point: 77 K, 305 K, and 480 K) which is available for platinum sensors
Best	Calibration	All sensors can be calibrated in the various pre-defined temperature ranges.  Lake Shore has defined calibration ranges available for each sensor type. The digits represent the lower range in kelvin, and the letter corresponds to high temperature limit, where:
		A = 6  K $B = 40  K$ $D = 100  K$ $L = 325  K$ $M = 420  K$ $H = 500  K$ $J = 800  K$
		For example: The calibration range "1.4L" would result in a sensor characterized from 1.4 K to 325 K



# **Sensor Characteristics**

#### Sensor packages and characteristics


	Sensor type/ packages		erature 1ge	Physical size <sup>1</sup>	Mass	Typical dimensionless sensitivity S <sub>p</sub>							
		low	high			1.4 K	4.2 K	20 K	77.4 K	295 K	475 K		
	DT-670-SD	1.4 K	500 K	1.08 mm high $\times$ 1.905 mm wide $\times$ 3.175 mm long	37 mg	-0.01	-0.08	-0.26	-0.13	-1.19	-7.5		
les	DT-670E-BR	30 K	500 K	0.178 mm × 0.432 mm × 0.406 mm	72.7 µg	-0.01	-0.08	-0.26	-0.13	-1.19	_		
Silicon diodes	DT-414	1.4 K	375 K	0.5 mm high $\times$ 0.635 mm $\times$ 1.524 mm long	3 mg	-0.01	-0.09	-0.29	-0.15	-1.3	_		
OI (	DT-421	1.4 K	325 K	$0.762 \text{ mm high} \times 1.27 \text{ mm dia.}$	23 mg	-0.01	-0.09	-0.29	-0.15	-1.3	_		
.≘	DT-470-SD	1.4 K	500 K	1.08 mm high × 1.905 mm wide × 3.175 mm long	37 mg	-0.01	-0.09	-0.20	-0.15	-1.3	-11.6		
S	DT-471-SD	10 K	500 K	1.08 mm high $\times$ 1.905 mm wide $\times$ 3.175 mm long	37 mg		_	-0.29	-0.15	-1.3	-11.6		
လ လ	TG-120-P	1.4 K	325 K	2.794 mm long × 3.048 mm dia.	79 mg	-0.03	-0.19	-0.77	-0.07	-0.9	-4.0		
GaAlAs	TG-120-PL	1.4 K	325 K	1.335±0.3175 mm long × 1.333±0.3175 mm thick	20 mg	-0.03	-0.19	-0.77	-0.07	-0.9	-4.0		
g g g	TG-120-SD	1.4 K	500 K	1.08 mm high $\times$ 1.905 mm wide $\times$ 3.175 mm long	38 mg	-0.03	-0.19	-0.77	-0.07	-0.9	-4.0		
	CX-1010-BC	0.1 K	325 K	3.175 mm × 8.89 mm × 7.620 mm	3.0 mg	-0.68	-0.49	-0.44	-0.56	-0.65	_		
	CX-1010-SD	0.1 K	325 K	1.08 mm high $\times$ 1.905 mm wide $\times$ 3.175 mm long	40 mg		—	_			_		
	CX-1010-AA	0.1 K	325 K	3.048 mm dia. × 8.509 mm long	400 mg		_	_	_	_	_		
	CX-1030-BC	0.30 K	325 K	3.175 mm × 8.89 mm × 7.620 mm	3.0 mg	-1.15	-0.71	-0.56	-0.63	-0.64	_		
	CX-1030-SD-HT	0.30 K	420 K	1.08 mm high $\times$ 1.905 mm wide $\times$ 3.175 mm long	40 mg	_	_	_	_	_	_		
	CX-1030-AA	0.30 K	325 K	3.048 mm dia. × 8.509 mm long	400 mg		_	_	_	_	_		
@_	CX-1050-BC	1.4 K	325 K	3.175 mm × 8.89 mm × 7.620 mm	3.0 mg	-2.5	-1.3	-0.9	-0.91	-0.87	_		
Cernox®	CX-1050-SD-HT	1.4 K	420 K	1.08 mm high $\times$ 1.905 mm wide $\times$ 3.175 mm long	40 mg	_	_	_	_	_	_		
Cer	CX-1050-AA	1.4 K	325 K	3.048 mm dia. × 8.509 mm long	400 mg		_	_		_	_		
	CX-1070-BC	4.2 K	325 K	3.175 mm × 8.89 mm × 7.620 mm	3.0 mg		-1.5	-1	-1	-0.9	_		
	CX-1070-SD-HT	4.2 K	420 K	1.08 mm high $\times$ 1.905 mm wide $\times$ 3.175 mm long	40 mg		_			_	_		
	CX-1070-AA	4.2 K	325 K	3.048 mm dia. × 8.509 mm long	400 mg		_	_	_	_	_		
	CX-1080-BC	20 K	325 K	3.175 mm × 8.89 mm × 7.620 mm	3.0 mg		_	-1.5	-1.4	-1.2	_		
	CX-1080-SD-HT	20 K	420 K	1.08 mm high $\times$ 1.905 mm wide $\times$ 3.175 mm long	40 mg		_	_	_	_	_		
	CX-1080-AA	20 K	325 K	3.048 mm dia. × 8.509 mm long	400 mg		_	_		_	_		
⊂	CGR-1-500	1.4 K	325 K	3.048 mm dia. × 8.509 mm long	330 mg	-6.9	-3.1	-0.98	-0.48	-0.33	_		
Carbon glass	CGR-1-1000	1.4 K	325 K	3.048 mm dia. × 8.509 mm long	330 mg	-7.8	-3.5	-1.1	-0.53	-0.35	_		
යු ල	CGR-1-2000	1.4 K	325 K	3.048 mm dia. × 8.509 mm long	330 mg	-8.4	-3.8	-1.2	-0.56	-0.37	_		
Ε	GR-50-AA	0.05 K	5 K	3.048 mm dia. × 8.509 mm long	355 mg	-0.74	-0.32	_	_	_	_		
ie	GR-300-AA	0.3 K	100 K	3.048 mm dia. × 8.509 mm long	355 mg	-1.8	-1.2	-0.93	-1.1	_	_		
Germanium	GR-1400-AA	1.4 K	100 K	3.048 mm dia. × 8.509 mm long	355 mg	-3.7	-2.1	-2.4	-1.1	_	_		
	RX-102A-BR	0.05 K	40 K	1.45 mm $\times$ 1.27 mm $\times$ 0.65 mm thick	2.8 mg	-0.47	-0.25	-0.07			_		
	RX-102A-AA	0.05 K	40 K	3.048 mm dia. × 8.509 mm long	350 mg	-0.47	-0.25	-0.07		_	_		
E.	RX-102B-CB	0.01 K <sup>2</sup>	40 K	14.6 mm high $\times$ 6.4 mm wide $\times$ 6.4 mm long	3.5 g	-0.16	-0.11	-0.12	_	_	_		
Rox™	RX-202A-AA	0.05 K	40 K	3.048 mm dia. × 8.509 mm long	350 mg	-0.34	-0.17	-0.10	_	_	_		
	RX-103A-BR	1.4 K	40 K	1.40 mm $\times$ 1.23 mm $\times$ 0.41 mm thick	3.7 mg	-0.62	-0.36	-0.17	_	_	_		
	RX-103A-AA	1.4 K	40 K	3.048 mm dia. × 8.509 mm long	350 mg	-0.62	-0.36	-0.17	_	_	_		
_	PT-102	14 K	873 K	2.007 mm dia. × 20.995 mm long	250 mg	_	_	+0.74	+1.6	+1.1	+1.0		
III I	PT-103	14 K	873 K	1.6 mm dia. × 12.192 mm long	120 mg		_	+0.74	+1.6	+1.1	+1.0		
Platinum	PT-111	14 K	673 K	1.8 mm dia. × 5 mm long	52 mg	_	_	+0.74	+1.6	+1.1	+1.0		
Ē	RF-100-BC	1.4 K	325 K	1.3 mm wide $\times$ 3.8 mm long $\times$ 0.38 mm	7 mg	+0.10	+0.21	+0.23	+1.1	+1.0	_		
hodium	RF-100-AA	1.4 K	325 K	3.048 mm dia. × 8.509 mm long	360 mg	+0.10	+0.21	+0.23	+1.1	+1.0	_		
뜶	RF-800-4	0.65 K	800 K	3.175 mm dia. × 20.32 mm long	735 mg	+0.16	+0.29	+0.29	+1.1	+1.0	_		
Thermocouples Capacitance	CS-501-GR	1.4 K	290 K	3.048 mm dia. × 8. 484 mm long	260 mg	+0.01	+0.02	+0.11	+0.46	-4.4	_		
es	Type K	3.2 K	1543 K	30 AWG (0.254 mm) & 36 AWG (0.127 mm)				1		1			
ldn	Type E	3.2 K	953 K	30 AWG (0.254 mm) & 36 AWG (0.127 mm)									
ermoco	Chromel-AuFe (0.07%)	1.2 K	610 K	30 AWG (0.254 mm) & 36 AWG (0.127 mm)				NA					

 $<sup>^{1} \ \</sup>text{Adapters will increase thermal response times} \underline{\hspace{0.5cm}} \text{see individual sensor specifications for thermal response times}$ 

<sup>&</sup>lt;sup>2</sup> Calibrations down to 20 mK available; 10 mK calibrations coming soon



#### Sensor package size versus temperature sensor characteristics

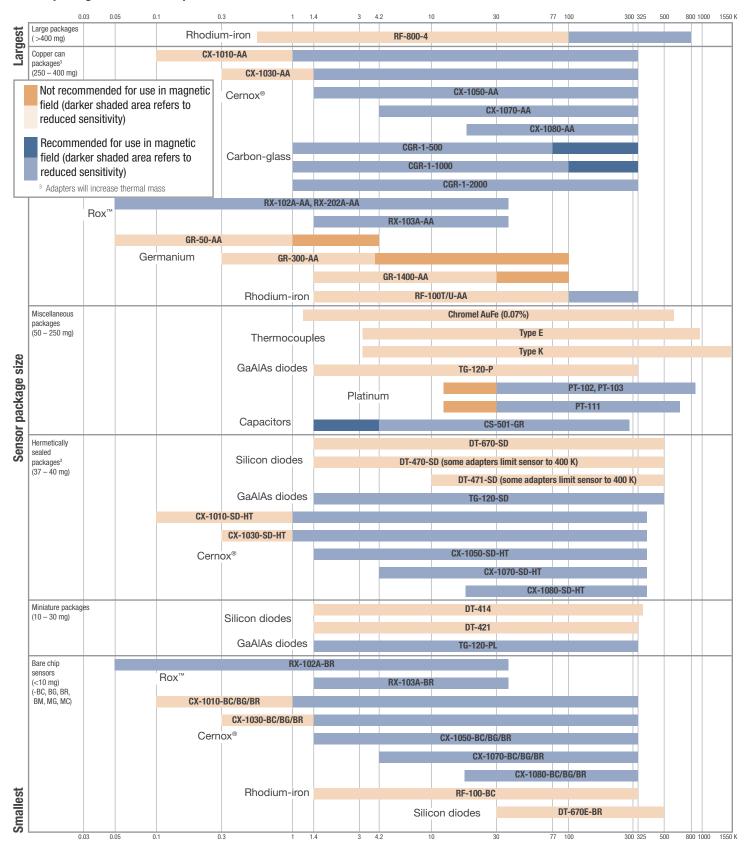



#### **Short and long term sensor characteristics**

	Interchangeability	Typical reproducibility at 4.2 K	Typical lo	ng-term stability
			Use to 305 K <sup>4</sup>	Use to 500 K <sup>5</sup>
Silicon diode	Yes—see page 18	±10 mK	4.2 K: ±10 mK/yr	4.2 K: ±40 mK/yr
			77 K: ±40 mK/yr	77 K: ±100 mK/yr
			305 K: ±25 mK/yr	305 K: ±50 mK/yr
				500 K: ±150 mK/yr
GaAlAs diode	No	±10 mK	4.2 K: ±15 mK/yr	
			77 K: ±15 mK/yr	
			330 K: ±50 mK/yr	
Cernox <sup>®</sup>	No	±3 mK	1 K to 100 K: ±25 m	√yr
			100 K to 300 K: 0.05	% of T
Germanium	No	±0.5 mK	4.2 K: ±1 mK/yr	
			77 K: ±10 mK/yr	
Rox™	Yes	±15 mK	4.2 K: ±15 to 50 mK/	yr (model dependent)
Platinum	Yes—see page 18	±5 mK <sup>6</sup>	77 K to 273 K: ±10 m	nK/yr
Rhodium-iron	No	±5 mK	1.4 K to 325 K: ±10 r	mK/yr
Capacitance	No	±0.01 K after cooling and stabilizing	±1.0 K/yr	
Thermocouples				
Type K	Yes—see ASTM standard	NA	NA	
Type E	Yes—see ASTM standard	NA	NA	
Type T	Yes—see ASTM standard	NA	NA	
Chromel-AuFe (0.07%)	Yes—see ASTM standard	NA	NA	

 $<sup>^{\</sup>rm 4}$  Long-term stability data is obtained by subjecting sensor to 200 thermal shocks from 305 K to 77 K

#### Sensor characteristics in various environments

		Use in vacuum		Use in radiation <sup>7</sup>	Use in magnetic fields <sup>7</sup>
	High	Very high	Ultra high		
	10 <sup>-1</sup> to 10 <sup>-4</sup> Pa	10-4 to 10-7 Pa	10 <sup>-7</sup> to 10 <sup>-10</sup> Pa		
Silicon diode	DT-421	_	DT-670-SD DT-414 DT-470-SD DT-471-SD	Not recommended	Not recommended for T<60 K, or for B>5 tesla above 60 K SD package has magnetic leads
GaAlAs diode	TG-120-P	_	TG-120-SD TG-120-PL	Not recommended	Relatively low field dependence DT/T(%) $\leq$ 4% for B<5 tesla and T $\geq$ 4.2 K; SD package with non-magnetic leads
Cernox <sup>®8</sup>	AA can	_	Bare chip SD	Recommended	Excellent for use in magnetic fields 1 K and up SD package with non-magnetic leads
Germanium <sup>8</sup>	AA can	_	Bare chip	Recommended	Not recommended for use except at low B due to large orientation-dependent magnetic field effect
Rox™	AA can	_	Bare chip	Recommended	Excellent for use in magnetic fields
Platinum	_	PT-111	_	Recommended	Moderately orientation dependent—suggested use only T $\geq$ 30 K
Rhodium-iron <sup>8</sup>	_	RF-800-4	_	Recommended	Not recommended below 77 K
Capacitance	CS-501	_	_	Not available	Recommended for control purposes
Thermocouples	Insulated wire	_	_	Recommended	Useful when T $\geq$ 10 K

<sup>&</sup>lt;sup>7</sup> See additional information in Appendix A: Overview of Thermometry

<sup>&</sup>lt;sup>5</sup> Based on 670 h of baking at 500 K

<sup>&</sup>lt;sup>6</sup> Platinum reproducibility tested at 77 K

<sup>8</sup> Adapters with epoxy are limited to a bakeout temperature of 127 °C



### Typical magnetic field-dependent temperature errors, ΔT/T (%), at B (magnetic induction)

	فيست		Magnetic	flux density B		
	T(K)	2.5 T	8 T	14 T	19 T	Notes
Cernox® 1050	2	1.3	3.1	3.9	5	Best sensor for use in magnetic field (T > 1 K)
(CX series)	4.2	0.1	-0.15	-0.85	-0.8	
(	10	0.04	-0.4	-1.1	-1.5	
	20	0.04	0.02	-0.16	-0.2	
	30	0.01	0.04	0.06	0.11	
	77	0.002	0.022	0.062	0.11	_
	300	0.003	0.004	0.004	0.006	
Rox <sup>™</sup> 102A	2	-1.4	-7.9	-13	-17	Recommended for use over the 0.05 K to 40 K temperature
	3	-1.5	-7	-14	-18	range. Consistent behavior between devices in magnetic
	4	-0.56	-6.7	-14	-18	fields.
	8	-1.3	-6.1	-13	-21	_
	16	-0.40	-3.4	-9.6	-16	_
DIM 100D	23	-0.31	-2.2	-6.2	-11	
Rox <sup>™</sup> 102B	2	3.29	13.82	22.53	27.95	
	3	3.96	14.68	23.12	29.12	_
	4	3.53	13.92	22.57	28.20	
	8 16	1.53 0.27	7.53 2.14	13.50	17.86 6.58	_
	23	0.27	0.79	4.66 2.01	3.11	
Rox™ 103A	23	0.58	1.5	2.01	2.6	Excellent for use in magnetic fields from 1.4 K to 40 K.
RUX TUSA	3	0.36	1.1	1.7	2.0	
	4	0.44	0.95	1.4	1.7	Predictable behavior.
	8	0.27	0.95	0.71	0.80	_
	16	0.018	0.49	0.089	0.040	_
	23	0.0051	0.0058	-0.0060	0.095	_
Rox™ 202A	2	-0.13	-2.2	-3.9	-5.2	Recommended for use over the 0.05 K to 40 K temperature
TION ZOZA	3	0.18	-0.68	-2.7	-3.7	·
	4	0.77	0.046	-1.8	-3.2	range. Consistent behavior between devices in magnetic
	8	-0.023	0.16	-0.65	-3.0	fields.
	16	0.03	0.16	-0.48	-1.5	
	23	-0.05	-0.08	-0.39	-0.92	_
Platinum Resistors	20	20	100	250		Recommended for use when $T \ge 40$ K.
(PT series)	40	0.5	3	6	8.8	_
(FT SCHOS)	87	0.04	0.4	1	1.7	_
	300	< 0.01	0.02	0.07	0.13	
Rhodium-iron	4.2	11	40	_	_	Not recommended for use below 77 K in
(RF series)	40	1.5	12	30	47	magnetic fields.
(III dollod)	87	0.2	1.5	4	6	_ magnotic notaci
	300	< 0.01	0.1	0.4	_	
Capacitance CS-501-GR series		$\Delta T/T(\%) < 0.01$	5 at 4.2 K and 18.7	7 tesla		Recommended for control purposes.
		$\Delta T/T(\%) < 0.05$	at 77 K and 305 K	and 18.7 tesla		Monotonic in C vs. T to nearly room temperature.
Germanium resistors	2.0	-8	-60	_	_	Not recommended except at low B owing to large,
(GR series)	4.2	-5 to -20	-30 to -55	-60 to -75	_	orientation-dependent temperature effect.
(dit series)	10	-4 to -15	-25 to -60	-60 to -75	_	onontation dependent temperature eneet.
	20	-3 to -20	-15 to -35	-50 to -80		
Chromel-AuFe (0.07%)	10	3	20	30		Data taken with entire thermocouple in field, cold junction a
,	45	1	5	7	_	4.2 K; errors in hot junction.
	100	0.1	0.8	_	_	, 0
	100	0.1	0.0			
Type E thermocouples	10	1	3	7	_	Useful when T ≥ 10 K.
Type E thermocouples (chromel-constantan)						Useful when T $\geq$ 10 K.  Refer to notes for Chromel-AuFe (0.07%).

	T(K)	1T	2 T	3 T	4 T	5 T	Notes
Silicon diodes	4.2	-200	-300	-350	-400	-500	Strongly orientation dependent.
Junction parallel to field	20	-10	-20	-25	-30	-40	
	40	-4	-6	-8	-10	-12	
(DT series)	60	-0.5	-1	-2	-3	-3.5	_
	80	<0.1	-0.5	-0.8	-1.1	-1.5	
	300	< 0.1	<-0.1	<-0.1	<-0.1	<-0.1	
Silicon diodes	4.2	-8	-9	-11	-15	-20	Strongly orientation dependent.
Junction perpendicular to field	20	-4	-5	-5	-5	-10	_
' '	40	-1.5	-3	-4	-5	-5.5	_
(DT series)	60	-0.5	-1	-2	-3	-3.5	_
	80	-0.1	-0.3	-0.5	-0.6	-0.7	_
	300	<0.1	0.2	0.5	0.6	0.6	
GaAlAs diodes	4.2	2.9	3.8	3.7	2.8	1	Shown with junction perpendicular (package base parallel)
(TG series)	30	0.2	0.2	0.3	0.3	0.2	to applied field B. When junction is parallel to B, induced
(1.0.001.00)	78	<0.1	<0.1	0.17	0.16	0.1	errors are typically less than or on the order of those shown
	300	-0.1	< 0.1	< 0.1	< 0.1	< 0.1	citors are typically less than of our tile order of those shown



#### Typical accuracy\* (interchangeability): uncalibrated sensors

	0.05 K	0.5 K	1.4 K	2 K	4.2 K	10 K	20 K	25 K	40 K	70 K	100 K	305 K	400 K	500 K	670 K
Silicon diode															
DT-470-SD, Band 11	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±0.5 K	±1.0 K	±1.0 K	_				
DT-470-SD, Band 11A	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±1% of temp	±1% of temp	±1% of temp					
DT-470-SD, Band 12	_	_	_	±0.5 K	±0.5 K	±0.5 K	±0.5 K	±1.0 K	±2.0 K	±2.0 K					
DT-470-SD, Band 12A	_	_	_	±0.5 K	±0.5 K	±0.5 K	±0.5 K	±1% of temp	±1% of temp	±1% of temp	_				
DT-470-SD, Band 13	_	_	_	±1.0 K	±1.0 K	±1.0 K	±1.0 K	±1% of temp	±1% of temp	±1% of temp					
DT-471-SD	_	_	_	_	_	±1.5 K	±1.5 K	±1.5 K	±1.5 K	±1.5 K	±1.5 K	±1.5% of temp	$\pm 1.5\%$ of temp	$\pm 1.5\%$ of temp	_
DT-414	_	_	_	±1.5 K	±1.5 K	±1.5 K	±1.5 K	$\pm 1.5\%$ of temp	_	_	_				
DT-421	_	_	_	_	_	_	±2.5 K	±2.5 K	±2.5 K	±2.5 K	±2.5 K	$\pm 1.5\%$ of temp	_	_	_
DT-670-SD, Band A	_	_		±0.25 K	±0.25 K	±0.25 K	±0.25 K	±0.5 K	±0.5 K	±0.5 K					
DT-670-SD, Band B	_	_	_	±0.5 K	±0.5 K	±0.5 K	±0.5 K	±0.5 K	±0.33% of temp	±0.33% of temp	_				
DT-670-SD, Band C	_	_		±1.0 K	±1.0 K	±1.0 K	±1.0 K	±1.0 K	$\pm 0.5\%$ of temp	$\pm 0.5\%$ of temp					
DT-670-SD, Band D	_	_	_	_	_	_	_	$\pm 0.25~\text{K}$	$\pm 0.25~\mathrm{K}$	$\pm 0.25~\mathrm{K}$	±0.25 K	±0.50 K	$\pm 0.2\%$ of temp	$\pm 0.2\%$ of temp	_
DT-670-SD, Band E	_	_	_	_	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±0.25% of temp	±0.25% of temp	±0.25% of temp	_
Platinum															
PT-102	_	_	_	_	_	_	_	_	_	±1.3 K	±1.2 K	±0.5 K	±0.9 K	±1.4 K	±2.3 K
PT-103	_	_	_		_					±1.3 K	±1.2 K	±0.5 K	±0.9 K	±1.4 K	±2.3 K
PT-111	_	_	_	_	_	_	_	_	_	±1.3 K	±1.2 K	±0.5 K	±0.9 K	±1.4 K	±2.3 K
Rox™															
RX-102A-AA	±10 mK	±25 mK	±50 mK	±75 mK	±125 mK	±300 mK	±1.25 K	±1.5 K	±4.0 K	_	_	_	_	_	_
RX-102A-AA-M	±5 mK	±20 mK	±25 mK	±40 mK	±75 mK	±200 mK	±500 mK	±750 mK	±1.5 K	_	_	_	_	_	
RX-202A-AA	±15 mK	±30 mK	±100 mK	±125 mK	±250 mK	±1 K	±2.5 K	±3 K	±5.0 K	_	_	_	_	_	_
RX-202A-AA-M	±10 mK	±25 mK	±50 mK	±75 mK	±150 mK	±500 mK	±1.0 K	±1.5 K	±2.0 K	_	_	_	_	_	
RX-103A-AA	_	_	±150 mK	±180 mK	±400 mK	±1 K	±2.0 K	±2.5 K	±4.0 K	_	_	_		_	_
RX-103A-AA-M	_	_	±50 mK	±75 mK	±100 mK	±300 mK	±700 mK	±1 K	±1.5 K	_	_	_	_	_	_

#### Typical accuracy\*: SoftCal™ (2-point and 3-point soft calibration sensors)

Silicon diode	2 K	4.2 K	10 K	30 K	70 K	305 K	400 K	475 K	FOO V	070.1/
Silicon diode					70 K	200 K	400 K	4/3 K	500 K	670 K
DT-470-SD-2S9 (Band 13) ±	1.0 K	±1.0 K	±1.0 K	±0.25 K	±0.15 K	±0.15 K	±1.0 K	±1.0 K	_	_
DT-471-SD-2S9 (Band 13)	_	_	±1.5 K	±0.25 K	±0.15 K	±0.15 K	±1.0 K	±1.0 K	_	_
DT-421-2S9 (Band 13)	_	_	±2.0 K	±0.25 K	±0.15 K	±0.15 K	_	_	_	_
DT-470-SD-3S10 (Band 13) ±	0.5 K	±0.5 K	±0.5 K	±0.25 K	±0.15 K	±0.15 K	±1.0 K	±1.0 K	_	_
DT-471-SD-3S10 (Band 13) ±	0.5 K	$\pm 0.5~\mathrm{K}$	$\pm 0.5~\mathrm{K}$	±0.25 K	±0.15 K	±0.15 K	±1.0 K	±1.0 K	_	_
Platinum										
PT-102-2S9	_	_	_	_	±0.25 K	±0.25 K	±0.9 K	±1.3 K	±1.4 K	±2.3 K
PT-103-2S9	_	_	_	_	±0.25 K	±0.25 K	±0.9 K	±1.3 K	±1.4 K	±2.3 K
PT-111-2S9	_	_	_	_	±0.25 K	±0.25 K	±0.9 K	±1.3 K	±1.4 K	±2.3 K
PT-102-3S11	_	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±1.4 K	±2.3 K
PT-103-3S11	_	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±1.4 K	±2.3 K
PT-111-3S11	_	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±1.4 K	±2.3 K

\*The use of the terms accuracy and uncertainty throughout this catalog are used in the more general and conventional sense as opposed to following the strict metrological definitions. For more information, see Appendix B: Accuracy versus Uncertainty, page 166.

<sup>&</sup>lt;sup>9</sup>2S (2-point at 77 K and 305 K)

 $<sup>^{10}3</sup>S$  (3-point at 4.2 K, 77 K, and 305 K)

<sup>&</sup>lt;sup>11</sup>3S (3-point at 77 K, 305 K, and 480 K)



#### Typical accuracy\*: calibrated sensors (in mK)<sup>12</sup>

Typical accuracy : calibrated Self	0.02 K		0.1 K	0.3 K	0.5 K	1 K	1.4 K	4.2 K	10 K	20 K	77 K	300 K	400 K	500 K
Silicon diode	0.02 K	0.00 K	0.1 IX	0.0 IX	0.5 K	110	1.410	7.2 IX	IUK	ZUK	// IX	300 K	400 K	300 K
							. 10	. 10	. 10	. 4.4	. 00	. 00	. 45	. 50
DT-670-SD/CO/CU-HT							±12	±12	±12	±14	±22	±32	±45	±50
DT-670-CU/CO/LR/CY/ET/B0							±12	±12	±12	±14	±22	±32		
DT-414								±12	±12	±14	±22	±32		
DT-421							±12	±12	±12	±14	±22	±32		
DT-470-SD/CO/CU-HT							±12	±12	±12	±14	±22	±32	±45	±50
DT-470-B0/BR/CU/CY/ET/LR/MT							±12	±12	±12	±14	±22	±32		
DT-471-SD/CO/CU-HT									±12	±14	±22	±32	±45	±50
DT-471-BO/BR/CU/CY/ET/LR/MT									±12	±14	±22	±32		
GaAlAs diode														
TG-120-P							±12	±12	±12	±14	±22	±32		
TG-120-PL				_			±12	±12	±12	±14	±22	±32		
TG-120-SD/C0							±12	±12	±12	±14	±22	±32	±45	±50
TG-120-CU							±12	±12	±12	±14	±22	±32		
TG-120-CU-HT	_	_	_	_	_	_	±12	±12	±12	±14	±22	±32	±45	±50
Cernox®														
CX-1010-AA/CD/C0/CU/LR/ET/MT/SD	_	_	±3	±3.5	±4.5	±5	±5	±5	±6	±9	±25	±75	_	_
CX-1010-BC	_	_	_		_	_	±5	±5	±6	±9	±25	±75		
CX-1030-AA/CD/C0/CU/LR/ET/MT/SD			_	±3	±4	±5	±5	±5	±6	±9	±25	±75		
CX-1030-BC		_	_				±5	±5	±6	±9	±25	±75		
CX-1050-AA/BC/CD/CO/CU/LR/ET/MT/SD			_				±5	±5	±6	±9	±16	±40		
CX-1070-AA/BC/CD/CO/CU/LR/ET/MT/SD	_	_	_	_		_	_	±5	±6	±9	±16	±40		_
CX-1080-AA/BC/CD/CO/CU/LR/ET/MT/SD			_	_			_		_	±9	±16	±40		
CX-1010-CO/SD/CU-HT	_	_	±3	±3.5	±4.5	±5	±5	±5	±6	±9	±25	±75	_	
CX-1030-C0/SD/CU-HT			_	±3	±4	±5	±5	±5	±6	±9	±16	±40	±65	
CX-1050-CO/SD/CU-HT	_	_	_	_	_	_	±5	±5	±6	±9	±16	±40	±65	
CX-1070-CO/SD/CU-HT								±5	±6	±9	±16	±40	±65	
CX-1080-CO/SD/CU-HT	_	_	_	_	_	_	_	_	_	±9	±16	±40	±65	_
Rox™														
RX-102A-AA/CD	_	±3	±3.5	±4	±4.5	±5.5	±5	±16	±18	±37	_	_	_	_
RX-102B-CB	±2	±4	±4.5	±5	±6	±9	±16	±16	±18	±39				
RX-103A-AA/CD							±5	±17	±22	±38				
RX-202A-AA/CD		±3	±3.5	±4	±4.5	±5.5	±5	±16	±18	±37				
Rhodium-iron		_0	_0.0		_ 1.0	_0.0		_10		_0,				
RF-100T-AA/CD/BC/MC							. 44	. 44	. 10	. 4.4	. 15	. 0.5		
							±11	±11	±12	±14	±15	±25		
RF-100U-AA/CD/BC							±11	±11	±12	±14	±15	±25	. 41	- 40
RF-800-4	_		_				±7	±7	±8	±10	±13	±23	±41	±46
Platinum														
PT-102										±10	±12	±23	±40	±46
PT-103										±10	±12	±23	±40	±46
PT-111										±10	±12	±23	±40	±46
Germanium														
GR-50-AA/CD	±5	±5	±5	±5	±6	±6	±6	_	_	_	_	_	_	_
GR-300-AA/CD	_	_	±4	±4	±4	±4	±4	±4	±8	±25	_	_	_	_
GR-1400-AA/CD	_	_	_	_	_	±4	±4	±4	±7	±15	_	_	_	

 $<sup>^{12}</sup>$ All accuracies are: 2  $\sigma$  figures; [(calibration uncertainty) $^2$  + (reproducibility) $^2$ ] $^{0.5}$ ; for additional information, please see Appendix D.




# Sensor Packages and Mounting Adapters


Temperature sensors are available in a variety of packages to facilitate mounting. Included are adapters that allow the sensor to be soldered in place, screwed on, bolted down, inserted into a hole, or inserted through a pressure seal in the form of a thermowell. Gold-plated copper bobbins are available for both diodes and resistors in order to heat sink leads. The chart below summarizes the standard Lake Shore sensor and packaging configurations. Appendix C: Sensor Packaging and Installation discusses techniques for the correct installation of temperature sensors. More specific installation notes are included for the bare chip sensors, the SD package, and the CU, DI, CY, and CD adapters. Special packaging is also available—consult Lake Shore for custom orders.

						Lak	e Sho	re sens	sors							
		S	ilicor	n dioc	le					P	latinu	ım				
Packa (see in	i <b>ging</b> dividual sensor pages for additional details)	DT-414/DT-670E-BR	DT-421	DT-670	DT-471	GaAIAs diode	Cernox®	Germanium	Rox™	PT-102	PT-103	PT-111	Rhodium-iron	Installation instructions		
Commo	n															
Bare ch	ip sensors Bare chip with 2 copper leads (42 AWG)						٠							Appendix C		
BG	Bare chip with 2 or 4 gold leads						-							Appendix C		
BR	Bare chip, no leads						-		-					Appendix C		
Hermeti SD	cally sealed package			•	٠	٠	•							Appendix C		
Mountin CO	g adapters for SD Clamp			•	•	•	-							Appendix C		
ET	Screw-in						-							Order from Lake Shore		
MT	Screw-in (metric)					•	-							Order from Lake Shore		
CU	Copper bobbin (small, 4-lead)			-		•	-							Appendix C		
DI	Copper bobbin (small, 2-lead)					•								Appendix C		
CY	Copper bobbin (large, 2-lead)													Appendix C		
LR	Half-rounded cylinder					•	-							Order from Lake Shore		
В0	Beryllium oxide heat sink block						-							Order from Lake Shore		
Platinun AL	n mounting adapters									•				Order from Lake Shore		
AM											-			Order from Lake Shore		
Copper	canister package						-	•	•					Appendix C		
CD	Copper bobbin						-	•	-					Appendix C		
Unique	packages															
See indi	vidual sensor specifications					•				•	-			Order from Lake Shore		



#### **Packages**

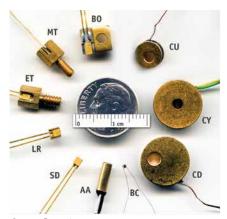

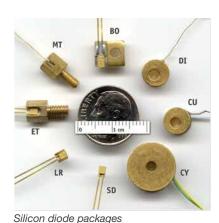



Germanium and Rox<sup>™</sup> packages (Rox -AA and -CD only)

Cernox® packages








**Unique packages** 

see individual sensor pages

TG-120-P DT-414
TG-120-PL DT-421-HR
PT-102 DT-670E-BR
PT-103 CX-10XX-BC
PT-111 CX-10XX-BG
CX-10XX-BR
RX-102B-CB
RF-800

The Lake Shore Hermetically Sealed SD Package

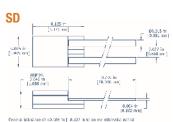

- Small package designed primarily for bonding or clamping to a flat surface
- Indium, silver epoxy, 2850 Stycast® epoxy, or a CO clamp may be used for mounting

Package material: Sapphire base with alumina body and lid.

Molybdenum/manganese metallization on base and lid top with nickel and gold plating. Gold tin solder as

hermetic lid seal.

Leads:

Lead material: Silicon diode: brazed Kovar

Cernox®: gold-plated copper soldered with

60/40 SnPb


GaAlAs: welded platinum

Mass: 0.03 g

Limitation: The useful upper temperature limit of this


configuration is 500 K

#### SD (Cernox®)

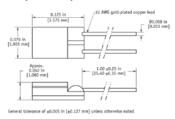

### The Lake Shore SD package—the most rugged, versatile package in the industry

The SD package, with its sapphire base, direct sensor-to-sapphire mounting, hermetic sealing, and brazed Kovar leads provides the industry's most rugged, versatile sensors with the best thermal connection between the sample and sensor chip. In addition, this

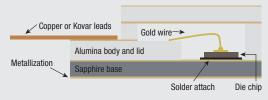


package is designed so heat coming down the leads bypasses the sensor chip. It can survive several thousand hours at 500 K and is compatible with most ultra high vacuum applications, and can be indium soldered to samples. The Lake Shore SD package is now available with Cernox® resistors and GaAlAs diodes, as well as silicon diodes. For the Cernox® resistors and GaAlAs diodes, the Kovar leads are replaced with nonmagnetic leads.



### Mounting adapters for SD package—CO, CU, DI, CY, LR, BO, ET, MT

#### CO

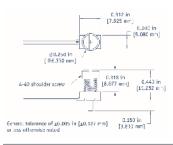

Spring-loaded clamp holds standard SD sensor in contact with the surface of the sample and allows the sensor to be easily changed or replaced

Extra clamps are available for frequent relocation of the sensor

4-40 stainless steel screw has a formed shoulder, thus applying correct pressure to the clamp

Package material: Adapter material:

See SD package Gold-plated copper (nickel strike);


spring is ASTM A313 302 Austenitic steel

See SD package Leads: Lead material: See SD package

1.8 g (including SD package and clamp) Mass: Limitation: The useful upper temperature limit of this

configuration is 500 K

#### CU/CU-HT & DI

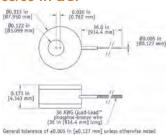

CU

SD packaged sensor indium-soldered into a flat copper bobbin with the leads thermally anchored to that same bobbin

HT (high temperature) version is soldered using high temperature (90% Pb, 10% Sn) solder

Can be mounted to any flat surface with a 4-40 screw

DI

2-lead version of the CU

Package material: See SD package

Adapter material: Gold-plated copper bobbin (SD indium-soldered

to adapter and wrapped in Stycast® epoxy); high temperature CU uses high temperature (90% Pb,

10% Sn) solder

CU leads: Four 0.91 m (36 in), 36 AWG, color-coded

Quad-Lead™

DI leads: 0.91 m (36 in), 36 AWG, color-coded, 2-lead

ribbon cable

Lead material: Phosphor bronze alloy

1.1 g (including SD package and bobbin, excluding Mass:


leads)

Limitation: The epoxy limits the upper useful temperature of

this configuration to 378 K (high temperature CU-HT upper temperature limit is 420 K with Cernox® and

500 K with silicon and GaAlAs diodes)

CY

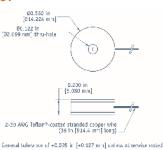

Similar to the DI package, except the bobbin is larger in diameter with a centered mounting hole

Relatively large-sized, robust

Package material: See SD package

Adapter material: Gold-plated copper bobbin (SD indium-soldered to

adapter and wrapped in Stycast® epoxy)

Leads: Two 0.91 m (36 in), 30 AWG Teflon®-coated leads

Lead material: Stranded copper

4.3 g (Including SD package and bobbin, excluding Mass:


leads)

The epoxy limits the upper useful temperature of Limitation:

this configuration to 400 K



LR

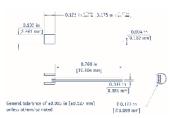

 With an SD packaged sensor mounted on a slightly-more-than half-rounded cylinder, this package is designed to be inserted into a 3.2 mm (1/8 in) diameter hole Package material: Adapter material: See SD package

Gold-plated flat cylindrical copper disk

(SD indium-soldered to adapter)


Leads: See SD package Lead material: See SD package

Mass: 0.2 g (Including SD package and disk)

Limitation: Indium solder limits the upper useful temperature of

this configuration to 420 K

**BO** 

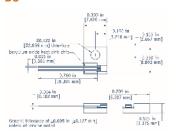

 SD package is soldered to a mounting block and the leads are thermally anchored (without epoxy) to the block via a beryllium oxide insert

Since leads can be a significant heat path to the sensing element and can lead to measurement errors when incorrectly anchored, this configuration helps maintain the leads at the same temperature as the sensor Package material: See SD package

Adapter material: Gold-plated bolt-on copper block with leads


thermally anchored to block (SD indium-soldered to adapter)

Leads: See SD package Lead material: See SD package

Mass: 1.5 g (including SD package and mounting block)
Limitation: Indium solder limits the upper useful temperature of

this configuration to 420 K

ET

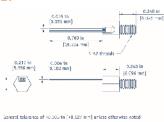

 Convenient screw-in package formed by indiumsoldering a basic SD configuration into a recess in one flat of a hexagonal screw head

The head terminates in a standard SAE 6-32 threaded stud allowing the sensor to be threaded into a mounting hole in the sample Package material: See SD package

Adapter material: ET: gold-plated copper SAE-threaded screw head

#6-32

MT: gold-plated copper metric threaded screw head

3 mm × 0.5 metric


Leads: See SD package Lead material: See SD package

Mass: 1.5 g (including SD package and screw head)

Limitation: Indium solder limits the upper useful temperature of

this configuration to 420 K

MT

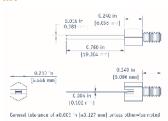


The MT package is similar to the ET version except the SD package is mounted in a slot in the center of the hexagonal head and the stud is a 3 mm × 0.5 metric thread

**Note:** A light coating of vacuum grease on the threads further enhances the thermal contact between the sensor package and the sample.



# Copper canister packages



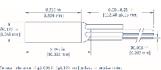

Used with Cernox®, germanium, and Rox™

sensors

Used only with germanium sensors

Adapter material: Gold-plated cylindrical copper canister, BeO header,

Stycast® epoxy

Leads: Four 32 AWG  $\times$  152 mm (6 in) long (Rox™: Two 32 AWG × 152 mm [6 in] long) Lead material:

Phosphor bronze insulated with polyimide (Rox™: copper insulated with Formvar®)

Mass: AA canister (empty): 0.091 g B canister (empty): 0.080 g

Once sensors are installed, total mass increases to 0.197 g to 0.416 g. Refer to individual sensor


specifications.

The epoxy limits the upper useful temperature of Limitation:

this configuration to 400 K

### Mounting adapter for AA canister package

#### CD

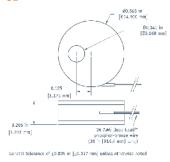

- AA canister sensor soldered into a flat, copper bobbin with the sensor leads thermally anchored
- Can be mounted to any flat surface with a 6-40 screw (not supplied)
- Used with Cernox®, Germanium, and Rox™ sensors

Adapter material:

Copper bobbin, gold-plated (AA canister epoxied to

bobbin with Stycast® epoxy)

Leads: 0.91 m (36 in), 36 AWG, color-coded, Quad-Lead™


Lead material: Phosphor bronze Grade A alloy

Limitation: The epoxy limits the upper useful temperature of

this configuration to 378 K

### Mounting adapters for platinum RTDs

#### PT-102-AL

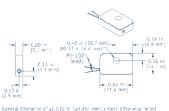

- PT-102 (AL) or PT-103 (AM) mounted into a flat aluminum block
- Can be mounted to any flat surface with a 6-32 or M3 screw (not included) and Inconel® Belleville washer (included)

Adapter material: 6061 Al block (PT mounted to adapter using Cotronics Durabond® 950 Al-based adhesive) AL leads:

Two 0.010-inch diameter;  $10.160 \pm 1.270 \text{ mm}$  $(0.400 \pm 0.050 \text{ in}) \text{ long}$ 

AM leads: Two 0.010-inch diameter; 15.240 ±1.270 mm  $(0.600 \pm 0.050 \text{ in}) \log$ 


Lead material: Platinum Mass: PT-102-AL: 3.8 g

PT-103-AM: 2.1 q

Limitation: The aluminum alloy limits the upper useful

temperature of these configurations to 800 K

PT-103-AM

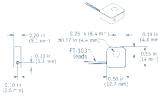

General colerance of +0.010 in (+0.254 mm) unless otherwise noted



### **Lead Extensions**

(formerly SMODs)

Adding extra wire to your sensor leads can be cumbersome and time consuming. Lake Shore offers this service for you at the time of order, allowing numerous options to best suit your application.

There are various options available when selecting a lead extension:

#### Number of Wires

4-wire: For accurate sensor measurements, 4-lead connections are by far the superior option when adding a lead extension to both diodes and resistive temperature sensors. See Appendix C and Appendix E for additional information.

2-wire: This option is useful if the number of electrical connections inside a system must be kept to a minimum. However, 2-lead connections add measureable resistance to sensor measurements as described in Appendix E. This additional resistance will cause a significant (but repeatable) shift on all sensors except diodes.

#### Wire Type

Phosphor bronze: This all-purpose cryogenic wire has a great balance of features.

- Low thermal conductivity minimizes heat leak (lower is generally better)
- Moderate electrical resistance (lower is generally better)
- Non-ferromagnetic and very low magnetoresistance, making this wire the best choice for applications where magnetic fields are present
- Available in several convenient configurations in addition to single strand, such as Quad-Lead™ and Quad-Twist™

Manganin: This wire has several interesting characteristics that make it useful in certain situations.

- Coefficient of thermal expansion very close to that of pure copper
- Very low thermal conductivity minimizes heat leak (lower is generally better)
- Somewhat high electrical resistance (lower is generally better)
- Heavy Formvar® insulation limits upper temperature of wire to 378 K
- Non-ferromagnetic
- Available as single strand wire only

#### Wire Gauge

Wire gauge (AWG)	Wire diameter (in)	Wire diameter (mm)
30	0.01	0.255
32	0.00795	0.202
36	0.0055	0.127
42	0.0025	0.0635

Various wire thicknesses are available, depending on the wire type selected. The wire gauge selection process usually involves a compromise between thermal conductivity and ease-of-use, with thinner wire being preferred to reduce thermal conductivity and thicker wire being easier to handle and work with. Lake Shore uses American wire gauge (AWG) for its wire. This conversion table is provided for your convenience.

32 AWG and 36 AWG are our preferred wire gauges to use with cryogenic sensors. By far they provide the best balance between reduced thermal conductivity and ease-ofuse.

Manganin is the only wire type available in 30 AWG as the extremely low thermal conductivity of the wire helps compensate for the "large" cross-sectional area associated with 30 AWG.

Phosphor bronze is the only wire type available in 42 AWG. This wire thickness reduces thermal conductivity substantially to the levels possible with manganin, with the same low magentoresistance of phosphor bronze. Unfortunately, this wire is extremely delicate and can break easily. Lake Shore suggests this wire be ordered only by users with extensive experience with system wiring.

#### Wire Length

Lead Extensions

Standard lengths of 2 m and 5 m are offered with all wire types and gauges. These lengths have been selected to suit a wide range of applications, most commonly wiring from a temperature sensor through the various stages of a cryostat, up to and terminating at an electrical feedthrough. Additional wire may be trimmed from both of these wire lengths if necessary. However, if a custom length is required, please contact Lake Shore to discuss custom wire lengths.

#### **Component temperature limits**

The lead extention components have different maximum temperatures. Use this chart to ensure the lead extensions you order are appropriate for your given application.

Lead extension component	Maximum temperature
Formvar	378 K (105 °C)
Bond Coat 999	433 K (160 °C)
Polyimide	500 K (227 °C)
63/37 Solder	450 K (177 °C)
90/10 Solder	548 K (275 °C)



#### **Recommended Standard Lead Extensions**

Lake Shore recommends selecting from one of these two configurations — our most popular configurations due to the wide range of applications they cover.

#### -QL

Quad-Lead<sup>™</sup> phosphor bronze, 32 AWG, 2 m For situations where ease-of-use and ruggedness is important.

- 32 AWG wire is easier to prepare and solder to that thinner gauges
- Quad-lead™ wire is easy to heat-sink around copper bobbins due to its ribbon structure
- Polyimide insulation is strong and is resistant to solvents, and also has a high temperature rating that protects it from heating that might be applied to help soften the bonding agent used to join the wires to one another

#### -QT

Quad-Twist<sup>™</sup> phosphor bronze, 36 AWG, 2 m For noisy environments where signal integrity must be protected.

- Quad-twist™ wire helps reject electromagnetic interference that may be present inside the measurement space
- 32 AWG wire is easier to prepare and solder to that thinner gauges
- Quad-twist<sup>™</sup> can be slightly more difficult to heat-sink, but the 36 AWG wire reduces thermal conductivity and therefore reduces heat-leak naturally
- Formvar® insulation has excellent mechanical properties such as abrasion resistance and flexibility, which is important when using 36 AWG wire. However, care should be taken as Formvar® can craze when exposed to solvents.

There are certain scenarios where these standard offerings are not adequate and alternative solutions should be selected. One such example is higher-temperature applications above 450 K where both Quad-Lead™ wire and Formvar® insulation become inappropriate. This application would require Quad-Twist, 32 AWG. In this scenario, please use the full part configurations to define the lead extension.

#### -XXYY-Z

**XX** = Wire type

YY = Wire gauge (AWG)

Z = Length in meters

#### **Method of Ordering**

When ordering a lead extension on the website, add the sensor to the shopping cart first, and then come to this page to add a lead extension.

If placing a purchase order, please append the lead extension part number to the sensor that requires the extension. Examples:

CX-1050-SD-HD-4L <b>-QL</b>	Quad-Lead™, 32 AWG, 2 m
DT-670-CU-HT-1.4L <b>-QT</b>	Quad-Twist <sup>™</sup> , 36 AWG, 2 m attached to 0.91 m of Quad-Twist <sup>™</sup> , 36 AWG wire that comes standard with the diode CU-HT package.
PT-102-14L-QT32 <b>-5</b>	Quad-Twist™, 32 AWG, 5 m
DT-670C-SD-DT32 <b>-2</b>	Duo-Twist™, 32 AWG, 2 m

Lead extensions are not available on devices with gold or no leads

For more information please visit www.lakeshore.com.



## **Temperature Probes**

#### **Temperature probe features**

- Stainless steel-encased probes that provide highly reliable sensor performance in a thermowell or direct cryogen contact
- Highly customizable to suit your particular application
- May be configured with many sensor types, including Cernox® for superior temperature performance from room temperature down to 4 K (-269.15 °C) and below
- Thin-walled probe tubing reduces thermal lag and heat leak from outside the measurement space
- Ideal for temperature measurements in fluid containers and tanks
- Full 3 year standard warranty



Lake Shore offers a variety of temperature sensors in packages that enable mounting in very tight areas. But for some applications (especially if the sensors have to be immersed in liquid) you need to do more to protect the sensor circuitry. For these applications, a cryogenic temperature probe is the optimum choice. Encased in one of these stainless steel thermowell fixtures, the sensor can perform as designed, unaffected by high pressure and sealed to keep electrical components and wiring protected from fluids and other elements.

#### **Typical applications**

Lake Shore temperature probes are ideal for thermometry applications where you need to measure inside:

- fluid containers, tanks, and pipes
- cryostats and cryogenic liquid flow meters
- other liquid storage systems.

#### **Highly customizable**

Lake Shore temperature probes are made-to-order with a wide range of configuration options available. These include:

- Multiple sensor types including our extremely popular Cernox® RTDs and DT-670 diodes
- Either 1/8 in or 1/4 in stem diameter in lengths up to 0.71 m (28 in) are standard
- Various mounting adapters suited for either positive or negative pressures, if required
- Numerous connectivity options including wire types and lengths as well as various terminating connectors for direct connection to Lake Shore temperature instruments or third party equipment

If you do not see an option available as part of our standard offerings, please contact Lake Shore to discuss further customization options.



#### **Specifications**

**Note:** These probes are not designed to be intrinsically safe. It is the responsibility of the user to operate these probes safely in explosive environments.

#### **Probe construction**

#### **Stem**

Material: 316 stainless steel (non-magnetic)<sup>1</sup>

	Wall thickness	Maximum length
1/4 in stem	0.028 in ±0.003 in	28 in*
1/8 in stem	0.010 in ±0.001 in	20 in

<sup>&</sup>lt;sup>2</sup>Not suitable for direct immersion in liquid oxygen or hydrogen environments.

#### **Internal components**


Internal atmosphere: Air

Internal atmosphere pressure: 98 kPa (14.2 psia)

**Internal sensor wire:** Quad-Twist™ 4-lead 36 AWG phosphor bronze wire with polyimide insulation


#### **Probe mount**

#### Swagelok® fittings

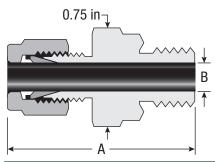

	1/4 in probe	1/8 in probe		
Swagelok® part number:	SS-400-1-4BT	SS-200-1-2BT		
Material	316 stainless steel			
Thread	0.25 in NPT male	0.125 in NPT male		
Α	1.59 in	1.5 in		
В	0.25 in	0.125 in		

#### **CF flange**

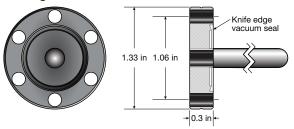

Material: 304L stainless steel Flange size: 11/3 in (DN16)

**Vacuum rating:**  $1 \times 10^{-13}$  torr (<1.3 × 10<sup>-13</sup> mbar)\*

#### **Connectors**

#### **BNC** connector

Standard male BNC connector. When ordering with 4-lead wire, two separate BNC connectors will be provided to maintain the 4-lead measurement.


#### **Configuration:**

	BN	C 1	BN	C 2
	Center pin	Shield	Center pin	Shield
2-lead cable	I/V+ (anode)	I/V- (cathode)	_	_
4-lead cable	l+	I-	V+	V-

#### 10-pin Detoronics® connector

The Detoronics connector is o-ring sealed to the temperature probe.

**Note:** This connector is mounted directly to the probe, meaning that no external cable can be selected with this option. It also eliminates the CF flange probe mount option.

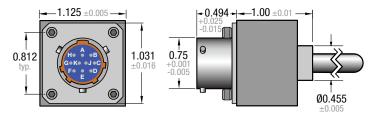

#### General specifications

Air leakage: 1  $\times$  10  $^{\text{-6}}$  cm  $^{3}$ /s at 15 psi Insulation resistance: 5,000 M $\Omega$  at 500 VDC

Operating temperature: -55 °C to +125 °C (-67 °F to +257 °F)

Finish is tin-plated shell and pins.

Materials

Shell, bayonet and flange: Carbon steel

Pins: 52 nickel alloy Insulator: Glass

#### 25-pin D-sub connector



The 25-pin D-sub is required to connect directly to particular Lake Shore temperature monitors.

#### Supported instruments:

- Model 211
- Model 218

#### 6-pin DIN connector



The 6-pin DIN is required to connect directly to particular Lake Shore temperature controllers and monitors.

#### Supported current instruments:

- Model 350
- Model 336
- Model 335
- Model 224

### Supported discontinued instruments:

- Model 340
- Model 331/332
- Model 330 (diodes only)
- Model 321 (silicon diodes only)

<sup>&</sup>lt;sup>2</sup>Longer lengths may be possible depending on the overall configuration. Please contact Lake Shore to discuss.

<sup>\*</sup>Requires the use of appropriate bolts, gasket and mating surface.



#### **Connector configurations**


Connector type	l+	V+	I-	V-	Shield*
2-lead BNC (1 connector)	Center pin		Outer cu	Not connected	
4-lead BNC (2 connectors)	Center pin of 'I' BNC	Center pin of 'V' BNC	Outer cup of 'I' BNC	Outer cup of 'V' BNC	Not connected
10-pin probe- mounted Detoronics connector®	Pin A	Pin C	Pin B	Pin D	NA
6-pin DIN	Pin 5	Pin 4	Pin 1	Pin 2	Pin 6
25-pin D-sub	Pin 3	Pin 4	Pin 15	Pin 16	Pin 2

<sup>\*</sup>Shield connection is only used in conjunction with external cable choices that include a braided shield (Cryocable™ and instrument cable)

#### Wire

#### Instrument cable

Robust 4-lead cable best for wiring to instrument where both the wire and instrument are at room temperature. The 30 AWG signal wires make these wires easier to work with than traditional cryogenic wire.

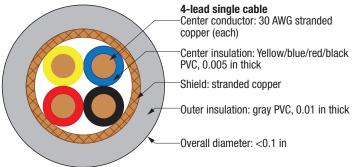

Rated temperature: -20 °C to 80 °C Thermal conductivity (300 K): 400 W/(m·K)

Resistance (300 K):  $0.32 \Omega/m$ 

Supported sensor types: Cernox® RTD, silicon diode, GaAlAs diode, platinum RTD

Maximum rated temperature: 378 K

#### Cryogenic wire


Phosphor-bronze wire combinations that limit heat transfer into the temperature probe and are themselves rated for use in cryogenic environments.

	Quad-Twist™ 36 AWG*	Quad-Twist™ 32 AWG	Quad-Lead™ 32 AWG	Duo-Twist™ 32 AWG
Configuration		4-lead	2-lead	
Wire		Phospho	or bronze	
Gauge	36 AWG		32 AWG	
Insulation	Formvar		Polyimide	
Structure	Two twis	ted pairs	One twisted pair	
Thermal conductivity (300 K)		48 W		
Resistance (300 K)	10.3 Ω/m			
Supported sensors	Cernox® RTD, sil	icon diode, GaAlAs RTD	Diodes only	

\*Also used for internal probe wiring. Ordering this cable will result in a continuous length of wire from the sensor through to the outside environment.

#### SS (stainless steel) coaxial cable

2-lead cabling solution that is extremely robust and limits heat transfer into the probe. Due to the 2-lead configuration, this cable is only compatible with diode sensors and will cause a predictable (potentially insignificant) offset in any temperature readings.

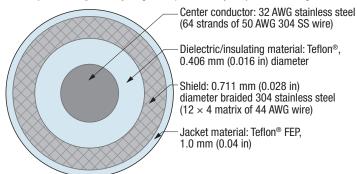

#### Electrical properties

Resistance—center conductor at 295 K (22 °C):  $23.62 \Omega/m$  (7.2  $\Omega/ft$ )

Resistance—shield at 295 K (22 °C):  $3.61~\Omega/m$  (1.1  $\Omega/ft$ )

Insulation temperature range: 10 mK to 473 K  $\,$ 

Supported sensor types: Silicon diode, GaAlAs diode, platinum RTD


#### Cryocable™

A robust, 4-wire cable for use in cryogenic environments to room temperature for the



#### **Cryocable**<sup>™</sup>

A robust, 4-wire cable for use in cryogenic environments to room temperature for the ultimate in thermal isolation from external heat sources. This cable is designed around 32 AWG (203  $\mu m)$  diameter superconductive wires consisting of a NbTi core (128  $\mu m$  diameter) and a Cu-10% Ni jacket. The wire is LTS, requiring very low temperatures for it to become superconducting.

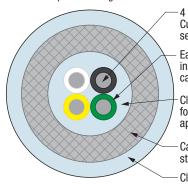

4 32 AWG wires: Nb-48wt%Ti core with Cu-10wt%Ni jacket, CuNi to NbTi cross sectional area ratio = 1.5:1

Each wire overcoated with Teflon® (PFA) insulation 0.003 in (75  $\mu$ m) thick; wires cabled with approx. 25 mm twist pitch

Clear Teflon® (PFA) extruded over the four-wire cable to an overall diameter of approx. 1.2 mm (0.048 in)

Cable overbraided with 304 stainless steel wire

Clear Teflon® (PFA) extruded over the entire cable; finished cable has an overall diameter 2.4 mm  $\pm 0.2$  mm (0.094 in  $\pm 0.008$  in)

 $\label{eq:minimum bend radius: 15 mm (0.6 in)} {\mbox{Superconducting critical temperature: } 9.8 \mbox{ K} \\ {\mbox{Superconducting critical magnetic field: } 10 \mbox{ T} \\ \mbox{}$ 

Supported sensor types: Cernox® RTD, silicon diode, GaAlAs diode, platinum RTD


Magnetic field	Critical current (per wire)
3 T	35 A
5 T	25 A
7 T	15 A
9 T	6 A

	Temperature (K)					
	295	77	4.2			
Wire resistance (Ω/m)	9.2	8.4	0*			
Overbraid resistance (Ω/m)	0.90	0.64	0.62			
Thermal conductivity—entire cable assembly (W/(m·K))	7.6	2.8	0.17			

<sup>\*</sup>Superconducting

#### Wire configurations

Wire type	l+	V+	I-	V-	Shield
Instrument cable	Black	Yellow	Red	Blue	Copper braid
Quad-Twist™ 36 AWG	Green (from red/green pair)	Green (from clear/green pair)	Red	Clear	None
Quad-Twist™ 32 AWG	Red	Black	Green	Clear	None
Quad-Lead™ 32 AWG	Clear	Black	Red	Green	None
Duo-Twist™ 32 AWG	CI	ear	Gree	en	None
Stainless steel coaxial	Center o	Center conductor		eld	None
Cryocable™	Black	Yellow	White	Green	Stainless steel braid

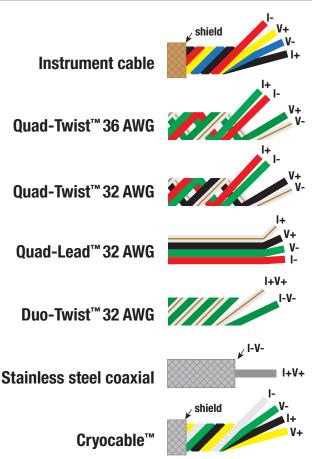

#### **Temperature sensors**

See the individual Cernox, DT-670, and platinum sensor pages for specifications:

Sensor type	Installed sensor package		
Cernox®	SD		
DT-670	SD		
Platinum	Standard PT-100 Series packages		

All temperature sensor calibrations are performed before the device is installed into the probe. At this time, Lake Shore does not perform recalibrations on finished probes.



#### Temperature probe ordering information

The easiest way to request a quote for a temperature probe is to use the online configurator at www.lakeshore.com. Otherwise contact our Sales department at sales@lakeshore.com and we can assist you.

# Specify TP-a-bcd-e-f-g, where:

a = probe length in inches—offered in whole inch increments from 1 to 28 inches

#### b = tube diameter<sup>1</sup>

2	1/8 in
4	1/4 in

<sup>&</sup>lt;sup>1</sup> Probes over 20 inches long are only available in 1/4-inch diameter

#### c = probe mount

N	no probe mount adapter
S	Swagelok® fitting²
F	CF™ flange mount <sup>3</sup>

<sup>&</sup>lt;sup>2</sup> For 1/8 in diameter probe, Swagelok® fitting uses a 1/8 in NPT male thread; for 1/4 in diameter probe, Swagelok® fitting uses a 1/4 in NPT male thread

#### d = external cable/wire type<sup>4</sup>

N	no external cable (usually used with Detoronics connector)			
S	S1 coaxial cable (2-lead)			
- 1	30 AWG instrument cable (4-lead)			
T	DT-32 (twisted pair of 32 AWG phosphor bronze wire)			
F	QT-32 (two twisted pairs of 32 AWG phosphor bronze wire)			
Q	QT-36 (two twisted pairs of 36 AWG phosphor bronze wire)			
L	QL-32 (four 32 AWG wires in a ribbon configuration)			
CryoCable™ (4-lead cryogenic coaxial cable)				
4 Lake Shore strongly recommends that all RTD temperature sensors use a 4-lead cable/wire type				

Lake Shore strongly recommends that all RTD temperature sensors use a 4-lead cable/wire type

#### e = terminator

N	no connector (leads stripped and tinned)
В	BNC connector
D	10-pin Detoronics connector <sup>5</sup>
Y	25-pin D-shell connector for temperature monitors
R	connector wired for temperature instruments (6-pin round)

<sup>&</sup>lt;sup>5</sup> Selecting a Detoronics connector limits the following selections: d = N and f = 0; the Detoronics connector is o-ring sealed to the probe

# **f = external cable length** — offered in whole meter increments from 1 to 10 m (enter '0' for no external cable)

# **g** = **temperature sensor type**<sup>6</sup>—specify sensor model number with calibration range, if applicable

#### **Ordering example**

**TP**- 06 - 2FS - B - 03 - S27

(6 in probe, 1/8 in diameter, flange, S1 coaxial cable, BNC connector, 3 m cable length, DT-670-SD calibrated 1.4 K to 325 K)

#### Calibration range suffix codes

Numeric figure is the low end of the calibration Letters represent the high end: B = 40 K, D = 100 K, L = 325 K, H = 500 K

#### Cernox® RTDs

Uncalibrated	C01	CX-1010-SD		
	C02	CX-1030-SD		
	C03	CX-1050-SD		
	C04	CX-1070-SD		
	C05	CX-1080-SD		
Calibrated	C07	CX-1010-SD-0.1L		
	C16	CX-1030-SD-0.3L		
	C25	CX-1050-SD-1.4L		
	C31	CX-1070-SD-4L		
	C32	CX-1080-SD-20L		
	C13	CX-1010-SD-1.4L		

#### **Platinum RTDs**

Uncalibrated	P01	PT-102		
	P02	PT-103		
	P03	PT-111		
Calibrated	P04	PT-102-2S		
	P05	PT-102-3S		
	P07	PT-102-14L		
	P08	PT-102-14H		
	P11	PT-103-2S		
	P12	PT-103-3S		
	P14	PT-103-14L		
	P15	PT-103-14H		
	P18	PT-111-2S		
	P19	PT-111-3S		
	P21	PT-111-14L		
	P22	PT-111-14H		

#### Silicon diodes

Uncalibrated	S07	DT-670A-SD		
	S08	DT-670B-SD		
	S09	DT-670C-SD		
	S10	DT-670D-SD		
	SOA	DT-670A1-SD		
	SOB	DT-670B1-SD		
Calibrated	S27	DT-670-SD-1.4L		
	S28	DT-670-SD-1.4H		
	S32	DT-670-SD-70L		
	S33	DT-670-SD-70H		

#### GaAlAs diodes

Uncalibrated	G01	TG-120-SD	
Calibrated	G04	TG-120-SD-1.4L	
	G05	TG-120-SD-1.4H	
	G10	TG-120-SD-70L	
	G11	TG-120-SD-70H	

 $<sup>^{\</sup>scriptscriptstyle 3}$  The  $\text{CF}^{\scriptscriptstyle \text{TM}}$  flange is welded to the probe

<sup>&</sup>lt;sup>6</sup> Due to indium solder use, all SD sensors have an upper temperature usage limit of 400 K



### Cernox® RTDs

#### Cernox® features

- Low magnetic field-induced errors
- Temperature range of 100 mK to 420 K (model dependent)
- High sensitivity at low temperatures and good sensitivity over a broad range
- Excellent resistance to ionizing radiation
- Bare die sensor with fast characteristic thermal response times: 1.5 ms at 4.2 K, 50 ms at 77 K
- Broad selection of models to meet your thermometry needs
- Excellent stability
- Variety of packaging options



CAUTION: These sensors are sensitive to electrostatic discharge (ESD). Use ESD precautionary procedures when handling, or making mechanical or electrical connections to these devices in order to avoid performance degradation or loss of functionality.

Cernox® thin film resistance temperature sensors offer significant advantages over comparable bulk or thick film resistance sensors. The smaller package size of these thin film sensors makes them useful in a broader range of experimental mounting schemes, and they are also available in a chip form. They are easily mounted in packages designed for excellent heat transfer, yielding a characteristic thermal response time much faster than possible with bulk devices requiring strain-free mounting. Additionally, they have been proven very stable over repeated thermal cycling and under extended exposure to ionizing radiation.

#### **Packaging options**

AA, BC, BG, BO, BR, CD, CO, CU, ET, LR, MT, SD



The CX-1010 is the first Cernox® designed to operate down to 100 mK, making it an ideal replacement for Germanium RTDs. Unlike Germanium, all Cernox models have the added advantage of being able to be used to room temperature. In addition, Cernox is offered in the incredibly robust Lake Shore SD package, giving researchers more flexibility in sensor mounting.








# The Lake Shore SD package — the most rugged, versatile package in the industry


CX-BR

The SD package, with direct sensor-to-sapphire base mounting, hermetic seal, and brazed Kovar leads, provides the industry's most rugged, versatile sensors with the best sample to chip connection. Designed so heat coming down the leads bypasses the chip, it can survive several thousand hours at 500 K (depending on model) and is compatible with most ultra high vacuum applications. It can be indium soldered to samples without shift in sensor calibration. If desired, the SD package is also available without Kovar leads.


#### Typical Cernox® resistance

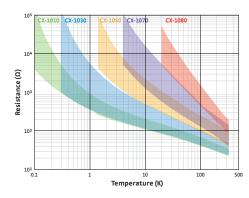

#### **Typical Cernox® sensitivity**

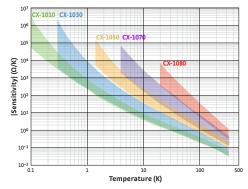

#### Typical Cernox® dimensionless sensitivity

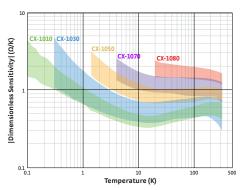



#### **Specifications**

#### Standard curve Not applicable

Recommended excitation<sup>1</sup> 20 μV (0.1 K to 0.5 K); 63 μV (0.5 K to 1 K); 10 mV or less for T > 1.2 K

Dissipation at recommended excitation Typical 10<sup>-5</sup> W at 300 K,  $10^{-7}$  W at 4.2 K,  $10^{-13}$  W at 0.3 K (model and temperature dependent)

Thermal response time BC, BR, BG: 1.5 ms at 4.2 K, 50 ms at 77 K, 135 ms at 273 K; SD: 15 ms at 4.2 K, 0.25 s at 77 K, 0.8 s at 273 K; AA: 0.4 s at 4.2 K, 2 s at 77 K, 1.0 s at 273 K

Use in radiation Recommended for use in radiation environments-see Appendix B

Use in magnetic field Recommended for use in magnetic fields at low temperatures. The magnetoresistance is typically negligibly small above 30 K and not significantly affected by orientation relative to the magnetic field—see Appendix B

Reproducibility<sup>2</sup> ±3 mK at 4.2 K

#### Soldering standard J-STD-001 Class 2

- <sup>1</sup> Recommended excitation for T < 1 K based on Lake Shore calibration procedures using an AC resistance bridge—for more information refer to Appendix D and Appendix E
- <sup>2</sup> Short-term reproducibility data is obtained by subjecting sensor to repeated thermal shocks from 305 K to 4.2 K

#### Range of use

	Minimum limit	
Cernox®	0.10 K <sup>3</sup>	420 K

<sup>3</sup> Model dependent

#### Calibrated accuracy4

	Typical sensor accuracy <sup>5</sup>	Long-term stability <sup>6</sup>
1.4 K	±5 mK	±3 mK
4.2 K	±5 mK	±3 mK
10 K	±6 mK	±6 mK
20 K	±9 mK	±12 mK
30 K	±10 mK	±18 mK
50 K	±13 mK	±30 mK
77 K	±16 mK	±46 mK
300 K	±40 mK	±180 mK
400 K	±65 mK	_

- <sup>4</sup> Bare chip sensors can only be calibrated after attaching gold wire leads—the user must remove the ball bonded leads if they are not desired (the bond pads are large enough for additional bonds)
- [(Calibration uncertainty)<sup>2</sup> + (reproducibility)<sup>2</sup>]<sup>0.5</sup> for more information see Appendices B, D, and E
- <sup>6</sup> Long-term stability data is obtained by subjecting sensor to 200 thermal shocks from 305 K to 77 K

Typical magnetic field-dependent temperature errors7 AT/T (%) at B (magnetic induction)

Cernox® RTDs

Cernox® 1050					
	2.5 T	8 T	14 T	19 T	
2 K	1.3	3.1	3.9	5	
4.2 K	0.1	-0.15	-0.85	-0.8	
10 K	0.04	-0.4	-1.1	-1.5	
20 K	0.04	0.02	-0.16	-0.2	
30 K	0.01	0.04	0.06	0.11	
77 K	0.002	0.022	0.062	0.11	
300 K	0.003	0.004	0.004	0.006	

Excellent for use in magnetic fields, depending on temperature range (>2 K)

#### Temperature response data table (typical)


	CX-1010			CX-1030			CX-1050			CX-1070			CX-1080		
	R <sup>8</sup> (Ω)	dR/dT (Ω/K)	(T/R)· (dR/dT)												
4.2	277.32	-32.209	-0.49	574.20	-97.344	-0.71	3507.2	-1120.8	-1.34	5979.4	-2225.3	-1.56	_	_	_
10	187.11	-8.063	-0.43	331.67	-19.042	-0.57	1313.5	-128.58	-0.98	1927.2	-214.11	-1.11	_	_	_
20	138.79	-3.057	-0.44	225.19	-6.258	-0.56	692.81	-30.871	-0.89	938.93	-46.553	-0.99	6157.5	-480.08	-1.56
30	115.38	-1.819	-0.47	179.12	-3.453	-0.58	482.88	-14.373	-0.89	629.90	-20.613	-0.98	3319.7	-165.61	-1.50
77.35	70.837	-0.510	-0.56	101.16	-0.820	-0.63	205.67	-2.412	-0.91	248.66	-3.150	-0.98	836.52	-15.398	-1.42
300	30.392	-0.065	-0.65	41.420	-0.088	-0.64	59.467	-0.173	-0.87	66.441	-0.201	-0.91	129.39	-0.545	-1.26
400 (HT)	_	_	_	34.779	-0.050	-0.57	46.782	-0.093	-0.79	51.815	-0.106	-0.81	91.463	-0.261	-1.14
420 (HT)	_			33.839	-0.045	-0.55	45.030	-0.089	-0.77	49.819	-0.094	-0.80	86.550	-0.231	-1.12

<sup>&</sup>lt;sup>7</sup> See Appendix G for expanded response table

<sup>8</sup> Cernox sensors do not follow a standard response curve — the listed resistance ranges are typical, but can vary widely; consult Lake Shore to choose a specific range


Cernox® RTDs

#### Magnetic field dependence data for sample CX RTDs

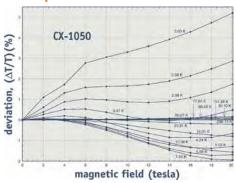


Typical temperature reading errors for operation of CX-1050 sensors in magnetic fields at temperatures from 2.03 K to 286 K. "Low temperature thermometry in high magnetic fields VII. Cernox® sensors to 32 T," B. L. Brandt, D. W. Liu and L. G. Rubin; Rev. Sci. Instrum., Vol. 70, No. 1, 1999, pp 104-110.

#### **Neutrons and gamma rays**

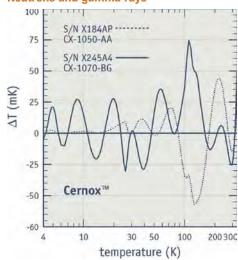

#### **Typical calibration shifts**

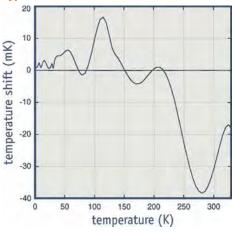


Typical calibration shift after 200 thermal shocks from 305 K to 77 K for a Model CX-1030 temperature sensor  $(\Delta T = 1 \text{ mK at } 4.2 \text{ K and } 10 \text{ mK at } 100 \text{ K}).$ 

# 0.005 in ±0.001 in [0.127 mm ±0.025 mm] **CX-BR** 0.008 in ±0.001 in [0.203 mm ±0.025 mm] 0.030 in [0.762 mm]

General tolerance of ±0.002 in [±0.051 mm] unless otherwise noted

# **CX-SD** 0.125 in [3.175 mm] Ø0.008 in [0.203 mm] Approx. 0.042 in [1.080 mm]

General tolerance of  $\pm 0.005$  in [ $\pm 0.127$  mm] unless otherwise noted

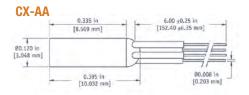


General tolerance of ±0.005 in [±0.127 mm] unless otherwise noted

#### **Physical specifications**

	Mass	Lead type	Internal atmosphere	Sensor materials used			
Bare chip (BC), (BG), (BR)	≤ 3.0 mg	BR: none BG: two 2 mil (44 AWG) bare gold 25 mm long wires BC: two 2.5 mil (42 AWG) bare copper 25 mm long wires	NA	Ceramic oxynitride, gold pads and sapphire substrate with Au Pt Mo back (chip in all models)			
Hermetic ceramic package (SD)	≈ 40 mg	2 gold-plated copper	Vacuum	Chip mounted on sapphire base with alumina body and lid, Mo/Mn with nickel and gold plating on base and lid, gold-tin solder as hermetic lid seal, 60/40 SnPb solder used to attach leads			
Copper canister ≈ 390 mg package (AA)		4 phosphor bronze with HML heavy build insulation attached with epoxy strain relief at sensor	Helium 4 ( <sup>4</sup> He) is standard	Chip mounted in a gold plated cylindrical copper can			

#### AA package

Wires with the same color code are connected to the same side of the sensor (looking at epoxy seal with leads toward user)

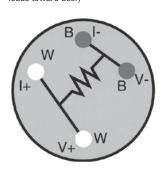



#### Ordering information

Uncalibrated sensor—Specify the model number in the left column only, for example CX-1050-CD.

Calibrated sensor—Add the calibration range suffix code to the end of the model number, for example CX-1050-CD-1.4L.



Cernox® RTD	Calibration range suffix codes Numeric figure is the low end of the calibration Letters represent the high end: L=325 K, M=420 K										
	Uncal	0.1L	0.1M	0.3L	0.3M	1.4L	1.4M	4L	4M	20L	20M
CX-1010-AA, -BC, -B0, -CD, -ET, -LR, -MT	•	-				-					
CX-1010-BG-HT, -BR-HT											
CX-1010-CO-HT, -CU-HT, SD-HT	•	-	•			-	•				
CX-1030-AA, -BC, -B0, -CD, -ET, -LR, -MT	•			-							
CX-1030-BG-HT, -BR-HT											
CX-1030-CO-HT, -CU-HT, -SD-HT	•			-	•		•				
CX-1050-AA, -BC, -B0, -CD, -ET, -LR, -MT	•										
CX-1050-BG-HT, -BR-HT											
CX-1050-CO-HT, -CU-HT, -SD-HT	•						•				
CX-1070-AA, -BC, -B0, -CD, -ET, -LR, -MT	•							-			
CX-1070-BG-HT, -BR-HT											
CX-1070-CO-HT, -CU-HT, -SD-HT	•							-	•		
CX-1080-AA, -BC, -B0, -CD, -ET, -LR, -MT	•									-	
CX-1080-BG-HT, -BR-HT											
CX-1080-CO-HT, -CU-HT, -SD-HT										-	-

ADD -P Add spot-welded platinum leads to the SD package for Cernox® sensors only

#### Accessories available for sensors

SN-CO-C1 CO style sensor clamps for SD package Expanded interpolation table **ECRIT** 8000 Calibration report on CD-ROM COC-SEN Certificate of conformance


## Accessories suggested for installation-

see Accessories section for full descriptions VGE-7031 varnish Stycast® epoxy Phosphor bronze wire Apiezon® grease 90% Pb, 10% Sn solder Manganin wire CryoCable<sup>™</sup> Indium solder

#### **Packaging options**

Cernox® RTDs

For more information on sensor packages and mounting adapters, see page 21.

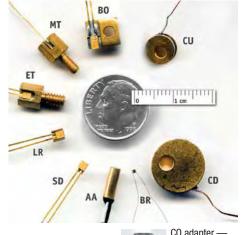



CO adapter spring loaded clamp for easy sensor interchangeability







#### DT-670 Silicon Diodes

#### **DT-670-SD features**

- Best accuracy across the widest useful temperature range—1.4 K to 500 K—of any silicon diode in the industry
- Tightest tolerances for 30 K to 500 K applications of any silicon diode to date
- Rugged, reliable Lake Shore SD package designed to withstand repeated thermal cycling and minimize sensor self-heating
- Conformance to standard DT-670 temperature response curve
- Variety of packaging options

#### **DT-670E-BR features**

- Temperature range: 1.4 K to 500 K
- Bare die sensors with the smallest size and fastest thermal response time of any silicon diode on the market today
- Non-magnetic sensor

#### **DT-621-HR features**

- Temperature range: 1.4 K to 325 K (uncalibrated down to 20 K)
- Non-magnetic package
- Exposed flat substrate for surface mounting



DT-670 Series silicon diodes offer better accuracy over a wider temperature range than any previously marketed silicon diodes. Conforming to the Curve DT-670 standard voltage versus temperature response curve, sensors within the DT-670 series are interchangeable, and for many applications do not require individual calibration. DT-670 sensors in the SD package are available in four tolerance bands—three for general cryogenic use across the 1.4 K to 500 K temperature range, and one that offers superior accuracy for applications from 30 K to room temperature.

DT-670-SD diodes are available with calibration across the full 1.4 K to 500 K temperature range.

The bare die sensor, the DT-670E-BR, provides the smallest physical size and fastest thermal response time of any silicon diode on the market today. This is an important advantage for applications where size and thermal response time are critical, including focal plane arrays and high temperature superconducting filters for cellular communication.

#### **Packaging options**

BO, BR, CO, CU, CY, DI, ET, LR, MT




CAUTION: These sensors are sensitive to electrostatic discharge (ESD). Use ESD precautionary procedures when handling, or making mechanical or electrical connections to these devices in order to avoid performance degradation or loss of functionality.


# The Lake Shore SD package — the most rugged, versatile package in the industry

The SD package, with direct sensor-to-sapphire base mounting, hermetic seal, and brazed Kovar leads, provides the industry's most rugged, versatile sensors with the best sample to chip connection. Designed so heat coming down the leads bypasses the chip, it can survive several thousand hours at 500 K (depending on model) and is compatible with most ultra high vacuum applications. It can be indium soldered to samples without shift in sensor calibration. If desired, the SD package is also available without Kovar leads.

#### Typical DT-670 diode voltage

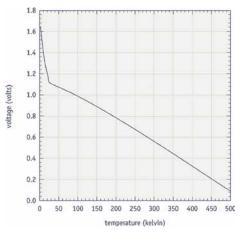

#### Typical DT-670 diode sensitivity

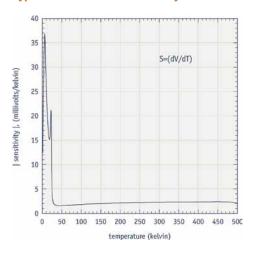



Standard curve Curve DT-670—see next page

Recommended excitation 10  $\mu A \pm 0.1\%$ 

Max reverse voltage 40 V

Max current before damage 1 mA continuous or 100 mA pulsed

Dissipation at recommended excitation 16  $\mu$ W at 4.2 K; 10  $\mu$ W at 77 K; 5  $\mu$ W at 300 K

Thermal response time SD: typical <10 ms at 4.2 K, 100 ms at 77 K, 200 ms at 305 K; BR: 1 ms at 4.2 K, 13 ms at 77 K, 20 ms at 305 K

**Use in radiation** Recommended for use only in low level radiation—see Appendix B

**Use in magnetic field** Not recommended for use in magnetic field applications below 60 K. Low magnetic field dependence when used in fields up to 5 tesla above 60 K—see Appendix B

Reproducibility¹ ±10 mK at 4.2 K

Soldering standard J-STD-001 Class 2

#### Range of use

Package	Minimum limit	Maximum limit
SD, CU-HT, BR	1.4 K	500 K
CU, LR, CY, ET, MT, BO, HR	1.4 K	420 K

#### DT-621-HR miniature silicon diode

The DT-621 miniature silicon diode temperature sensor is configured for installation on flat surfaces. Due to the absence of magnetic materials in its construction, this package is suited for applications where minimal interaction between the diode and sample space magnetic field is desired. The DT- 621 sensor package exhibits precise, monotonic temperature response over its useful range. The sensor chip is in direct contact with the epoxy dome, which causes increased voltage below 20 K and prevents full range Curve DT-670 conformity. For use below 20 K, calibration is required.



#### **Calibrated accuracy**

	Typical sensor accuracy <sup>2</sup>
1.4 K	±12 mK
4.2 K	±12 mK
10 K	±12 mK
77 K	±22 mK
300 K	±32 mK
500 K	±50 mK

<sup>&</sup>lt;sup>2</sup> [(Calibration uncertainty)<sup>2</sup> +(reproducibility)2]<sup>0.5</sup> for more information see Appendices B, D, and E

#### **Temperature response data table (typical)**

	DT-	670	DT-62	1-HR
	V (volts) dV/dT (mV/K)		V (volts)	dV/dT (mV/K)
1.4 K	1.64	-12.5	_	_
4.2 K	1.58	-31.6	1.678	-35
10 K	1.38	-26.8	_	_
77 K	1.03	-1.73	1.03	1.73
305 K	0.560	-2.30	0.560	-2.3

See Appendix G for expanded response table

#### Long-term stability

4.2 K     ±10 mK     ±40 mK       77 K     ±40 mK     ±100 mK       305 K     ±25 mK     ±50 mK		Use to 305 K <sup>3</sup>	Use to 500 K <sup>4</sup>
210 1111	4.2 K	±10 mK	±40 mK
305 K ±25 mK ±50 mK	77 K	±40 mK	±100 mK
	305 K	±25 mK	±50 mK
500 K — ±150 mK	500 K	_	±150 mK

<sup>&</sup>lt;sup>3</sup> Long-term stability data is obtained by subjecting sensor to 200 thermal shocks from 305 K to 77 K

#### Standard curve DT-670 tolerance bands

	2 K to 100 K	100 K to 305 K	305 K to 500 K
Band A	±0.25 K	±0.5 K	±0.5 K
Band A1	±0.25 K	$\pm 1.5\%$ of temp	$\pm 1.5\%$ of temp
Band B	±0.5 K	±0.5 K	±0.33% of temp
Band B1	±0.5 K	$\pm 1.5\%$ of temp	±1.5% of temp
Band C	±1 K	±1 K	±0.50% of temp

	30 K to 100 K	100 K to 305 K	305 K to 500 K
Band D <sup>5</sup>	±0.25 K	±0.50 K	$\pm 0.20\%$ of temp

 $<sup>^{5}~</sup>$  For T  $<30~K~\pm1.5~K$ 

	2 K to 100 K	100 K to 500 K
DT-670E-BR	±1.5 K typical	±1.5% of temp typical


	20 K to 325 K
DT-621-HR	$\pm 2.5$ K or $\pm 1.5\%$ of temperature,
	whichever is greater

#### **Physical specifications**

	Mass	Lead type	Lead polarity	Sensor materials used
DT-670-SD	37 mg	2—nickel and gold plated Kovar	Positive lead on right with package lid up and leads towards user	Sapphire base with alumina body & lid. Molybdenum/manganese metallization on base and lid top with nickel and gold plating. Gold tin solder as hermetic seal.
DT-670E-BR (bare die)	72.7 µg	None	Positive connection made through bottom of chip; negative connection made on base pad on top of chip	Silicon chip with aluminum metallization on chip contacts.
DT-621-HR	23 mg	2—platinum ribbon with tinned 60/40 SnPb solder	Positive lead is right-hand ribbon with platinum disk down and leads towards user	Sensing element is mounted to a platinum disk and covered with a dome of Stycast® 2850 epoxy

Short-term reproducibility data is obtained by subjecting sensor to repeated thermal shocks from 305 K to 4.2 K

<sup>&</sup>lt;sup>4</sup> Based on 670 h of baking at 500 K

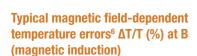


Package base parallel to field B							
	1T 2T 3T 4T 5T						
4.2 K	-200	-300	-350	-400	-500		
20 K	-10	-20	-25	-30	-40		
40 K	-4	-6	-8	-10	-12		
60 K	-0.5	-1	-2	-3	-3.5		
80 K	< 0.1	-0.5	-0.8	-1.1	-1.5		
300 K	<-0.1	<-0.1	<-0.1	<-0.1	<-0.1		

Package base perpendicular to field B							
	1T	2 T	3 T	4 T	5 T		
4.2 K	-8	-9	-11	-15	-20		
20 K	-4	-5	-5	-5	-10		
40 K	-1.5	-3	-4	-5	-5.5		
60 K	-0.5	-0.7	-0.8	-1	-1.1		
80 K	-0.1	-0.3	-0.5	-0.6	-0.7		
300 K	< 0.1	0.2	0.5	0.6	0.6		

To minimize magnetic field-induced temperature errors, the sensor should be oriented so that the package base is perpendicular to the magnetic field flux lines this results in the diode current being parallel to the magnetic field


#### **DT-670-SD**

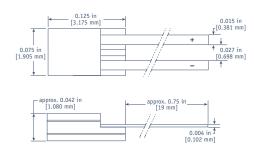




**CAUTION:** (+) lead connected electrically to external braze ring—take care not to cause a short

#### DT-670E-BR

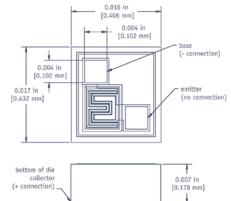


#### **DT-621-HR**

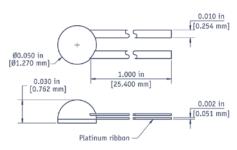


General tolerance of  $\pm 0.005$  in [ $\pm 0.127$  mm] unless otherwise noted

#### **DT-670 temperature response curve**

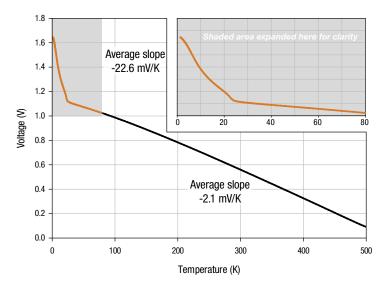

#### **Curve DT-670 tolerance bands**

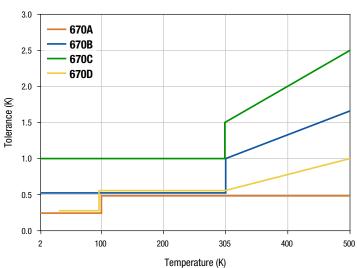



#### DT-670 Series expanded temperature response data table

T (K)	Voltage (V)	dV/dT (mV/K)									
1.4	1.644290	-12.5	6.0	1.51541	-36.7	28.0	1.110421	-2.25	160.0	0.868518	-2.07
1.5	1.642990	-13.6	6.5	1.49698	-36.9	29.0	1.108261	-2.08	170.0	0.847659	-2.10
1.6	1.641570	-14.8	7.0	1.47868	-36.2	30.0	1.106244	-1.96	180.0	0.826560	-2.12
1.7	1.640030	-16.0	7.5	1.46086	-35.0	31.0	1.104324	-1.88	190.0	0.805242	-2.14
1.8	1.638370	-17.1	8.0	1.44374	-33.4	32.0	1.102476	-1.82	200.0	0.783720	-2.16
1.9	1.636600	-18.3	8.5	1.42747	-31.7	33.0	1.100681	-1.77	210.0	0.762007	-2.18
2.0	1.634720	-19.3	9.0	1.41207	-29.9	34.0	1.098930	-1.73	220.0	0.740115	-2.20
2.1	1.632740	-20.3	9.5	1.39751	-28.3	35.0	1.097216	-1.70	230.0	0.718054	-2.21
2.2	1.630670	-21.1	10.0	1.38373	-26.8	36.0	1.095534	-1.69	240.0	0.695834	-2.23
2.3	1.628520	-21.9	10.5	1.37065	-25.5	37.0	1.093878	-1.64	250.0	0.673462	-2.24
2.4	1.626290	-22.6	11.0	1.35820	-24.3	38.0	1.092244	-1.62	260.0	0.650949	-2.26
2.5	1.624000	-23.2	11.5	1.34632	-23.2	39.0	1.090627	-1.61	270.0	0.628302	-2.27
2.6	1.621660	-23.6	12.0	1.33499	-22.1	40.0	1.089024	-1.60	273.0	0.621141	-2.28
2.7	1.619280	-24.0	12.5	1.32416	-21.2	42.0	1.085842	-1.59	280.0	0.605528	-2.28
2.8	1.616870	-24.2	13.0	1.31381	-20.3	44.0	1.082669	-1.59	290.0	0.582637	-2.29
2.9	1.614450	-24.4	13.5	1.30390	-19.4	46.0	1.079492	-1.59	300.0	0.559639	-2.30
3.0	1.612000	-24.7	14.0	1.29439	-18.6	48.0	1.076303	-1.60	310.0	0.536542	-2.31
3.1	1.609510	-25.1	14.5	1.28526	-17.9	50.0	1.073099	-1.61	320.0	0.513361	-2.32
3.2	1.606970	-25.6	15.0	1.27645	-17.3	52.0	1.069881	-1.61	330.0	0.490106	-2.33
3.3	1.604380	-26.2	15.5	1.26794	-16.8	54.0	1.066650	-1.62	340.0	0.466760	-2.34
3.4	1.601730	-26.8	16.0	1.25967	-16.3	56.0	1.063403	-1.63	350.0	0.443371	-2.34
3.5	1.599020	-27.4	16.5	1.25161	-15.9	58.0	1.060141	-1.64	360.0	0.419960	-2.34
3.6	1.596260	-27.9	17.0	1.24372	-15.6	60.0	1.056862	-1.64	370.0	0.396503	-2.35
3.7	1.59344	-28.4	17.5	1.23596	-15.4	65.0	1.048584	-1.67	380.0	0.373002	-2.35
3.8	1.59057	-29.0	18.0	1.22830	-15.3	70.0	1.040183	-1.69	390.0	0.349453	-2.36
3.9	1.58764	-29.6	18.5	1.22070	-15.2	75.0	1.031651	-1.72	400.0	0.325839	-2.36
4.0	1.58465	-30.2	19.0	1.21311	-15.2	77.35	1.027594	-1.73	410.0	0.302161	-2.37
4.2	1.57848	-31.6	19.5	1.20548	-15.3	80.0	1.022984	-1.75	420.0	0.278416	-2.38
4.4	1.57202	-32.9	20.0	1.197748	-15.6	85.0	1.014181	-1.77	430.0	0.254592	-2.39
4.6	1.56533	-34.0	21.0	1.181548	-17.0	90.0	1.005244	-1.80	440.0	0.230697	-2.39
4.8	1.55845	-34.7	22.0	1.162797	-21.1	100.0	0.986974	-1.85	450.0	0.206758	-2.39
5.0	1.55145	-35.2	23.0	1.140817	-20.8	110.0	0.968209	-1.90	460.0	0.182832	-2.39
5.2	1.54436	-35.6	24.0	1.125923	-9.42	120.0	0.949000	-1.94	470.0	0.159010	-2.37
5.4	1.53721	-35.9	25.0	1.119448	-4.60	130.0	0.929390	-1.98	480.0	0.135480	-2.33
5.6	1.53000	-36.2	26.0	1.115658	-3.19	140.0	0.909416	-2.01	490.0	0.112553	-2.25
5.8	1.52273	-36.5	27.0	1.112810	-2.58	150.0	0.889114	-2.05	500.0	0.090681	-2.12



#### Ordering information

#### **Uncalibrated sensor**

Step 1: Choose diode series, for example DT-670.

Step 2: Choose tolerance band (if applicable), for example DT-670A.

**Step 3:** Choose package or mounting adapter—if ordering adapter, substitute the adapter suffix for the SD suffix, for example DT-670A-CU.

#### **Calibrated sensor**

Step 1: Choose diode series, for example DT-670.

**Step 2:** Choose package or mounting adapter—if ordering adapter, substitute the adapter suffix for the SD suffix, for example DT-670-CU.

**Step 3:** Specify the calibration range suffix code after the model number and package suffix, for example DT-670-CU-1.4L.

DT-670	Numeric figure	Calibration range suffix codes Numeric figure is the low end of the calibration Letters represent the high end: L=325 K, H=500 K					
Model number	Uncal	1.4L	1.4H	70L	70H		
DT-621-HR							
DT-670A-SD							
DT-670A1-SD	-						
DT-670B-SD							
DT-670B1-SD							
DT-670C-SD							
DT-670D-SD	-						
DT-670-SD				-	-		
Mounting adapters are replace SD suffix with n			nge—				
CO	-	•	-	-	-		
CU, LR, CY, ET, BO, MT	-	-		-			
CU-HT	-	•		-			
DI							
DT-670E-BR-10	-	bare chip silic	con diode sensor, o	quantity 10			

Note: upper temperature limit package dependent—see Sensor Packages section Other packaging available by special order—please consult Lake Shore

#### Accessories available for sensors

SN-CO-C1 CO style sensor clamps for SD package

ECRIT Expanded interpolation table
8000 Calibration report on CD-ROM
COC-SEN Certificate of conformance












Accessories suggested for installation—
see Accessories section for full descriptions

Stycast® epoxy
Apiezon® grease
90% Pb, 10% Sn solder
Indium solder
VGE-7031 varnish
Phosphor bronze wire
Manganin wire

#### **Packaging options**

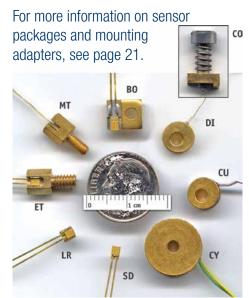

CO-spring loaded clamp for easy sensor interchangeability

#### **Upgrade conversion chart**

	From:	То:
Sensor	DT-470	DT-670
Band	11	Α
	11A	A1
	12	В
	12A	B1
	13	С



# See the appendices for a detailed description of: Installation

Uncalibrated sensors SoftCal™ Calibrated sensors CalCurve™ Sensor packages






#### **GaAlAs Diodes**

#### **TG-120-SD features**

- Monotonic temperature response from 1.4 K to 500 K
- Excellent sensitivity (dV/dT) at temperatures below 50 K
- Relatively low magnetic field-induced errors
- Rugged, reliable Lake Shore SD package designed to withstand repeated thermal cycling and minimize sensor self-heating
- Variety of packaging options

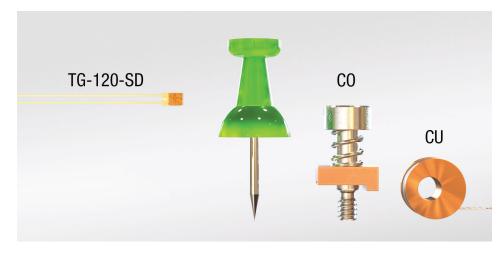

The TG-120 gallium-aluminum-arsenide (GaAlAs) diode temperature sensors are particularly well suited for low to moderate magnetic field applications at low temperatures. The GaAlAs sensing element exhibits high sensitivity (dV/dT) at low temperatures. Voltage-temperature characteristics are monotonic over the sensor's useful range from 1.4 K to 500 K (see data plots below).

Gallium-aluminum-arsenide diodes are direct band-gap, single junction devices that produce small output variances in the presence of magnetic fields. Consequently, their low magnetic field dependence makes them ideally suited for applications in moderate magnetic fields up to five tesla.


#### **Packaging options**

P, PL, SD, CO, CU


#### The Lake Shore SD package the most rugged, versatile package in the industry

The SD package, with direct sensor-to-sapphire base mounting, hermetic seal, and brazed Kovar leads, provides the industry's most rugged, versatile sensors with the best sample to chip connection. Designed so heat coming down the leads bypasses the chip, it can survive several thousand hours at 500 K (depending on model) and is compatible with most ultra high vacuum applications. It can be indium soldered to samples without shift in sensor calibration. If desired, the SD package is also available without Kovar leads.

#### Typical GaAlAs diode sensitivity

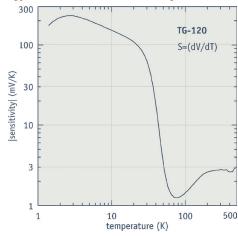

#### Typical GaAlAs diode voltage

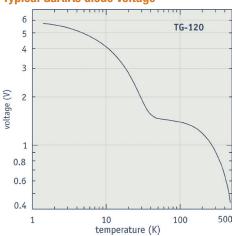



Standard curve Not applicable

Recommended excitation 10  $\mu A \pm 0.1\%$ 

Maximum reverse voltage (diode) 2 V

Maximum forward current (diode) 500 mA

Dissipation at recommended excitation Typical 50  $\mu W$  max at 4.2 K, 14  $\mu W$  at 77 K, 10  $\mu W$  at 300 K

**Thermal response time** (typical) SD: <10 ms at 4.2 K; P and PL: 100 ms at 4.2 K, 250 ms at 77 K, 3 s at 305 K;

**Use in radiation** Recommended for use only in low level radiation—see Appendix B

**Use in magnetic field** Low magnetic field dependence when used in fields up to 5 T above 60 K—see Appendix B

Reproducibility<sup>1</sup> ±10 mK at 4.2 K

Soldering standard J-STD-001 Class 2

Short-term reproducibility data is obtained by subjecting sensor to repeated thermal shocks from 305 K to 4.2 K

#### Range of use

	Minimum Iimit	Maximum limit
TG-120-SD, CU-HT	1.4 K	500 K
TG-120-P, PL (discontinued)	1.4 K	325 K

#### **Calibrated accuracy**

	Typical sensor accuracy <sup>2</sup>	Long-term stability <sup>3</sup>
1.4 K <sup>4</sup>	±12 mK	±25 mK
4.2 K <sup>4</sup>	±12 mK	±15 mK
10 K	±12 mK	±25 mK
77 K	±22 mK	±15 mK
300 K	±32 mK	±50 mK
500 K	±50 mK	_

- <sup>2</sup> [(Calibration uncertainty)<sup>2</sup> + (reproducibility)<sup>2</sup>]<sup>0.5</sup> for more information see Appendices B, D, and E
- Long-term stability data is obtained by subjecting sensor to 200 thermal shocks from 305 K to 77 K
- <sup>4</sup> Under 10 K calibration valid in vacuum only

#### **Temperature response data table (typical)**

	TG-120				
	V (volts)	dV/dT (mV/K)			
1.4 K	5.3909	-97.5			
4.2 K	4.7651	-214			
20 K	2.5341	-97.5			
77 K	1.4222	-1.24			
300 K	0.8978	-2.85			
500 K	0.3778	-3.15			

See Appendix G for expanded response table

# Typical magnetic field-dependent temperature errors<sup>5</sup> ΔT/T (%) at B (magnetic induction)


	Package base parallel to field B						
	1T	2 T	3 T	4 T	5 T		
4.2 K	2.9	3.8	3.7	2.8	1		
30 K	0.2	0.2	0.3	0.3	0.2		
78 K	<0.1	< 0.1	0.17	0.16	0.1		
300 K	0.1	0.1	0.1	0.1	0.1		


To minimize magnetic field-induced temperature errors, the sensor should be oriented so that the package base is perpendicular to the magnetic field flux lines this results in the diode current being parallel to the magnetic field

#### Physical specifications


	Mass	Lead type	Internal atmosphere	Lead polarity	Sensor materials
TG-120-SD	38 mg	2 platinum, welded to package; CAUTION: leads are delicate	Hermetically sealed in vacuum	Positive lead on right with package lid up and leads toward user	Chip mounted on sapphirebase with alumina body andlid, Mo/ Mn metallization onbase & lid top with nickel and gold plating
TG-120-P (discontinued)	79 mg	2 phosphor bronze, insulated with heavy build Polyimide	Air	Short (+) Long (–)	BeO ceramic header set into a gold plated copper cylinder
TG-120-PL (discontinued)	20 mg	2 platinum	Solid epoxy	Short (+) Long (–)	Constructed with platinum, Stycast® epoxy, and alumina

#### **TG-120-SD**

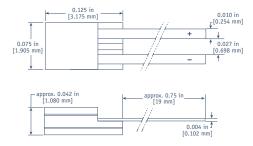

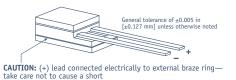


#### TG-120-P

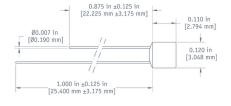

General tolerance of ±0.005 in [±0.127 mm] unless otherwise noted

#### **TG-120-PL**

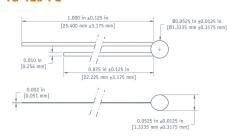

General tolerance of  $\pm 0.005$  in [ $\pm 0.127$  mm] unless otherwise noted



#### Ordering information

#### **Uncalibrated sensor**

Specify the model number in the left column only, for example TG-120-SD.

#### **Calibrated sensor**

Add calibration range suffix code to the end of the model number, for example TG-120-SD-1.4L.

GaAlAs diode	Numeric figure is	Calibration range suffix codes Numeric figure is the low end of the calibration Letters represent the high end: L=325 K, H=500 K					
Part number	Uncal	1.4L	1.4H	70L	70H		
TG-120-SD							
TG-120-C0		-	•	-	•		
TG-120-CU		-		-			
TG-120-CU-HT			•				

\*Below 10 K, calibration is valid in vaccuum only

Other packaging available by special order—please consult Lake Shore

#### Accessories available for sensors

Expanded interpolation table **ECRIT** 8000 Calibration report on CD-ROM COC-SEN Certificate of conformance









#### Accessories suggested for installationsee Accessories section for full descriptions

Stycast® epoxy Apiezon® grease 90% Pb, 10% Sn solder Indium solder VGE-7031 varnish Phosphor bronze wire Manganin wire CryoCable™

#### **Packaging options**

**GaAlAs Diodes** 

For more information on sensor packages and mounting adapters, see page 21.



CO adapter spring loaded clamp for easy sensor interchangeability



SoftCal<sup>™</sup> Calibrated sensors CalCurve™ Sensor packages



#### Germanium RTDs

#### **Germanium features**

- Recognized as a "Secondary Standard Thermometer"
- High sensitivity provides submillikelvin control at 4.2 K and below
- Excellent reproducibility better than ±0.5 mK at 4.2 K
- Various models for use from 0.05 K to 100 K
- Excellent resistance to ionizing radiation

Lake Shore germanium resistance temperature sensors are recognized as "Secondary Standard Thermometers" and have been employed in the measurement of temperature from 0.05 K to 30 K for more than 40 years.

Germanium sensors have a useful temperature range of about two orders of magnitude. The exact range depends upon the doping of the germanium element. Sensors with ranges from below 0.05 K to 100 K are available. Between 100 K and 300 K, dR/dT changes sign and dR/dT above 100 K is very small for all models. Sensor resistance varies from several ohms at its upper useful temperature to several tens of kilohms at its lower temperature. Because device sensitivity increases rapidly with decreasing temperature, a high degree of resolution is achieved at lower temperatures, making these resistors very useful for submillikelvin control at 4.2 K and below.




The sensors offer excellent stability, and  $\pm 0.5$  mK reproducibility at 4.2 K. The germanium resistor is usually the best choice for high-accuracy work below 30 K. Use in a magnetic field is not recommended.


#### **Packaging options**

AA,CD


#### **Typical germanium resistance**

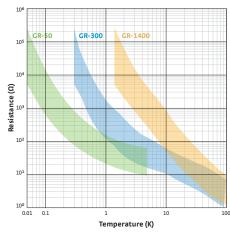

#### Typical germanium sensitivity

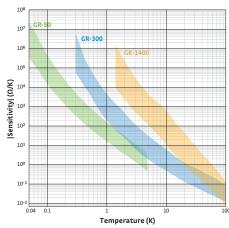

#### Typical germanium dimensionless sensitivity

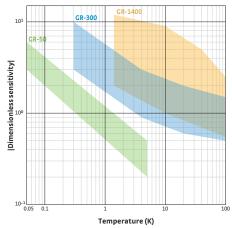



#### Standard curve Not applicable

**Recommended excitation**<sup>1</sup> 20  $\mu$ V (0.05 K to 0.1 K); 63  $\mu$ V (0.1 K to 1 K); 10 mV or less for T > 1 K

Dissipation at recommended excitation  $10^{-13}$  W at 0.05 K,  $10^{-7}$  W at 4.2 K (temperature and model dependent)

Thermal response time 200 ms at 4.2 K, 3 s at 77 K

**Use in radiation** Recommended for use in ionizing radiation environments—see Appendix B

**Use in magnetic field** Because of their strong magnetoresistance and associated orientation effect, germanium sensors are of very limited use in magnetic fields—see Appendix B

Soldering standard J-STD-001 Class 2

#### Reproducibility

	Short term <sup>2</sup>	Long term <sup>3</sup>
4.2 K	±0.5 mK	±1 mK/yr
77 K	_	±10 mK/yr

- Recommended excitation for T < 1 K based on Lake Shore calibration procedures using an AC resistance bridge—for more information refer to Appendix D and Appendix E
- <sup>2</sup> Short-term reproducibility data is obtained by subjecting sensor to repeated thermal shocks from 305 K to 4.2 K
- <sup>3</sup> Long-term stability data is obtained by subjecting sensor to 200 thermal shocks from 305 K to 77 K

#### Range of use


	Minimum limit	Maximum limit
GR-50-AA	<0.05 K	5 K
GR-300-AA	0.3 K	100 K
GR-1400-AA	1.4 K	100 K

#### Calibrated accuracy<sup>4</sup>

	Typic GR-50	al sensor accu GR-300	ıracy⁴   GR-1400
0.05 K	±5 mK	_	_
0.3 K	±5 mK	±4 mK	_
0.5 K	±5 mK	±4 mK	_
1.4 K	±6 mK	±4 mK	±4 mK
4.2 K	±6 mK	±4 mK	±4 mK
77 K	_	±25 mK	±15 mK
100 K	_	±32 mK	±18 mK

<sup>&</sup>lt;sup>4</sup> [(Calibration uncertainty)2 + (reproducibility)2]0.5 for more information see Appendices B, D, and E

#### AA package

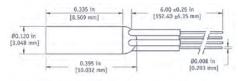

General tolerance of  $\pm 0.005$  in [ $\pm 0.127$  mm] unless otherwise noted

# Typical magnetic field-dependent temperature errors<sup>5</sup> ΔT/T (%) at B (magnetic induction)

Germanium						
	2.5 T	8 T	14 T			
2.0 K	-8	-60	_			
4.2 K	-5 to -20	-30 to -55	-60 to -75			
10 K	-4 to -15	-25 to -60	-60 to -75			
20 K	-3 to -20	-15 to -35	-50 to -80			

<sup>&</sup>lt;sup>5</sup> Long axis of thermometer parallel to applied field

#### **Typical resistance values**

GR-AA	Typical resistance at 4.2 K	Typical resistance range at 4.2 K
50	30 Ω	9 Ω to 65 Ω
300	95 Ω	15 $\Omega$ to 155 $\Omega$
1400	1750 Ω	350 $\Omega$ to 6500 $\Omega$

#### Temperature response data table (typical)—see Appendix G for expanded response table

		GR-50-AA			GR-300-AA			GR-1400-AA	
	R <sup>8</sup> (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	R <sup>8</sup> (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	R <sup>8</sup> (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)
0.05 K	35000	-3642000	-5.2	_	_	_	_	_	_
0.1 K	2320	-71860	-3.1	_	_	_	_	_	_
0.2 K	364.6	-4043	-2.2	_	_	_	_	_	_
0.3 K	164.0	-964.0	-1.8	35180	-512200	-4.4	_	_	_
0.5 K	73.75	-202.9	-1.4	5443	-34800	-3.2	_	_	_
1.0 K	33.55	-31.33	-0.93	875.7	-1901	-2.2	_	_	_
1.4 K	24.73	-13.15	-0.74	448.6	-581.3	-1.8	35890	-94790	-3.7
2.0 K	19.32	-6.167	-0.64	248.8	-187.4	-1.5	11040	-16670	-3.0
4.2 K	13.66	-1.036	-0.32	94.46	-26.56	-1.2	1689	-861.9	-2.1
10 K	_	_	_	33.20	-3.97	-1.2	252.8	-61.95	-2.5
40 K	_	_	_	7.79	-0.235	-1.2	9.57	-0.449	-1.9
77.4 K	_	_	_	3.50	-0.050	-1.1	3.55	-0.050	-1.1
100 K	_	_	_	2.72	-0.024	-0.88	2.80	-0.021	-0.74



# Proper selection of germanium sensors for use below 1 K

Germanium resistance thermometers are often classified according to their 4.2 K resistance value. However, for devices to be used below 1 K, there is no close correlation between the 4.2 K resistance and the suitability of the device as a thermometer. As a result, the Lake Shore low resistance germanium sensors (GR-50-AA and GR-300-AA) are classified according to their lowest useful temperatures, not their 4.2 K resistance values.

The resistance vs. temperature behavior for these devices is typical of all the germanium sensors. As the temperature is lowered, both the resistance and sensitivity (dR/dT) increase logarithmically. The lowest useful temperature is generally limited by the rapidly increasing resistance and the difficulties encountered in measuring high resistance values.

The following recommendations are made concerning the optimum temperature range for using these devices:

GR-50-AA 0.05 K to 1.0 K GR-300-AA 0.3 K to 100 K

Increasingly better temperature resolution is achievable at lower temperatures.

In general, it is recommended you do not purchase a device which has a lower temperature limit than required, since some sensitivity (dR/dT) will be sacrificed at the higher temperatures. For example, a GR-300-AA will have more sensitivity at 1 K than a GR-50-AA.


# GD package 36-inch long Quad-Lead<sup>®</sup> 36 AWG phosphor bronze wire Sensor leads are anchored by a Stycast® coating

#### **Physical specifications**

	Mass	Lead type	Internal atmosphere	Sensor materials used
GR-50-AA GR-300-AA GR-1400-AA	395 mg	4 color coded phosphor bronze with heavy build polyimide, attached with epoxy strain relief at sensor	Helium 4 (4He) at ≥500 Ω, air at <500 Ω	Doped germanium chip mounted strain-free in a gold plated cylindrical copper can


#### **Germanium series construction detail**

The epoxy holding the chip to the header is omitted for germanium devices designed for

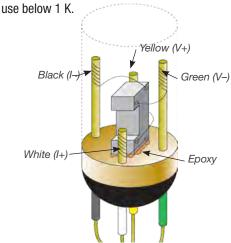

Looking at the wiring end with leads toward user


	Key	Lead Col
N	/+	White
G	V-	Green
Y	V+	Yellow
В	<i>I</i> -	Black

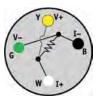

#### **Packaging options**

For more information on sensor packages and mounting adapters, see page 21.



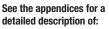

Installation Uncalibrated sensors SoftCal™

Calibrated sensors
CalCurve™
Sensor packages



#### Ordering information

Uncalibrated sensor—Specify the model number in the left column only, for example GR-50-AA.

**Calibrated sensor**—Add the calibration range suffix code to the end of the model number, for example GR-50-AA-0.05A.

Germanium RTD	Calibration range suffix codes  Numeric figure is the low end of the calibration  Letters represent the high end: A=5 K, D=100 K					
Part number	Uncal	0.05A	0.3D	1.4D		
GR-50-AA						
GR-300-AA						
GR-1400-AA						
GR-50-CD						
GR-300-CD			•	•		
GR-1400-CD						

\*NOTE: The GR-50-AA calibration is not useful above 5 K Other packaging available through special order—consult Lake Shore

#### Accessories available for sensors

ECRIT Expanded interpolation table COC-SEN Certificate of conformance











# Ruthenium Oxide (Rox<sup>®</sup>) RTDs

#### **RX-102A features**

- Standard curve interchangeable
- Good radiation resistance
- Useful down to 50 mK
- Low magnetic field-induced errors

#### **RX-102B features**

- Useful below 10 mK; calibrations down to 10 mK available
- Monotonic from 10 mK to 300 K

#### **RX-202A features**

- Standard curve interchangeable
- Good radiation resistance
- Monotonic from 50 mK to 300 K
- 4× improvement in magnetic fieldinduced errors over other ruthenium oxides

#### **RX-103A features**

- Standard curve interchangeable
- Good radiation resistance
- Best choice for interchangeability from 1.4 K to 40 K
- Low magnetic field-induced errors

Ruthenium oxide temperature sensors are thick-film resistors used in applications involving magnetic fields. These composite sensors consist of bismuth ruthenate, ruthenium oxides, binders, and other compounds that allow them to obtain the necessary temperature and resistance characteristics. Each Lake Shore Rox™ model adheres to a single resistance versus temperature curve.

#### **RX-102A**

The RX-102A (1000  $\Omega$  at room temperature) is useful down to 50 mK and has better interchangeability than the RX-202A as well as low magnetic field-induced errors below 1 K.

#### **RX-102B-RS**

The RX-102B-CB (1000  $\Omega$  at room temperature) is useful below 10 mK (calibrations available down to 10 mK) and monotonic from 10 mK to 300 K. The unique package design maximizes thermal connection and minimizes heat capacity at ultra low temperatures. The RX-102B-CB is not interchangeable to a standard curve and not recommended for use in magnetic fields.

#### **RX-202A**

The RX-202A (2000  $\Omega$  at room temperature) has a 4× improvement in magnetic field-induced errors over other commercially available ruthenium oxide temperature sensors with similar resistances and sensitivities. Most ruthenium oxide sensors have a maximum useful temperature limit well below room temperature, where the sensitivity changes from negative to positive. The RX-202A however, is designed to have a monotonic response from 0.05 K up to 300 K.

#### **RX-103A**


The RX-103A (10,000  $\Omega$  at room temperature) has a unique resistance and temperature response curve combined with low magnetic field-induced errors, and is the best choice for interchangeability from 1.4 K to 40 K.

#### **Packaging options**


AA, BR




#### **Typical Rox™ resistance**

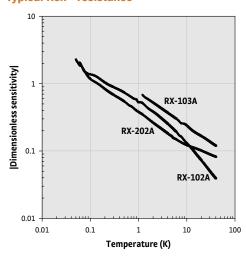

Typical Rox™ sensitivity

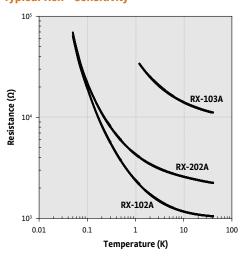

Typical Rox™ dimensionless sensitivity

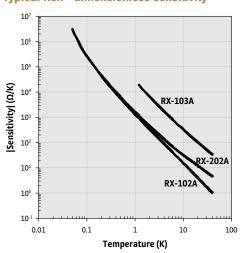

Lake Shore Cryotronics, Inc. | t. 614.891.2244 | f. 614.818.1600 | e. info@lakeshore.com | www.lakeshore.com



**Standard Curve**<sup>1</sup> 102 and 202: 0.05 K to 40 K; 103: 1.4 K to 40 K

**Recommended excitation**<sup>2</sup> RX-102 and RX-202: 20  $\mu$ V (0.05 K to 0.1 K); 63  $\mu$ V (0.1 K to 1.2 K); 10 mV or less for T > 1 K. RX-103: 10 mV or less for T > 1 K

Dissipation at recommended excitation 102 and 202:  $7.5\times10^{-8}$  W at 4.2 K; 103:  $3.2\times10^{-9}$  W at 1.4 K,  $5.5\times10^{-9}$  W at 4.2 K,  $9.6\times10^{-9}$  W at 77 K

Thermal response time  $0.5\ s$  at  $4.2\ K$ ,  $2.5\ s$  at  $77\ K$ 

Use in radiation Recommended—see Appendix B

Use in magnetic field<sup>3</sup> Recommended—see Appendix B

Reproducibility<sup>4</sup> ±15 mK

Soldering standard J-STD-001 Class 2

- 1 102B does not follow a standard curve
- Recommended excitation for T < 1 K based on Lake Shore calibration procedures using an AC resistance bridge—for more information refer to Appendix D and Appendix E
- 3 102B not recommended for use in magnetic fields
- Short-term reproducibility data is obtained by subjecting sensor to repeated thermal shocks from 305 K to 4.2 K

#### Range of use

	Minimum limit	Maximum limit
RX-102A-AA	0.05 K	40 K
RX-102B-CB	0.01 K <sup>5</sup>	40 K
RX-202A-AA	0.05 K	40 K
RX-103A-AA	1.4 K	40 K

<sup>5</sup> Calibrations down to 20 mK available; 10 mK calibrations coming soon

#### **Long-term stability**

	RX-	RX-	RX-	RX-
	102A-AA	102B-RS	202A-AA	103A-AA
4.2 K	±30 mK	±30 mK	±50 mK	±15 mK

#### Calibrated accuracy<sup>6</sup>

	RX-102A- AA	RX-102B- CB	RX-202A- AA	RX-103A- AA
20 mK	_	±2 mK	_	_
50 mK	_	±4 mK	_	_
1.4 K	±16 mK	±16 mK	±16 mK	±16 mK
4.2 K	±16 mK	±16 mK	±16 mK	±17 mK
10 K	±18 mK	±18 mK	±18 mK	±22 mK

<sup>&</sup>lt;sup>6</sup> [(Calibration uncertainty)<sup>2</sup> + (reproducibility)<sup>2</sup>]<sup>0.5</sup> for more information see Appendices B, D, and E

#### **Accuracy: interchangeability**

	RX-102A-AA-M matched	RX-102A-AA unmatched	RX-202A-AA-M matched	RX-202A-AA unmatched	RX-103A-AA-M matched	RX-103A-AA unmatched
0.05 K	±5 mK	±10 mK	±10 mK	±15 mK	_	_
0.3 K	±15 mK	±20 mK	±20 mK	±25 mK	_	_
0.5 K	±20 mK	±25 mK	±25 mK	±30 mK	_	_
1.4 K	±25 mK	±50 mK	±50 mK	±100 mK	±50 mK	±150 mK
4.2 K	±75 mK	±125 mK	±150 mK	±250 mK	±100 mK	±400 mK
20 K	±500 mK	±1.25 K	±1 K	±2.5 K	±700 mK	±2 K
40 K	±1.5 K	±4 K	±2 K	±5 K	±1.5 K	±4 K

#### **Temperature response data table (typical)**—See Appendix G for expanded response table

		102A			202A			103A	
	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)
0.05 K	63765	-2888654	-2.27	69191	-3186379	-2.3	_	_	_
0.1 K	19400	-266199	-1.37	21927	-256913	-1.17	_	_	_
0.3 K	5615	-16647	-0.89	8079	-18420	-0.68	_	_	_
1.4 K	2005	-667	-0.47	3820	-879	-0.32	30745	-13571	-0.62
4.2 K	1370	-80.4	-0.25	2929	-124	-0.18	18149	-1559	-0.36
10 K	1167	-15.2	-0.13	2582	-31.6	-0.12	14083	-337	-0.24
20 K	1089	-3.96	-0.07	2389	-12.1	-0.1	12289	-102	-0.17
40 K	1049	-1.03	-0.04	2243	-4.6	-0.08	11137	-33	-0.12
40 K	1049	-1.06	-0.04	2244	-4.58	-0.08	11150	-21.7	-0.08



# Typical magnetic field-dependent temperature errors $\Delta T/T$ (%) at B (magnetic induction)


Rox™ 102A						
	2.5 T	8 T	14 T	19 T		
2 K	-1.4	-7.9	-13	-17		
3 K	-1.5	-7	-14	-18		
4 K	-0.56	-6.7	-14	-18		
8 K	-1.3	-6.1	-13	-21		
16 K	-0.40	-3.4	-9.6	-16		
23 K	-0.31	-2.2	-6.2	-11		


Rox™ 102B						
	2.5 T	8 T	14 T	19 T		
2 K	3.29	13.82	22.53	27.95		
3 K	3.96	14.68	23.12	29.12		
4 K	3.53	13.92	22.57	28.20		
8 K	1.53	7.53	13.50	17.86		
16 K	0.27	2.14	4.66	6.58		
23 K	0.06	0.79	2.01	3.11		


	Rox™ 202A					
	2.5 T	8 T	14 T	19 T		
2 K	-0.13	-2.2	-3.9	-5.2		
3 K	0.18	-0.68	-2.7	-3.7		
4 K	0.77	0.046	-1.8	-3.2		
8 K	-0.023	0.16	-0.65	-3.0		
16 K	0.03	0.16	-0.48	-1.5		
23 K	-0.05	-0.08	-0.39	-0.92		


	Rox™ 103A					
	2.5 T	8 T	14 T	19 T		
2 K	0.58	1.5	2.2	2.6		
3 K	0.44	1.1	1.7	2.0		
4 K	0.27	0.95	1.4	1.7		
8 K	0.11	0.49	0.71	0.80		
16 K	0.018	0.076	0.089	0.040		
23 K	0.0051	0.0058	-0.0060	-0.095		

#### Magnetic field dependance data for sample Rox™ RTDs

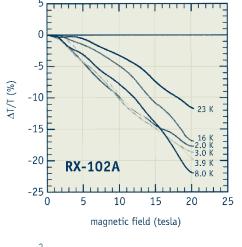


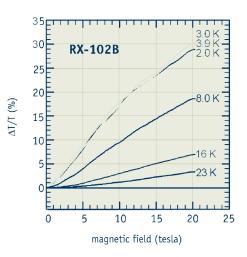

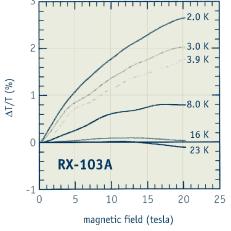

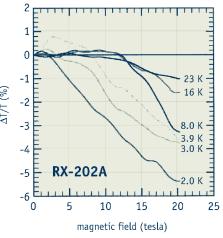

#### **RX-AA**

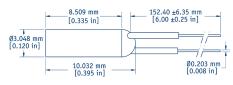


General tolerance of  $\pm 0.127$  mm [ $\pm 0.005$  in] unless otherwise noted

# RX-102B-RS 3.2 mm (0.125 in) 3.2 mm (0.125 in) 2.5 mm (0.100 in) (0.200 in) 1.9 mm (0.075 in) 1.9 mm (0.075 in) 1.4.6 mm (0.375 in) 1.4.6 mm (0.355 in) 3.2 mm (0.125 in) 3.2 mm (0.125 in) 3.2 mm (0.125 in) 3.2 mm (0.125 in)

General tolerance of  $\pm 0.127$  mm (0.005 in) on X.XXX and  $\pm 0.254$  mm (0.01 in) on X.XX unless otherwise noted Mount using a #6 or M3 screw

#### Bare chip (see table on page 50)



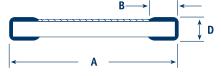



#### **Packaging options**

For more information on sensor packages and mounting adapters, see page 21.



See the appendices for a detailed description of: Installation
Uncalibrated sensors
SoftCal™
Calibrated sensors
CalCurve™
Sensor packages



#### **Packaging**

The  $\text{Rox}^{\text{\tiny{TM}}}$  202A, 102A, and 103A sensors are available in the Lake Shore standard copper AA canister and the 102B is available in the RS copper block package. Two are available as bare chips for applications requiring a smaller sensor or a faster thermal response time. The RX-102A-BR is a bare chip version of RX-102A. This bare chip features wrap-around noble metal contacts that can be soldered to using standard lead/tin solder. The RX-103A-BR is a bare chip version of the RX-103A. This bare chip has wrap-around pretinned contacts that can be soldered to using standard lead/tin solder. The pretinned contacts increase the sensor thickness from 0.25 mm to 0.41 mm. Leads are not attached to these models, so they are not available as matched or calibrated.

See the Physical Specifications for details and individual dimensions.

#### **Physical specifications**

	Mass	Lead type	Internal atmosphere	Materials used	
RX-102A-AA	3.3 g	Two 6 in 32 AWG copper	Air	Thick ruthenium dioxide and bismuth ruthenate	
RX-202A-AA	3.28 g	leads with heavy build Formvar® attached with		films with palladium silver contacts, indium solder, aluminum oxide substrate, sapphire	
RX-103A-AA	3.36 g	epoxy strain relief at sensor–user should branch to 4 (no polarity)		header and copper canister with epoxy sea	
RX-102B-RS	3.5 g	Two 36 AWG copper leads with heavy build polyimide insulation, 15 cm length	NA	Thick ruthenium dioxide and bismuth ruthenate films on aluminum dioxide substrate with palladium silver contacts; epoxy attachment to OFHC adapter; copper leads indium soldered to chip and heat sunk to copper adapter using varnish, covered with copper radiation shield	

Bare chip	A (chip length)	B (pad width)	C (chip width)	D (thickness)	Materials used
RX-102A-BR	1.45 mm (0.057 in)	0.30 mm (0.012 in)	1.27 mm (0.050 in)	0.65 mm (0.022 in)	Thick ruthenium dioxide and bismuth ruthenate
RX-103A-BR	1.40 mm (0.070 in)	0.21 mm (0.010 in)	1.23 mm (0.060 in)	0.41 mm (0.016 in)	films with palladium silver contacts

#### Ordering information

Rox™ RTD	Numeric fig Letters repr	Calibration range suffix codes Numeric figure is the low end of the calibration Letters represent the high end: C=1 K, B=40 K, M = matched (calibration of matched sensors is available—consult Lake Shore)				
Part number	Uncal	0.02C	0.02B	0.05B	0.3B	1.4B
RX-102B-RS						
RX-202A-AA, CD						•
RX-202A-AA-M						
RX-102A-AA, CD						
RX-102A-AA-M						
RX-102A-BR						
RX-103A-AA, CD						
RX-103A-AA-M						
RX-103A-BR	-					

**Note:** the RX-102B-CB is not interchangeable to a standard curve and is not available as matched. Other packaging available through special order—consult Lake Shore

#### Accessories available for sensors

ECRIT Expanded interpolation table 8000 Calibration report on CD-ROM COC-SEN Certificate of conformance









Accessories suggested for installation—
see Accessories section for full descriptions
Styrast® enorgy VGE-7031 varnish

Stycast® epoxy
Apiezon® grease
90% Pb, 10% Sn solder
Indium solder

VGE-7031 varnish Phosphor bronze wire Manganin wire



#### PT-100 Series Platinum RTDs

#### **PT-100 Series features**

- Temperature range: 14 K to 873 K (model dependant)
- Conforms to IEC 751 standards down to 70 K
- High reproducibility: ±5 mK at 77 K
- Low magnetic field dependence above 40 K
- Excellent for use in ionizing radiation
- SoftCal™ calibration available
- Non-magnetic packages available (all PT-102 and PT-103 variants)

#### **Matching**

If your application requires more than one platinum resistor, up to five platinum resistors can be matched to one another to within  $\pm 0.1$  K at liquid nitrogen temperature with the purchase of only one calibration.

PT-100 platinum resistance thermometers (PRTs) are an excellent choice for use as cryogenic temperature sensing and control elements in the range from 30 K to 873 K (-243 °C to 600 °C). Over this temperature span, PRTs offer high repeatability and nearly constant sensitivity (dR/dT). Platinum resistors are also useful as control elements in magnetic field environments where errors approaching one degree can be tolerated. PRTs are interchangeable above 70 K. The use of controlled-purity platinum assures uniformity from one device to another.


PRTs experience rapidly decreasing sensitivity below approximately 30 K. They should be calibrated in order to achieve maximum accuracy for use below 100 K. The plot illustrates platinum sensor conformance to the IEC 751 curve.

#### **Packaging options**


AL, AM




#### **Typical platinum resistance**

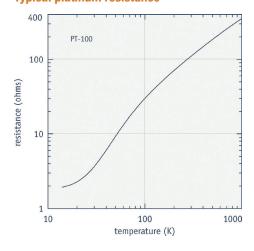

#### Typical platinum sensitivity

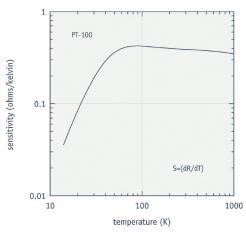

#### Typical platinum dimensionless sensitivity

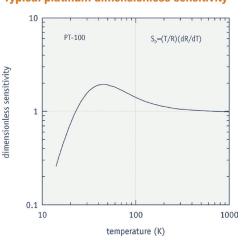



Standard curve IEC 751

Recommended excitation 1 mA

Dissipation at recommended excitation 100  $\mu W$  at 273 K

**Thermal response time** PT-102 & PT-103: 1.75 s at 77 K, 12.5 s at 273 K; PT-111: 2.5 s at 77 K, 20 s at 273 K

**Use in radiation** Recommended for use in ionizing radiation environments—see Appendix B

Use in magnetic field Because of their relatively low magnetic field dependence above 30 K, platinum sensors are useful as control elements in magnetic field applications when some error can be tolerated—see Appendix B

Reproducibility<sup>1</sup> ±5 mK at 77 K

Soldering standard J-STD-001 Class 2

Short-term reproducibility data is obtained by subjecting sensor to repeated thermal shocks from 305 K to 77 K

#### Range of use

	Minimum limit	Maximum limit
PT-102	14 K	873 K
PT-103	14 K	873 K
PT-111	14 K	673 K

#### SoftCal™ accuracy

	30 K to 305 K	305 K to 400 K	400 K to 475 K	475 K to 500 K	500 K to 670 K
2S	±0.25 K	±0.9 K	±1.3 K	±1.4 K	±2.3 K
3S	±0.25 K	±0.25 K	±0.25 K	±1.4 K	±2.3 K

2S: 77 K and 305 K 3S: 77 K, 305 K and 480 K

#### Calibrated accuracy

	Typical sens	Long-term	
	Calibrations to 800 K	All other calibrations	stability <sup>3</sup>
30 K	±10 mK	±10 mK	_
77 K	±12 mK	±12 mK	±10 mK
305 K	±23 mK	±23 mK	_
400 K	±210 mK	±41 mK	_
500 K	±210 mK	±46 mK	_
800 K	±310 mK	_	_

<sup>2</sup> [(Calibration uncertainty)<sup>2</sup> + (reproducibility)<sup>2</sup>]<sup>0.5</sup> for more information see Appendices B, D, and E

# Typical magnetic field-dependent temperature errors<sup>4</sup> $\Delta T/T$ (%) at B (magnetic field)


Package parallel to field B					
	2.5 T	5 T	8 T	14 T	19 T
20 K	20	_	100	250	_
40 K	0.5	1.5	3	6	8.8
87 K	0.04	0.14	0.4	1	1.7
300 K	0.01	0.001	0.02	0.07	0.13

 $<sup>^4</sup>$  Recommended for use when T  $\geq$  30 K

#### **Temperature response data table (typical)** See Appendix G for expanded response table

		PT-100	
	R (Ω)	S (Ω/K)	(T/R)(dR/dT)
20 K	2.2913	0.085	0.74
50 K	9.3865	0.360	1.90
77 K	20.380	0.423	1.60
150 K	50.788	0.409	1.20
300 K	110.354	0.387	1.10
600 K	221.535	0.372	1.00
800 K	289.789	0.360	1.00

#### PT-100 Series interchangeability

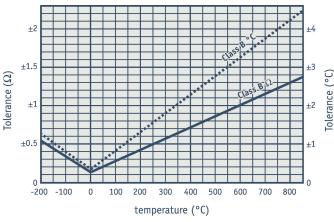


#### **Physical specifications**

	Mass	Lead type	Internal atmosphere	Materials used	UHV compatible	Non- magnetic package
PT-102	250 mg	2, platinum	Partially filled powder	Platinum winding partially supported by a high temperature alumina powder inside a ceramic tube, platinum lead wires	No	Yes
PT-103	120 mg	2, platinum	Fully filled powder	Platinum winding fully supported by a high temperature alumina powder inside a ceramic tube, platinum/rhodium lead wires	Somewhat	Yes
PT-111	52 mg	2, platinum- coated nickel	Solid glass	One platinum band wound onto a glass tube which is protected from the environment by a layer of glaze, platinum coated nickel lead wires	Yes	No

If not heated above 475 K—long-term stability data is obtained by subjecting sensor to 200 thermal shocks from 305 K to 77 K




#### PT-102

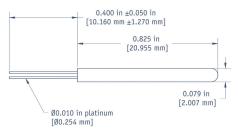


General tolerance of  $\pm 0.010$  in  $[\pm 0.254$  mm] unless otherwise noted

#### PT-103

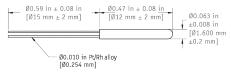


General tolerance of  $\pm 0.010$  in  $[\pm 0.254$  mm] unless otherwise noted

#### PT-111

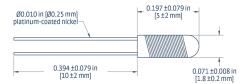


#### PT-102-AL

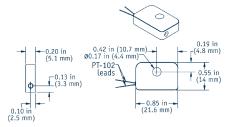

General tolerance of ±0.010 in (±0.254 mm) unless otherwise noted

#### PT-103-AM

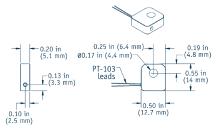

General tolerance of  $\pm 0.010$  in ( $\pm 0.254$  mm) unless otherwise noted

#### Ordering information

**Uncalibrated sensor**—Specify the model number in the left column only, for example PT-103. **Calibrated sensor**—Add the calibration range suffix code to the end of the model number, for example PT-103-14L.

Platinum RTD	Calibration range suffix codes Numeric figure is the low end of the calibration Letters represent the high end: S=SoftCal™, L=325 K, H = 500 K, J = 800 K					
Part number	Uncal	2S	3S	14L	14H	14J
PT-102	•					
PT-102-AL	•	•	•	•	•	•
PT-103						
PT-103-AM	•	•				•
PT-111	•	•	•	•	•	

#### ADD -LN\* Matching PT sensors to ±0.1 K at 77 K

\*MUST be purchased with all matching sensors, as well as with the sensor to be matched

one other uncalibrated PT-102 to within  $\pm 0.1$  K at liquid nitrogen temperature.

#### Motes

- 1. Upper temperature of AL and AM packages is limited to 800 K.
- 2. If your application requires more than one platinum resistor, up to five platinum resistors can be matched with one another to within ±0.1 K at liquid nitrogen temperature with the purchase of only one calibration. If absolute accuracy is required, one of these matched RTDs can be calibrated. For larger quantities, or for different requirements, consult Lake Shore. At the time of order, add -LN to the model number.
  Example: PT-102-14L-LN is a PT-102-LN RTD with a calibration range of 14 K to 325 K that is matched with at least
- 3. For metrological applications below 30 K, use a germanium RTD. PT-100 sensors are not useful below 14 K for metrology and are of limited use below 30 K for temperature control, due to rapid decline in sensitivity.
- 4. For use above 500 K, anneal at  $T_{\text{max}}$  +10 °C for 4 h.

#### Accessories available for sensors

ECRIT Expanded interpolation table 8000 Calibration report on CD-ROM COC-SEN Certificate of conformance









#### Accessories suggested for installation see Accessories section for full descriptions

Stycast® epoxy Apiezon® grease 90% Pb, 10% Sn solder Indium solder VGE-7031 varnish Phosphor bronze wire Manganin wire CryoCable™


#### **Packaging options**

For more information on sensor packages and mounting adapters, see page 21.



See the appendices for a detailed description of: Installation

Uncalibrated sensors SoftCal™ Calibrated sensors CalCurve™ Sensor packages To add length to sensor leads (SMOD), see page 25.

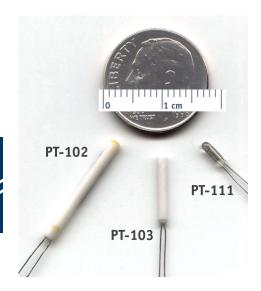



#### Rhodium-Iron RTDs

#### **Rhodium-iron features**

- Good long term stability: ±10 mK from 1.4 K to 325 K
- RF-800 offers a wide temperature range from 0.65 K to 500 K
- Linear response above 100 K
- Excellent resistance to ionizing radiation

Rhodium-iron temperature sensors offer a positive temperature coefficient, monotonic response over a wide temperature range, and high resistance to ionizing radiation.


The RF-800 rhodium-iron resistance sensor features monotonically decreasing resistivity from 500 K to 0.65 K, although sensitivity (dR/dT) falls off in the region of 30 K. From 100 K to 273 K the resistance changes linearly with temperature to within 1 K. RF-800-4 sensors also exhibit monotonic response at higher temperatures, hence their adaptability for use over the broad range from 1.4 K to 500 K.

#### **Packaging options**


RF-800-4



#### **RF-800**

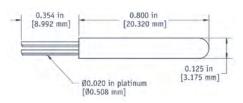


General tolerance of ±0.005 in [±0.127 mm] unless otherwise noted

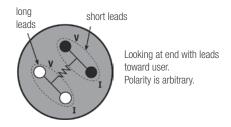


#### **Typical rhodium-iron resistance**

# 100 RF-800-4 1 1 10 100 500 temperature (K)

#### Typical rhodium-iron sensitivity

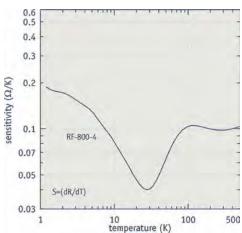

#### Typical rhodium-iron dimensionless sensitivity

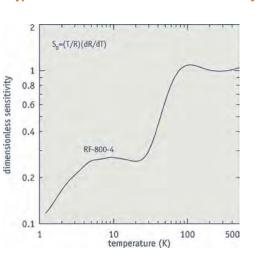



Standard curve Not applicable

Recommended excitation 1 mA

Dissipation at recommended excitation 10  $\mu W$  at 4.2 K, 250  $\mu W$  at 273 K

Thermal response time 10 s at 273 K

**Use in radiation** Recommended for use in ionizing radiation environments—see Appendix B

**Use in magnetic field** Not recommended for use in magnetic fields below 77 K—see Appendix B

Reproducibility<sup>1</sup> ±5 mK at 4.2 K

Soldering standard J-STD-001 Class 2

Short-term reproducibility data is obtained by subjecting sensor to repeated thermal shocks from 305 K to 4.2 K

#### Range of use

	Minimum limit	Maximum limit
RF-800-4	0.65 K	500 K <sup>2</sup>

<sup>&</sup>lt;sup>2</sup> Usable to 800 K, but large and erratic temperature shifts can occur at lower temperatures without proper thermal conditioning

#### **Calibrated accuracy**

	Typical sensor accuracy <sup>3</sup>	Long-term stability <sup>4</sup>
1.4 K	±7 mK	±10 mK
4.2 K	±7 mK	±10 mK
10 K	±8 mK	±10 mK
77 K	±13 mK	±10 mK
305 K	±23 mK	±10 mK
400 K	±41 mK	_
500 K	±42 mK	_

<sup>&</sup>lt;sup>3</sup> [(Calibration uncertainty)<sup>2</sup> + (reproducibility)<sup>2</sup>]<sup>0.5</sup> for more information see Appendices B, D, and E

#### **Physical specifications**

	Size	Mass	Lead type	Internal atmosphere	Materials used
RF-800-4	0.51 mm × 9 mm long	735 mg	4 platinum wire	solid	Alumina and glass cylindrical case—rhodium-iron alloy wire encapsulated in ceramic

#### Typical magnetic field-dependent temperature errors<sup>5</sup> ΔT/T (%) at B (magnetic induction)

Package parallel to field B					
	2.5 T	8 T	14 T	19 T	
4.2 K	11	40	_	_	
20 K	4	_	_	_	
40 K	1.5	12	30	47	
66 K	0.3	2.5	6	9	
87 K	0.2	1.5	4	6	
110 K	0.1	0.9	2.4	_	
190 K	0.03	0.3	0.9	_	
300 K	-0.01	0.1	0.4	_	

 $<sup>^{\</sup>rm 5}~$  Not recommended for use in magnetic fields below 77 K

#### **Temperature response data table (typical)** See Appendix G for expanded response table

		RF-800-4	
	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)
1.4 K	1.5204	0.178	0.16
4.2 K	1.9577	0.135	0.29
20 K	3.1632	0.0461	0.29
77 K	6.8341	0.0959	1.1
150 K	14.463	0.105	1.1
300 K	29.697	0.101	1.0
400 K	39.824	0.103	1.0

#### **Packaging options**

For information on mounting adapters available for use with the SD package, see page 25.





See the appendices for a detailed description of: Installation Uncalibrated sensors SoftCal™ Calibrated sensors CalCurve™ Sensor packages

To add length to sensor leads (SMOD), see page 25.

#### Ordering information

Uncalibrated sensor—Specify the part number in the left column only, for example RF-800-4.

**Calibrated sensor**—Add the calibration range suffix code to the end of the model number, for example RF-800-4-1.4L.

Rhodium-iron RTD	Calibration range suffix codes Numeric figure is the low end of the calibration Letters represent the high end: B=40 K, D=100 K, L=325 K, H=500 K									
Part number	Uncal	1.4B	1.4D	1.4L	1.4H	4B	4D	4L	4H	70L
RF-800-4				-	•			-	•	

Other packaging available through special order—consult Lake Shore

Accessories suggested for installation see Accessories section for full descriptions Stycast® epoxy VGE-7031 varnish

CryoCable™
Apiezon® grease
Manganin wire

VGE-7031 varnish Indium solder 90% Pb, 10% Sn solder Phosphor bronze wire



<sup>&</sup>lt;sup>4</sup> Long-term stability data is obtained by subjecting sensor to 200 thermal shocks from 305 K to 77 K



# Capacitance Temperature Sensors\*

#### **Capacitance features**

- Virtually no magnetic field-induced errors
- Capable of mK control stability in the presence of strong magnetic fields
- Monotonic in C versus T to nearly room temperature
- \* Patent #3,649,891, exclusively assigned to Lake Shore Cryotronics, Inc.


#### Temperature reproducibility

Over a period of days, thermal cycling of capacitance sensors can provide variations in their capacitance/ temperature values equivalent to several tenths of a degree at 4.2 K, 77 K, and room temperature. Over longer periods of time, variations can reach one degree or more. However, any reduced capacitance, C(T)/C(4.2 K), is generally stable to within ±0.5 K. These variations, or shifts, in the temperature response curve have no effect on the sensor's stability when operating at a given temperature and, therefore, do not impair the sensor's primary function as a control element.

Capacitance sensors (CS) are ideally suited for use as temperature control sensors in strong magnetic fields because they exhibit virtually no magnetic field dependence. Displacement current is not affected by magnetic fields. Consequently, temperature control fluctuations are kept to a minimum when sweeping magnetic field or when changing field values under constant temperature operation.

Because small variations in the capacitance/ temperature curves occur upon thermal cycling, calibrations must be transferred to the capacitor from another sensor after cooling for the best accuracy. It is recommended that temperature in zero field be measured with another temperature sensor and that the capacitance sensor be employed as a control element only.

#### **CS-501GR**

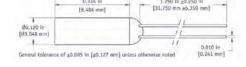




#### Temperature stability/ temperature transfer accuracy

Capacitance sensors will provide very stable control conditions for long periods of time at operating temperature, but because an operational "aging" phenomenon exists, care must be taken to account for this occurrence in their use.


The variation in capacitance/temperature characteristics is likely the result of the time dependence of the dielectric constant and the dielectric loss, or "aging", that all ferroelectric dielectrics exhibit. This time dependence, which occurs as a short term drift (minutes to hours) in capacitance/temperature value, is initiated by disturbing the sensor thermally or by changing the voltage or frequency of excitation. To compensate for this, the sensor should be stabilized for one hour after initial cool-down to desired operating temperature and whenever significant adjustments in control temperature are made.

After the one hour stabilization, this short-term drift is on the order of a few tenths of a millikelvin per minute at 4.2 K, and several millikelvin per minute at 305 K. The drift is always in the direction of decreasing capacitance; consequently, it corresponds to decreasing temperature below 290 K.


#### Typical CS capacitance

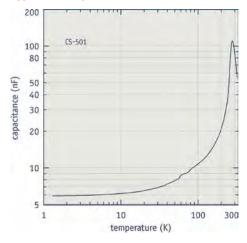

#### Typical CS sensitivity

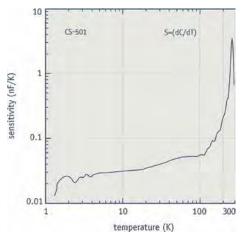

#### Typical CS dimensionless sensitivity

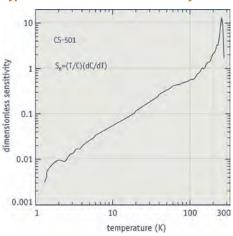



Standard curve Not applicable

Nominal capacitance 6.1 nF

Nominal sensitivity 26 pF/K

Accuracy (interchangeability) Not applicable

Accuracy (calibrated) Calibration should be performed in situ

**Recommended excitation** 1 to 5 kHz, 0 to 7 V peak to peak or any other acceptable capacitance measuring method

**Dissipation at recommended excitation** Not applicable

Expected long-term stability ±1.0 K/yr

Thermal response time Minutes, dominated by electronic setting time

Radiation effects Not available

Magnetic fields See table on right

**Reproducibility** See shaded box on previous page for detailed discussion

Soldering standard J-STD-001 Class 2

#### Range of use

	Minimum limit	Maximum limit
CS-501GR	1.4 K	290 K

# Typical magnetic field-dependent temperature errors¹ ΔT/T (%) at B (magnetic induction)

Package parallel to field B				
18.7 T				
4.2 K	-0.15			
77 K <0.05				

Recommended for control purposes; monotonic in C vs T to nearly room temperature; frequency dependent

#### **Packaging options**

For more information on sensor packages and mounting adapters, see page 21.





# See the appendices for a detailed description of:

Installation
Uncalibrated sensors
SoftCal™
Calibrated sensors
CalCurve™
Sensor packages



#### **Physical specifications**

	Size	Mass	Lead Type	Internal atmosphere
CS-501	3.0 mm × 8.5 mm long	260 mg	2 phosphor bronze with heavy build polyimide attached with epoxy strain relief at sensor	Air





## Thermocouple Wire

#### **Thermocouple features**

- Chromel-Au/Fe (0.07%) consists of a gold (Au)-0.07 at % iron (Fe) as the negative thermoelement and a Ni-Cr alloy (chromel) as the positive thermoelement. This thermocouple is more widely used because of its relatively high thermoelectric sensitivity (>15 µV/K above 10 K).
- Type E (chromel-constantan) has the highest sensitivity among the three standard thermocouple types typically used at low temperatures (types E, K, and T). The best choice for temperatures down to 40 K.
- Type K (chromel-alumel) Recommended for continuous use in inert atmospheres. Has a sensitivity of 4.1 mV/K at 20 K (about ½ of Type E).

Thermocouples are used in a variety of cryogenic applications, but special techniques must be employed to approach temperature accuracies of 1% of temperature, even without consideration for the effects of high magnetic fields or high radiation fluxes. The problems are further complicated by exposure to variable gradient conditions at cryogenic temperatures.

Many Lake Shore temperature controllers offer inputs that accommodate most common types of cryogenic thermocouples in use.

#### Range of use

	Minimum Iimit	Maximum limit
Chromel-AuFe (0.07%)	1.2 K	610 K
Type E	3.15 K	953 K
Type K	3.15 K	1543 K

Upper limit dependent on wire size; to achieve higher than 473 K. insulation must be removed

#### Note:

Heat conduction down the thermocouple wire is the same as with lead wire going to any other sensing device. Refer to Appendix C: Conduction (Lead Attachment) for more detailed information.

See Appendix G for thermocouple curve data



# Typical magnetic field-dependent temperature errors¹ ΔT/T (%) at B (magnetic induction)

Chromel-AuFe (0.07%) <sup>1</sup>					
	2.5 T	8 T	14 T		
4.2 K	_	5	_		
10 K	3	20	30		
45 K	1	5	7		
100 K	0.1	8.0	_		

Data taken with entire thermocouple in field, cold junction at 4.2 K, errors in hot junction temperature

	Type E thermocouple <sup>2</sup>									
	2.5 T	8 T	14 T							
10 K	1	3	7							
20 K	<1	2	4							
45 K	<1	<1	2							

 $<sup>^2</sup>$  Useful when T  $\geq$  10 K. Refer to comments for chromel-AuFe (0.07%)

#### Part number Explanation

**TC** = Thermocouple

 $\mathbf{Y}$  = Wire type, E or K

**ZZ** = Wire diameter excluding insulation

**XX** = Wire length in meters

	Wire (	gauge
	30 AWG	36 AWG
Type E	TC-E-30-XX	TC-E-36-XX
Type K	TC-K-30-XX	TC-K-36-XX

#### Ordering information

Thermocouple wire	36 AWG = 0.005 in (0.127 mm) diameter wire, excluding insulation 30 AWG = 0.010 in (0.254 mm) diameter wire, excluding insulation All thermocouple wire is Teflon® insulated—76.2 μm wall
Part number	Description
C-Y-ZZ-03	Thermocouple wire — 3 m
C-Y-ZZ-06	Thermocouple wire — 6 m
TC-Y-ZZ-10	Thermocouple wire — 10 m
TC-Y-ZZ-20	Thermocouple wire — 20 m
TC-Y-ZZ-50	Thermocouple wire — 50 m




# Cryogenic Hall Sensors

#### **Cryogenic Hall sensor features**

- Low temperature dependence
- Low resistance, low power dissipation
- Low linearity error:-150 kG to +150 kG
- Axial and transverse configurations available
- Small active area

# Attaching Hall sensors to the Model 425, 455, 460, and 475 gaussmeters

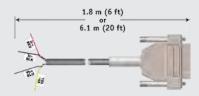

The MCBL-6 cable allows discrete Hall sensors to be mated to the Model 425, 455, 460, and 475 gaussmeters. The cable is shipped with a CD-ROM containing the **Hallcal.exe** file to program the cable PROM through the gaussmeter RS-232C port. Because of the many intricacies involved with proper calibration, the user is responsible for the measurement accuracy.

# Certain Hall sensor sensitivity constraints are applicable:

Sensitivities between 5.6 mV/kG and 10.4 mV/kG at 100 mA current.

Sensitivities between 0.56 mV/kG and 1.04 mV/kG at 100 mA current.


#### **System requirements**

- Lake Shore gaussmeter connected via RS-232C to the PC
- Hall sensor meeting the sensitivities given above
- Calibration or sensitivity constant and serial number of the Hall sensor

#### Hall sensor theory

A Hall sensor is a solid state sensor which provides an output voltage proportional to magnetic flux density. As implied by its name, this device relies on the Hall effect principle. The Hall effect principle is the development of a voltage across a sheet of conductor when current is flowing and the conductor is placed in a magnetic field.

Electrons (the majority carrier most often used in practice) "drift" in the conductor when under the influence of an external driving electric field. When exposed to a magnetic field, these moving charged

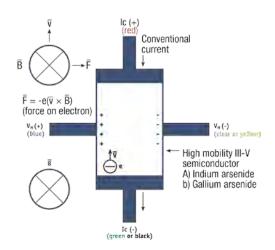

particles experience a force perpendicular to both the velocity and magnetic field vectors. This force causes the charging of the edges of the conductor, one side positive with respect to the other. This edge charging sets up an electric field which exerts a force on the moving electrons equal and opposite to that caused by the magnetic-field-related Lorentz force. The voltage potential across the width of the conductor is called the Hall voltage. This Hall voltage can be utilized in practice by attaching two electrical contacts to the sides of the conductor.

The Hall voltage can be given by the expression:  $V_{\mu} = \gamma_{R} B \sin \phi$ 

where  $V_H$  = Hall voltage (mV)

 $\gamma_{\rm B}$  = Magnetic sensitivity (mV/kG) at a fixed current

B = Magnetic field flux density (kG)

Φ = Angle between magnetic flux vector and the plane of Hall sensor

As can be seen from the formula above, the Hall voltage varies with the angle of the sensed magnetic field, reaching a maximum when the field is perpendicular to the plane of the Hall sensor.


#### **Using a Hall sensor**

A Hall sensor is a 4-lead device. The control current ( $I_c$ ) leads are normally attached to a current source such as the Lake Shore Model 121. The Model 121 provides several fixed current values compatible with various Hall sensors.

The Hall voltage leads may be connected directly to a readout instrument, such as a high impedance voltmeter, or can be attached to electronic circuitry for amplification or conditioning. Device signal levels will be in the range of microvolts to hundreds of millivolts.

The Hall sensor input is not isolated from its output. In fact, impedance levels on the order of the input resistance are all that generally exist between the two ports. To prevent erroneous current paths which can cause large error voltages the current supply must be isolated from the output display or the downstream electronics.

#### A typical Hall effect measurement scheme

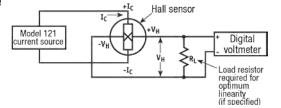



#### **Configurations**


Hall sensors come in two main configurations, axial and transverse. **Transverse** devices are generally thin and rectangular in shape. They are applied successfully in magnetic circuit gaps, surface measurements, and general open field measurements. **Axial** sensors are mostly cylindrical in shape. Their applications include ring magnet center bore measurements, solenoids, surface field detection, and general field sensing.

#### **Active area**


The Hall sensor assembly contains the sheet of semiconductor material to which the four contacts are made. This is normally called a "Hall plate." The Hall plate is, in its simplest form, a rectangular shape of fixed length, width, and thickness. Due to the shorting effect of the current supply contacts, most of the sensitivity to magnetic fields is contained in an area approximated by a circle, centered in the Hall plate, with a diameter equal to the plate width. Thus, when the active area is given, the circle as described above is the common estimation.

Specifications	HGCA-3020	HGCI-3020
Description	Cryogenic axial; phenolic package	Cryogenic transverse; ceramic package
Active area (approximate)	0.030 in (0.762 mm) diameter	0.040 in (1.016 mm) diameter
Input resistance (approximate)	1 Ω	1 Ω
Output resistance (approximate)	1 Ω	1 Ω
Nominal control current (I <sub>CN</sub> )	100 mA	100 mA
Maximum continuous current (non-heat sunk)	300 mA	300 mA
Magnetic sensitivity at I <sub>CN</sub>	0.55 mV/kG to 1.05 mV/kG	0.55 mV/kG to 1.05 mV/kG
Magnetic sensitivity change with temperature	+0.7% at 200 K; +0.8% at 100 K; +1.0% at 3.8 K	+0.7% at 200 K; +0.8% at 100 K; +1.0% at 3.8 k
Maximum linearity error (sensitivity versus field)	±1.0% RDG (-30 kG to +30 kG) ±2.0% RDG (-150 kG to +150 kG)	$\pm 1.0\%$ RDG (–30 kG to +30 kG) $\pm 2.0\%$ RDG (–150 kG to +150 kG)
Zero field offset voltage (maximum)(I <sub>c</sub> = nominal control current)	±200 μV	±200 μV
Operating temperature range	1.5 K to 375 K	1.5 K to 375 K
Mean temperature coefficient of magnetic sensitivity (approximate)	±0.01%/K	±0.01%/K
Mean temperature coefficient of offset (maximum) ( $I_{\text{c}} = \text{nominal control current})$	±0.4 μV/K	±0.4 μV/K
Mean temperature coefficient of resistance (maximum)	±0.6%/K	±0.6%/K
Leads	34 AWG copper with Teflon® insulation	34 AWG copper with Teflon® insulation
Data	Room temp; 30 kG data supplied	Room temp; 30 kG data supplied

#### **HGCA-3020**

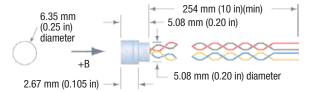

#### **HGCT-3020**

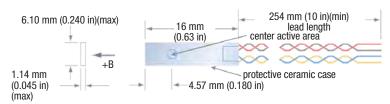

#### **Gaussmeters**

Lake Shore gaussmeters offer a straightforward and cost effective solution to measure magnetic fields. Hall sensors or factory calibrated probes connect to the gaussmeter rear panel and the sensor data is automatically uploaded into the instrument. Lake Shore gaussmeters offer easy-to-make flux density measurements with high accuracy, resolution, and stability, and are available with RS-232C and IEEE interfaces, analog outputs, relays, and alarms.



For more information call, or visit www.lakeshore.com.



# Cryogenic Hall Probes

Lake Shore offers cryogenic Hall sensors mounted into gaussmeter probes, which work in a variety of magnetic measurement applications. Our probes are factory calibrated for accuracy and interchangeability. Pre-calibrated probes feature a programmable read-only memory (PROM) in the probe connector so that Hall sensor calibration data can be read automatically by a Lake Shore gaussmeter.

Lake Shore also offers a complete line of axial, transverse, flexible, tangential, gamma, brass stem, and multi-axis Hall probes. For more information call us, or visit www.lakeshore.com.

#### **Axial probes for cryogenic applications**



	L mm (in)	D mm (in)	A mm (in)	Active area mm (in)	Stem material	Frequency range	Usable full scale ranges	Corrected accuracy (% rdg at 25 °C)	Operating temp range	Temp coefficient (max) zero	(app	error prox) ration	Contains temp sensor
for Models 4	75, 455, a	nd 425									300 K	ref	
HMCA-2560-WN	1524 ±12.7 (60 ±0.50)	6.35 dia ±0.15 (0.25 dia ±0.006)	0.64 ±0.13 (0.025 ±0.005)	Approx 0.76 dia (0.030 dia)	Stainless steel	DC	HST-3 35 G, 350 G, 3.5 kG, 35 kG, 350 kG	±2% to 100 kG	1.5 K to 350 K	±0.13 G/°C		+0.05% -0.04% -0.09%	No
for Models 4	60, 450, a	nd 421									20 K	-0.40%	
MCA-2560-WN	1524 ±12.7 (60 ±0.50)	6.35 dia ±0.15 (0.25 dia ±0.006)	0.64 ±0.13 (0.025 ±0.005)	Approx 0.76 dia (0.030 dia)	Stainless steel	DC	HST-1 300 G, 3 kG, 30 kG, 300 kG	±2% to 100 kG	1.5 K to 350 K	±0.13 G/°C	4 K 1.5 K	-0.70% -1.05%	No

#### Transverse probes for cryogenic applications



	L mm (in)	D mm (in)	A mm (in)	Active area mm (in)	Stem material	Frequency range	Usable full scale ranges	Corrected accuracy (% rdg at 25 °C)	Operating temp range	Temp coefficient (max) zero	(ap	error orox) ration	Contains temp sensor
for Models 4	or Models 475, 455, and 425											ref	
	1549	6.35 dia	5.33 ±1.27	Approx 1			HST-3				200 K	+0.05%	
HMCT-3160-WN	±25.4 (61	±0.25 (0.25	0.25 (0.25 (0.210	dia (0.040	Stainless steel	DC to 800 Hz	35 G, 350 G, 3.5 kG,	±2% to 100 kG	1.5 K to 350 K	±0.13 G/°C	100 K	-0.04%	No
	±1)	dia ±0.010)	±0.050)	dia)	01001	000112	35 kG, 350 kG	, Ku	00010		80 K	-0.09%	
for Models 4	<del>1</del> 60, 450,	and 421									20 K	-0.40%	
MOT O4 CO MAI	1549	6.35 dia	5.33 ±1.27	Approx 1	Stainless	DC and	HST-1	±2% to	1.5 K to	0.40.0/00	4 K	-0.70%	
MCT-3160-WN	±25.4 (61 ±1)	±0.25 (0.25 dia ±0.010)	(0.210 ±0.050)	dia (0.040 dia)	steel	10 Hz to 400 Hz	300 G, 3 kG, 30 kG, 300 kG	100 kG	350 K	±0.13 G/°C	1.5 K	-1.05%	No

#### Ordering information Other Hall sensors and probes are available — please consult Lake Shore Part number Description **Cryogenic Hall sensors** Cryogenic axial Hall sensor HGCA-3020 Cryogenic transverse Hall sensor HGCT-3020 Cryogenic Hall probes HMCA-2560-WN Cryogenic axial gaussmeter probe for Models 475, 455, and 425 MCA-2560-WN Cryogenic axial gaussmeter probe for Models 460, 450, and 421 HMCT-3160-WN Cryogenic transverse gaussmeter probe for Models 475, 455, and 425 MCT-3160-WN Cryogenic transverse gaussmeter probe for Models 460, 450, and 421






# Instruments

- 64 Instrument Selection Guide
- 69 Model 372 AC Resistance Bridge
- **83 Model 350 Temperature Controller**
- 91 Model 336 Temperature Controller
- 99 Model 335 Temperature Controller
- **106 Model 325 Temperature Controller**
- 112 Model 224 Temperature Monitor
- 118 Model 218 Temperature Monitor
- 123 Model 211 Temperature Monitor
- 127 230 Series Temperature Transmitters/Monitors
- 132 Model 121 DC Current Source
- 136 Model 625 Superconducting Magnet Power Supply





#### Instrument Selection Guide

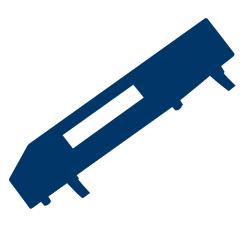

# How to select a temperature instrument for your application

Lake Shore offers the most comprehensive line of cryogenic temperature instruments in the world. The instruments described in this section are designed and manufactured for both general and specific temperature research applications in mind. For much of its 35-year history, Lake Shore has focused on instrumentation used for the precise measurement of temperatures from near absolute zero to well above room temperature.

Unfortunately, you can't have it all in one instrument. The most precise and accurate temperature instruments optimized for operation below 100 mK work with fewer sensors and provide lower heater power. The stable and high-resolution instruments designed for general cryogenic use work well for nearly any application, but can have limitations in rare circumstances. Choosing the appropriate instrument for a particular application necessitates prioritizing the requirements for that application.

Any one or several of the following factors may be important to you in selecting an instrument, whether temperature control or temperature monitoring is required:

- Operating temperature range
- Number of sensor inputs
- Sensor type compatibility
- Sensor input resistance and voltage ranges
- Current excitation ranges and methods
- High measurement resolution
- High electronic accuracy
- Control stability
- Number of reading displays
- Interfaces

Ethernet

USB

**IEEE-488** 

RS-232C

Alarms

Relays

Analog outputs


Digital I/O

Data logging


- Number of control loops, control type, and operating parameters
- Heater power and ranges
- Low cost

The tables on the following pages are designed to compare the instruments more easily with regard to sensor compatibility, operating temperature range, control capability, display features, and interface flexibility.

Our experienced sales staff is here to answer your questions. If you already know what your needs are, please inform us. Otherwise we ask a lot of questions to inform, educate, and to assist you in selecting the correct instrument.

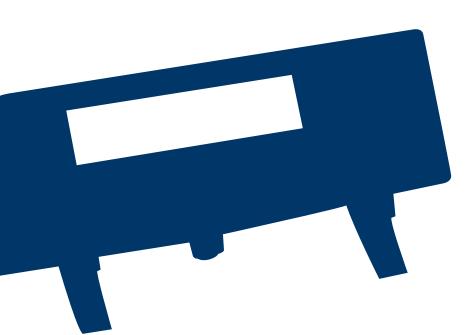

#### **Temperature controller temperature ranges**

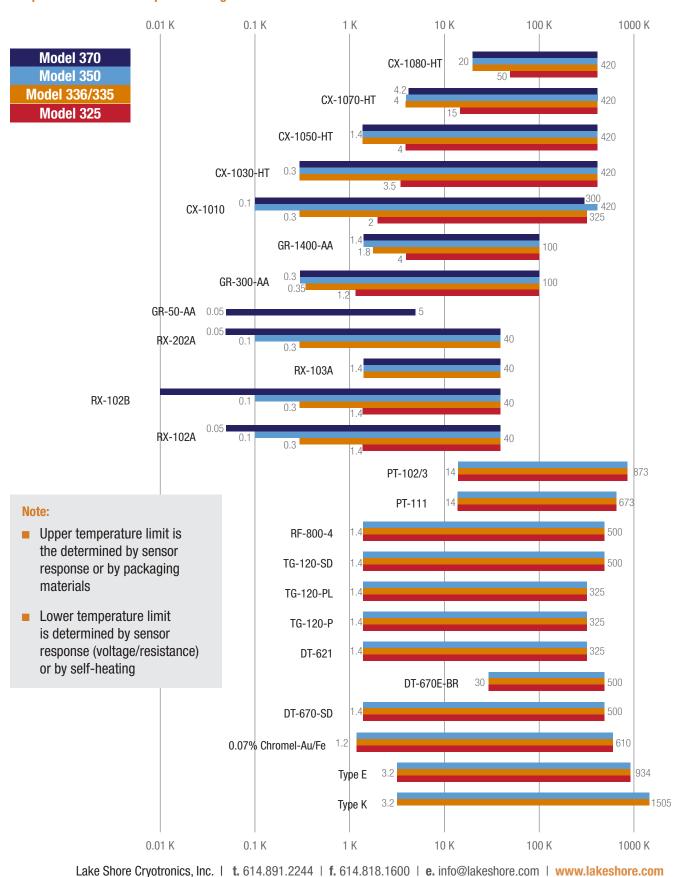




#### **Current excitation ranges**

	3.16 pA	10 pA	31.6 pA	100 pA	316 pA	1.0 nA	3.16 nA	10 nA	31.6 nA	100 nA	316 nA	1 μΑ	3.16 μΑ	10 μΑ	31.6 μΑ	100 μΑ	316 µА	1 mA	3.16 mA	10 mA	31.6 mA
370																					
350																					
336																					
335																					
325																					

	AC bridge		Contr	ollers	
	370	350	336	335	325
Number of sensor inputs	1 to 16¹	8 <sup>1</sup>	4	2	2
Number of user curves	20	39	39	39	15
Minimum operating temperature	<20 mK	100 mK	300 mK	300 mK	1.2 K
Maximum operating temperature	420 K	1505 K	1505 K	1505 K	1505 K
Current reversal	Yes	Yes	Yes	Yes	Yes
Current excitation autoranging	Yes	Yes	Yes	Yes	_
Number of reading displays	1 to 8	1 to 8	1 to 8	1 to 4	1 to 4
Interfaces					
Ethernet	_	Yes	Yes	_	_
USB	_	Yes	Yes	Yes	_
IEEE-488.2	Yes	Yes	Yes	Yes	Yes
RS-232C	Yes	_	_	_	Yes
Number of alarms	32	8	4	2	_
Number of relays	2	2	2	2	_
Analog voltage output	2 at ±10 V	2 at ±10 V	2 at ±10 V	1 at ±10 V	0 to 10 V
Number of autotuning control loops	1	1	2	2	2
Maximum heater output power					
Control loop 1	1 W	75 W	100 W	75 W <sup>2</sup>	25 W
Control loop 2	_	50 W	50 W	25 W	2 W
Number of heater ranges	8	5	3	3	2

Optional input card or scanner 75 W only available when output 2 is in voltage mode; maximum in other modes 50 W

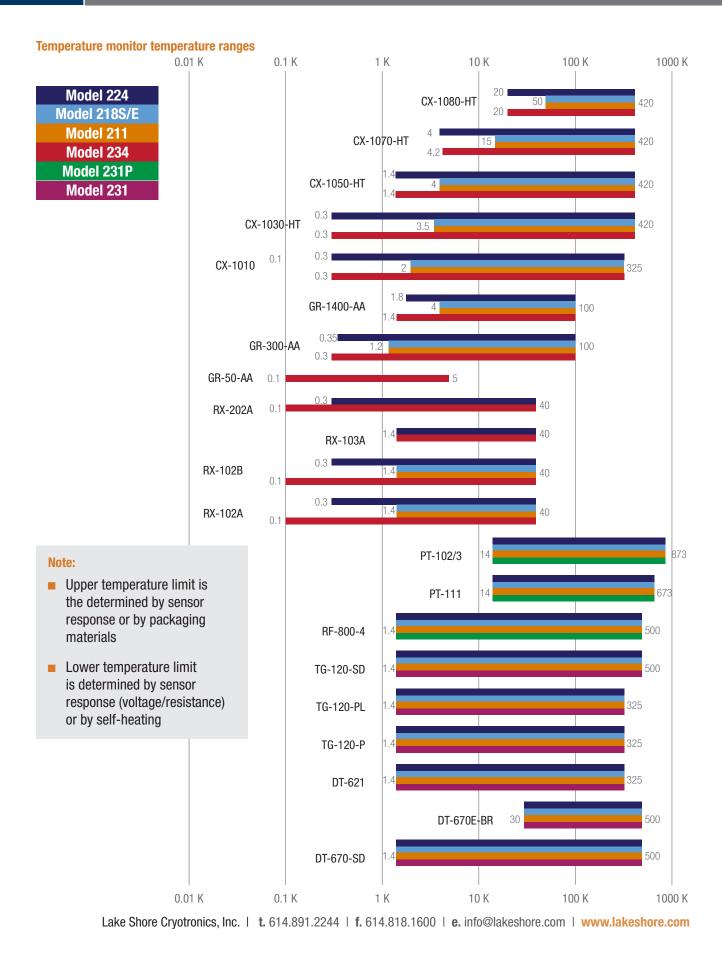



#### **Current excitation ranges**

	31.6 nA	100 nA	316 nA	1 μΑ	3.16 μΑ	10 μΑ	31.6 µА	100 μΑ	316 μΑ	500 μΑ	1 mA
224											
218S/E											
211											
234											
231P											
231											

		Moni	tors			Transmitters	
	224	218\$	218E	211	234	231P	231
Number of sensor inputs	12	8	8	1	1	1	1
Number of user curves	39	8	8	1	1	1	1
Minimum operating temperature	0.3 K	1.2 K	1.2 K	1.2 K	100 mK	1.4 K	1.4 K
Maximum operating temperature	873 K	800 K	800 K	800 K	420 K	800 K	500 K
Current reversal	Yes	_	_	_	Yes	_	_
Current excitation autoranging	Yes	<del>_</del>	_	_	Yes	<del>_</del>	_
Number of reading displays	1 to 16	1 to 8	1 to 8	1	1	<del>-</del>	_
Interfaces							
IEEE-488.2	Yes	Yes	_	_	_	_	_
USB	Yes	_	_	_	_	_	_
Ethernet	Yes	_	_	_	_	_	_
RS-232C	_	Yes	Yes	Yes	Yes	_	_
Number of alarms	12	16	16	2	_	_	_
Number of relays	2	8	_	2	_	_	_
Analog voltage output	_	2 at ±10 V	_	0 – 10 V	0 – 10 V	0 – 10 V	0 – 10 V
4 – 20 mA output	_	_	_	Yes	Yes	Yes	Yes
Data logging	Yes	Yes	Yes	_	_	_	_

<sup>&</sup>lt;sup>3</sup> Uses 5 mV or 10 mV constant voltage



# Model 372 AC Resistance Bridge

and Temperature Controller



- Patented noise rejection technology
- Highly versatile and reliable measurement input
- Ability to increase the number of measurement channels to a maximum of 16 with optional 3726 scanner
- Dedicated input for ultra-low temperature control
- Powerful impedance measurement capabilities such as quadrature measurements
- Multiple PID controllable outputs with up to 10 W of heater power available
- Latest generation front panel for ease of use
- 3-year standard warranty





#### Introduction

The Model 372 AC resistance bridge and temperature controller builds on the solid foundation provided by the original Lake Shore AC resistance bridge. The Model 372 provides the best possible temperature measurement and control capabilities for dilution refrigerators (DRs) that are intended to be operated below 100 mK. The Model 372 makes it easy to perform multiple tasks that were once very difficult to perform reliably at ultra-low temperatures:

- Temperature measurement
- Automatic or manual temperature control
- Device or sample impedance measurements

#### Targeted applications

#### **Ultra-low temperature measurement**

Making measurements below 100 mK is far from a trivial exercise, with even the smallest amounts of added energy leading to self-heating and unwanted temperature shifts. Every design decision made on the Model 372 aims to minimize the amount of energy needed to take measurements.

#### **Noise rejection**

Externally generated electronic noise can be a major cause of self-heating if it is allowed to couple into the device under test. Thankfully, multiple noise-rejection strategies have been implemented to reduce this effect substantially:

- Our patented\* balanced noise-rejecting current source ensures that external signals have no path to ground through the measurement circuit, effectively making the Model 372 unaltered by these noise sources.
- The measurement signal cables use a driven guard that reduces parasitic capacitance in the cables that connect a scanner to the Model 372. This helps to further balance the measurement network and bolster the integrity of the noise rejection circuitry.
- All measurement circuitry is isolated from other instrument components, limiting the impact of any small electrical disturbances.
- The AC frequency options used for the measurement signal are selected to be naturally resilient to line voltage frequencies (50 and 60 Hz).

#### **AC** measurement signals

By using alternating current (AC) measurement in tandem with a specially designed internal lock-in amplifier, the Model 372 is able to extract very small measurement signals from background noise. This allows for much lower excitation levels to be used when compared to traditional direct current (DC) systems, minimizing the amount of energy that is dissipated into the device under test.


These AC excitation levels can be set to as low as 10 pA, while still maintaining accuracy of better than 1% over quite a wide range of resistances. This enables impedance and temperature measurements to be made while adding power levels so small that they are measured in the attowatt range (10<sup>-18</sup> W). These features are vital in allowing accurate measurement to be made while minimizing the negative effects of self-heating.

#### Low noise signal recovery

Due to the very low excitation level used for measurement, the resulting voltage levels must first be boosted to allow those signals to be measured. The internal lock-in amplifier in the Model 372 has been specifically designed to minimize the amount of noise added to the signal. This results in an input noise figure that is less than 10 nV/\JHz, thereby increasing the resolution of measurements and limiting the amount of post-measurement filtering that needs to be applied.

#### **Temperature measurement**

Extremely accurate and reliable ultra-low temperature measurements can be achieved by combining the Model 372 with a negative temperature coefficient (NTC) resistive temperature device (RTD), such as the Lake Shore Cernox®, Rox<sup>TM</sup> or germanium temperature sensors. Multiple calibration curves can easily be uploaded to the Model 372, allowing highly accurate conversion of sensor resistance to equivalent temperature using cubic spline interpolation (an improved interpolation technique compared to older instruments).

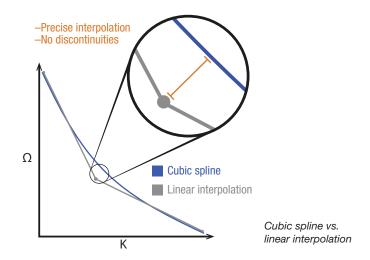

### Model 372 rear panel



Model 372 AC Resistance Bridge

- 1 Scanner control input (DA-15)
- 2 Sensor voltage/current input (6-pin DIN)
- 3 Secondary control input (6-pin DIN)
- 4 Monitor output (BNC)

- 5 Reference output (BNC)
- 6 Sample heater output, warm-up heater output, and still heater output (terminal block)
- 7 Relay 1 and 2 (terminal block)

- B Ethernet interface (RJ-45)
- 9 USB interface (USB Type A)
- 10 IEEE-488.2 interface
- 11 Line power/fuse assembly

User-generated calibration curves can also be created and loaded into the Model 372, allowing great flexibility in the type of resistive sensors that are used. A maximum of 39 calibration curves can be stored on the instrument, and when used with a 3726 scanner, up to 17 sensors can be connected simultaneously, each with their own curve.

#### Measure a wide range of resistive devices

With up to 22 different current (I) excitation levels available, the Model 372 is able to perform accurate impedance measurements from several microohms ( $10^{-6}~\Omega$ ) to many megohms ( $10^{6}~\Omega$ ), all while keeping power dissipation levels to an absolute minimum.

The addition of full quadrature measurements means that both the resistive and reactive components of an impedance can now be measured. This enables much better characterization of the device under test by allowing capacitive or inductive components to be measured.




The new 3726 scanner option

#### **Expandability**

For situations where temperature measurements must be taken at multiple locations, the 3726 scanner and preamp can be paired with the Model 372 to provide up to 16 connections for 4-wire resistance measurements. The Model 372 can switch measurement to any one of these connections as required, removing the need to physically switch cables on the instrument to look at different sensors. The measurement signal is also boosted by a pre-amp circuit in the 3726, preserving the signal-to-noise ratio between the sensor and measurement circuitry of the Model 372. This allows connection cables of up to 10 m to be used between the Model 372 and the 3726.

In cases where measurements are required at multiple locations simultaneously within an experiment space, additional Model 372 units may be used together. Five different AC excitation frequencies are available for this purpose, ensuring that up to five simultaneous measurements can be performed without the risk of co-channel interference.

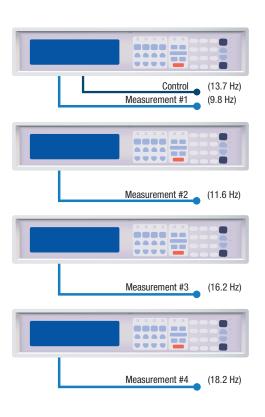



### Dilution Refrigerator Temperature Control

Accurate measurement at ultra-low temperatures are no easy feat, especially when working in the ranges seen by modern dilution refrigerators. The Model 372 has many features specifically developed for dilution refrigerator applications.

#### **Dedicated temperature control input**

Taking measurements at ultra-low temperatures deserves uninterrupted attention from measurement devices. The Model 372 uses a dedicated temperature control input that is designed specifically for connection to a negative temperature coefficient resistive sensor. This input is designed to continuously monitor the temperature of the dilution refrigerator sample holder, while the measurement input scans through the multiple other temperature sensors placed throughout the dilution refrigerator.




The dedicated control input ensures uninterrupted dilution refrigerator temperature control


#### **Multiple heater options**

Three separate heater outputs are available on the Model 372:

- Sample heater—for fine control of the sample stage at ultra-low temperatures with up to 1 W of power available.
- Warm-up heater—supplying up to 10 W of power and featuring a warm-up mode specifically for the purpose of bringing the system temperature up to allow work to be performed on the sample stage.
- Still heater—an additional 1 W heater is available for the purpose of controlling the temperature of a dilution refrigerator's still. Alternatively, this output can provide an analog out signal to other devices if required.

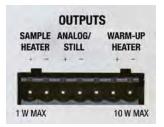

The sample and warm-up heaters have many powerful control options, including PID control that allows both the setting of fixed temperature setpoints as well as ramp rates.

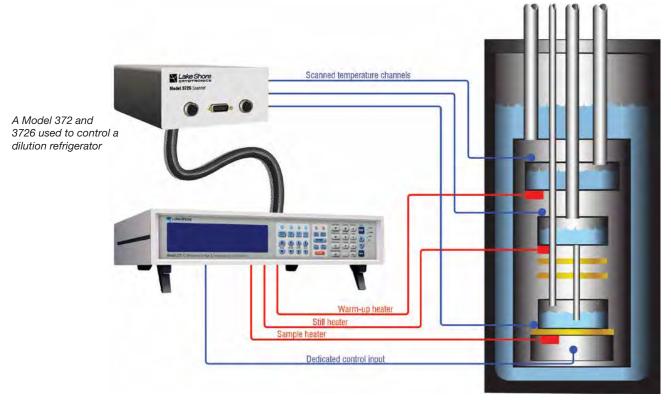



#### Stable temperature control

When operating at ultra-low temperatures, even small amounts of added energy can cause unwanted spikes in system temperature. The Model 372 heater outputs implement several protection mechanisms to reduce or eliminate this potential:

- The circuitry for the sample and still heaters are electrically isolated from other instrument sections
- Multiple power range settings allow extremely fine or coarse power transitions, depending on the need
- Heater outputs are shunted during power up and power range changes, eliminating the potential for unwanted power surges
- Terminal connections allow twisted pair cabling to be easily used for heater wiring; additional shielding of these wires can also be added to further reduce the potential of injecting noise into a system via the heater cabling

#### **Temperature zone control**

Thermal response characteristics of a dilution refrigerator system can change quite dramatically over the useful range of operation, particularly down towards the lower temperature limit of a system, where cooling power is reduced. To accommodate these system variations, different PID values can be set for different temperature ranges (zones). This allows for more aggressive transition settings to be used at higher temperatures where system response is faster, and less reactive settings at low temperatures when temperature overshoots result in long recovery times.

#### **Heater fail-safes**

The Model 372 has several features that will protect your system and experiment from accidental deviations in planned temperature settings:

- Temperature thresholds can be set for all heater outputs, meaning the heaters will automatically shut down if it is detected that the system is being overheated.
- An easy-to-hit "ALL OFF" button is provided that shuts all heaters down instantly. This eliminates the terrible experience of having to hurriedly search through menu options while your experiment continues to heat.



# Low-Power Impedance Characterization—the 3708 Scanner

Model 372 AC Resistance Bridge



Many material characterization experiments require measurements to be performed at cryogenic temperatures. This can be because the material behavior changes in interesting ways at these temperatures, or because background thermal noise must be minimized for useful measurement data to be extracted. The standard inputs of the Model 372 accurately measure higher-impedance devices such as temperature sensors, but begin to lose resolution and accuracy when extremely low impedances are encountered such as in Hall effect or superconducting material measurements. However, by adding a 3708 preamp and scanner to the Model 372, these materials can be characterized with the same accuracy and stability as when measuring higher-impedance devices.

To accomplish this, the 3708 produces higher levels of DC bias current than both the Model 372 and the 3726 scanner and preamp. This means the 3708 would cause self-heating in a temperature sensor used at ultra-low temperatures. The new dedicated control input resolves this issue by providing the ability to make highly reliable measurements of a temperature control sensor.

#### Lower input voltage noise

The limiting factor for making extremely low-impedance measurements directly with the Model 372 is the input voltage noise figure of 10 nV/\/Hz. The preamp in the 3708 reduces this by a factor of 5 to an impressive 2 nV/\/Hz. By reducing the amount of input noise, even smaller return signals can be recovered with excellent accuracy. When combined with the ability of the Model 372 to smooth measurement values with user-settable filters ranging from 1 to 200 s, the 3708 preamp and scanner provides the best solution to measuring low-impedance devices at cryogenic temperatures.



#### **Multiple simultaneous connections**

The 3708 scanner and preamp allows up to eight simultaneous connections to be made, with the scanner feature enabling measurement to be switched between those connections. Unlike the 3726 scanner, all connections that are not actively being measured are left open, allowing the 3708 to be connected to Hall bar devices.

#### **Overcoming cable length**

With such small resultant voltages needing to be measured, it can be very helpful to have these signals amplified slightly as close as possible to the source of these signals. The compact size of the 3708 scanner and preamp allows it be mounted close to the device or sample being measured, thereby maintaining signal-to-noise ratio for the measurement signal between the sample and the Model 372 that will ultimately perform the measurements. Cable lengths of up to 10 m are supported by the 3708, allowing the Model 372 to be located away from the experiment area if needed.


### Connectivity and Usability—Communication Options

#### **Physical connectivity**

Various methods for communicating with the Model 372 are made available:

- Ethernet: allows full control and reporting throughout an IP network.
- USB: provides direct serial communication by emulating a standard RS-232 connection.
- IEE-488.2: allows connection to GPIB systems.



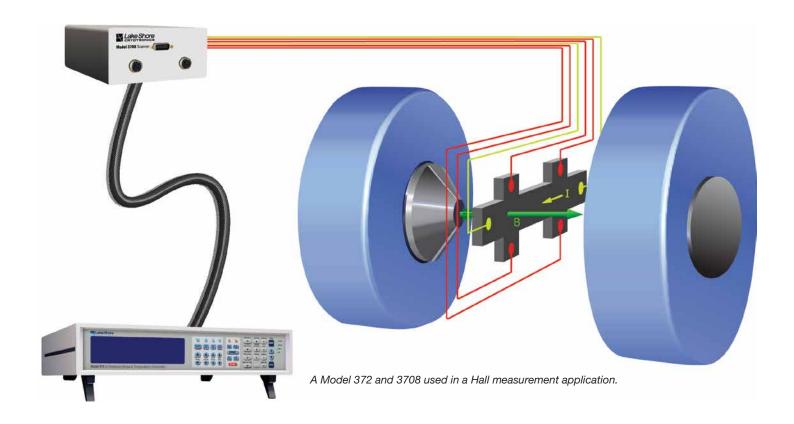



#### **Available functions**

Multiple actions can be performed when connected to the Model 372 through one of its various remote access options:

- Send any command to the instrument that could be entered via the front panel
- Read and store measurement data that is generated by the instrument
- Live graphical viewing of data using the Lake Shore Cryotronics Chart Recorder software
- Load new calibration curves for use with new temperature sensors
- Upload new firmware if required

#### **Backwards compatibility**

The Model 372 is designed for trouble-free integration with existing equipment and software that has been built around the previous generation Model 370.


- Emulation mode on the Model 372 is designed to imitate all important communication functions of the Model 370. In most cases, programming that was previously written for the Model 370 can be used to interact with the Model 372.
- A convenient heater connector adapter (117-071) can also be purchased. This adapter replicates the BNC heater connections that were available on the Model 370, allowing connection swapping between the Model 372 and Model 370 without the need to rewire experiment cabling.

#### The Rox™ RX-102B-CB

Model 372 AC Resistance Bridge

The RX-102B-CB (1000  $\Omega$  at room temperature) is useful down to 10 mK (calibrations available down to 20 mK) and monotonic from 10 mK to 300 K. The unique package design maximizes thermal connection and minimizes heat capacity at ultra low temperatures. The RX-102B-CB is not interchangeable to a standard curve and not recommended for use in magnetic fields.



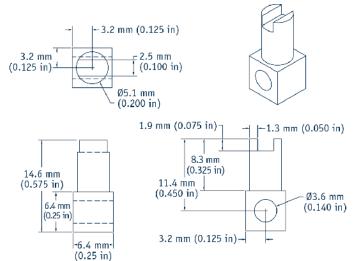

General tolerance of ±0.127 mm (0.005 in) on X.XXX and ±0.254 mm (0.01 in) on X.XX unless otherwise noted Mount using a #6 or M3 screw



#### Sensor performance

Excitation ranges in sensor tables were selected to minimize sensor self-heating. Excitation power = actual current<sup>2</sup> × example resistance

**Measurement resolution** comes from electronic instrumentation and sensor thermal noises. Measurement resolution is given by:  $Resolution (\Omega) = ((instrument \ noise)^2 + (sensor \ thermal \ noise)^2)^{0.5}$ 

Electronic instrumentation noise is taken at ambient temperature, while sensor thermal noise is taken at the temperature specified in the following tables.

Resolution (K) = 
$$\frac{Resolution(\Omega)}{(dR/dT)}$$

**Electronic accuracy** is influenced by the measurement range used and sensor resistance value. Electronic accuracy is given by:  $Electronic accuracy(\Omega) = Accuracy(\%) \times example resistance + 0.005\% of resistance range$  Where: Accuracy (%) is given in the instrument performance table (pages 10–11) at the selected current and voltage range

Electronic accuracy (K) = 
$$\frac{Electronic \ accuracy(\Omega)}{(dR/dT)}$$

**Self-heating errors** are measurement errors due to power dissipation in the sensor causing unwanted temperature rises.

Self-heating error is given by:

Self heating error=thermal resistance×power

Thermal resistances specified are typical values resulting from minimal heat sinking.

Improved values can be achieved with permanent installation.

Calibration accuracies are based on Lake Shore sensor calibration uncertainty and repeatability values—see Appendices B, D & E of the Temperature Measurement and Control Catalog for more information.

**Interpolation errors** are due to the linear interpolation method used by the Model 372 to convert resistance values to temperatures when using a temperature sensor. These errors are not present when resistance is measured directly.

Overall accuracy is a combination of all listed sources of potential error and is given by:

Overall accuracy = (measurement resolution<sup>2</sup> + electronic accuracy<sup>2</sup> + self heating errors<sup>2</sup> + calibration accuracy<sup>2</sup> + interpolation error<sup>2</sup>)<sup>0.5</sup>

#### Lake Shore Rox™ RX-102B-CB with 0.02 to 40 K calibration

Values given are for measurement input. If the value is different for the control input, it is shown in blue.

	Sensor p	roperties		Excitat	tion and ir	strument	ation	Instrument pe	erformance		Overall pe	rformance	
Temperature	Nominal resistance	Typical sensor sensitivity	Thermal resistance	Resistance range	Excitation voltage limit	Excitation current	Power	Measurement resolution	Electronic accuracy	Calibration accuracy	Self-heating errors	Interpolation error	Overall accuracy
20 mK	7.3 kΩ	-171 kΩ/K	17.2 K/nW	20 kΩ 632 kΩ	6.32 μV 200 μV	316 pA	730 aW	7.3 Ω (42.7 μK) 33.9 Ω (198 μK)	8.3 Ω (48.5 μK) 35.3 Ω (206 μK)	±2 mK	12.6 µK	±0.2 mK	2 mK
30 mK	6.0 kΩ	-100 kΩ/K	8.2 K/nW	6.32 kΩ 200 kΩ	6.32 μV 200 μV	1 nA	6 fW	485 m $\Omega$ (4.9 $\mu$ K) 7.3 $\Omega$ (73 $\mu$ K)	6.3 Ω (63 μK) 13.0 Ω (130 μK)	±4 mK	49.2 μK	±0.2 mK	4 mK
40 mK	5.2 kΩ	-62 kΩ/K	635.8 mK/nW	6.32 kΩ 63.2 kΩ	20 μV 200 μV	3.16 nA	52 fW	502 mΩ (8.1 μK) 1.5 Ω (24.2 μK)	2.9 Ω (46.8 μK) 4.7 Ω (75.8 μK)	±4 mK	33.1 µK	±0.2 mK	4 mK
50 mK	4.7 kΩ	-41 kΩ/K	415.1 mK/nW	6.32 kΩ 63.2 kΩ	20 μV 200 μV	3.16 nA	47 fW	502 mΩ (12.2 μK) 1.5 Ω (36.6 μK)	2.7 Ω (65.9 μK) 4.6 Ω (112 μK)	±4 mK	19.5 µK	±0.2 mK	4 mK
100 mK	3.5 kΩ	-13 kΩ/K	33.2 mK/nW	6.32 kΩ 20 kΩ	63.2 μV 200 μV	10 nA	350 fW	48.6 mΩ (3.7 μK) 338 mΩ (26 μK)	2.1 Ω (162 μK)	±4 mK	11.6 µK	±0.2 mK	4 mK
300 mK	2.5 kΩ	-2.4 kΩ/K	2.8 mK/nW	6.32 kΩ	200 μV	31.6 nA	2.5 pW	50.2 mΩ (20.9 μK) 87 mΩ (36.3 μK)	1.1 Ω (458 μΚ)	±4 mK	7.0 µK	±0.2 mK	4 mK
1 K	1.9 kΩ	-351 Ω/Κ	609.6 μK/nW	6.32 kΩ	200 μV	31.6 nA	1.9 pW	50.2 mΩ (143 μK) 87 mΩ (248 μK)	0.9 Ω (2.6 mK)	±4 mK	1.2 μΚ	±0.2 mK	4.7 mK

#### Lake Shore GR-50-AA with 0.05 to 6 K calibration

Values given are for measurement input. If the value is different for the control input, it is shown in blue.

	Sensor p	roperties		Excitat	ion and ir	strument	ation	Instrument p	erformance		curacy errors error acc 4 mK 7.0 μK ±0.2 mK 4 4 mK 4.6 μK ±0.2 mK 4 4 mK 66 μK ±0.2 mK 4		
Temperature	Nominal resistance	Typical sensor sensitivity	Thermal resistance	Resistance range	Excitation voltage limit	Excitation current	Power	Measurement resolution	Electronic accuracy	Calibration accuracy			Overall accuracy
50 mK	35 kΩ	-3.6 MΩ/K	200 mK/nW	63.2 kΩ 200 kΩ	63.2 μV 200 μV	1 nA	35 fW	$3.4~\Omega$ (944 nK) $7.3~\Omega$ (2 $\mu$ K)	20.7 Ω (5.8 μK) 27.5 Ω (7.6 μK)	±4 mK	7.0 µK	±0.2 mK	4 mK
100 mK	2317 Ω	-72 kΩ/K	20 mK/nW	6.32 kΩ 20 kΩ	63.2 μV 200 μV	10 nA	232 fW	$48.5 \text{ m}\Omega$ (674 nK) $338 \text{ m}\Omega$ (4.7 μK)	1.5 Ω (20.8 μK) 1.7 Ω (23.6 μK)	±4 mK	4.6 μΚ	±0.2 mK	4 mK
300 mK	164 Ω	-964 Ω/K	4 mK/nW	632 Ω 2 kΩ	200 μV	316 nA 100 nA	16 pW 1.6 pW	$3.6~\text{m}\Omega$ (3.7 µK) $29~\text{m}\Omega$ (30.1 µK)	81 mΩ (84 μK) 149 mΩ (155 μK)	±4 mK	66 μK 6.6 μK	±0.2 mK	4 mK
500 mK	73.8 Ω	-202.9 Ω/K	1.2 mK/nW	632 Ω 2 kΩ	200 μV	316 nA 100 nA	7.4 pW 738 fW	3.6 mΩ (17.7 μK) 29 mΩ (143 μK)	54 mΩ (266 μK) 122 mΩ (601 μK)	±4 mK	8.9 μK 886 nK	±0.2 mK	4 mK 4.1 mK
1 K	34 Ω	-31 Ω/Κ	100 μK/nW	200 Ω 2 kΩ	200 μV	1 μA 100 nA	34 pW 340 fW	1.2 mΩ (38.7 μK) 29 mΩ (935 μK)	20 mΩ (645 μK) 110 mΩ (3.5 mK)	±4 mK	3.4 μK 34 nK	±0.2 mK	4.1 mK 5.4 mK
1.4 K	24.7 Ω	-13.15 Ω/K	75 μK/nW	200 Ω 2 kΩ	200 μV	1 μA 100 nA	25 pW 247 fW	1.2 mΩ (91.3 μK) 29 mΩ (2.2 mK)	17 mΩ (1.3 mK) 107 mΩ (8.1 mK)	±5 mK	1.9 μK 19 nK	±0.2 mK	5.2 mK 9.8 mK
4.2 K	13.7 Ω	-1.036 Ω/K	25 μK/nW	20 Ω 2 kΩ	200 μV	10 μA 100 nA	1.4 nW 137 fW	120 μΩ (116 μK) 29 mΩ (28 mK)	5.1 mΩ (4.9 mK) 104 mΩ (100 mK)	±5 mK	3.5 μK 3.4 nK	±0.2 mK	7 mK 104 mK

Model 372 AC Resistance Bridge

#### Lake Shore CX-1010-SD with 0.1 to 325 K calibration

Values given are for measurement input. If the value is different for the control input, it is shown in blue.

	Sensor p	roperties		Excitat	tion and ir	strument	ation	Instrument p	erformance		Overall pe	errors error accu 30 μK ±0.2 mK 4 n 52 μK ±0.2 mK 4 n	
Temperature	Nominal resistance	Typical sensor sensitivity	Thermal resistance	Resistance range	Excitation voltage limit	Excitation current	Power	Measurement resolution	Electronic accuracy	Calibration accuracy	Self-heating errors	•	Overall accuracy
100 mK	21.389 kΩ	-558 kΩ/K	1.4 K/nW	63.2 kΩ 200 kΩ	63.2 μV 200 μV	1 nA	21 fW	3.4 Ω (6.1 μK) 7.4 Ω (13.3 μK)	13.9 Ω (24.9 μK) 20.7 Ω (37.1 μK)	±4 mK	30 μK	±0.2 mK	4 mK
300 mK	2.3224 kΩ	-10.8 kΩ/K	26.8 mK/nW	6.32 kΩ	200 μV	31.6 nA	2.3 pW	$50.2~\text{m}\Omega$ (4.6 µK) $87.0~\text{m}\Omega$ (8.1 µK)	1.0 Ω (92.6 μK) 1.0 Ω (93.8 μK)	±4 mK	62 μK	±0.2 mK	4 mK
500 mK	1.2482 kΩ	-2.7 kΩ/K	4.3 mK/nW	2 kΩ	200 μV	100 nA	12.5 pW	$\begin{array}{c} \text{14.5 m}\Omega\text{ (5.4 }\mu\text{K)} \\ \text{29.2 m}\Omega\text{ (10.8 }\mu\text{K)} \end{array}$	475 mΩ (176 μK) 474 mΩ (176 μK)	±4 mK	54 μK	±0.2 mK	4 mK
4.2 K	277.32 Ω	-32.2 Ω/K	2 μK/nW	632 Ω 2 kΩ	6.32 mV 200 μV	10 μA 100 nA	28 nW 2.8 pW	1.3 mΩ (40.4 μK) 29.2 mΩ (907 μK)	115 mΩ (3.6 mK) 183 mΩ (5.7 mK)	±4 mK	56 μK 5.6 nK	±0.2 mK	5.4 mK 7 mK
300 K	30.392 Ω	-65.4 mΩ/K	426 fK/nW	63.2 Ω 2 kΩ	6.32 mV 200 μV	100 μA 100 nA	304 nW 304 fW	130 μ $\Omega$ (2.0 mK) 29.2 m $\Omega$ (446 mK)	12.3 mΩ (188 mK) 109 mΩ (1.7 K)	±78 mK	130 pK 129 aK	±0.2 mK	203 mK 1.7 K

#### 372/3726 performance specification table

The values below apply to the measurement input. The control input operates over a reduced range indicated by the black-bordered cells. These cells contain bracketed numbers to indicate the resolution that applies to the control input.

						Voltage	e range					
	632 mV	200 mV	63.2 mV	20 mV	6.32 mV	2 mV	632 μV	200 μV	63.2 μV	20 μV	6.32 μV	2 μV
31.6 mA	20 Ω 20 μΩ 10 mW	6.32 Ω 6.3 μΩ 3.2 mW	2 Ω 2 μΩ 1 mW	632 mΩ 1.3 μΩ 320 μW	200 mΩ 400 nΩ 100 μW	63.2 mΩ 95 nΩ 32 μW	20 mΩ 36 nΩ 10 μW	6.32 mΩ 35 nΩ 3.2 μW	2 mΩ 40 nΩ 1 μW	* *	* *	* *
10 mA	63.2 Ω 63 μΩ 3.2 mW	20 Ω 20 μΩ 1 mW	6.32 Ω 6.3 μΩ 320 μW	2 Ω 4 μΩ 100 μW	632 mΩ 1.3 μΩ 32 μW	200 mΩ 300 nΩ 10 μW	63.2 mΩ 120 nΩ 3.2 μW	20 mΩ 120 nΩ 1 μW	6.32 mΩ 130 nΩ 320 nW	2 mΩ 120 nΩ 100 nW	* *	* *
3.16 mA	200 Ω 200 μΩ 1 mW	63.2 Ω 63 μΩ 320 μW	20 Ω 20 μΩ 100 μW	6.32 Ω 13 μΩ 32 μW	2 Ω 4 μΩ 10 μW	632 mΩ 950 nΩ 3.2 μW	200 mΩ 390 nΩ 1 μW	63.2 mΩ 370 nΩ 320 nW	20 mΩ 400 nΩ 100 nW	6.32 mΩ 380 nΩ 32 nW	2 mΩ 400 nΩ 10 nW	* *
1 mA	632 Ω 630 μΩ 3.2E-04	200 Ω 200 μΩ 100 μW	63.2 Ω 63 μΩ 32 μW	20 Ω 40 μΩ 10 μW	6.32 Ω 13 μΩ 3.2 μW	2 Ω 3 μΩ 1 μW	632 mΩ 1 μΩ 320 nW	200 mΩ 1 μΩ 100 nW	63.2 mΩ 1.3 μΩ 32 nW	20 mΩ 1.2 μΩ 10 nW	6.32 mΩ 1.3 μΩ 3.2 nW	2 mΩ 1 μΩ 1 nW
316 μΑ	2 kΩ 2 mΩ 100 μW	632 Ω 630 μΩ 32 μW	200 Ω 200 μΩ 10 μW	63.2 Ω 130 μΩ 3.2 μW	20 Ω 40 μΩ 1 μW	6.32 Ω 9.5 μΩ 320 nW	2 Ω 3.8 μΩ 100 nW	632 mΩ 3.7 μΩ 32 nW	200 mΩ 4 μΩ 10 nW	63.2 mΩ 3.8 μΩ 3.2 nW	20 mΩ 4 μΩ 1 nW	6.32 mΩ 3.7 μΩ 320 pW
100 μΑ	6.32 kΩ 6.3 mΩ 32 μW	2 kΩ 2 mΩ 10 μW	632 Ω 630 μΩ 3.2 μW	200 Ω 400 μΩ 1 μW	63.2 Ω 130 μΩ 320 nW	20 Ω 30 μΩ 100 nW	6.32 Ω 12 μΩ 32 nW	2 Ω 12 μΩ 10 nW	632 mΩ 13 μΩ 3.2 nW	200 mΩ 12 μΩ 1 nW	63.2 mΩ 13 μΩ 320 pW	20 mΩ 12 μΩ 100 pW
31.6 μΑ	20 kΩ 20 mΩ 10 μW	6.32 kΩ 6.3 mΩ 3.2 μW	2 kΩ 2 mΩ 1 μW	632 Ω 1.3 mΩ 320 nW	200 Ω 400 μΩ 100 nW	63.2 Ω 95 μΩ 32 nW	20 Ω 37 μΩ 10 nW	6.32 Ω 37 μΩ 3.2 nW	2 Ω 40 μΩ 1 nW	632 mΩ 38 μΩ 320 pW	200 mΩ 40 μΩ 100 pW	63.2 mΩ 37 μΩ 32 pW
10 μΑ	63.2 kΩ 63 mΩ 3.2 μW	20 kΩ 20 mΩ 1 μW	6.32 kΩ 6.3 mΩ 320 nW	2 kΩ 4 mΩ 100 nW	632 Ω 1.3 mΩ 32 nW	200 Ω 300 μΩ 10 nW	63.2 Ω 120 μΩ 3.2 nW	20 Ω 120 μΩ 1 nW	6.32 Ω 130 μΩ 320 pW	2 Ω 120 μΩ 100 pW	632 mΩ 130 μΩ 32 pW	200 mΩ 120 μΩ 10 pW
3.16 μΑ	200 kΩ 200 mΩ 1 μW	63.2 kΩ 63 mΩ 320 nW	20 kΩ 20 mΩ 100 nW	6.32 kΩ 13 mΩ 32 nW	2 kΩ 4 mΩ 10 nW	632 Ω 950 μΩ 3.2 nW	200 Ω 370 μΩ 1 nW	63.2 Ω 370 μΩ 320 pW	20 Ω 400 μΩ 100 pW	6.32 Ω 380 μΩ 32 pW	2 Ω 400 μΩ 10 pW	632 mΩ 370 μΩ 3.2 pW
1 μΑ	632 kΩ 630 mΩ 320 nW	200 kΩ 200 mΩ 100 nW	63.2 kΩ 63 mΩ 32 nW	20 kΩ 40 mΩ 10 nW	6.32 kΩ 13 mΩ 3.2 nW	2 kΩ 3 mΩ 1 nW	632 Ω 1.2 mΩ 320 pW	200 Ω 1.2 mΩ 100 pW	63.2 Ω 1.3 mΩ 32 pW	20 Ω 1.2 mΩ 10 pW	6.32 Ω 1.3 mΩ 3.2 pW	2 Ω 1.2 mΩ 1 pW
316 nA	2 MΩ 2 Ω 100 nW	632 kΩ 630 mΩ 32 nW	200 kΩ 200 mΩ 10 nW	63.2 kΩ 130 mΩ 3.2 nW	20 kΩ 40 mΩ 1 nW	6.32 kΩ 13 mΩ 320 pW	2 kΩ 4 mΩ 100 pW	632 Ω 3.8 mΩ 32 pW	200 Ω 4 mΩ 10 pW	63.2 Ω 3.8 mΩ 3.2 pW	20 Ω 4 mΩ 1 pW	6.32 Ω 3.8 mΩ 320 fW
316 nA	6.32 MΩ ** 32 nW	2 MΩ 2 Ω 10 nW	632 kΩ 630 mΩ 3.2 nW	200 kΩ 400 mΩ 1 nW	63.2 kΩ 130 mΩ 320 pW	20 kΩ 30 mΩ 100 pW	6.32 kΩ 13 mΩ 32 pW	2 kΩ 16 [30] mΩ 10 pW	632 Ω 13 mΩ 3.2 pW	200 Ω 12 mΩ 1 pW	63.2 Ω 13 mΩ 320 fW	20 Ω 12 mΩ 100 fW
31.6 nA	20 MΩ ** 10 nW	6.32 MΩ ** 3.2 nW	2 MΩ 2 Ω 1 nW	632 kΩ 1.3 Ω 320 pW	200 kΩ 300 mΩ 100 pW	63.2 kΩ 160 mΩ 32 pW	20 kΩ 100 mΩ 10 pW	6.32 kΩ 63 [95] mΩ 3.2 pW	2 kΩ 40 mΩ 1 pW	632 Ω 38 mΩ 320 fW	200 Ω 40 mΩ 100 fW	63.2 Ω 38 mΩ 32 fW
10 nA	63.2 MΩ **	20 MΩ **	6.32 MΩ ** 320 pW	2 MΩ 3 Ω 100 pW	632 kΩ 1.6 Ω 32 pW	200 kΩ 600 mΩ 10 pW	63.2 kΩ 470 mΩ 3.2 pW	20 kΩ 300 [400] mΩ 1 pW	6.32 kΩ 130 mΩ 320 fW	2 kΩ 160 mΩ 100 fW	632 Ω 130 mΩ 32 fW	200 Ω 120 mΩ 10 fW
3.16 nA	* *	63.2 MΩ ** 320 pW	20 MΩ ** 100 pW	6.32 MΩ ** 32 pW	2 MΩ 9 Ω 10 pW	632 kΩ 4.7 Ω 3.2 pW	200 kΩ 3 Ω 1 pW	63.2 kΩ 1.6 [1.9] Ω 320 fW	20 kΩ 1 Ω 100 fW	6.32 kΩ 630 mΩ 32 fW	2 kΩ 500 mΩ 10 fW	632 Ω 380 mΩ 3.2 fW
1 nA	* *	* *	63.2 MΩ ** 32 pW	20 MΩ ** 10 pW	6.32 MΩ ** 3.2 pW	2 MΩ 30 Ω 1 pW	632 kΩ 16 Ω 320 fW	200 kΩ 6 [10] Ω 100 fW	63.2 kΩ 5.1 Ω 32 fW	20 kΩ 3 Ω 10 fW	6.32 kΩ 1.3 Ω 3.2 fW	2 kΩ 1.6 Ω 1 fW
316 pA		* resistance		63.2 MΩ ** 3.2 pW	20 MΩ ** 1 pW	6.32 MΩ ** 320 fW	2 MΩ 90 Ω 100 fW	632 kΩ 47 [51] Ω 32 fW	200 kΩ 30 Ω 10 fW	63.2 kΩ 16 Ω 3.2 fW	20 kΩ 10 Ω 1 fW	6.32 kΩ 6.3 Ω 320 aW
100 pA	1.0 fW	100 $\Omega$ [150 $\Omega$ ] — measurement resolution [control resolution] 1.0 fW — power		* *	63.2 MΩ ** 320 fW	20 MΩ ** 100 fW	6.32 MΩ ** 32 fW	2 MΩ 300 Ω 10 fW	632 kΩ 160 Ω 3.2 fW	200 kΩ 100 Ω 1 fW	63.2 kΩ 51 Ω 320 aW	20 kΩ 30 Ω 100 aW
31.6 pA	Resistance range: Full scale resistance range, nominal 20% over range.  Resolution: RMS noise with 18 s		20% Acc	wracy	*	63.2 MΩ ** 32 fW	20 MΩ ** 10 fW	6.32 MΩ ** 3.2 fW	2 MΩ 900 Ω 1 fW	632 kΩ 470 Ω 320 aW	200 kΩ 300 Ω 100 aW	63.2 kΩ 160 Ω 32 aW
10 pA	filter settlin 3 s analog	g time (approxii time constant). t ½ full scale re	nates Noise sistance	±0.05% + 0.0	005% of range 008% of range 15% of range	. *	63.2 MΩ ** 3.2 fW	20 MΩ ** 1 fW	6.32 MΩ ** 320 aW	2 MΩ 3 kΩ 100 aW	632 kΩ 1.6 kΩ 32 aW	200 kΩ 1 kΩ 10 aW
3.16 pA	Power: Exc one-half fu	riperature. citation power a Il scale resistan Dominated by	t e	±0.3% + 0.05 ±0.5% + 0.05	5% of range	* * *	* *	63.2 MΩ ** 320 aW	20 MΩ ** 100 aW	6.32 MΩ ** 32 aW	2 MΩ 9 kΩ 10 aW	632 kΩ 4.7 kΩ 3.2 aW
1 pA	measurem coefficient	ent temperature (±0.0015% of roof range)/°C	eading * F	±1.0% + 0.15		* * *	* *	* *	63.2 MΩ ** 32 aW	20 MΩ ** 10 aW	6.32 MΩ ** 3.2 aW	2 MΩ 30 kΩ 1 aW




#### 372/3708 performance specification table

					Voltage rang	е		
	6.32 mV	2.0 mV	632 μV	200 μV	63.2 μV	20 μV	6.32 μV	2.0 μV
31.6 mA	200 mΩ 200 nΩ 100 μW	63.2 mΩ 63 nΩ 32 μW	20 mΩ 40 nΩ 10 μW	6.32 mΩ 13 nΩ 3.2 μW	2.0 mΩ 10 nΩ 1.0 μW	* *	* *	* *
10 mA	632 mΩ 630 nΩ 32 μW	200 mΩ 200 nΩ 10 μW	63.2 mΩ 130 nΩ 3.2 μW	20 mΩ 40 nΩ 1.0 μW	6.32 mΩ 32 nΩ 320 nW	2.0 mΩ 32 nΩ 100 nW	* *	* *
3.16 mA	2.0 Ω 2.0 μΩ 10 μW	632 mΩ 630 nΩ 3.2 μW	200 mΩ 400 nΩ 1.0 μW	63.2 mΩ 130 nΩ 320 nW	20 mΩ 100 nΩ 100 nW	6.32 mΩ 100 nΩ 32 nW	2.0 mΩ 100 nΩ 10 nW	* *
1 mA	6.32 Ω 6.3 μΩ 3.2 μW	2.0 Ω 2.0 μΩ 1.0 μW	632 mΩ 1.3 μΩ 320 nW	200 mΩ 400 nΩ 100 nW	63.2 mΩ 320 nΩ 32 nW	20 mΩ 320 nΩ 10 nW	6.32 mΩ 320 nΩ 3.2 nW	2.0 mΩ 320 nΩ 1.0 nW
316 μΑ	20 Ω 20 μΩ 1.0 μW	6.32 Ω 6.3 μΩ 320 nW	2.0 Ω 4.0 μΩ 100 nW	632 mΩ 1.3 μΩ 32 nW	200 mΩ 1.0 μΩ 10 nW	63.2 mΩ 1.0 μΩ 3.2 nW	20 mΩ 1.0 μΩ 1.0 nW	6.32 mΩ 1.0 μΩ 320 pW
100 μΑ	63.2 Ω 63 μΩ 320 nW	20 Ω 20 μΩ 100 nW	6.32 Ω 13 μΩ 32 nW	2.0 Ω 4.0 μΩ 10 nW	632 mΩ 3.2 μΩ 3.2 nW	200 mΩ 3.2 μΩ 1.0 nW	63.2 mΩ 3.2 mΩ 320 pW	20 mΩ 3.2 μΩ 100 pW
31.6 µА	200 Ω 200 μΩ 100 nW	63.2 Ω 63 μΩ 32 nW	20 Ω 40 μΩ 10 nW	6.32 Ω 13 μΩ 3.2 nW	2.0 Ω 10 μΩ 1.0 nW	632 mΩ 10 μΩ 320 pW	200 mΩ 10 μΩ 100 pW	63.2 mΩ 10 μΩ 32 pW
10 μΑ	632 Ω 630 μΩ 32 nW	200 Ω 200 μΩ 10 nW	63.2 Ω 130 μΩ 3.2 nW	20 Ω 40 μΩ 1.0 nW	6.32 Ω 32 μΩ 320 pW	2.0 Ω 32 μΩ 100 pW	632 mΩ 32 μΩ 32 pW	200 mΩ 32 μΩ 10 pW
3.16 μΑ	2.0 kΩ 2.0 mΩ 10 nW	632 Ω 630 μΩ 3.2 nW	200 Ω 400 μΩ 1.0 nW	63.2 Ω 130 μΩ 320 pW	20 Ω 100 μΩ 100 pW	6.32 Ω 100 μΩ 32 pW	2.0 Ω 100 μΩ 10 pW	632 mΩ 100 μΩ 3.2 pW
πtion ation	$\begin{array}{c} 6.32 \text{ k}\Omega \\ 6.3 \text{ m}\Omega \\ 3.2 \text{ nW} \end{array}$	2.0 kΩ 2.0 mΩ 1.0 nW	$\begin{array}{c} 632~\Omega \\ 1.3~\text{m}\Omega \\ 320~\text{pW} \end{array}$	200 Ω 400 μΩ 100 pW	63.2 Ω 320 μΩ 32 pW	20 Ω 320 μΩ 10 pW	6.32 Ω 320 μΩ 3.2 pW	2.0 Ω 320 μΩ 1.0 pW
Current excitation	20 kΩ 20 mΩ 1.0 nW	$\begin{array}{c} 6.32 \text{ k}\Omega \\ 6.3 \text{ m}\Omega \\ 320 \text{ pW} \end{array}$	2.0 kΩ 4.0 mΩ 100 pW	632 Ω 1.3 mΩ 32 pW	200 Ω 1.0 mΩ 10 pW	63.2 Ω 1.0 mΩ 3.2 pW	20 Ω 1.0 mΩ 1.0 pW	6.32 Ω 1.0 mΩ 320 fW
100 nA	63.2 kΩ 63 mΩ 320 pW	20 kΩ 40 mΩ 100 pW	6.32 kΩ 13 mΩ 32 pW	2.0 kΩ 6.0 mΩ 10 pW	632 Ω 3.2 mΩ 3.2 pW	200 Ω 3.2 mΩ 1.0 pW	63.2 Ω 3.2 mΩ 320 fW	$\begin{array}{c} 20~\Omega\\ 3.2~\text{m}\Omega\\ 100~\text{fW} \end{array}$
31.6 nA	200 kΩ 400 mΩ 100 pW	63.2 kΩ 130 mΩ 32 pW	20 kΩ 60 mΩ 10 pW	6.32 kΩ 20 mΩ 3.2 pW	2.0 kΩ 20 mΩ 1.0 pW	$\begin{array}{c} 632~\Omega \\ 10~\text{m}\Omega \\ 320~\text{fW} \end{array}$	200 Ω 10 mΩ 100 fW	63.2 Ω 10 mΩ 32 fW
10 nA	632 kΩ 1.9 Ω 32 pW	200 kΩ 600 mΩ 10 pW	63.2 kΩ 200 mΩ 3.2 pW	20 kΩ 200 mΩ 1.0 pW	$\begin{array}{c} 6.32 \text{ k}\Omega \\ 63 \text{ m}\Omega \\ 320 \text{ fW} \end{array}$	2.0 kΩ 63 mΩ 100 fW	632 Ω 32 Ω 32 fW	200 Ω 32 mΩ 10 fW
3.16 nA	2.0 MΩ 6.0 Ω 10 pW	632 kΩ 2.0 Ω 3.2 pW	200 kΩ 2.0 Ω 1.0 pW	63.2 kΩ 630 mΩ 320 fW	20 kΩ 600 mΩ 100 fW	$\begin{array}{c} 6.32 \text{ k}\Omega \\ 200 \text{ m}\Omega \\ 32 \text{ fW} \end{array}$	2.0 kΩ 200 mΩ 10 fW	$\begin{array}{c} 632~\Omega \\ 100~\text{m}\Omega \\ 3.2~\text{fW} \end{array}$
1 nA	6.32 MΩ ** 3.2 pW	2.0 MΩ 20 Ω 1.0 pW	632 kΩ 6.3 Ω 320 fW	200 kΩ 6.0 Ω 100 fW	63.2 kΩ 3.2 Ω 32 fW	20 kΩ 2.0 Ω 10 fW	6.32 kΩ 630 mΩ 3.2 fW	2.0 kΩ 1.0 Ω 1.0 fW
316 pA	* *	6.32 MΩ ** 320 fW	2.0 MΩ 60 Ω 100 fW	632 kΩ 19 Ω 32 fW	200 kΩ 20 Ω 10 fW	$63.2 \text{ k}\Omega$ $6.3 \Omega$ $3.2 \text{ fW}$	20 kΩ 3.0 Ω 1.0 fW	6.32 kΩ 3.2 Ω 320 aW
100 pA	* *	* *	6.32 MΩ ** 32 fW	2.0 MΩ 200 Ω 10 fW	632 kΩ 63 Ω 3.2 fW	200 kΩ 60 Ω 1.0 fW	63.2 kΩ 32 Ω 320 aW	20 kΩ 20 Ω 100 aW
31.6 pA	* *	* *	* *	6.32 MΩ ** 3.2 fW	2.0 MΩ 600 Ω 1.0 fW	632 kΩ 190 Ω 320 aW	200 kΩ 200 Ω 100 aW	63.2 kΩ 63 Ω 32 aW
10 pA	* *	* *	* *	* *	6.32 MΩ ** 320 aW	2.0 MΩ 2.0 kΩ 100 aW	632 kΩ 630 Ω 32 aW	200 kΩ 600 Ω 10 aW
3.16 pA	* *	* *	* *	* *	* *	6.32 MΩ ** 32 aW	2.0 MΩ 6.0 kΩ 10 aW	632 kΩ 1.9 kΩ 3.2 aW

$\begin{array}{cccc} 200 \ k\Omega & & & \\ 100 \ \Omega & & & \\ & & \\ 1.0 \ fW & & \\ & $
Resistance range: Full scale
resistance range, nominal 20% over
range.
Resolution: RMS noise with 18 s filter
settling time (approximates 3 s analog
time constant). Noise specified at ½ full
scale resistance at room temperature.
Power: Excitation power at
one-half full scale resistance.

Precision: Dominated by measurement temperature coefficient (±0.0015% of reading ±0.0002% of range)/°C.

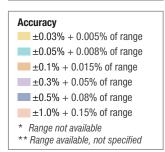

#### **Specifications**

#### Measurement input

Temperature coefficient

Input type AC, four-lead differential, resistance

Number of inputs

Maximum channels 16 (with optional scanner) Measurement units  $\Omega$ , K (with temperature curve)

22 ranges from 2 m $\Omega$  to 63.2 M $\Omega$  (excitation dependent) Resistance ranges

Maximum update rate 10 rdg/s (single range and input)

Range change settling 3 s + filter settling Channel change (scan) settling 3 s + filter settling

Sensor and range dependent, refer to Measurement Input Resolution

Specifications table

Sensor and range dependent, refer to Measurement Input Accuracy

> Specifications table ±0.0015%/°C of rdg

Maximum lead resistance  $100 \Omega + 10\%$  of resistance range per lead for current

 $\leq$ 3.16 mA; 10  $\Omega$  + 10% of resistance range per lead for

current ≥10 mA

Isolation Isolated from chassis and heater grounds Lead connections V+, V-, I+, I-, V shield, I shield, individual guards V+, V-, I+, I-, for each sensor, shield common to all Scanner lead connections Matched impedance voltage input and current output, Common mode rejection

active CMR

Excitation Sinusoidal AC current source

**Excitation frequency** 9.8 Hz, 11.6 Hz, 13.7 Hz (default), 16.2 Hz, or 18.2 Hz

**Excitation currents** 22 ranges from 1 pA to 31.6 mA RMS

**Excitation accuracy** ±2% of nominal

Minimum excitation power  $10^{-18}$  W into a 100 k $\Omega$  (see Measurement Input

Specifications table for other ranges)

Typical DC bias current 2 pA + 1% of excitation current (4.0  $\times$   $10^{\text{-19}}\,\text{W}$  into

100 kΩ)

4 pA + 1% of excitation current (1.6  $\times$  10  $^{\!-18}$  W into Maximum DC bias current

 $100 \text{ k}\Omega$ )

Power up current protection Current output shunted on power up Voltage input ranges 12 ranges from 2 µV to 632 mV RMS

Voltage input over-range Voltage input impedance

 ${>}5\times10^{\scriptscriptstyle13}\,\Omega$ Maximum input voltage noise 10 nV/J/Hz at 10 Hz

Range selection modes

Manual, voltage excitation, current excitation, autorange Scanner modes Manual or autoscan

1 s to 200 s settling time, 1% to 80% filter window Additional software features

Min/Max reading capture, pause (3 s to 60 s) on range and/or channel change, scanner dwell time (1 s to 200 s)

Supported temperature sensors NTC resistive sensors including germanium, Cernox®,

Rox<sup>™</sup>, PTC resistive sensors including rhodium-iron

Quadrature display Real and Imaginary

Connectors 6-pin DIN (current out), 6-pin DIN (voltage in), and DA-15

(scanner control)

Supported scanners Lake Shore 3726 and 3708

#### **Control** input

Input type AC, four-lead differential, resistance

Number of inputs

Measurement units Ω, K (with temperature curve)

6 ranges from 2 k $\Omega$  to 632 k $\Omega$  (excitation dependent) Resistance ranges

10 rdg/s (single range) Maximum update rate Range change settling 3 s + filter settling

Resolution Sensor and range dependent, refer to Control Input

Specifications table

Sensor and range dependent, refer to Control Input Accuracy

Specifications table

Temperature coefficient ±0.0015%/°C of reading

Maximum lead resistance  $100 \Omega + 10\%$  of resistance range per lead

Isolation Isolated from chassis, common to measurement input

V+, V-, I+, I-, shield Lead connections

Common mode rejection Matched impedance voltage input and current output

Excitation Sinusoidal AC current source

**Excitation frequency** 9.8 Hz, 11.6 Hz, 13.7 Hz, 16.2 Hz (default), or 18.2 Hz

**Excitation currents** 6 ranges from 316 pA to 100 nA RMS **Excitation accuracy** ±8% of nominal for 316 pA and 1 nA ranges; ±2% of nominal for the other ranges

Power up current protection Current output shunted on power up Voltage input range 200 µV

Voltage input over-range 20%

Maximum input voltage noise 20 nV/√Hz at 10 Hz

Manual, standard autorange, and Rox™ RX-102B-CB Range selection modes

optimized autorange

Filter 1 s to 200 s settling time, 1% to 80% filter window

Min/Max reading capture Additional software features

Supported sensors NTC resistive sensors (optimized for Rox™ RX-102B-CB

Down to 10 mK using a Rox™ RX-102B-CB sensor in a Minimum temperature

well-designed system

Connector 6-pin DIN

#### Temperature conversion

Sensor temperature coefficient Negative or positive

User curves Up to 39 CalCurves™ or user curves (200-point)

Curve entry Via front panel or computer interface

 $\Omega/K$ , Log  $\Omega/K$ Curve format Curve interpolation Cubic spline, linear

#### Sample heater output

Туре	Variable DC current source
Control modes	Closed loop PID, PID zones, open loop
Setpoint units	Ω, K (with temperature curve)
D/A resolution	16-bit
Ranges	100 mA, 31.6 mA, 10 mA, 3.16 mA, 1 mA, 316 μA,
	100 μΑ, 31.6 μΑ
Output compliance voltage (min)	±10 V
Maximum power of output ranges	1 W, 100 mW, 10 mW, 1 mW, 100 μW, 10 μW, 1 μW,
, , , ,	0.1 μW
Resistance range	1 $\Omega$ to 2 k $\Omega$ , 100 $\Omega$ for maximum power
Heater offset (at 0%)	±0.02% of range
Heater gain accuracy	±1% of setting
Heater noise	<0.005% of range
Isolation	Isolated from chassis ground, measurement and
	control inputs; shared ground with analog/still
	output
Heater connector	Detachable terminal block
Safety limits	Curve temperature, power up heater off, shunted
	with a relay on power up, short-circuit protection.
	compliance voltage limit detection, input
	temperature limit
Additional software features	Heater power display based on user entered
	resistance
Heater offset (at 0%) Heater gain accuracy Heater noise Isolation  Heater connector Safety limits	±1% of setting <0.005% of range Isolated from chassis ground, measurement and control inputs; shared ground with analog/still output Detachable terminal block Curve temperature, power up heater off, shunter with a relay on power up, short-circuit protection compliance voltage limit detection, input temperature limit Heater power display based on user entered

#### Warm-un heater outnut

warm-up neater output					
Туре	Variable DC current source	9			
Control modes	Closed loop PID, PID zones, open loop, warm-up				
	mode				
Setpoint units	Ω, K (with temperature cu	rve)			
D/A resolution	16-bit				
	25 Ω setting	50 Ω setting			
Maximum power	10 W	10 W			
Maximum current	0.63 A	0.45 A			
Voltage compliance (min)	+15.8 V	+22.4 V			
Heater load for maximum power	25 Ω	50 Ω			
Resistance range	10 Ω to 100 Ω				
Isolation	Chassis ground reference				
Heater connector	Detachable terminal block				
Safety limits	Curve temperature, power	up heater off, shunted			
	with a relay on power up,	short-circuit protection,			
	compliance voltage limit, i	relay disconnects output			
	when off, input temperatu				



#### Analog/still output

Variable DC voltage source Type Control modes Open loop, still heater, monitor output

Isolation Isolated from chassis ground, measurement and control

inputs; shared ground with sample heater

Output voltage range ±10 V Maximum current 100 mA Maximum power 1 W into 100  $\Omega$ 

100  $\Omega$  (short-circuit protected) Minimum load resistance

Accuracy ±2.5 mV Noise (resolution) <0.003% of range

Monitor output settings

Scale User selected

Data source Temperature or sensor units

Settings Input, source, top of scale, and bottom of scale

Connector Detachable terminal block

#### **Heater control**

Number of control loops 2 (sample heater, warm-up heater)

Update rate

Manual PID, zone Tuning

PID control settings

Proportional (gain) 0.0 to 1,000 Integral (reset) 0 to 10,000 s Derivative (rate) 0 to 2,500 s

Manual output 0 to 100% with 0.01% setting resolution

10 temperature zones with P, I, D, manual heater out, Zone control

heater range, setpoint, relays, and analog output (still)

Setpoint ramping 0.001 K/min to 100 K/min

Scanner support Control with scanned channel (reduced stability) Control stability Below 10 µK peak-to-peak at 50 mK (system dependent)

Warm-up heater mode settings

Warm-up percentage 0 to 100% with 1% resolution Warm-up mode Continuous control or auto-off

#### Front panel

Display 8-line by 40-character (256 × 64 pixel) graphic VF display

module 1 to 8

Number of reading displays

Display units mK, K, mΩ, Ω, kΩ, MΩ

Reading source Resistance, temperature, max, min

Display update rate 2 rdg/s

Other displays Input name, channel number, resistance range, excitation

voltage, excitation current, excitation power, control setpoint, PID, heater range, heater output, and quadrature

Same as display resolution (sensor-dependent) Setpoint setting resolution Heater output display Numeric display in percent of full scale for power or

current

Display annunciators Control input and alarm

LED annunciators Autorange, excitation mode, autoscan, control outputs,

remote, Ethernet status, alarm, still output

Kevpad 34-key silicone elastomer keypad Front panel features Front panel curve entry, and keypad lock-out

#### Interface

IEEE-488.2

Capabilities SH1, AH1, T5, L4, SR1, RL1, PP0, DC1, DT0, C0, E1

Update rate To 10 rdg/s on each input

Software support LabVIEW™ driver (see www.lakeshore.com)

**USB** 

Function Emulates a standard RS-232 serial port

Baud rate 57.600

Connector B-type USB connector Update rate To 10 rdg/s on each input

LabVIEW™ driver (see www.lakeshore.com) Software support

Ethernet

Function TCP/IP, web interface, curve handler, configuration backup.

chart recorder

Connector R.I-45 Update rate To 10 rdg/s on each input

Software support LabVIEW™ driver (see www.lakeshore.com) Special interface feature Model 370 command emulation mode

300, 1,200, 9,600, 57,600

Alarms

Number 34, high and low for each measurement channel and the

control input

Data source Temperature or sensor units

Source, high setpoint, low setpoint, deadband, latching or Settings

non-latching, audible on/off, visible on/off Display annunciator, beeper, and relays

Actuators Relays

> Number 2

Available baud rates

Normally open (NO), normally closed (NC), and common (C) Contacts

30 VDC at 2 A

Activate relays on high, low, or both alarms for any Operation

measurement channel or control input, manual mode, or

zone control mode

Connector Diagnostic monitor output

Contact rating

Operation

User selects one of several analog voltage diagnostic

Detachable terminal block monitor output

points (must remain isolated)

1. AC voltage driving positive/negative side of current Available signals

source programming resistor

2. AC voltage present on the positive/negative side of the

differential input amplifier

3. AC voltage present on the output of the differential input

amplifier

4. AC voltage into the measurement channel or control

input AD converter

Reference output

Connector

Phase-sensitive detector reference (must remain isolated) Signal type

Amplitude 0 to +5 V nominal Waveform Square wave Connector **BNC** 

#### General

Ambient temperature 15 °C to 35 °C at rated accuracy; 5 °C to 40 °C at reduced

Power requirement 100, 120, 220, 240 VAC, ±10%, 50 or 60 Hz, 90 VA 435 mm W  $\times$  89 mm H  $\times$  368 mm D (17 in  $\times$  3.5 in  $\times$ Size

14.5 in), full rack 6.8 kg (15 lb)

Weight Approval CE mark, RoHS

135 mm W  $\times$  66 mm H  $\times$  157 mm D (5.3 in  $\times$  2.6 in  $\times$ Scanner size 6.2 in), plus connector clearance of 125 mm (5 in)



## Ordering information

#### **Description** Part number

372N AC resistance bridge and temperature controller with no connection cable 372U AC resistance bridge with 3708 scanner and standard 3 m (10 ft) connection cable AC resistance bridge with 3708 scanner and 6 m (20 ft) connection cable 372U-6 372U-10 AC resistance bridge with 3708 scanner and 10 m (33 ft) connection cable **372S** AC resistance bridge with 3726 scanner and standard 3 m (10 ft) connection cable AC resistance bridge with 3726 scanner and 6 m (20 ft) connection cable 372S-6 AC resistance bridge with 3726 scanner and 10 m (33 ft) connection cable

#### Please indicate your power/cord configuration:

**1** 100 V—U.S. cord (NEMA 5-15) 2 120 V—U.S. cord (NEMA 5-15) 3 220 V—Euro cord (CEE 7/7) 4 240 V—Euro cord (CEE 7/7)

5 240 V—U.K. cord (BS 1363) 6 240 V—Swiss cord (SEV 1011) **7** 220 V—China cord (GB 1002)

#### **Scanners**

372S-10

3708 Ultra-low resistance 8-channel scanner with standard 3 m (10 ft) connection cable—

includes one scanner cable and bracket kit (117-035)

3708-6 Ultra-low resistance 8-channel scanner with 6 m (20 ft) connection cable Ultra-low resistance 8-channel scanner with 10 m (33 ft) connection cable 3708-10 16-channel scanner with standard 3 m (10 ft) connection cable (Model 372 only) 3726 16-channel scanner with 6 m (20 ft) connection cable (Model 372 only) 3726-6 16-channel scanner with 6 m (20 ft) connection cable (Model 372 only) 3726-10

#### **Accessories/options**

G-112-374 3 m (10 ft) AC resistance bridge cable G-112-375 6 m (20 ft) AC resistance bridge cable G-112-376 10 m (33 ft) AC resistance bridge cable

117-017 1 m (3.3 ft long) IEEE-488 (GPIB) computer interface cable assembly

RM-1 Kit for mounting one full rack instrument in a 482.6 mm (19 in) rack mount cabinet

G-106-233 Sensor input mating connector (6-pin DIN plug)

117-071 Model 372 heater adapter cable CAL-372-CERT Instrument recalibration with certificate Instrument recalibration with certificate and data CAL-372-DATA

119-372 Model 372 user manual

All specifications are subject to change without notice











- platforms down to 100 mK
- Optimized performance with Cernox® RTDs
- Patented low-noise input circuitry enables super low excitation power for minimal self-heating and high resolution measurement
- 4 independent control loops and a broad range of I/O configurations can eliminate need for additional instrumentation
- 4 PID-controlled outputs: 75 W warm-up heater, 1 W sample heater, and 2 auxiliary 1 W ±10 V outputs
- Proven, intuitive interface
- Performance assurance even at the extremes, with verifiable product specifications
- CE certification
- Full 3 year standard warranty





#### A powerful ultra-low temperature physics tool

The Model 350 is designed for the demands of pumped He-3 refrigerators and other ultra-low and low temperature platforms. It provides excellent measurement performance, superior control accuracy, and convenient operation in a wide range of advanced research applications. Whether the need is for high accuracy with minimal thermal impact, or precise temperature control in high magnetic fields, or dependable measurement in radiation environments, the new Model 350 controller matched with Lake Shore's industry-leading Cernox® sensors provides a cryogenic solution that's demonstrably best-in-class.



The patented noise reduction input circuitry of the Model 350 is just one reason why this controller works so well for ultra-low temperature (ULT) applications, all the way down to 100 mK. When combined with precision Cernox sensors, this performance-optimized design allows as little as 10 nA of excitation current to be used, minimizing self-heating effects, and ensures best possible measurement accuracy throughout the entire temperature range.

This single instrument offers extraordinary capability and flexibility, often eliminating the need for additional instrumentation in a refrigeration control system. Its four input channels and four independent control outputs are configurable to support a broad range of I/O requirements, including the heaters and auxiliary devices typical of ULT refrigeration systems, as well as other cryogenic sensor types like ruthenium oxide and platinum RTDs. Standard computer interfaces enable remote communications, control and coordination with other systems.

In short, the Model 350 cryogenic temperature controller brings a new level of power, precision, and performance to critical low temperature physics research. It is ideal for use with He-3 systems, adiabatic demagnetization refrigerators (ADRs), certain dilution refrigerators, and many other applications demanding low thermal power and high measurement precision.

#### 4 standard sensor input channels

The Model 350 comes with four standard sensor inputs supporting Cernox®, ruthenium oxide, platinum RTDs, and other NTC RTD sensors. Inputs can be configured to accept any of the supported input types. Each sensor input channel has its own current source, providing fast settling times. The four sensor inputs are optically isolated from other circuits to reduce noise and to provide repeatable sensor measurements. Current reversal eliminates thermal electromotive force (EMF) errors in resistance sensors. Nine excitation currents facilitate temperature measurement and control down to 100 mK, with the nominal temperature range (using Cernox® sensors) spanning to 420 K. The instrument automatically selects the optimal current and gain levels for you once the sensor type is selected, and automatically scales current to minimize self-heating effects at low temperatures. The patented input circuitry eliminates any errors associated with grounding inconsistencies, making it easier to achieve reliable measurements at ultra-low temperatures. With the ability to label each sensor input channel with a customized name, it's also easy to identify the measured values being displayed.

#### **Application versatility**

Designed to support a broad range of sensor types, the Model 350 is performance-optimized for use over the entire temperature range of Cernox® sensors, making it the instrument of choice for ULT environments as well as other cryogenic systems where errors due to magneto-resistive or radiation effects need to be minimized.



# 3 option cards for more inputs and a wider range of applications

Field installable input option cards can expand vour sensor selection to include silicon diodes (like DT-670), capacitance sensors or thermocouples. Once installed, the option input can be selected and named from the front panel like any other input type. These option cards further expand the application versatility of the Model 350 temperature controller by allowing specialized sensors to be switched in and out to achieve specific measurement objectives. For example, addition of the thermocouple input option enables continuous measurement to 1000 K and above. Alternatively, the capacitance sensor option card enables a magneticsimpervious capacitance temperature sensor to be temporarily switched in for elimination of magneto-resistive effects while taking low temperature sample measurements under high or changing fields. Diode sensor support is provided by the 4-channel expansion card, which also enables use of additional Cernox® sensors for supplemental monitoring.



#### 4 PID controlled outputs

For convenient integration into a wide range of systems, the Model 350 offers four PID-controlled outputs. Variable DC current source outputs include a 75 W output for direct control of the typical main warm-up heater, and a 1 W output for fine control of the sample heater. Two additional 1 W variable DC voltage source outputs can be used to power auxiliary devices like a still heater in a dilution refrigerator, or to control a magnet power supply driving an ADR. The ability to dynamically select an input to associate with the controlled output provides additional flexibility in setting up the control scheme.

#### **Precision temperature control**

The Model 350 calculates the precise control output based on your temperature setpoint and feedback from the control sensor. You can manually set the PID values for fine control, or the temperature control loop autotuning feature can automate the tuning process for you. The setpoint ramp feature provides smooth, continuous setpoint changes and predictable setpoint approaches without the worry of overshoot or excessive settling times. When combined with the zone setting feature, which enables automatic switching of sensor inputs and scales current excitation through ten different preloaded temperature zones, the Model 350 provides continuous measurement and control over the entire temperature range required.

#### Simple and increased productivity

With remote control and automated features. the Model 350 will simplify your temperature control processes and increase your productivity in the laboratory.

Model 350 Temperature Controller

#### 3 interfaces for remote control

The Model 350 temperature controller includes Ethernet, USB, and IEEE-488 interfaces. In addition to gathering data, nearly every function of the instrument can be controlled through a computer interface. Ethernet provides the ability to access and monitor instrument activities via the internet from anywhere in the world, allowing distributed sharing of the controller and the controlled system. You can download the Lake Shore curve handler software to your computer to easily enter and manipulate sensor calibration curves for storage in the instrument's non-volatile flash memory.

#### Simple automation

Each sensor input has a high and low alarm that offer latching and non-latching operation. The two relays can be used in conjunction with the alarms to alert you of a fault condition and perform simple on/off control. Relays can be assigned to any alarm or operated manually. Choosing appropriate PID control settings for a closed loop system can be tedious, but the Model 350 provides the temperature control loop autotuning feature to simplify the process. It's an automated


process that measures system characteristics and computes setting values for P, I, and D for you. Once PID tuning parameters are chosen for a given setpoint, the zone tuning feature automatically switches sensor inputs for new setpoints, enabling you to control temperatures from 100 mK to over 1000 K without interrupting your experiment.

#### Performance you can count on

As with all Lake Shore products, the Model 350 product specifications are documented and verifiable in keeping with Lake Shore's tradition of performance assurance even at application extremes. The product is supported by a 3-year standard warranty, our confirmation of quality and commitment for the long term. Choosing the Model 350 for your ultra-low temperature application means you'll have the ultimate confidence in meeting your integration, measurement and control needs, now and into the future.

#### Use additional input types with option cards

The field installable input option cards add additional input types. The Model 3060 adds thermocouple capability. The Model 3061 adds capacitance sensor inputs. The Model 3062 adds 4 Cernox®/diode inputs. While the option cards can be easily removed, it is not necessary as the standard inputs remain functional when the options are not being used.

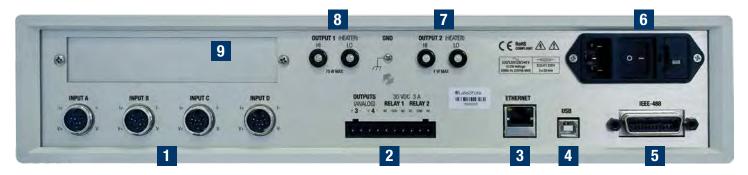

#### Model 350 rear panel

- Sensor inputs
- Terminal block (analog output and relays)
- Ethernet interface

- 4 LISR interface
- 5 IEEE-488 interface
- Line input assembly
- Output 2 heater
- Output 1 heater
- Option card slot



#### **Configurable display**

The Model 350 offers a bright, graphic liquid crystal display with an LED backlight that simultaneously displays up to eight readings. You can show all four loops, all inputs, or if you need to monitor one input, you can display just that one in greater detail. Or you can custom configure each display location to suit your experiment. Data from any input can be assigned to any of the locations, and your choice of temperature or sensor units can be displayed. For added convenience, you can also custom label each sensor input, eliminating the guesswork in remembering or determining the location to which a sensor input is associated.



Four input/output display with labels Standard display option featuring all four inputs and associated outputs.



Intuitive menu structure
Logical navigation allows you to spend more time on research and less time on setup.

# # Sample Space # Rad Shield 27.8645 K Li # Sample Space P 50.0 MOut Setp 273.000 K I 20.0 MOut Heat 0.00% of Low I 0.0 0.0 0.00%

Two input/output display with labels Reading locations can be user configured to meet application needs. Here, the input name is shown above each measurement reading along with the designated input letter.

#### Sensor selection

#### Sensor temperature range (sensors sold separately)

		Model	Useful range	Magnetic field use
Negative	Cernox®	CX-1010-HT	0.1 K to 420 K <sup>1, 2</sup>	$T > 2 K \& B \le 19 T$
temperature	Cernox®	CX-1030-HT	0.3 K to 420 K <sup>1, 2</sup>	$T > 2 K \& B \le 19 T$
coefficient RTDs	Cernox®	CX-1050-HT	1.4 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1070-HT	4 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1080-HT	20 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
	Germanium	GR-300-AA	0.3 K to 100 K	Not recommended
	Germanium	GR-1400-AA	1.4 K to 100 K	Not recommended
	Rox™	RX-102B	0.1 K to 40 K <sup>2</sup>	$T > 2 K \& B \le 10 T$
	Rox™	RX-103	1.4 K to 40 K	$T > 2 K \& B \le 10 T$
	Rox™	RX-202	0.1 K to 40 K <sup>2</sup>	$T > 2 K \& B \le 10 T$
Positive temperature	100 Ω platinum	PT-102/3	14 K to 873 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
coefficient RTDs	100 Ω platinum	PT-111	14 K to 673 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
	Rhodium-iron	RF-800-4	1.4 K to 500 K	$T > 77 \text{ K \& B} \le 8 \text{ T}$
Diodes	Silicon diode	DT-670-SD	1.4 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
Option - 3062	Silicon diode	DT-670E-BR	30 K to 500 K	T≥60 K & B≤3 T
	Silicon diode	DT-414	1.4 K to 375 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-421	1.4 K to 325 K	T≥60 K & B≤3 T
	Silicon diode	DT-470-SD	1.4 K to 500 K	T≥60 K & B≤3 T
	Silicon diode	DT-471-SD	10 K to 500 K	T ≥ 60 K & B ≤ 3 T
	GaAlAs diode	TG-120-P	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-PL	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-SD	1.4 K to 500 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
Capacitance Option—3061		CS-501	1.4 K to 290 K	$T > 4.2 \text{ K \& B} \le 18.7 \text{ T}$
Thermocouples	Type K	9006-006	3.2 K to 1505 K	Not recommended
Option-3060	Type E	9006-004	3.2 K to 934 K	Not recommended
	Chromel- AuFe 0.07%	9006-002	1.2 K to 610 K	Not recommended

<sup>&</sup>lt;sup>1</sup> Non-HT version maximum temperature: 325 K

**Cernox®** thin-film RTDs offer high sensitivity and low magnetic field-induced errors over the 0.1 K to 420 K temperature range. Cernox sensors require calibration.

**Platinum RTDs** offer high uniform sensitivity from 30 K to over 800 K. With excellent reproducibility, they are useful as thermometry standards. They follow a standard curve above 70 K and are interchangeable in many applications.

**Silicon diodes** are the best choice for general cryogenic use from 1.4 K to above room temperature. Silicon diodes are economical to use because they follow a standard curve and are interchangeable in many applications. They are not suitable for use in ionizing radiation or magnetic fields.

**Capacitance** sensors are ideally suited for use in strong magnetic fields because they exhibit virtually no magnetic field dependence. They can be used from 1.4 K to 290 K.

Low temperature specified with self-heating error:  $\leq 5$  mK



#### Typical sensor performance—see Appendix F for sample calculations of typical sensor performance

	Example Lake Shore sensor	Temperature (K)	Nominal resistance/ voltage	Typical sensor sensitivity <sup>3</sup>	Measurement resolution: temperature equivalents	Electronic accuracy: temperature equivalents	Temperature accuracy including electronic accuracy, CalCurve™, and calibrated sensor	Electronic control stability <sup>4</sup> : temperature equivalents
Cernox® (1 mV)	CX-1010-SD with 0.1L calibration	0.1 0.3 0.5 4.2 300	21389 Ω 2322.4 Ω 1248.2 Ω 277.32 Ω 30.392 Ω	-558110 Ω/K -10785 Ω/K -2665.2 Ω/K -32.209 Ω/K -0.0654 Ω/K	5.4 μK 28 μK 113 μK 931 μK 153 mK	±69 µK ±272 µK ±938 µK ±6.5 mK ±1.7 K	±3.1 mK ±3.8 mK ±5.4 mK ±11.5 mK ±1.8 K	±10.8 µK ±56.0 µK ±225 µK ±1.9 mK ±306 mK
Cernox® (10 mV)	CX-1050-SD-HT <sup>5</sup> with 1.4M calibration	1.4 4.2 77 420	26566 Ω 3507.2 Ω 205.67 Ω 45.03 Ω	-48449 Ω/K -1120.8 Ω/K -2.4116 Ω/K -0.0829 Ω/K	6.2 μK 89 μK 1.2 mK 12 mK	±261 μK ±2.1 mK ±38 mK ±338 mK	±5.3 mK ±7.1 mK ±54 mK ±378 mK	±12.4 μK ±178 μK ±2.4 mK ±24 mK
Germanium (1 mV)	GR-50-AA with 0.05A calibration	0.1 0.3 0.5 1.4 4.2	2317 Ω 164 Ω 73.8 Ω 24.7 Ω 13.7 Ω	-71858 Ω/K -964 Ω/K -202.9 Ω/K -13.15 Ω/K -1.036 Ω/K	4.2 μK 31.1 μK 49.3 μK 228 μK 2.9 mK	±14 µK ±78 µK ±195 µK ±904 µK ±7.2 mK	±3.2 mK ±3.8 mK ±4.5 mK ±4.9 mK ±11 mK	±8.4 μK ±62.2 μK ±98.6 μK ±456 μK ±5.8 mK
Germanium (10 mV)	GR-300-AA with 0.3D calibration	0.3 1.4 4.2 100	35180 Ω 448.6 Ω 94.46 Ω 2.72 Ω	-512200 Ω/K -581.3 Ω/K -26.56 Ω/K -0.024 Ω/K	2 μK 17 μK 38 μK 4.2 mK	±47 μK ±481 μK ±1.8 mK ±151 mK	±3.7 mK ±4.5 mK ±5.8 mK ±181 mK	±4.0 μK ±34 μK ±76 μK ±8.4 mK
Germanium (10 mV)	GR-1400-AA with 1.4D calibration	1.4 4.2 77 100	35890 Ω 1689 Ω 3.55 Ω 2.8 Ω	-94790 Ω/K -861.9 Ω/K -0.05 Ω/K -0.021 Ω/K	11 μK 35 μK 2 mK 4.8 mK	±257 μK ±900 μK ±83 mK ±175 mK	±4.3 mK ±4.9 mK ±99 mK ±191 mK	±21.1 μK ±69.6 μK ±4 mK ±9.5 mK
Rox <sup>™</sup> (1 mV)	RX-102B-CB with 0.02C calibration	0.1 0.5 1.4 4.2 40	3549 Ω 2188 Ω 1779 Ω 1546 Ω 1199 Ω	-12578 Ω/K -1056 Ω/K -198 Ω/K -40.0 Ω/K -3.41 Ω/K	79.5 μK 284 μK 1.5 mK 7.5 mK 88 mK	±908 µK ±2.7 mK ±13.7 mK ±65.4 mK ±727 mK	±3.8 mK ±5.7 mK ±18.7 mK ±81.4 mK ±764 mK	±159 µK ±568 µK ±3.0 mK ±15.0 mK ±176 mK
Platinuim RTD 500 Ω full scale	PT-103 with 14J calibration	30 77 300 500	3.66 Ω 20.38 Ω 110.35 Ω 185.668 Ω	0.191 Ω/K 0.423 Ω/K 0.387 Ω/K 0.378 Ω/K	0.5 mK 0.7 mK 7.8 mK 7.9 mK	±22 mK ±34 mK ±140 mK ±223 mK	$\pm 32 \text{ mK}  \pm 46 \text{ mK}  \pm 163 \text{ mK}  \pm 269 \text{ mK}$	±1.0 mK ±1.4 mK ±15.6 mK ±15.8 mK
Silicon diode	DT-670-CO-13 with 1.4H calibration	1.4 77 300 500	1.664 V 1.028 V 0.5596 V 0.0907 V	-12.49 mV/K -1.73 mV/K -2.3 mV/K -2.12 mV/K	0.8 mK 5.8 mK 4.3 mK 4.7 mK	±13 mK ±76 mK ±47 mK ±40 mK	±25 mK ±98 mK ±79 mK ±90 mK	±1.6 mK ±11.6 mK ±8.7 mK ±9.4 mK
Silicon diode	DT-470-SD-13 with 1.4H calibration	1.4 77 300 475	1.6981 V 1.0203 V 0.5189 V 0.0906 V	-13.1 mV/K -1.92 mV/K -2.4 mV/K -2.22 mV/K	0.8 mK 5.2 mK 4.2 mK 4.5 mK	±13 mK ±68 mK ±44 mK ±38 mK	±25 mK ±90 mK ±76 mK ±88 mK	±1.6 mK ±10.4 mK ±8.4 mK ±9.0 mK
GaAIAs diode	TG-120-SD with 1.4H calibration	1.4 77 300 475	5.3909 V 1.4222 V 0.8978 V 0.3778 V	-97.5 mV/K -1.24 mV/K -2.85 mV/K -3.15 mV/K	0.21 mK 16 mK 7 mK 6.3 mK	±8.8 mK ±373 mK ±144 mK ±114 mK	±21 mK ±395 mK ±176 mK ±164 mK	±410 µK ±32.3 mK ±14.0 mK ±12.6 mK
Thermocouple 50 mV Option — 3060	Type K	75 300 600 1505	-5862.9 μV 1075.3 μV 13325 μV 49998.3 μV	15.6 μV/K 40.6 μV/K 41.7 μV/K 36.0 μV/K	26 mK 9.9 mK 9.6 mK 11 mK	±252 mK <sup>6</sup> ±38 mK <sup>6</sup> ±184 mK <sup>6</sup> ±730 mK <sup>6</sup>	Calibration not available from Lake Shore	±52 mK ±19.6 mK ±19.2 mK ±22.2 mK
Capacitance Option — 3061	CS-501	4.2 77 200	6.0 nF 9.1 nF 19.2 nF	27 pF/K 52 pF/K 174 pF/K	1.9 mK 1.0 mK 2.9 mK	Not applicable	Calibration not available from Lake Shore	±3.8 mK ±2.0 mK ±5.8 mK

Model 350 Temperature Controller

<sup>&</sup>lt;sup>3</sup> Typical sensor sensitivities were taken from representative calibrations for the sensor listed

<sup>&</sup>lt;sup>4</sup> Control stability of the electronics only, in an ideal thermal system

<sup>&</sup>lt;sup>5</sup> Non-HT version maximum temperature: 325 K

<sup>&</sup>lt;sup>6</sup> Accuracy specification does not include errors from room temperature compensation

#### **Input specifications**

Standard inputs	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution <sup>7</sup>	Electronic accuracy (at 25 °C)	Measurement temperature coefficient	Electronic control stability <sup>8</sup>
NTC RTD/	Negative/	0 $\Omega$ to 10 $\Omega$	1 mA <sup>10</sup>	$0.1~\text{m}\Omega$	0.1 mΩ	$\pm 0.002~\Omega~\pm 0.06\%$ of rdg	$(0.01 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	$\pm 0.2~\text{m}\Omega$
PTC RTD	Positive	0 Ω to 30 Ω	300 μA <sup>10</sup>	0.1 mΩ	0.3 mΩ	$\pm 0.002~\Omega~\pm 0.06\%$ of rdg	$(0.03 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±0.6 mΩ
10 mV		0 Ω to 100 Ω	100 μA <sup>10</sup>	1 mΩ	1 mΩ	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	$(0.1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±2 mΩ
		0 Ω to 300 Ω	30 μA <sup>10</sup>	1 mΩ	3 mΩ	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	$(0.3 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±6 mΩ
		0 Ω to 1 kΩ	10 μA <sup>10</sup>	10 mΩ	10 mΩ	$\pm 0.1 \Omega \pm 0.04\%$ of rdg	$(1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±20 mΩ
		0 Ω to 3 kΩ	3 μA <sup>10</sup>	10 mΩ	30 mΩ	$\pm 0.1 \Omega \pm 0.04\%$ of rdg	$(3 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±60 mΩ
		$0~\Omega$ to $10~k\Omega$	1 μA <sup>10</sup>	$100~\text{m}\Omega$	100 mΩ	$\pm 1.0 \Omega \pm 0.04\%$ of rdg	$(10 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±200 mΩ
		0 Ω to 30 kΩ	300 nA <sup>10</sup>	100 mΩ	300 mΩ	$\pm 2.0 \Omega \pm 0.04\%$ of rdg	$(30 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±600 mΩ
		0 Ω to 100 kΩ	100 nA <sup>10</sup>	1 Ω	1 Ω	$\pm 10.0 \Omega \pm 0.04\%$ of rdg	$(100 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±2 Ω
		0 $\Omega$ to 300 $k\Omega$	30 nA <sup>10</sup>	1 Ω	3 Ω	$\pm 30 \Omega \pm 0.04\%$ of rdg	$(300 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±6 Ω
NTC RTD	Negative	0 Ω to 10 Ω	100 μA <sup>10</sup>	$0.1~\text{m}\Omega$	1 mΩ	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	$(0.1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±2 mΩ
1 mV		0 Ω to 30 Ω	30 μA <sup>10</sup>	$0.1~\text{m}\Omega$	3 mΩ	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	$(0.3 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±6 mΩ
		0 $\Omega$ to 100 $\Omega$	10 μA <sup>10</sup>	1 mΩ	10 mΩ	$\pm 0.1 \Omega \pm 0.04\%$ of rdg	$(1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	$\pm 20~\text{m}\Omega$
		0 Ω to 300 Ω	3 μA <sup>10</sup>	1 mΩ	30 mΩ	$\pm 0.1 \Omega \pm 0.04\%$ of rdg	$(3 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±60 mΩ
		0 Ω to 1 kΩ	1 μA <sup>10</sup>	10 mΩ	100 mΩ	$\pm 1.0 \Omega \pm 0.04\%$ of rdg	$(10 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±200 mΩ
		0 Ω to 3 kΩ	300 nA <sup>10</sup>	10 mΩ	300 mΩ	$\pm 2.0 \Omega \pm 0.04\%$ of rdg	$(30 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±600 mΩ
		0 Ω to 10 kΩ	100 nA <sup>10</sup>	100 mΩ	1 Ω	$\pm 10.0 \Omega \pm 0.04\%$ of rdg	$(100 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±2 Ω
		0 $\Omega$ to 30 k $\Omega$	30 nA <sup>10</sup>	100 mΩ	3 Ω	$\pm 30~\Omega~\pm 0.04\%$ of rdg	$(300 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±6 Ω
		0 $\Omega$ to 100 $k\Omega$	10 nA <sup>10</sup>	1 Ω	10 Ω	$\pm 100~\Omega~\pm 0.04\%$ of rdg	$(1 \Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±20 Ω

Scanner option Model 3062	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy (at 25 °C)	Measurement temperature coefficient	Electronic control stability <sup>a</sup>
Diode	Negative	0 V to 2.5 V	10 μA ±0.05% <sup>9</sup>	10 μV	10 μV	±80 μV ±0.005% of rdg	(10 μV + 0.0005% of rdg)/°C	±20 μV
	Negative	0 V to 10 V	10 μA ±0.05% <sup>9</sup>	100 μV	20 μV	±160 μV ±0.01% of rdg	(20 μV + 0.0005% of rdg)/°C	±40 μV
PTC RTD	Positive	0 $\Omega$ to 10 $\Omega$	1 mA <sup>10</sup>	$0.1~\text{m}\Omega$	$0.2~\text{m}\Omega$	$\pm 0.002 \Omega \pm 0.01\%$ of rdg	$(0.01 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	$\pm 0.2~\text{m}\Omega$
		0 Ω to 30 Ω	1 mA <sup>10</sup>	0.1 mΩ	0.2 mΩ	$\pm 0.002 \Omega \pm 0.01\%$ of rdg	$(0.03 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±0.4 mΩ
		0 Ω to 100 Ω	1 mA <sup>10</sup>	1 mΩ	2 mΩ	$\pm 0.004 \Omega \pm 0.01\%$ of rdg	$(0.1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±4 mΩ
		0 Ω to 300 Ω	1 mA <sup>10</sup>	1 mΩ	2 mΩ	±0.004 Ω ±0.01% of rdg	$(0.3 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±4 mΩ
		0 Ω to 1 kΩ	1 mA <sup>10</sup>	10 mΩ	20 mΩ	±0.04 Ω ±0.02% of rdg	$(1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±40 mΩ
		0 Ω to 3 kΩ	1 mA <sup>10</sup>	10 mΩ	20 mΩ	±0.04 Ω ±0.02% of rdg	$(3 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±40 mΩ
		0 Ω to 10 kΩ	1 mA <sup>10</sup>	100 mΩ	200 mΩ	$\pm 0.4 \Omega \pm 0.02\%$ of rdg	$(10 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±400 mΩ
NTC RTD	Negative	0 $\Omega$ to 10 $\Omega$	1 mA <sup>10</sup>	$0.1~\text{m}\Omega$	0.15 mΩ	$\pm 0.002~\Omega~\pm 0.06\%$ of rdg	$(0.01 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	$\pm 0.3~\text{m}\Omega$
10 mV		0 Ω to 30 Ω	300 μA <sup>10</sup>	0.1 mΩ	0.45 mΩ	$\pm 0.002~\Omega~\pm 0.06\%$ of rdg	$(0.03 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$	±0.9 mΩ
		0 $\Omega$ to 100 $\Omega$	100 μA <sup>10</sup>	1 mΩ	1.5 mΩ	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	$(0.1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±3 mΩ
		0 $\Omega$ to 300 $\Omega$	30 μA <sup>10</sup>	1 mΩ	4.5 mΩ	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	$(0.3 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$	±9 mΩ
		$0~\Omega$ to $1~k\Omega$	10 μA <sup>10</sup>	10 mΩ	15 m $\Omega$ +0.002% of rdg	$\pm 0.1~\Omega~\pm 0.04\%$ of rdg	$(1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	$\pm 30$ m $\Omega$ $\pm 0.004\%$ of rdg
		0 Ω to 3 kΩ	3 μA <sup>10</sup>	10 mΩ	45 mΩ+0.002% of rdg	$\pm 0.1 \Omega \pm 0.04\%$ of rdg	$(3 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$	$\pm 90 \text{ m}\Omega \pm 0.004\%$ of rdg
		0 Ω to 10 kΩ	1 µA <sup>10</sup>	100 mΩ	150 mΩ+0.002% of rdg	±1.0 Ω ±0.04% of rdg	$(10 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±300 mΩ ±0.004% of rdg
		0 Ω to 30 kΩ	300 nA <sup>10</sup>	100 mΩ	450 mΩ+0.002% of rdg	±2.0 Ω ±0.04% of rdg	$(30 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±900 mΩ ±0.004% of rdg
		0 $\Omega$ to 100 $k\Omega$	100 nA <sup>10</sup>	1 Ω	1.5 Ω+0.005% of rdg	$\pm 10.0~\Omega~\pm 0.04\%$ of rdg	$(100 \text{ m}\Omega + 0.002\% \text{ of rdg})/^{\circ}\text{C}$	$\pm 3~\Omega \pm 0.01\%$ of rdg
		1	1			- · · · · · · · · · · · · · · · · · · ·		

Thermocouple option Model 3060	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy (at 25 °C)	Measurement temperature coefficient	Electronic control stability <sup>8</sup>
Thermocouple	Positive	±50 mV	NA	0.1 μV	0.4 μV	$\pm 1~\mu V~\pm 0.05\%~of~rdg^{11}$	$(0.1  \mu V + 0.001\% \text{ of rdg})/^{\circ}C$	±0.8 μV

Capacitance option Model 3061	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy (at 25 °C)	Measurement temperature coefficient	Electronic control stability <sup>8</sup>
Capacitance	Positive or	0.1 to 15 nF	3.496 kHz 1 mA square wave	0.1 pF	0.05 pF	$\pm 50$ pF $\pm 0.4\%$ of rdg	2.5 pF/°C	0.1 pF
	Negative	1 to 150 nF	3.496 kHz 10 mA square wave	1 pF	0.5 pF	$\pm 50 \text{ pF} \pm 0.4\% \text{ of rda}$	5 pF/°C	1 pF

<sup>&</sup>lt;sup>7</sup> Measurement resolution measured at 4.2 K to remove the thermal noise of the resistor

<sup>&</sup>lt;sup>8</sup> Control stability of the electronics only, in ideal thermal system

<sup>&</sup>lt;sup>9</sup> Current source error has negligible effect on measurement accuracy

<sup>&</sup>lt;sup>10</sup> Current source error is removed during calibration

<sup>&</sup>lt;sup>11</sup> Accuracy specification does not include errors from room temperature compensation

#### **Thermometry**

Number of inputs 4 (8 with scanner option)

Input configuration Inputs can be configured from the front panel to accept any of the supported input types. Thermocouple, capacitance and diode inputs require an optional input card that can be installed in the field.

Isolation Sensor inputs optically isolated from other circuits but not each other  ${\bf A/D}$  resolution 24-bit

Input accuracy Sensor dependent, refer to Input Specifications table

Measurement resolution Sensor dependent, refer to Input Specifications table

Maximum update rate 10 rdg/s on each non-scanned input

**Maximum update rate (scanner)** The maximum update rate for a scanned input is 10 rdg/s distributed among the enabled channels. Any channel configured as 100 k $\Omega$  RTD with reversal on changes the update rate for the channel to 5 rdg/s.

Scanner channels enabled*	Update rate
1	10 rdg/s (100 ms/rdg)
2	5 rdg/s (200 ms/rdg)
3	31/3 rdg/s (300 ms/rdg)
4	2½ rdg/s (400 ms/rdg)
5	2 rda/s (500 ms/rda)

<sup>\*</sup> No channels configured for 100 k $\Omega$  NTC RTD

Autorange Automatically selects appropriate NTC RTD or PTC RTD range

User curves Room for 39 200-point CalCurves™ or user curves

SoftCal™ Improves accuracy of DT-470 diode to ±0.25 K from 30 K to 375 K; improves accuracy of platinum RTDs to ±0.25 K from 70 K to 325 K; stored as user curves

Math Maximum and minimum

Filter Averages 2 to 64 input readings

#### **Control**

#### **Control outputs** 4

#### Heater outputs (Outputs 1 & 2)

Control type Closed loop digital PID with manual heater output or open loop  ${\bf Update\ rate\ 10/s}$ 

Tuning Autotune (one loop at a time), PID, PID zones

Control stability Sensor dependent, see Input Specifications table

#### PID control settings

Proportional (gain) 0 to 9999 with 0.1 setting resolution

Integral (reset) 1 to 1000 (1000/s) with 0.1 setting resolution

Derivative (rate) 1 to 200% with 1% resolution

Manual output 0 to 100% with 0.01% setting resolution

Zone control 10 temperature zones with P, I, D, manual heater out, heater range, control channel, ramp rate

Setpoint ramping 0.001 K/min to 100 K/min

#### Analog outputs (Outputs 3 & 4)

Control type Closed loop PID, PID zones, warm up heater mode, still heater, manual output, or monitor output

#### Warm up heater mode settings

Warm up percentage 0 to 100% with 1% resolution Warm up mode Continuous control or auto-off

#### **Monitor output settings**

Scale User selected

Data source Temperature or sensor units

Settings Input, source, top of scale, bottom of scale, or manual

Type Variable DC voltage source

Update rate 10/s Range ±10 V

Resolution 16-bit, 0.3 mV

Accuracy ±2.5 mV

Noise 0.3 mV RMS

Maximum current 100 mA

Maximum power 1 W (into 100  $\Omega$ )

Minimum load resistance 100  $\Omega$  (short-circuit protected)

**Connector** Detachable terminal block

#### **Output 1**

	25 Ω setting	50 Ω setting				
Туре	Variable DC c	urrent source				
D/A resolution	16-	bit				
Max power	75 W	50 W				
Max current	1.732 A	1 A				
Voltage compliance (min)	50 V	50 V				
Heater load for max power	25 Ω	50 Ω				
Heater load range	10 Ω to 100 Ω					
Ranges	5 (decade ste	eps in power)				
Heater noise	1.2 µA RMS (dominated by line frequency and its harmonics)					
Grounding	Output referenced to chassis ground					
Heater connector	Dual b	anana				
Safety limits	Curve temperature, power up he	eater off, short circuit protection				

#### **Output 2**

Туре	Variable DC current source			
D/A resolution	16-bit			
Max power	1 W			
Max current	100 mA			
Voltage compliance	10 V			
(min)	10 ν			
Heater load for max	100 Ω			
power	100 Ω			
Heater load range	25 Ω to 2 kΩ			
Ranges (100 Ω load)	1 W, 100 mW, 10 mW, 1 mW, 100 μW			
Heater noise	<0.005% of range			
Grounding	Output referenced to measurement common			
Heater connector	Dual banana			
Safety limits	Curve temperature, power up heater off, short circuit protection			

#### **Sensor input configuration**

	RTD	Diode (option)	Thermocouple (option)	Capacitance (option)
Measurement	4-lead	4-lead	2-lead differential,	4-lead
type	differential	differential	room temperature	differential,
			compensated	variable duty
				cycle
Excitation	Constant current	10 μΑ	N/A	Constant
	with current	constant		current,
	reversal	current		3.496 kHz
				square wave
Supported	100 Ω Platinum,	Silicon,	Most thermocouple	CS-501GR
sensors	1000 Ω	GaAlAs	types	
	Platinum,			
	Germanium,			
	Carbon-Glass,			
	Cernox®, and			
	Rox™			
Standard	PT-100,	DT-470,	Type E, Type K, Type T,	N/A
curves	PT-1000,	DT-670,	AuFe 0.07% vs. Cr,	
	RX-102A,	DT-500-D,	AuFe 0.03% vs Cr	
	RX-202A	DT-500-E1		
Input connector	6-pin DIN	6-pin DIN	Screw terminals in a	6-pin DIN
			ceramic isothermal block	



#### **Front panel**

Display 8-line by 40-character (240 × 64 pixel) graphic LCD display module with LED backlight

Number of reading displays 1 to 8 Display units K, °C, V, mV, Ω, nF

Reading source Temperature, sensor units, max, and min

Display update rate 2 rdg/s

Temperature display resolution 0.00001° from 0° to 9.99999°, 0.0001° from 10° to

99.9999°, 0.001° from 100° to 999.999°, 0.01° above 1000° Sensor units display resolution Sensor dependent, to 6 digits

Other displays Input name, setpoint, heater range, heater output, and PID Setpoint setting resolution Same as display resolution (actual resolution is sensor dependent)

Heater output display Numeric display in percent of full scale for power or current Heater output resolution 0.01%

Display annunciators Control input, alarm, tuning

**LED annunciators** Remote, Ethernet status, alarm, control outputs

Keypad 27-key silicone elastomer keypad

Front panel features Front panel curve entry, display contrast control, and keypad lock-out

#### **Interface**

#### IEEE-488.2

SH1, AH1, T5, L4, SR1, RL1, PP0, DC1, DT0, C0, E1 Capabilities

Reading rate To 10 rdg/s on each input

Software support LabVIEW™ driver (see www.lakeshore.com)

**USB** 

Emulates a standard RS-232 serial port Function

Baud rate 57,600

Connector B-type USB connector To 10 rdg/s on each input Reading rate

Software support LabVIEW™ driver (see www.lakeshore.com)

**Ethernet** 

Function TCP/IP, web interface, curve handler, configuration backup,

chart recorder **RJ-45** Connector

Reading rate To 10 rdg/s on each input

LabVIEW<sup>™</sup> driver (see www.lakeshore.com) Software support

**Alarms** 

Number 4 (8 with scanner option), high and low for each input

Data source Temperature or sensor units

Settings Source, high setpoint, low setpoint, deadband,

latching or non-latching, audible on/off, and visible on/off

Display annunciator, beeper, and relays Actuators

Relays

Number

Contacts Normally open (NO), normally closed (NC), and common (C)

Contact rating 30 VDC at 3 A

Activate relays on high, low, or both alarms for any input, Operation

or manual mode

Connector Detachable terminal block

#### General

Ambient temperature 15 °C to 35 °C at rated accuracy; 5 °C to 40 °C at reduced accuracy

Power requirement 100, 120, 220, 240 VAC, ±10%, 50 or 60 Hz, 220 VA Size 435 mm W  $\times$  89 mm H  $\times$  368 mm D (17 in  $\times$  3.5 in  $\times$  14.5 in), full rack

Weight 7.6 kg (16.8 lb) Approval CE mark, RoHS

## Ordering information

#### Part number **Description**

350 2 diode/resistor inputs temperature controller, includes one

dual banana jack heater output connector, four 6-pin DIN plug sensor input mating connectors, one 10-pin terminal block, a

calibration certificate and a user's manual Model 350 with a 3060 option card installed Model 350 with a 3061 option card installed

350-3062 Model 350 with a 3062 option card installed 3060 2-thermocouple input option for 350/336, field-installable 3061 Capacitance input option for 350/336, field-installable 4-channel scanner option for diodes and RTD sensors for

350/336, field-installable

#### Please indicate your power/cord configuration:

1 100 V-U.S. cord (NEMA 5-15)

2 120 V-U.S. cord (NEMA 5-15) 3 220 V—Euro cord (CEE 7/7)

4 240 V—Euro cord (CEE 7/7)

5 240 V—U.K. cord (BS 1363)

240 V-Swiss cord (SEV 1011)

**7** 220 V—China cord (GB 1002)

#### **Accessories**

350-3060

350-3061

112-177 Temperature controller cable, 3 m (10 ft)—IN STOCK

112-178 Temperature controller cable, 6 m (20 ft) 112-180 Temperature controller cable, 10 m (33 ft)

117-017 1 m (3.3 ft long) IEEE-488 (GPIB) computer interface cable

assembly

RM-1 Rack mount kit for mounting one full rack temperature

instrument

Sensor input mating connector (6-pin DIN plug) G-106-233

G-106-755 Terminal block, 10-pin 106-009 Banana plug, dual

CAL-350-CERT Instrument calibration with certificate

CAL-350-DATA Instrument recalibration with certificate and data 119-057 Model 350 temperature controller manual

All specifications are subject to change without notice











## Model 336 Temperature Controller



Model 336 Temperature Controller

#### **Model 336 features**

- Operates down to 300 mK with appropriate NTC RTD sensors
- Four sensor inputs and four independent control outputs
- Two PID control loops: 100 W and 50 W into a 50  $\Omega$  or 25  $\Omega$  load
- Autotuning automatically calculates PID parameters
- Automatically switch sensor inputs using zones to allow continuous measurement and control from 300 mK to 1505 K
- Custom display setup allows you to label each sensor input
- Ethernet, USB and IEEE-488 interfaces
- Supports diode, RTD, and thermocouple temperature sensors
- Sensor excitation current reversal eliminates thermal EMF errors for resistance sensors
- ±10 V analog voltage outputs, alarms, and relays
- CE certification
- Full 3 year standard warranty





#### Introduction

The first of a new generation of innovative temperature measurement and control solutions by Lake Shore, the Model 336 temperature controller comes standard equipped with many advanced features promised to deliver the functionality and reliable service you've come to expect from the world leader in cryogenic thermometry. The Model 336 is the only temperature controller available with four sensor inputs, four control outputs and 150 W of low noise heater power. Two independent heater outputs providing 100 W and 50 W can be associated with any of the four sensor inputs and programmed for closed loop temperature control in proportional-integral-derivative (PID) mode. The improved autotuning feature of the Model 336 can be used to automatically calculate PID parameters, so you spend less time tuning your controller and more time conducting experiments.



The Model 336 supports the industry's most advanced line of cryogenic temperature sensors as manufactured by Lake Shore, including diodes, resistance temperature detectors (RTDs) and thermocouples. The controller's zone tuning feature allows you to measure and control temperatures seamlessly from 300 mK to over 1,500 K by automatically switching temperature sensor inputs when your temperature range goes beyond the usable range of a given sensor. You'll never again have to be concerned with temperature sensor over or under errors and measurement continuity issues. Alarms, relays, and  $\pm 10$  V analog voltage outputs are available to help automate secondary control functions.

Another innovative first from Lake Shore, the ability to custom label sensor inputs eliminates the guesswork in remembering or determining the location to which a sensor input is associated. As we strive to maintain increasingly demanding workloads, ease of use and the ability to stay connected from anywhere in the world are critical attributes. With standard Ethernet, USB, and IEEE-488 interfaces and an intuitive menu structure and logic, the Model 336 was designed with efficiency, reliable connectivity, and ease of use in mind. While you may need to leave your lab, Ethernet ensures you'll always be connected to your experiments. The new intuitive front panel layout and keypad logic, bright graphic display, and LED indicators enhance the user friendly front panel interface of the Model 336.

In many applications, the unparalleled feature set of the Model 336 allows you to replace several instruments with one, saving time, money and valuable laboratory space. Delivering more feedback, tighter control, and faster cycle times, the Model 336 keeps up with increasingly complex temperature measurement and control applications. It is the ideal solution for general purpose to advanced laboratory applications. Put the Model 336 temperature controller to use in your lab and let it take control of your measurement environment.

#### **Sensor inputs**

The Model 336 offers four standard sensor inputs that are compatible with diode and RTD temperature sensors. The field installable Model 3060 thermocouple input option provides support for up to two thermocouple inputs by adding thermocouple functionality to inputs C and D.

Sensor inputs feature a high-resolution 24-bit analog-to-digital converter; each input has its own current source, providing fast settling times. All four sensor inputs are optically isolated from other circuits to reduce noise and to provide repeatable sensor measurements. Current reversal eliminates thermal electromotive force (EMF) errors in resistance sensors. Nine excitation currents facilitate temperature measurement and control down to 300 mK using appropriate negative temperature coefficient (NTC) RTDs. Autorange mode automatically scales excitation current in NTC RTDs to reduce self heating at low temperatures as sensor resistance changes by many orders of magnitude. Temperatures down to 1.4 K can be measured and controlled using silicon or GaAlAs diodes. Software selects the appropriate excitation current and signal gain levels when the sensor type is entered via the instrument front panel. The unique zone setting feature automatically switches sensor inputs, enabling you to measure temperatures from 300 mK to over 1,500 K without interrupting your experiment.

The Model 336 includes standard temperature sensor response curves for silicon diodes, platinum RTDs, ruthenium oxide RTDs, and thermocouples. Non-volatile memory can also store up to 39 200-point CalCurves for Lake Shore calibrated temperature sensors or user curves. A built-in SoftCal algorithm can be used to generate curves for silicon diodes and platinum RTDs that can be stored as user curves. Temperature sensor calibration data can be easily uploaded and manipulated using the Lake Shore curve handler software.



#### **Temperature control**

Providing a total of 150 W of heater power, the Model 336 is the most powerful temperature controller available. Delivering very clean heater power, it precisely controls temperature throughout the full scale temperature range for excellent measurement reliability, efficiency, and throughput. Two independent PID control outputs supplying 100 W and 50 W of heater power can be associated with any of the four standard sensor inputs. Precise control output is calculated based on your temperature setpoint and feedback from the control sensor. Wide tuning parameters accommodate most cryogenic cooling systems and many high-temperature ovens commonly used in laboratories. PID values can be manually set for fine control, or the improved autotuning feature can automate the tuning process. Autotune calculates PID parameters and provides information to help build zone tables. The setpoint ramp feature provides smooth, continuous setpoint changes and predictable setpoint approaches without the worry of overshoot or excessive settling times. When combined with the zone setting feature, which enables automatic switching of sensor inputs and scales current excitation through ten different preloaded temperature zones, the Model 336 provides continuous measurement and control from 300 mK to 1505 K.

Control outputs 1 and 2 are variable DC current sources referenced to chassis ground. Output 1 can provide 100 W of continuous power to a 25  $\Omega$  load or 50 W to a 50  $\Omega$  or 25  $\Omega$  load. Output 2 provides 50 W to 25  $\Omega$ or 50  $\Omega$  heater loads. Outputs 3 and 4 are variable DC voltage source outputs providing two ±10 V analog outputs. When not in use to extend the temperature controller heater power, these outputs can function as manually controlled voltage sources.

Temperature limit settings for inputs are provided as a safeguard against system damage. Each input is assigned a temperature limit, and if any input exceeds that limit, all control channels are automatically disabled.

#### Interface

The Model 336 is standard equipped with Ethernet, universal serial bus (USB) and parallel (IEEE-488) interfaces. In addition to gathering data, nearly every function of the instrument can be controlled through a computer interface. You can download the Lake Shore curve handler software to your computer to easily enter and manipulate sensor calibration curves for storage in the instruments non-volatile memory.

Model 336 Temperature Controller



Four input/output display with labels

Standard display option featuring all four inputs and associated outputs.



#### Two input/output display with labels

Reading locations can be user configured to meet application needs. Here, the input name is shown above each measurement reading along with the designated input letter.



#### Intuitive menu structure

Logical navigation allows you to spend more time on research and less time on setup.

Ethernet provides the ability to access and monitor instrument activities via the internet from anywhere in the world. The USB interface emulates an RS-232C serial port at a fixed 57,600 baud rate, but with the physical connections of a USB. It also allows you to download firmware upgrades, ensuring the most current firmware version is loaded into your instrument without having to physically change anything.

Each sensor input has a high and low alarm that offer latching and non-latching operation. The two relays can be used in conjunction with the alarms to alert you of a fault condition and perform simple on/off control. Relays can be assigned to any alarm or operated manually.

The  $\pm 10$  V analog voltage outputs on outputs 3 and 4 can be configured to send a voltage proportional to temperature to a strip chart recorder or data acquisition system. You may select the scale and data sent to the output, including temperature or sensor units.

#### Model 336 rear panel



- Sensor input connectors
- Terminal block (analog outputs and relays)
- 3 Ethernet interface
- 4 LISR interface
- 5 IEEE-488 interface

- 6 Line input assembly
- 7 Output 2 heater
- 8 Output 1 heater
- 9 Thermocouple option inputs



#### **Configurable display**

The Model 336 offers a bright, graphic liquid crystal display with an LED backlight that simultaneously displays up to eight readings. You can show all four loops, or If you need to monitor one input, you can display just that one in greater detail. Or you can custom configure each display location to suit your experiment. Data from any input can be assigned to any of the locations, and your choice of temperature or sensor units can be displayed. For added convenience, you can also custom label each sensor input, eliminating the guesswork in remembering or determining the location to which a sensor input is associated.

#### Model 3060 thermocouple input option

The field installable Model 3060 thermocouple input option adds thermocouple functionality to inputs C and D. While the option can be easily removed, this is not necessary as the standard inputs remain fully functional when they are not being used to measure thermocouple temperature sensors. Calibration for the option is stored on the card so it can be installed in the field and used with multiple Model 336 temperature controllers without recalibration.

#### Sensor selection

#### Sensor temperature range (sensors sold separately)

		Model	Useful range	Magnetic field use
Diodes	Silicon diode	DT-670-SD	1.4 K to 500 K	T≥60 K & B < 3 T
Diougs	Silicon diode	DT-670E-BR	30 K to 500 K	T260 K & B ≤ 3 T
	Silicon diode	DT-414	1.4 K to 375 K	T260 K & B ≤ 3 T
	Silicon diode	DT-414	1.4 K to 375 K	
				T≥60 K & B ≤ 3 T
	Silicon diode	DT-470-SD	1.4 K to 500 K	T≥60 K & B ≤ 3 T
	Silicon diode	DT-471-SD	10 K to 500 K	T≥60 K & B ≤ 3 T
	GaAlAs diode	TG-120-P	1.4 K to 325 K	T > 4.2 K & B ≤ 5 T
	GaAlAs diode	TG-120-PL	1.4 K to 325 K	$T > 4.2 \text{ K & B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-SD	1.4 K to 500 K	$T > 4.2 \text{ K & B} \le 5 \text{ T}$
Positive temperature	100 Ω platinum	PT-102/3	14 K to 873 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
coefficient RTDs	100 Ω platinum	PT-111	14 K to 673 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
	Rhodium-iron	RF-800-4	1.4 K to 500 K	$T > 77 \text{ K \& B} \le 8 \text{ T}$
	Rhodium-iron	RF-100T/U	1.4 K to 325 K	$T > 77 K \& B \le 8 T$
Negative	Cernox®	CX-1010	0.3 K to 325 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
temperature	Cernox®	CX-1030-HT	0.3 K to 420 K <sup>1, 3</sup>	$T > 2 K \& B \le 19 T$
coefficient RTDs	Cernox®	CX-1050-HT	1.4 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1070-HT	4 K to 420 K <sup>1</sup>	T > 2 K & B ≤ 19 T
	Cernox®	CX-1080-HT	20 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
	Germanium	GR-300-AA	0.35 K to 100 K <sup>3</sup>	Not recommended
	Germanium	GR-1400-AA	1.8 K to 100 K <sup>3</sup>	Not recommended
	Carbon-glass	CGR-1-500	1.4 K to 325 K	T > 2 K & B ≤ 19 T
	Carbon-glass	CGR-1-1000	1.7 K to 325 K <sup>2</sup>	T > 2 K & B ≤ 19 T
	Carbon-glass	CGR-1-2000	2 K to 325 K <sup>2</sup>	T > 2 K & B ≤ 19 T
	Rox™	RX-102	0.3 K to 40 K <sup>3</sup>	T > 2 K & B ≤ 10 T
	Rox™	RX-103	1.4 K to 40 K	T > 2 K & B ≤ 10 T
	Rox™	RX-202	0.3 K to 40 K <sup>3</sup>	T > 2 K & B ≤ 10 T
Thermocouples	Type K	9006-006	3.2 K to 1505 K	Not recommended
Option—3060	Type E	9006-004	3.2 K to 934 K	Not recommended
,	Chromel- AuFe 0.07%	9006-002	1.2 K to 610 K	Not recommended

<sup>&</sup>lt;sup>1</sup> Non-HT version maximum temperature: 325 K

Silicon diodes are the best choice for general cryogenic use from 1.4 K to above room temperature. Silicon diodes are economical to use because they follow a standard curve and are interchangeable in many applications. They are not suitable for use in ionizing radiation or magnetic fields.

**Cernox®** thin-film RTDs offer high sensitivity and low magnetic field-induced errors over the 0.3 K to 420 K temperature range. Cernox sensors require calibration.

**Platinum RTDs** offer high uniform sensitivity from 30 K to over 800 K. With excellent reproducibility, they are useful as thermometry standards. They follow a standard curve above 70 K and are interchangeable in many applications.

<sup>&</sup>lt;sup>2</sup> Low temperature limited by input resistance range

<sup>&</sup>lt;sup>3</sup> Low temperature specified with self-heating error: ≤ 5 mK



#### **Typical sensor performance**

	Example Lake Shore sensor	Temperature	Nominal resistance/ voltage	Typical sensor sensitivity⁴	Measurement resolution: temperature equivalents	Electronic accuracy: temperature equivalents	Temperature accuracy including electronic accuracy, Calcurve™, and calibrated sensor	Electronic control stability <sup>5</sup> : temperature equivalents
Silicon diode	DT-670-C0-13	1.4 K	1.664 V	-12.49 mV/K	0.8 mK	±13 mK	±25 mK	±1.6 mK
	with 1.4H	77 K	1.028 V	-1.73 mV/K	5.8 mK	±76 mK	±98 mK	±11.6 mK
	calibration	300 K	0.5596 V	-2.3 mV/K	4.3 mK	±47 mK	±79 mK	±8.7 mK
		500 K	0.0907 V	-2.12 mV/K	4.7 mK	±40 mK	±90 mK	±9.4 mK
Silicon diode	DT-470-SD-13	1.4 K	1.6981 V	-13.1 mV/K	0.8 mK	±13 mK	±25 mK	±1.6 mK
	with 1.4H	77 K	1.0203 V	-1.92 mV/K	5.2 mK	±68 mK	±90 mK	±10.4 mK
	calibration	300 K	0.5189 V	-2.4 mV/K	4.2 mK	±44 mK	±76 mK	±8.4 mK
		475 K	0.0906 V	-2.22 mV/K	4.5 mK	±38 mK	±88 mK	±9 mK
GaAlAs diode	TG-120-SD	1.4 K	5.391 V	-97.5 mV/K	0.2 mK	±8.8 mK	±21 mK	±0.4 mK
	with 1.4H	77 K	1.422 V	-1.24 mV/K	16 mK	±373 mK	±395 mK	±32 mK
	calibration	300 K	0.8978 V	-2.85 mV/K	7 mK	±144 mK	±176 mK	±14 mK
		475 K	0.3778 V	-3.15 mV/K	6.4 mK	±114 mK	±164 mK	±12.6 mK
100 Ω platinum RTD	PT-103 with	30 K	3.660 Ω	0.191 Ω/Κ	1.1 mK	±13 mK	±23 mK	±2.2 mK
500 Ω full scale	14J calibration	77 K	20.38 Ω	0.423 Ω/K	0.5 mK	±10 mK	±22 mK	±1.0 mK
		300 K	110.35 Ω	0.387 Ω/K	5.2 mK	±39 mK	±62 mK	±10.4 mK
		500 K	185.668 Ω	0.378 Ω/Κ	5.3 mK	±60 mK	±106 mK	±10.6 mK
Cernox®	CX-1010-SD	0.3 K	2322.4 Ω	-10785 Ω/K	8.5 µK	±0.1 mK	±3.6 mK	±17 μK
	with 0.3L	0.5 K	1248.2 Ω	-2665.2 Ω/K	26 μK	±0.2 mK	±4.7 mK	±52 μK
	calibration	4.2 K	277.32 Ω	-32.209 Ω/K	140 µK	±3.8 mK	±8.8 mK	±280 μK
		300 K	30.392 Ω	-0.0654 Ω/K	23 mK	±339 mK	±414 mK	±46 mK
Cernox®	CX-1050-SD-HT <sup>6</sup>	1.4 K	26566 Ω	-48449 Ω/K	20 μK	±0.3 mK	±5.3 mK	±40 μK
COLLICK	with 1.4M	4.2 K	3507.2 Ω	-1120.8 Ω/K	196 µK	±2.1 mK	±7.1 mK	±392 μK
	calibration	77 K	205.67 Ω	-2.4116 Ω/K	1.9 mK	±38 mK	±54 mK	±3.8 mK
		420 K	45.03 Ω	-0.0829 Ω/K	1.9 mK	±338 mK	±403 mK	±36 mK
Germanium	GR-300-AA	0.35 K	45.03 Ω	-0.0829 Ω/K -193453 Ω/K	4 μK	±330 IIIK ±48 µK	±4.2 mK	±30 IIIK ±8 µK
ucilialiulii	with 0.3D	1.4 K	449 Ω	-581 Ω/K	41 μK	±481 μK	±4.7 mK	±82 μK
	calibration							
		4.2 K	94 Ω	-26.6 Ω/K	56 μK	±1.8 mK	±6.8 mK	±112 μK
Germanium	GR-1400-AA	100 K 1.8 K	2.7 Ω 15288 Ω	-0.024 Ω/K -26868 Ω/K	6.3 mK 28 μK	±152 mK	±175 mK ±4.5 mK	±12.6 mK
dermanium	with 1.4D					±302 μK		±56 μK
	calibration	4.2 K	1689 Ω	-862 Ω/K	91 μK	±900 μK	±5.1 mK	±182 μK
		10 K	253 Ω	-62.0 Ω/K	73 μK	±1.8 mK	±6.8 mK	±146 μK
0	000 4 500	100 K	2.8 Ω	-0.021 Ω/K	7.1 mK	±177 mK	±200 mK	±14.2 mK
Carbon-glass	CGR-1-500 with 1.4L	1.4 K	103900 Ω	-520000 Ω/K	13 µK	±0.1 mK	±4.1 mK	±26 μK
	calibration	4.2 K	584.6 Ω	-422.3 Ω/K	63 μK	±0.8 mK	±4.8 mK	±126 μK
	ouoruuru	77 K	14.33 Ω	-0.098 Ω/K	4.6 mK	±108 mK	±133 mK	±9.2 mK
- TM	B)/ / 200 1 1 1	300 K	8.55 Ω	-0.0094 Ω/K	16 mK	±760 mK	±865 mK	±32 mK
Rox™	RX-102A-AA	0.5 K	3701 Ω	-5478 Ω/K	41 μK	±0.5 mK	±5 mK	±82 μK
	with 0.3B calibration	1.4 K	2005 Ω	-667 Ω/K	128 µK	±1.4 mK	±6.4 mK	±256 μK
	Ganbration	4.2 K	1370 Ω	-80.3 Ω/K	902 μΚ	±8 mK	±24 mK	±1.8 mK
_		40 K	1049 Ω	-1.06 Ω/K	62 mK	±500 mK	±537 mK	±124 mK
Thermocouple	Type K	75 K	-5862.9 μV	15.6 μV/K	26 mK	±0.25 K <sup>7</sup>	Calibration not available	±52 mK
<b>50 mV</b> Option—3060		300 K	1075.3 μV	40.6 μV/K	10 mK	±0.038 K <sup>7</sup>	from Lake Shore	±19.6 mK
Option—3000		600 K	13325 μV	41.7 μV/K	10 mK	±0.184 K <sup>7</sup>		±20 mK
		1505 K	49998.3 μV	36.0 μV/K	11 mK	±0.73 K <sup>7</sup>		±22 mK
Capacitance	CS-501	4.2 K	6.0 nF	27 pF/K	1.9 mK	NA	Calibration not available	±3.8 mK
Option—3061		77 K	9.1 nF	52 pF/K	1.0 mK		from Lake Shore	±2.0 mK
		200 K	19.2 nF	174 pF/K	2.9 mK			±5.8 mK

Model 336 Temperature Controller

 $<sup>^{\</sup>rm 4}\,$  Typical sensor sensitivities were taken from representative calibrations for the sensor listed

 $<sup>^{\</sup>mbox{\tiny 5}}$  Control stability of the electronics only, in an ideal thermal system

<sup>&</sup>lt;sup>6</sup> Non-HT version maximum temperature: 325 K

 $<sup>^{\</sup>rm 7}\,$  Accuracy specification does not include errors from room temperature compensation



## Model 336 Specifications

#### **Input specifications**

Standard inputs and scanner option Model 3062	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy (at 25 °C)	Measurement temperature coefficient	Electronic control stability <sup>8</sup>
Diode	Negative	0 V to 2.5 V	10 μA ±0.05% <sup>9,10</sup>	10 μV	10 μV	±80 μV ±0.005% of rdg	(10 μV + 0.0005% of rdg)/°C	±20 μV
Diode	Negative	0 V to 10 V	10 μA ±0.05% <sup>9,10</sup>	100 μV	20 μV	±160 μV ±0.01% of rdg	(20 $\mu$ V + 0.0005% of rdg)/°C	±40 μV
		0 Ω to 10 Ω	1 mA <sup>11</sup>	0.1 mΩ	0.2 mΩ	$\pm 0.002~\Omega$ $\pm 0.01\%$ of rdg	(0.01 m $\Omega$ + 0.001% of rdg)/°C	±0.4 mΩ
		0 Ω to 30 Ω	1 mA <sup>11</sup>	0.1 mΩ	0.2 mΩ	$\pm 0.002~\Omega$ $\pm 0.01\%$ of rdg	(0.03 m $\Omega$ + 0.001% of rdg)/°C	±0.4 mΩ
		0 Ω to 100 Ω	1 mA <sup>11</sup>	1 mΩ	2 mΩ	$\pm 0.004~\Omega$ $\pm 0.01\%$ of rdg	(0.1 m $\Omega$ + 0.001% of rdg)/°C	±4 mΩ
PTC RTD	Positive	0 Ω to 300 Ω	1 mA <sup>11</sup>	1 mΩ	2 mΩ	$\pm 0.004~\Omega$ $\pm 0.01\%$ of rdg	(0.3 m $\Omega$ + 0.001% of rdg)/°C	±4 mΩ
		0 Ω to 1 kΩ	1 mA <sup>11</sup>	10 mΩ	20 mΩ	±0.04 Ω ±0.02% of rdg	(1 m $\Omega$ + 0.001% of rdg)/°C	±40 mΩ
		0 Ω to 3 kΩ	1 mA <sup>11</sup>	10 mΩ	20 mΩ	$\pm 0.04~\Omega~\pm 0.02\%$ of rdg	(3 m $\Omega$ + 0.001% of rdg)/°C	±40 mΩ
		0 Ω to 10 kΩ	1 mA <sup>11</sup>	100 mΩ	200 mΩ	$\pm 0.4~\Omega~\pm 0.02\%$ of rdg	(10 m $\Omega$ + 0.001% of rdg)/°C	±40 mΩ
		0 Ω to 10 Ω	1 mA <sup>11</sup>	0.1 mΩ	0.15 mΩ	$\pm 0.002~\Omega$ $\pm 0.06\%$ of rdg	(0.01 m $\Omega$ + 0.001% of rdg)/°C	±0.3 mΩ
		0 Ω to 30 Ω	300 μA <sup>11</sup>	0.1 mΩ	0.45 mΩ	$\pm 0.002~\Omega$ $\pm 0.06\%$ of rdg	(0.03 m $\Omega$ + 0.0015% of rdg)/°C	±0.9 mΩ
		0 Ω to 100 Ω	100 μA <sup>11</sup>	1 mΩ	1.5 mΩ	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	(0.1 m $\Omega$ + 0.001% of rdg)/°C	±3 mΩ
		0 Ω to 300 Ω	30 μA <sup>11</sup>	1 mΩ	4.5 mΩ	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	$(0.3 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$	±9 mΩ
NTC RTD 10 mV	Negative	0 Ω to 1 kΩ	10 μA <sup>11</sup>	10 mΩ	15 mΩ +0.002% of rdg	$\pm 0.1~\Omega~\pm 0.04\%$ of rdg	(1 m $\Omega$ + 0.001% of rdg)/°C	$\pm 30~\text{m}\Omega~\pm 0.004\%$ of rdg
		0 Ω to 3 kΩ	3 μA <sup>11</sup>	10 mΩ	45 mΩ +0.002% of rdg	$\pm 0.1~\Omega~\pm 0.04\%$ of rdg	$(3 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$	$\pm 90~\text{m}\Omega~\pm 0.004\%$ of rdg
		0 Ω to 10 kΩ	1 μA <sup>11</sup>	100 mΩ	$\begin{array}{c} 150 \text{ m}\Omega + 0.002\% \\ \text{of rdg} \end{array}$	$\pm 1.0~\Omega~\pm 0.04\%$ of rdg	(10 m $\Omega$ + 0.001% of rdg)/°C	$\pm 300~\text{m}\Omega~\pm 0.004\%$ of rdg
		0 Ω to 30 kΩ	300 nA <sup>11</sup>	100 mΩ	450 mΩ +0.002% of rdg	$\pm 2.0~\Omega~\pm 0.04\%$ of rdg	(30 m $\Omega$ + 0.001% of rdg)/°C	$\pm 900~\text{m}\Omega~\pm 0.004\%$ of rdg
Thermocouple option Model 3060	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy (at 25 °C)	Measurement temperature coefficient	Electronic control stability <sup>8</sup>
Thermocouple 3060	Positive	±50 mV	NA	0.1 μV	0.4 μV	±1 µV ±0.05% of rdg <sup>12</sup>	$(0.1 \ \mu V + 0.001\% \ of \ rdg)/^{\circ}C$	±0.8 μV

Capacitance option Model 3061	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy (at 25 °C)	Measurement temperature coefficient	Electronic control stability <sup>8</sup>
Capacitance Positive or 3061 negative	0.1 nF to 15 nF	3.496 kHz 1 mA square wave	0.1 pF	0.05 pF	±50 pF ±0.1% of rdg	2.5 pF/°C	0.1 pF	
	1 nF to 150 nF	3.496 kHz 10 mA square wave	1 pF	0.5 pF	±50 pF ±0.1% of rdg	5 pF/°C	1 pF	

 $<sup>^{\</sup>rm 8}$  Control stability of the electronics only, in ideal thermal system

<sup>&</sup>lt;sup>9</sup> Current source error has negligible effect on measurement accuracy

<sup>&</sup>lt;sup>10</sup> Diode input excitation can be set to 1 mA

<sup>&</sup>lt;sup>11</sup> Current source error is removed during calibration

 $<sup>^{\</sup>rm 12}$  Accuracy specification does not include errors from room temperature compensation



#### **Sensor input configuration**

	Diode/RTD	Thermocouple
Measurement type	4-lead differential	2-lead differential, room
		temperature compensated
Excitation	Constant current with current	NA
	reversal for RTDs	
Supported sensors	Diodes: Silicon, GaAlAs	Most thermocouple types
	<b>RTDs:</b> 100 Ω Platinum,	
	1000 Ω Platinum,	
	Germanium, Carbon-Glass,	
	Cernox®, and Rox™	
Standard curves	DT-470, DT-670, DT-500-D,	Type E, Type K, Type T,
	DT-500-E1, PT-100, PT-1000,	AuFe 0.07% vs. Cr,
	RX-102A, RX-202A	AuFe 0.03% vs. Cr
Input connector	6-pin DIN	Screw terminals in a ceramic
		isothermal block

#### **Thermometry**

Number of inputs 4 (8 with scanner option)

Input configuration Inputs can be configured from the front panel to accept any of the supported input types. Thermocouple and capacitance inputs require an optional input card that can be installed in the field.

**Supported option cards** Thermocouple (3060), capacitance (3061), or scanner (3062) Option slots 1

**Isolation** Sensor inputs optically isolated from other circuits but not each other A/D resolution 24-bit

**Input accuracy** Sensor dependent, refer to Input Specifications table

Measurement resolution Sensor dependent, refer to Input Specifications table

Maximum update rate 10 rdg/s on each input, 5 rdg/s when configured as 100 kΩ NTC RTD with reversal on, 2 rdg/s on each scanned input (scanner option only)

Autorange Automatically selects appropriate NTC RTD or PTC RTD range

**User curves** Room for 39 200-point CalCurves<sup>™</sup> or user curves

SoftCal<sup>™</sup> Improves accuracy of DT-470 diode to ±0.25 K from 30 K to 375 K; improves accuracy of platinum RTDs to  $\pm 0.25$  K from 70 K to 325 K; stored as user curves

Math Maximum and minimum

Filter Averages 2 to 64 input readings

#### Control

**Control outputs** 4

#### **Heater outputs (Outputs 1 & 2)**

Control type Closed loop digital PID with manual heater output or open loop Update rate 10/s

Tuning Autotune (one loop at a time), PID, PID zones

Control stability Sensor dependent, see Input Specifications table

PID control settings

Proportional (gain) 0 to 1000 with 0.1 setting resolution

Integral (reset) 1 to 1000 (1000/s) with 0.1 setting resolution

Derivative (rate) 1 to 200% with 1% resolution

Manual output 0 to 100% with 0.01% setting resolution

Zone control 10 temperature zones with P, I, D, manual heater out, heater range, control channel, ramp rate

Setpoint ramping 0.1 K/min to 100 K/min

#### **Output 1**

Model 336 Temperature Controller

	25 Ω setting	50 Ω setting					
Туре	Variable DC current source						
D/A resolution	16-	-bit					
Max power	100 W	50 W					
Max current	2 A	1 A					
Voltage compliance	50 V	50 V					
Heater load for max	25 Ω	50 Ω					
power							
Heater load range	10 Ω to	100 Ω					
Ranges	3 (decade ste	eps in power)					
Heater noise	0.12 µA RMS (dominated by line frequency and its harmonics)						
Grounding	Output referenced to chassis ground						
Heater connector	Dual banana						
Safety limits	Curve temperature, power up he	eater off, short circuit protection					

#### **Output 2**

	25 Ω setting	50 Ω setting			
Туре	Variable DC current source				
D/A resolution	16-bit				
Max power	50 W	50 W			
Max current	1.41 A	1 A			
Voltage compliance	35.4 V	50 V			
Heater load for max	25 Ω	50 Ω			
power					
Heater load range	10 Ω to 100 Ω				
Ranges	3 (decade steps in power)				
Heater noise	0.12 µA RMS (dominated by line frequency and its harmonics)				
Grounding	Output referenced to chassis ground				
Heater connector	Dual banana				
Safety limits	Curve temperature, power up heater off, short circuit protection				

#### Unpowered analog outputs (Outputs 3 & 4)

Control type Closed loop PID, PID zones, warm up heater mode, manual output, or monitor output

Tuning Autotune (one loop at a time), PID, PID zones

Control stability Sensor dependent, see Input Specifications table

PID control settings

Proportional (gain) 0 to 1000 with 0.1 setting resolution

Integral (reset) 1 to 1000 (1000/s) with 0.1 setting resolution

Derivative (rate) 1 to 200% with 1% resolution

0 to 100% with 0.01% setting resolution Manual output

Zone control 10 temperature zones with P, I, D, manual heater out, heater range, control channel, ramp rate

Setpoint ramping 0.1 K/min to 100 K/min

Warm up heater mode settings

Warm up percentage 0 to 100% with 1% resolution Warm up mode Continuous control or auto-off

**Monitor output settings** 

User selected Scale

Data source Temperature or sensor units

Settings Input, source, top of scale, bottom of scale, or manual

Type Variable DC voltage source

Update rate 10/s Range ±10 V

Resolution 16-bit, 0.3 mV Accuracy ±2.5 mV Noise 0.3 mV RMS

**Minimum load resistance** 1 k $\Omega$  (short-circuit protected)

**Connector** Detachable terminal block



#### Front panel

**Display** 8-line by 40-character (240 × 64 pixel) graphic LCD display module with LED backlight

Number of reading displays 1 to 8

Display units K,  $^{\circ}$ C, V, mV,  $\Omega$ 

Reading source Temperature, sensor units, max, and min

Display update rate 2 rdg/s

Temperature display resolution 0.0001° from 0° to 99.9999°, 0.001° from 100° to 999.999°, 0.01° above 1000°

Sensor units display resolution Sensor dependent, to 6 digits

Other displays Input name, setpoint, heater range, heater output, and PID Setpoint setting resolution Same as display resolution (actual resolution is sensor dependent)

Heater output display Numeric display in percent of full scale for power or current

Heater output resolution 0.01%

Display annunciators Control input, alarm, tuning

LED annunciators Remote, Ethernet status, alarm, control outputs

Keypad 27-key silicone elastomer keypad

Front panel features Front panel curve entry, display contrast control, and keypad lock-out

#### Interface

#### IEEE-488.2

**Capabilities** SH1, AH1, T5, L4, SR1, RL1, PP0, DC1, DT0, C0, E1

Reading rate To 10 rdg/s on each input

Software support LabVIEW™ driver (see www.lakeshore.com)

USB

Emulates a standard RS-232 serial port Function

**Baud rate** 57.600

Connector B-type USB connector Reading rate To 10 rdg/s on each input

Software support LabVIEW™ driver (see www.lakeshore.com)

**Ethernet** 

**Function** TCP/IP, web interface, curve handler, configuration backup,

chart recorder

Connector **RJ-45** 

Reading rate To 10 rdg/s on each input

Software support LabVIEW™ driver (see www.lakeshore.com)

**Alarms** 

Number 4, high and low for each input **Data source** Temperature or sensor units

Source, high setpoint, low setpoint, deadband, latching or non-Settings

> latching, audible on/off, and visible on/off Display annunciator, beeper, and relays

Relays

Number

**Actuators** 

Contacts Normally open (NO), normally closed (NC), and common (C)

**Contact rating** 30 VDC at 3 A

Operation Activate relays on high, low, or both alarms for any input, or

manual mode

Connector Detachable terminal block

#### General

Ambient temperature 15 °C to 35 °C at rated accuracy; 5 °C to 40 °C at reduced accuracy

Power requirement 100, 120, 220, 240 VAC, ±10%, 50 or 60 Hz, 250 VA Size 435 mm W  $\times$  89 mm H  $\times$  368 mm D (17 in  $\times$  3.5 in  $\times$  14.5 in), full rack

Weight 7.6 kg (16.8 lb) Approval CE mark, RoHS

#### Ordering information

#### Part number **Description**

4 diode/RTD inputs and 4 control outputs, including one dual 336 banana jack heater input connector (106-009), four 6-pin DIN

plug sensor input mating connectors (G-106-233), one 10-pin terminal block (G-106-750), a calibration certificate and a

user's manual

336-3060 Model 336 with a 3060 option card installed 336-3061 Model 336 with a 3061 option card installed 336-3062 Model 336 with a 3062 option card installed 3060 2-thermocouple input option, uninstalled 3061 Capacitance input option for 350/336, uninstalled 3062 4-channel scanner option for diodes and RTD sensors for

350/336, uninstalled

#### Please indicate your power/cord configuration:

1 100 V-U.S. cord (NEMA 5-15)

2 120 V-U.S. cord (NEMA 5-15)

3 220 V—Euro cord (CEE 7/7)

4 240 V-Euro cord (CEE 7/7)

5 240 V-U.K. cord (BS 1363) **6** 240 V—Swiss cord (SEV 1011)

**7** 220 V—China cord (GB 1002)

#### **Accessories**

Temperature controller cable, 3 m (10 ft)-IN STOCK 112-177

112-178 Temperature controller cable, 6 m (20 ft) 112-180 Temperature controller cable, 10 m (33 ft)

117-017 1 m (3.3 ft long) IEEE-488 (GPIB) computer interface cable

assembly

RM-1 Rack mount kit for mounting one full rack temperature

instrument

G-106-233 Sensor input mating connector (6-pin DIN plug)

G-106-755 Terminal block, 10-pin 106-009 Banana plug, dual

CAL-336-CERT Instrument recalibration with certificate

CAL-336-DATA Instrument recalibration with certificate and data CAL-3060-CERT Thermocouple board recalibration with certificate 119-048 Model 336 temperature controller manual

All specifications are subject to change without notice





Model 335 Temperature Controller

#### **Model 335 features**

- Operates down to 300 mK with appropriate NTC RTD sensors
- Two sensor inputs
- Two configurable PID control loops providing 50 W and 25 W or 75 W and 1 W
- Autotuning automatically calculates PID parameters
- Automatically switch sensor inputs using zones to allow continuous measurement and control from 300 mK to 1505 K
- Custom display set-up allows you to label each sensor input
- USB and IEEE-488 interfaces
- Supports diode, RTD, and thermocouple temperature sensors
- Sensor excitation current reversal eliminates thermal EMF errors for resistance sensors
- ±10 V analog voltage output, alarms, and relays
- CE certification
- Full 3 year standard warranty





#### Introduction

Designed with the user and ease of use in mind, the Model 335 temperature controller offers many user-configurable features and advanced functions that until now have been reserved for more expensive, high-end temperature controllers. The Model 335 is the first two-channel temperature controller available with user configurable heater outputs delivering a total of 75 W of low noise heater power—50 W and 25 W, or 75 W and 1 W. With that much heater power packed into an affordable half-rack sized instrument, the Model 335 gives you more power and control than ever.

The Model 335 supports the industry's most advanced line of cryogenic temperature sensors as manufactured by Lake Shore, including diodes, resistance temperature detectors (RTDs), and thermocouples. The controller's zone tuning feature allows you to measure and control temperatures seamlessly from 300 mK to over 1,500 K. This feature automatically switches temperature sensor inputs when your temperature range goes beyond the useable range of a given sensor. You'll never again have to be concerned with temperature sensor over or under errors and measurement continuity issues.

As a replacement to our popular Model 331 and 332 temperature controllers, the Model 335 offers software emulation modes for literal drop-in compatibility. The commands you are accustomed to sending to the Model 331 and 332 will either be interpreted directly or translated to the most appropriate Model 335 setting. The Model 335 comes standard-equipped with all of the functionality of the controllers it replaces, but offers additional features that save you time and money.

With the Model 335, you get a temperature controller you control from the world leader in cryogenic thermometry.



Control outputs are equipped with both hardware and software features allowing you, and not your temperature controller, to easily control your experiments. Output one functions as a current output while output two can be configured in either current or voltage mode. With output two in voltage mode, it functions as a ±10 V analog output while still providing 1 W of heater power and full closed loop proportional-integral-derivative (PID) control capability. Alarms and relays are included to help automate secondary control functions. The improved autotuning feature of the Model 335 can be used to automatically calculate PID control parameters, so you spend less time tuning your controller and more time conducting experiments.

The intuitive front panel layout and keypad logic, bright vacuum fluorescent display, and LED indicators enhance the user-friendly front panel interface of the Model 335. Four standard display modes are offered to accommodate different instrument configurations and user preferences. Say goodbye to sticky notes and hand written labels, as the ability to custom label sensor inputs eliminates the guesswork in remembering or determining the location to which a sensor input is associated. These features, combined with USB and IEEE-488 interfaces and intuitive menu structure and logic supports efficiency and ease of use.

#### **Sensor inputs**

The Model 335 offers two standard sensor inputs that are compatible with diode and RTD temperature sensors. The field-installable Model 3060 option adds thermocouple functionality to both inputs.

Sensor inputs feature a high-resolution 24-bit analog-to-digital converter and each of the two powered outputs function as separate current sources. Both sensor inputs are optically isolated from other circuits to reduce noise and to deliver repeatable sensor measurements. Current reversal eliminates thermal electromagnetic field (EMF) errors in resistance sensors. Ten excitation currents facilitate temperature measurement and control down to 300 mK using appropriate negative temperature coefficient (NTC) RTDs. Autorange mode automatically scales excitation current in NTC RTDs to reduce self heating at low temperatures as sensor resistance changes by many orders of magnitude. Temperatures down to 1.4 K can be measured and controlled using silicon or GaAlAs diodes. Software selects the appropriate excitation current and signal gain levels when the sensor type is entered via the instrument front panel. To increase your productivity, the unique zone setting feature automatically switches sensor inputs, enabling you to measure temperatures from 300 mK to over 1,500 K without interrupting your experiment.

The Model 335 includes standard temperature sensor response curves for silicon diodes, platinum RTDs, ruthenium oxide RTDs, and thermocouples. Non-volatile memory can also store up to 39 200-point CalCurves for Lake Shore calibrated temperature sensors or user curves. A built-in SoftCal algorithm can be used to generate curves for silicon diodes and platinum RTDs that can be stored as user curves. Temperature sensor calibration data can be easily loaded into the Model 335 temperature controller and manipulated using the Lake Shore curve handler software program.

#### **Temperature control**

Providing a total of 75 W of heater power, the Model 335 is the most powerful half rack temperature controller available. Designed to deliver very clean heater power, precise temperature control is ensured throughout your full scale temperature range for excellent measurement reliability, efficiency and throughput. Two independent PID control outputs can be configured to supply 50 W and 25 W or 75 W and 1 W of heater power. Precise control output is calculated based on your temperature setpoint and feedback from the control sensor. Wide tuning parameters accommodate most cryogenic cooling systems and many high-temperature ovens commonly used in laboratories. PID values can be manually set for fine control or the improved autotuning feature can automate the tuning process.

The Model 335 autotuning method calculates PID parameters and provides feedback to help build zone tables. The setpoint ramp feature provides smooth, continuous setpoint changes and predictable approaches to setpoint without the worry of overshoot or excessive settling times. The instrument's zone tuning feature automatically switches temperature sensor inputs when your temperature range goes beyond the useable range of a given sensor. This feature combined with the instrument's ability to scale the sensor excitation through ten pre-loaded current settings allows the Model 335 to provide continuous measurement and control from 300 mK to 1505 K.

Model 335 Temperature Controller

Both control outputs are variable DC current sources referenced to chassis ground. As a factory default, outputs 1 and 2 provide 50 W and 25 W of continuous power respectively, both to a 50  $\Omega$  or 25  $\Omega$  load. For increased functionality, output 2 can also be set to voltage mode. When set to voltage mode, it functions as a ±10 V analog output while still providing 1 W of heater power and full closed loop PID control capability. While in this mode, output 1 can provide up to 75 W of heater power to a 25  $\Omega$  load.

Temperature limit settings for inputs are provided as a safeguard against system damage. Each input is assigned a temperature limit, and if any input exceeds that limit, both control channels are automatically disabled.

#### Interface

The Model 335 is standard equipped with universal serial bus (USB) and parallel (IEEE-488) interfaces. In addition to gathering data, nearly every function of the instrument can be controlled via computer interface. You can download the Lake Shore curve handler software program to your computer to easily enter and manipulate sensor calibration curves for storage in the instrument's nonvolatile memory.

The USB interface emulates an RS-232C serial port at a fixed 57,600 baud rate, but with the physical plug-ins of a USB. It also allows you to download firmware upgrades, ensuring the most current firmware version is loaded into your instrument without having to physically change your instrument.

Both sensor inputs are equipped with a high and low alarm which offers latching and nonlatching operation. The two relays can be used in conjunction with the alarms to alert you of a fault condition and perform simple on-off control. Relays can be assigned to any alarm or operated manually.

The ±10 V analog voltage output can be configured to send a voltage proportional to temperature to a strip chart recorder or data acquisition system. You may select the scale and data sent to the output, including temperature or sensor units.

- 1 Sensor input connectors
- 2 Terminal block (analog outputs/ relays)
- 3 USB interface
- 4 IEEE-488 interface
- 5 Line input assembly
- 6 Output 2 heater
- 7 Output 1 heater
- 8 Thermocouple option inputs

#### Model 335 rear panel





#### **Configurable display**

The Model 335 offers a bright, vacuum fluorescent display that simultaneously displays up to four readings. You can display both control loops, or if you need to monitor just one input, you can display just that one in greater detail. Or you can custom configure each display location to suit your experiment. Data from any input can be assigned to any of the locations, and your choice of temperature sensor units can be displayed. For added convenience, you can also custom label each senor input, eliminating the guesswork in remembering or determining the location to which a sensor input is associated.



Two input/one loop display with labels
Standard display option featuring two inputs and associated outputs.

			Shield
H:	5.202K	800	34.635к

#### **Custom display with labels**

Reading locations can be user configured to accommodate application needs. Here, the input names are shown above the measurement readings along with the designated input letters.



#### Intuitive menu structure

Logical navigation allows you to spend more time on research and less time on setup.

#### Model 3060 thermocouple input option

The field installable Model 3060 thermocouple input option adds thermocouple functionality to both inputs. While the option can be easily removed, this is not necessary as the standard inputs remain fully functional when they are not being used to measure thermocouple temperature sensors. Calibration for the option is stored on the card so it can be installed in the field and used with multiple Model 335 temperature controllers without recalibration.

#### Sensor selection

#### Sensor temperature range (sensors sold separately)

		Model	Useful range	Magnetic field use
Diodes	Silicon diode	DT-670-SD	1.4 K to 500 K	T≥60 K & B≤3 T
	Silicon diode	DT-670E-BR	30 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-414	1.4 K to 375 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-421	1.4 K to 325 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-470-SD	1.4 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-471-SD	10 K to 500 K	T≥60 K & B≤3 T
	GaAlAs diode	TG-120-P	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-PL	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-SD	1.4 K to 500 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
Positive temperature	100 Ω platinum	PT-102/3	14 K to 873 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
coefficient RTDs	100 Ω platinum	PT-111	14 K to 673 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
	Rhodium-iron	RF-800-4	1.4 K to 500 K	$T > 77 \text{ K \& B} \le 8 \text{ T}$
	Rhodium-iron	RF-100T/U	1.4 K to 325 K	$T > 77 \text{ K \& B} \le 8 \text{ T}$
Negative	Cernox®	CX-1010	0.3 K to 325 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
temperature coefficient RTDs	Cernox®	CX-1030-HT	0.3 K to 420 K <sup>1, 3</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1050-HT	1.4 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1070-HT	4 K to 420 K <sup>1</sup>	T > 2 K & B ≤ 19 T
	Cernox®	CX-1080-HT	20 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
	Germanium	GR-300-AA	0.35 K to 100 K <sup>3</sup>	Not recommended
	Germanium	GR-1400-AA	1.8 K to 100 K <sup>3</sup>	Not recommended
	Carbon-glass	CGR-1-500	1.4 K to 325 K	$T > 2 K \& B \le 19 T$
	Carbon-glass	CGR-1-1000	1.7 K to 325 K <sup>2</sup>	$T > 2 K \& B \le 19 T$
	Carbon-glass	CGR-1-2000	2 K to 325 K <sup>2</sup>	$T > 2 K \& B \le 19 T$
	Rox™	RX-102	0.3 K to 40 K <sup>3</sup>	$T > 2 K \& B \le 10 T$
	Rox™	RX-103	1.4 K to 40 K	$T > 2 K \& B \le 10 T$
	Rox™	RX-202	0.3 K to 40 K <sup>3</sup>	$T > 2 K \& B \le 10 T$
Thermocouples	Type K	9006-006	3.2 K to 1505 K	Not recommended
Option—3060	Type E	9006-004	3.2 K to 934 K	Not recommended
	Chromel-	9006-002	1.2 K to 610 K	Not recommended
	AuFe 0.07%			

- <sup>1</sup> Non-HT version maximum temperature: 325 K
- <sup>2</sup> Low temperature limited by input resistance range
- <sup>3</sup> Low temperature specified with self-heating error: ≤ 5 mK

**Silicon diodes** are the best choice for general cryogenic use from 1.4 K to above room temperature. Silicon diodes are economical to use because they follow a standard curve and are interchangeable in many applications. They are not suitable for use in ionizing radiation or magnetic fields.

**Cernox**® thin-film RTDs offer high sensitivity and low magnetic field-induced errors over the 0.3 K to 420 K temperature range. Cernox sensors require calibration.

**Platinum RTDs** offer high uniform sensitivity from 30 K to over 800 K. With excellent reproducibility, they are useful as thermometry standards. They follow a standard curve above 70 K and are interchangeable in many applications.



#### **Typical sensor performance**

	Example Lake Shore sensor	Temperature	Nominal resistance/ voltage	Typical sensor sensitivity⁴	Measurement resolution: temperature equivalents	Electronic accuracy: temperature equivalents	Temperature accuracy including electronic accuracy, CalCurve™, and calibrated sensor	Electronic control stability <sup>5</sup> : temperature equivalents
Silicon diode	DT-670-CO-13	1.4 K	1.664 V	-12.49 mV/K	0.8 mK	±13 mK	±25 mK	±1.6 mK
	with 1.4H	77 K	1.028 V	-1.73 mV/K	5.8 mK	±76 mK	±98 mK	±11.6 mK
	calibration	300 K	0.5596 V	-2.3 mV/K	4.3 mK	±47 mK	±79 mK	±8.7 mK
		500 K	0.0907 V	-2.12 mV/K	4.7 mK	±40 mK	±90 mK	±9.4 mK
Silicon diode	DT-470-SD-13	1.4 K	1.6981 V	-13.1 mV/K	0.8 mK	±13 mK	±25 mK	±1.6 mK
	with 1.4H	77 K	1.0203 V	-1.92 mV/K	5.2 mK	±68 mK	±90 mK	±10.4 mK
	calibration	300 K	0.5189 V	-2.4 mV/K	4.2 mK	±44 mK	±76 mK	±8.4 mK
		475 K	0.0906 V	-2.22 mV/K	4.5 mK	±38 mK	±88 mK	±9 mK
GaAlAs diode	TG-120-SD	1.4 K	5.3909 V	-97.5 mV/K	0.2 mK	±8.8 mK	±21 mK	±0.4 mK
da iii lo diodo	with 1.4H	77 K	1.4222 V	-1.24 mV/K	16 mK	±373 mK	±395 mK	±32 mK
	calibration	300 K	0.8978 V	-2.85 mV/K	7 mK	±144 mK	±176 mK	±14 mK
100 Ω platinum RTD	PT-103 with	475 K 30 K	0.3778 V 3.660 Ω	-3.15 mV/K 0.191 Ω/K	6.4 mK 1.1 mK	±114 mK ±13 mK	±164 mK ±23 mK	±13 mK ±2.2 mK
$500 \Omega$ full scale	14J calibration	77 K	20.38 Ω		0.5 mK		±23 IIIK ±22 mK	±1.0 mK
ooo ir iun ooulo	1 10 odiibration			0.423 Ω/Κ		±10 mK		
		300 K	110.35 Ω	0.387 Ω/Κ	5.2 mK	±39 mK	±62 mK	±10.4 mK
Camana	CV 1010 CD	500 K	185.668 Ω	0.378 Ω/K	5.3 mK	±60 mK	±106 mK	±10.6 mK
Cernox®	CX-1010-SD with 0.3L	0.3 K	2322.4 Ω	-10785 Ω/K	8.5 µK	±0.1 mK	±3.6 mK	±17 μK
	calibration	0.5 K	1248.2 Ω	-2665.2 Ω/K	26 μΚ	±0.2 mK	±4.7 mK	±52 μK
	Gailbration	4.2 K	277.32 Ω	-32.209 Ω/K	140 µK	±3.8 mK	±8.8 mK	±280 μK
		300 K	30.392 Ω	-0.0654 Ω/K	23 mK	±339 mK	±414 mK	±46 mK
Cernox®	CX-1050-SD-HT <sup>6</sup>	1.4 K	26566 Ω	-48449 Ω/K	20 μΚ	±0.3 mK	±5.3 mK	±40 μK
	with 1.4M calibration	4.2 K	3507.2 Ω	-1120.8 Ω/K	196 µK	±2.1 mK	±7.1 mK	±392 μK
	Gailbration	77 K	205.67 Ω	-2.4116 Ω/K	1.9 mK	±38 mK	±54 mK	±3.8 mK
		420 K	45.03 Ω	-0.0829 Ω/K	18 mK	±338 mK	±403 mK	±36 mK
Germanium	GR-300-AA	0.35 K	18225 Ω	-193453 Ω/K	4 μΚ	±48 μK	±4.2 mK	±8 μK
	with 0.3D	1.4 K	449 Ω	-581 Ω/K	41 μK	±481 μK	±4.7 mK	±82 μK
	calibration	4.2 K	94 Ω	-26.6 Ω/K	56 μK	±1.8 mK	±6.8 mK	±112 μK
		100 K	2.7 Ω	-0.024 Ω/K	6.3 mK	±152 mK	±175 mK	±12.6 mK
Germanium	GR-1400-AA	1.8 K	15288 Ω	-26868 Ω/K	28 μΚ	±302 μK	±4.5 mK	±56 μK
	with 1.4D	4.2 K	1689 Ω	-862 Ω/K	91 μK	±900 μK	±5.1 mK	±182 μK
	calibration	10 K	253 Ω	-62.0 Ω/K	73 μK	±1.8 mK	±6.8 mK	±146 μK
		100 K	2.8 Ω	-0.021 Ω/K	7.1 mK	±177 mK	±200 mK	±14.2 mK
Carbon-glass	CGR-1-500	1.4 K	103900 Ω	-520000 Ω/K	13 μK	±0.1 mK	±4.1 mK	±26 μK
	with 1.4L	4.2 K	584.6 Ω	-422.3 Ω/K	63 μK	±0.8 mK	±4.8 mK	±126 μK
	calibration	77 K	14.33 Ω	-0.098 Ω/K	4.6 mK	±108 mK	±133 mK	±9.2 mK
		300 K	8.55 Ω	-0.0094 Ω/K	16 mK	±760 mK	±865 mK	±32 mK
Rox™	RX-102A-AA	0.5 K	3701 Ω	-5478 Ω/K	41 µK	±0.5 mK	±5 mK	±82 μK
	with 0.3B	1.4 K	2005 Ω	-667 Ω/K	128 µK	±1.4 mK	±6.4 mK	±256 μK
	calibration	4.2 K	1370 Ω	-80.3 Ω/K	902 μK	±8 mK	±24 mK	±1.8 mK
		40 K	1049 Ω	-1.06 Ω/K	62 mK	±500 mK	±537 mK	±124 mK
Thermocouple	Type K	75 K	-5862.9 μV	15.6 μV/K	26 mK	±0.25 K <sup>7</sup>	Calibration not available	±52 mK
50 mV	130011	300 K	1075.3 μV	40.6 μV/K	10 mK	±0.038 K <sup>7</sup>	from Lake Shore	±20 mK
Option—3060		600 K	1073.3 μV 13325 μV	40.0 μV/K 41.7 μV/K	10 mK	±0.036 K ±0.184 K <sup>7</sup>		±20 mK
		1505 K	49998.3 μV	36.006 μV/K	11 mK	±0.73 K <sup>7</sup>		±22 mK

Model 335 Temperature Controller

 $<sup>^{\</sup>rm 4}\,$  Typical sensor sensitivities were taken from representative calibrations for the sensor listed

<sup>&</sup>lt;sup>5</sup> Control stability of the electronics only, in an ideal thermal system

<sup>&</sup>lt;sup>6</sup> Non-HT version maximum temperature: 325 K

<sup>&</sup>lt;sup>7</sup> Accuracy specification does not include errors from room temperature compensation



#### Model 335 Specifications

#### Input specifications

	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy <sup>6</sup>	Measurement temperature coefficient	Electronic control stability <sup>1</sup>
Diode	Negative	0 V to 2.5 V	10 μA ±0.05% <sup>2,3</sup>	100 μV	10 μV	±80 µV ±0.005% of rda	(10 µV + 0.0005% of rdg)/°C	±20 μV
		0 V to 10 V	10 μA ±0.05% <sup>2,3</sup>	1 mV	20 μV	±320 µV ±0.01% of rdg	(20 µV + 0.0005% of rdg)/°C	±40 μV
PTC RTD	Positive	0 $\Omega$ to 10 $\Omega$	1 mA <sup>4</sup>	1 mΩ	0.2 mΩ	±0.002 Ω ±0.01% of rda	(0.01 m $\Omega$ + 0.001% of rdg)/°C	±0.4 mΩ
		0 Ω to 30 Ω	1 mA <sup>4</sup>	1 mΩ	0.2 mΩ	$\pm 0.002 \Omega$ $\pm 0.01\%$ of rdg	$(0.03 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±0.4 mΩ
		0 Ω to 100 Ω	1 mA <sup>4</sup>	10 mΩ	2 mΩ	$\pm 0.004 \Omega$ $\pm 0.01\%$ of rdg	(0.1 m $\Omega$ + 0.001% of rdg)/°C	±4 mΩ
		0 Ω to 300 Ω	1 mA <sup>4</sup>	10 mΩ	2 mΩ	±0.004 Ω ±0.01% of rda	(0.3 m $\Omega$ + 0.001% of rdg)/°C	±4 mΩ
		0 Ω to 1 kΩ	1 mA <sup>4</sup>	100 mΩ	20 mΩ	±0.04 Ω ±0.02% of rda	$(1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±40 mΩ
		0 Ω to 3 kΩ	1 mA <sup>4</sup>	100 mΩ	20 mΩ	±0.04 Ω ±0.02% of rda	$(3 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±40 mΩ
		0 Ω to 10 kΩ	1 mA <sup>4</sup>	1 Ω	200 mΩ	±0.4 Ω ±0.02% of rda	(10 m $\Omega$ + 0.001% of rdg)/°C	±400 mΩ
NTC RTD 10 mV	Negative	0 $\Omega$ to 10 $\Omega$	1 mA <sup>4</sup>	1 mΩ	0.15 mΩ	±0.002 Ω ±0.06% of rda	(0.01 m $\Omega$ + 0.001% of rdg)/°C	±0.3 mΩ
101114		0 Ω to 30 Ω	300 μA⁴	1 mΩ	0.45 mΩ	±0.002 Ω ±0.06% of rda	$(0.03 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$	±0.9 mΩ
		0 Ω to 100 Ω	100 μA <sup>4</sup>	10 mΩ	1.5 mΩ	$\pm 0.01 \Omega$ $\pm 0.04\%$ of rdg	$(0.1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±3 mΩ
		0 Ω to 300 Ω	30 μA <sup>4</sup>	10 mΩ	4.5 mΩ	±0.01 Ω ±0.04% of rda	$(0.3 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$	±9 mΩ
		0 Ω to 1 kΩ	10 μA <sup>4</sup>	100 mΩ	15 mΩ +0.002% of rda	±0.1 Ω ±0.04% of rda	$(1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$	±30 mΩ ±0.004% of rdg
		0 Ω to 3 kΩ	3 μA <sup>4</sup>	100 mΩ	of rdg 45 mΩ +0.002% of rda	±0.1 Ω ±0.04% of rda	$(3 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$	±90 mΩ ±0.004% of rda
		0 Ω to 10 kΩ	1 μA <sup>4</sup>	1 Ω	of rdg 150 mΩ +0.002% of rda	of rdg $\pm 1.0 \Omega \pm 0.04\%$ of rdg	(10 mΩ + 0.001% of rdg)/°C	±300 mΩ ±0.004% of rdg
		0 Ω to 30 kΩ	300 nA <sup>4</sup>	1 Ω	of rdg 450 mΩ +0.002% of rda	$\begin{array}{c} \text{ of rdg} \\ \pm 2.0 \ \Omega \ \pm 0.04\% \\ \text{ of rdg} \end{array}$	$(30 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$	±900 mΩ ±0.004% of rda
		0 Ω to 100 kΩ	100 nA⁴	10 Ω	of rdg 1.5 Ω +0.005% of rdg	±10.0 Ω ±0.04% of rdg	$(100 \text{ m}\Omega + 0.002\% \text{ of rdg})/^{\circ}\text{C}$	$\pm 3 \Omega \pm 0.01\%$
Thermocouple Option—3060	Positive	±50 mV	NA	1 μV	0.4 μV	±1 µV ±0.05% of rdg <sup>5</sup>	(0.1 $\mu$ V + 0.001% of rdg)/°C	±0.8 μV

- Control stability of the electronics only, in ideal thermal system
- <sup>2</sup> Current source error has negligible effect on measurement accuracy
- <sup>3</sup> Diode input excitation can be set to 1 mA

- Current source error is removed during calibration
- <sup>5</sup> Accuracy specification does not include errors from room temperature compensation
- <sup>6</sup> Accuracy at T<sub>cal</sub>, typically 23.5 °C ±1.5 °C

#### **Sensor input configuration**

	Diode/RTD	Thermocouple
Measurement type	4-lead differential	2-lead differential, room
		temperature compensated
Excitation	Constant current with current	NA
	reversal for RTDs	
Supported sensors	Diodes: Silicon, GaAlAs	Most thermocouple types
	<b>RTDs:</b> 100 Ω Platinum,	
	1000 Ω Platinum,	
	Germanium, Carbon-Glass,	
	Cernox®, and Rox™	
Standard curves	DT-470, DT-670, DT-500-D,	Type E, Type K, Type T,
	DT-500-E1, PT-100, PT-1000,	AuFe 0.07% vs. Cr,
	RX-102A, RX-202A	AuFe 0.03% vs. Cr
Input connector	6-pin DIN	Screw terminals in a ceramic
		isothermal block

#### **Thermometry**

Number of inputs 2

Input configuration Inputs can be configured from the front panel to accept any of the supported input types. Thermocouple inputs require an optional input card that can be installed in the field. Once installed the thermocouple input can be selected from the front panel like any other input type.

Isolation Sensor inputs optically isolated from other circuits but not each other A/D resolution 24-bit

Input accuracy Sensor dependent, refer to Input Specifications table

Measurement resolution Sensor dependent, refer to Input Specifications table

Maximum update rate 10 rdg/s on each input, 5 rdg/s when configured as 100 kΩ NTC RTD with reversal on

Autorange Automatically selects appropriate NTC RTD or PTC RTD range

User curves Room for 39 200-point CalCurves™ or user curves

SoftCal™ Improves accuracy of DT-470 diode to ±0.25 K from 30 K to 375 K; improves accuracy of platinum RTDs to  $\pm 0.25$  K from 70 K to 325 K; stored as user curves Math Maximum and minimum

Filter Averages 2 to 64 input readings

#### **Control**

**Control outputs** 2

#### **Heater outputs**

Control type Closed loop digital PID with manual heater output or open loop; warm up mode (output 2 only)

Update rate 10/s

Tuning Autotune (one loop at a time), PID, PID zones

Control stability Sensor dependent, see Input Specifications table

PID control settings

0 to 1000 with 0.1 setting resolution Proportional (gain) Integral (reset) 1 to 1000 (1000/s) with 0.1 setting resolution 1 to 200% with 1% resolution Derivative (rate) Manual output 0 to 100% with 0.01% setting resolution

Zone control 10 temperature zones with P, I, D, manual heater out, heater range, control channel, ramp rate

Setpoint ramping 0.1 K/min to 100 K/min
Warm up heater mode settings (output 2 only)

Warm up percentage 0 to 100% with 1% resolution Warm up mode Continuous control or auto-off

Monitor output settings (output 2 voltage only)

Scale User selected

Data source Temperature or sensor units

Settings Input, source, top of scale, bottom of scale, or manual

#### **Output 1**

Type	Variable DC current source					
Control modes	Closed loop digital PID with manual output or open loop					
D/A resolution	16-bit					
	25 Ω s	etting	50 Ω setting			
Max power	75 W*	50 W	50 W			
Max current	1.73 A	1.41 A	1 A			
Voltage compliance (min)	43.3 V	35.4 V	50 V			
Heater load for max	25 Ω 25 Ω		50 Ω			
power						
Heater load range	10 Ω to 100 Ω					
Ranges	3 (decade steps in power)					
Heater noise	0.12 µA RMS (dominated by line frequency and its harmonics)					
Heater connector	Dual banana					
Grounding	Output referenced to chassis ground					
Safety limits	Curve temperature, power up heater off, short circuit protection					

<sup>\*75</sup> W only available when output 2 is in voltage mode

#### **Output 2**

Type	Variable DC current source or voltage source				
	Curren	t mode	Voltage mode		
Control modes	Closed loop d	igital PID with	Closed loop digital PID with		
	manual output,	zone, open loop	manual output, zone, open		
			loop, warm up, monitor out		
D/A resolution	15-	-bit	16-bit (bipolar)/15-bit		
			(unipolar)		
	25 Ω setting	50 Ω setting	N/A		
Max power	25 W 25 W		1 W		
Max current	1 A 0.71 A		100 mA		
Voltage compliance (min)	25 V	35.4 V	±10 V		
Heater load for max	25 Ω 50 Ω		100 Ω		
power					
Heater load range	10 Ω to	100 Ω	100 Ω min (short circuit		
			protected)		
Ranges	3 (decade steps in power)		N/A		
Heater noise	0.12 μA RMS		0.3 mV RMS		
Heater connector	Dual banana		Detachable terminal block		
Grounding		put referenced to			
Safety limits	Curve temperatu	re, power up heat	er off, short circuit protection		

Update rate 10/s Range ±10 V Resolution 16-bit, 0.3 mV Accuracy ±2.5 mV Noise 0.3 mV RMS

Minimum load resistance 100 Ω (short-circuit protected)

Connector Detachable terminal block

#### Front panel

Display 2-line by 20-character, 9 mm character height, vacuum fluorescent display Number of reading displays 1 to 4

Display units K, °C, V, mV, Ω

Reading source Temperature, sensor units, max, and min

Display update rate 2 rdg/s

Temperature display resolution 0.001° from 0° to 99.999°, 0.01° from 100° to 999.99°,

0.1° above 1000°

Sensor units display resolution Sensor dependent, to 5 digits

Other displays Sensor name, setpoint, heater range, heater output, and PID Setpoint setting resolution Same as display resolution (actual resolution is

sensor dependent)

Heater output display Numeric display in percent of full scale for power or current

Heater output resolution 1%

Display annunciators Control input, alarm, tuning LED annunciators Remote, alarm, control outputs

Keypad 25-key silicone elastomer keypad

Front panel features Front panel curve entry, display brightness control, and keypad lock-out

#### Interface

IEEE-488.2

Capabilities SH1, AH1, T5, L4, SR1, RL1, PP0, DC1, DT0, C0, E1

Reading rate To 10 rdg/s on each input

Software support LabVIEW<sup>™</sup> driver (see www.lakeshore.com)

IISR

Model 335 Temperature Controller

**Function** Emulates a standard RS-232 serial port

**Baud rate** 57,600

Connector B-type USB connector

Reading rate
Software support
To 10 rdg/s on each input
LabVIEW™ driver (see www.lakeshore.com) Special interface features Model 331/332 command emulation mode

Number 2, high and low for each input Data source Temperature or sensor units

Settings Source, high setpoint, low setpoint, deadband, latching or non-

latching, audible on/off, and visible on/off **Actuators** Display annunciator, beeper, and relays

Relays

Number

Normally open (NO), normally closed (NC), and common (C) Contacts

**Contact rating** 

Operation Activate relays on high, low, or both alarms for any input, or manual

Connector Detachable terminal block

#### General

Ambient temperature 15 °C to 35 °C at rated specifications; 5 °C to 40 °C at reduced specifications

Power requirement 100, 120, 220, 240 VAC, ±10%, 50 or 60 Hz, 210 VA Size 217 mm W  $\times$  90 mm H  $\times$  317 mm D (8.5 in  $\times$  3.5 in  $\times$  14.5 in), half rack

Weight 5.1 kg (11.3 lb) Approval CE mark, RoHS

#### Ordering information

#### Part number **Description**

2 diode/RTD inputs and 2 control outputs temperature 335

controller-includes one dual banana jack heater output connector (106-009), two 6-pin DIN plug sensor input mating connectors (G-106-233), one 8-pin terminal block (G-107-773),

a calibration certificate and user manual Model 335 with 3060 option card installed

335-3060 2-thermocouple input option for Model 335, uninstalled

#### Please indicate your power/cord configuration:

100 V-U.S. cord (NEMA 5-15)

120 V—U.S. cord (NEMA 5-15)

220 V—Euro cord (CEE 7/7)

240 V—Euro cord (CEE 7/7)

240 V-U.K. cord (BS 1363)

240 V-Swiss cord (SEV 1011)

220 V—China cord (GB 1002)

#### **Accessories**

112-177 Temperature controller cable, 3 m (10 ft)—IN STOCK

112-178 Temperature controller cable, 6 m (20 ft) 112-180 117-017 Temperature controller cable, 10 m (33 ft)

1 m (3.3 ft long) IEEE-488 (GPIB) computer interface

cable assembly

RM-2 Kit for mounting two 1/2-rack temperature instruments in a 483

mm (19 in) rack

RM-1/2 Kit for mounting one 1/2-rack temperature instrument in a 483

mm (19 in) rack

G-106-233 Sensor input mating connector (6-pin DIN plug)

G-106-773 Terminal block, 8-pin

106-009 Banana plug, dual CAL-335-CERT Instrument recalibration with certificate

Instrument recalibration with certificate and data CAL-335-DATA Model 335 temperature controller manual

All specifications are subject to change without notice





# Model 325 Temperature Controller



- Operates down to 1.2 K with appropriate sensor
- Two sensor inputs
- Supports diode, RTD, and thermocouple sensors
- Sensor excitation current reversal eliminates thermal EMF errors in resistance sensors
- Two autotuning control loops: 25 W and 2 W maximum
- Control loop 2: variable DC voltage source from 0 to 10 V maximum
- IEEE-488 and RS-232C interfaces
- **CE** certification
- Full 3 year standard warranty





#### Introduction

The Model 325 dual-channel temperature controller is capable of supporting nearly any diode, RTD, or thermocouple temperature sensor. Two independent PID control loops with heater outputs of 25 W and 2 W are configured to drive either a 50  $\Omega$  or 25  $\Omega$ load for optimal cryocooler control flexibility. Designed with ease of use, functionality, and value in mind, the Model 325 is ideal for general-purpose laboratory and industrial temperature measurement and control applications.

#### **Sensor inputs**

The Model 325 temperature controller features two inputs with a high-resolution 24-bit analog-to-digital converter and separate current sources for each input. Constant current excitation allows temperature to be measured and controlled down to 2.0 K using appropriate Cernox® RTDs or down to 1.4 K using silicon diodes. Thermocouples allow for temperature measurement and control above 1.500 K. Sensors are optically isolated from other instrument functions for quiet and repeatable sensor measurements. The Model 325 also uses current reversal to eliminate thermal EMF errors in resistance sensors. Sensor data from each input is updated up to ten times per second, with display outputs twice each second.

Standard temperature response curves for silicon diodes, platinum RTDs, ruthenium oxide RTDs, and many thermocouples are included. Up to fifteen 200-point CalCurves™ (for Lake Shore calibrated temperature sensors) or user curves can be stored into non-volatile memory. A built-in SoftCal™1 algorithm can be used to generate curves for silicon diodes and platinum RTDs for storage as user curves. The Lake Shore curve handler software program allows sensor curves to be easily loaded and manipulated.

Model 325 Temperature Controller

Sensor inputs for the Model 325 are factory configured and compatible with either diodes/ RTDs or thermocouple sensors. Your choice of two diode/RTD inputs, one diode/RTD input and one thermocouple input, or two thermocouple inputs must be specified at time of order and cannot be reconfigured in the field. Software selects appropriate excitation current and signal gain levels when the sensor type is entered via the instrument front panel.

#### Temperature control

The Model 325 temperature controller offers two independent proportionalintegral-derivative (PID) control loops. A PID algorithm calculates control output based on temperature setpoint and feedback from the control sensor. Wide tuning parameters accommodate most cryogenic cooling systems and many small high-temperature ovens. A high-resolution digital-to-analog converter generates a smooth control

output. The user can set the PID values or the Autotuning feature of the Model 325 can automate the tuning process.

Control loop 1 heater output for the Model 325 is a well-regulated variable DC current source. The output can provide up to 25 W of continuous power to a 50  $\Omega$  or 25  $\Omega$  heater load, and includes a lower range for systems with less cooling power. Control loop 2 heater output is a single-range, variable DC voltage source. The output can source up to 0.2 A, providing 2 W of heater power at the 50  $\Omega$ setting or 1 W at the 25  $\Omega$  setting. When not being used for temperature control, the loop 2 heater output can be used as a manually controlled voltage source. The output voltage can vary from 0 to 10 V on the 50  $\Omega$  setting, or 0 to 5 V on the 25  $\Omega$  setting. Both heater outputs are referenced to chassis ground.

The setpoint ramp feature allows smooth continuous setpoint changes and can also make the approach to setpoint more predictable. The zone feature can automatically change control parameter values for operation over a large temperature range. Ten different temperature zones can be loaded into the instrument, which will select the next appropriate value on setpoint change.

Temperature limit settings for inputs are provided as a safeguard against system damage1. Each input is assigned a temperature limit, and if any input exceeds that limit, all control channels are automatically disabled.

1 Firmware version 1.5 and later



#### Model 325 rear panel

- Loop 1 heater output
- 2 Serial (RS-232C) I/O (DTE)
- Line input assembly

- Loop 2 heater output
- 5 Sensor input connectors
- IEEE-488 interface



#### Interface

The Model 325 includes both parallel (IEEE-488) and serial (RS-232C) computer interfaces. In addition to data gathering, nearly every function of the instrument can be controlled via computer interface. Sensor curves can also be entered and manipulated through either interface using the Lake Shore curve handler software program.



#### Normal (default) display configuration

The display provides four reading locations. Readings from each input and the control setpoint can be expressed in any combination of temperature or sensor units, with heater output expressed as a percent of full scale current or power.

#### Flexible configuration

Reading locations can be configured by the user to meet application needs. The character preceding the reading indicates input A or B or setpoint S. The character following the reading indicates measurement units.

	mer ca	and the same of th	Section 2	and an analysis of
1:	54.	NORU	28.	0000k
2:	38.	2000	30.	0000k

#### **Curve entry**

The Model 325 display offers the flexibility to support curve, SoftCal™, and zone entry. Curve entry may be performed accurately and to full resolution via the display and keypad as well as computer interface.

#### **Configurable display**

The Model 325 offers a bright, easy to read LCD display that simultaneously displays up to four readings. Display data includes input and source annunciators for each reading. All four display locations can be configured by the user. Data from either input can be assigned to any of the four locations, and the user's choice of temperature or sensor units can be displayed. Heater range and control output as current or power can be continuously displayed for immediate feedback on control operation. The channel A or B indicator is underlined to indicate which channel is being controlled by the displayed control loop.

#### Sensor Selection

#### Sensor temperature range (sensors sold separately)

		Model	Useful range	Magnetic field use
Diodes	Silicon diode	DT-670-SD	1.4 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-670E-BR	30 K to 500 K	$T \geq 60~K~\&~B \leq 3~T$
	Silicon diode	DT-414	1.4 K to 375 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-421	1.4 K to 325 K	$T \geq 60~K~\&~B \leq 3~T$
	Silicon diode	DT-470-SD	1.4 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-471-SD	10 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	GaAlAs diode	TG-120-P	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-PL	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-SD	1.4 K to 500 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
Positive temperature	100 Ω platinum	PT-102/3	14 K to 873 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
coefficient RTDs	100 Ω platinum	PT-111	14 K to 673 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
	Rhodium-iron	RF-800-4	1.4 K to 500 K	$T > 77 \text{ K \& B} \le 8 \text{ T}$
	Rhodium-iron	RF-100T/U	1.4 K to 325 K	$T > 77~K~\&~B \leq 8~T$
Negative	Cernox®	CX-1010	2 K to 325 K <sup>5</sup>	$T>2~K~\&~B\leq 19~T$
temperature	Cernox®	CX-1030-HT	3.5 K to 420 K <sup>3,6</sup>	$T>2~K~\&~B\leq19~T$
coefficient RTDs <sup>2</sup>	Cernox®	CX-1050-HT	4 K to 420 K <sup>3, 6</sup>	$T>2~K~\&~B\leq19~T$
	Cernox®	CX-1070-HT	15 K to 420 K <sup>3</sup>	$T>2~K~\&~B\leq19~T$
	Cernox®	CX-1080-HT	50 K to 420 K <sup>3</sup>	$T>2~K~\&~B\leq19~T$
	Germanium	GR-300-AA	1.2 K to 100 K <sup>4</sup>	Not recommended
	Germanium	GR-1400-AA	4 K to 100 K <sup>4</sup>	Not recommended
	Carbon-glass	CGR-1-500	4 K to 325 K <sup>5</sup>	$T > 2 K \& B \le 19 T$
	Carbon-glass	CGR-1-1000	5 K to 325 K <sup>5</sup>	$T>2~K~\&~B\leq 19~T$
	Carbon-glass	CGR-1-2000	6 K to 325 K <sup>5</sup>	$T>2~K~\&~B\leq19~T$
	Rox™	RX-102A	1.4 K to 40 K <sup>5</sup>	$T > 2 K \& B \le 10 T$
Thermocouples	Type K	9006-006	3.2 K to 1505 K	Not recommended
	Type E	9006-004	3.2 K to 934 K	Not recommended
	Chromel-AuFe	9006-002	1.2 K to 610 K	Not recommended
	0.07%			

<sup>&</sup>lt;sup>2</sup> Single excitation current may limit the low temperature range of NTC resistors

**Silicon diodes** are the best choice for general cryogenic use from 1.4 K to above room temperature. Diodes are economical to use because they follow a standard curve and are interchangeable in many applications. They are not suitable for use in ionizing radiation or magnetic fields.

**Cernox®** thin-film RTDs offer high sensitivity and low magnetic field-induced errors over the 2 K to 420 K temperature range. Cernox sensors require calibration.

**Platinum RTDs** offer high uniform sensitivity from 30 K to over 800 K. With excellent reproducibility, they are useful as thermometry standards. They follow a standard curve above 70 K and are interchangeable in many applications.

<sup>&</sup>lt;sup>3</sup> Non-HT version maximum temperature: 325 K

<sup>&</sup>lt;sup>4</sup> Low temperature limited by input resistance range

<sup>&</sup>lt;sup>5</sup> Low temperature specified with self-heating error: ≤5 mK

<sup>&</sup>lt;sup>6</sup> Low temperature specified with self-heating error: ≤12 mK



	Example Lake Shore sensor	Temperature	Nominal resistance/ voltage	Typical sensor sensitivity <sup>7</sup>	Measurement resolution: temperature equivalents	Electronic accuracy: temperature equivalents	Temperature accuracy including electronic accuracy,	Electronic control stability8: temperature
					<b>0 1</b>		CalCurve™, and calibrated sensor	equivalents
Silicon diode	DT-670-SD-13	1.4 K	1.664 V	-12.49 mV/K	0.8 mK	±13 mK	±25 mK	±1.6 mK
	with 1.4H	77 K	1.028 V	-1.73 mV/K	5.8 mK	±76 mK	±98 mK	±11.6 mK
	calibration	300 K	0.5596 V	-2.3 mV/K	4.3 mK	±47 mK	±79 mK	±8.7 mK
		500 K	0.0907 V	-2.12 mV/K	4.8 mK	±40 mK	±90 mK	±9.6 mK
Silicon diode	DT-470-SD-13	1.4 K	1.6981 V	-13.1 mV/K	0.8 mK	±13 mK	±25 mK	±1.6 mK
	with 1.4H	77 K	1.0203 V	-1.92 mV/K	5.2 mK	±68 mK	±90 mK	±10.4 mK
	calibration	300 K	0.5189 V	-2.4 mV/K	4.2 mK	±44 mK	±76 mK	±8.4 mK
		475 K	0.0906 V	-2.22 mV/K	4.6 mK	±39 mK	±89 mK	±9.2 mK
GaAlAs diode	TG-120-SD	1.4 K	5.391 V	-97.5 mV/K	0.2 mK	±8.8 mK	±21 mK	±0.4 mK
	with 1.4H	77 K	1.422 V	-1.24 mV/K	16.2 mK	±373 mK	±395 mK	±32.4 mK
	calibration	300 K	0.8978 V	-2.85 mV/K	7 mK	±144 mK	±176 mK	±14 mK
		475 K	0.3778 V	-3.15 mV/K	6.4 mK	±114 mK	±164 mK	±12.8 mK
00 Ω platinum RTD	PT-103 with	30 K	3.660 Ω	0.191 Ω/Κ	10.5 mK	±23 mK	±33 mK	±21 mK
500 Ω full scale	1.4J calibration	77 K	20.38 Ω	0.423 Ω/K	4.8 mK	±15 mK	±27 mK	±9.6 mK
		300 K	110.35 Ω	0.387 Ω/Κ	5.2 mK	±39 mK	±62 mK	±10.4 mK
		500 K	185.668 Ω	0.378 Ω/Κ	5.3 mK	±60 mK	±106 mK	±10.6 mK
Cernox®	CX-1050-SD-	4.2 K	3507.2 Ω	-1120.8 Ω/K	36 μK	±1.4 mK	±6.4 mK	±72 μK
	HT <sup>9</sup> with 4M	77 K	205.67 Ω	-2.4116 Ω/K	16.6 mK	±76 mK	±92 mK	±33.2 mK
	calibration	300 K	59.467 Ω	-0.1727 Ω/K	232 mK	±717 mK	±757 mK	±464 mK
		420 K	45.030 Ω	-0.0829 Ω/K	483 mK	±1.42 K	±1.49 K	±966 mK
Germanium	GR-300-AA	1.2 K	600 Ω	-987 Ω/K	51 μK	±345 μK	±4.5 mK	±101 μK
	with 0.3D	1.4 K	449 Ω	-581 Ω/K	86 μK	±481 μK	±4.7 mK	±172 μK
	calibration	4.2 K	94 Ω	-27 Ω/K	1.9 mK	±5.19 mK	±10.2 mK	±3.8 mK
		100 K	2.72 Ω	-0.024 Ω/K	2.1 K	±4.25 K	±4.27 K	±4.20 K
Germanium	GR-1400-AA	4 K	1873 Ω	-1008 Ω/K	50 μK	±842 μK	±5.0 mK	±99 μK
	with 1.4D	4.2 K	1689 Ω	-862 Ω/K	58 μK	±900 μK	±5.1 mK	±116 μK
	calibration	10 K	253 Ω	-62 Ω/K	807 μK	±3.2 mK	±8.2 mK	±1.6 mK
		100 K	2.80 Ω	-0.021 Ω/K	2.4 K	±4.86 K	±4.884 K	±4.81 K
Carbon-glass	CGR-1-2000	4.2 K	2260 Ω	-2060 Ω/K	20 μK	±0.5 mK	±4.5 mK	±40 μK
- <b>3</b>	with 4L	77 K	21.65 Ω	-0.157 Ω/K	255 mK	±692 mK	±717 mK	±510 mK
	calibration	300 K	11.99 Ω	-0.015 Ω/K	2.667 K	±7 K	±7.1 K	±5.334 K
Thermocouple	Type K	75 K	-5862.9 μV	15.6 µV/K	26 mK	±0.25 K <sup>10</sup>	Calibration not	±52 mK
50 mV	., 60	300 K	1075.3 μV	40.6 μV/K	10 mK	±0.038 K <sup>10</sup>	available from Lake	±20 mK
1				41.7 μV/K				
		600 K	13325 μV	41./ IIV/K	10 mK	±0.184 K <sup>10</sup>	Shore	±20 mK

Model 325 Temperature Controller

 $<sup>^{7}</sup>$  Typical sensor sensitivities were taken from representative calibrations for the sensor listed

<sup>&</sup>lt;sup>8</sup> Control stability of the electronics only, in an ideal thermal system

<sup>&</sup>lt;sup>9</sup> Non-HT version maximum temperature: 325 K

<sup>&</sup>lt;sup>10</sup> Accuracy specification does not include errors from room temperature compensation



#### Model 325 Specifications

#### **Input specifications**

	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy	Electronic control stability <sup>11</sup>
Diode	negative	0 V to 2.5 V	10 $\mu$ A $\pm 0.05\%^{12,13}$	100 μV	10 μV	$\pm 80~\mu V~\pm 0.005\%$ of rdg	±20 μV
_	negative	0 V to 7.5 V	10 μA ±0.05% <sup>12, 13</sup>	100 μV	20 μV	±320 μV ±0.01% of rdg	±40 μV
PTC RTD	positive	0 $\Omega$ to 500 $\Omega$	1 mA <sup>14</sup>	10 mΩ	2 mΩ	$\pm 0.004~\Omega~\pm 0.01\%$ of rdg	±4 mΩ
	positive	0 Ω to 5000 Ω	1 mA <sup>14</sup>	100 mΩ	20 mΩ	$\pm 0.04~\Omega~\pm 0.02\%$ of rdg	±40 mΩ
NTC RTD	negative	0 $\Omega$ to 7500 $\Omega$	$10 \ \mu A \pm 0.05\%^{14}$	100 mΩ	40 mΩ	$\pm 0.1~\Omega~\pm 0.04\%$ of rdg	±80 mΩ
Thermocouple	positive	±25 mV	NA	1 μV	0.4 μV	$\pm 1~\mu V~\pm 0.05\%$ of rdg $^{15}$	±0.8 μV
	positive	±50 mV	NA	1 μV	0.4 μV	$\pm 1~\mu V~\pm 0.05\%$ of rdg $^{15}$	±0.8 μV

<sup>&</sup>lt;sup>11</sup> Control stability of the electronics only, in an ideal thermal system

#### **Thermometry**

Number of inputs 2

**Input configuration** Each input is factory configured for either diode/RTD or thermocouple

**Isolation** Sensor inputs optically isolated from other circuits but not each other **A/D resolution** 24-bit

Input accuracy Sensor dependent—refer to Input Specifications table

 $\textbf{Measurement resolution} \ \textbf{Sensor dependent} - \textbf{refer to Input Specifications table}$ 

Maximum update rate 10 rdg/s on each input (except 5 rdg/s on input A when configured as thermocouple)

User curves Room for 15 200-point CalCurves  $^{\text{\tiny TM}}$  or user curves

SoftCal™ Improves accuracy of DT-470 diode to ±0.25 K

#### **Sensor input configuration**

	Diode/RTD	Thermocouple
Measurement type	4-lead differential	2-lead differential, room
		temperature compensated
Excitation	Constant current with current	N/A
	reversal for RTDs	
Supported sensors	Diodes: Silicon, GaAlAs	Most thermocouple types
	RTDs: 100 $\Omega$ Platinum, 1000 $\Omega$	
	Platinum, Germanium, Carbon-	
	Glass, Cernox®, and Rox™	
Standard curves	DT-470, DT-500D, DT-670,	Type E, Type K, Type T,
	PT-100, PT-1000, RX-102A,	AuFe 0.07% vs. Cr,
	RX-202A	AuFe 0.03% vs Cr
Input connector	6-pin DIN	Ceramic isothermal block

from 30 K to 375 K; improves accuracy of platinum RTDs to  $\pm 0.25$  K from 70 K to 325 K; stored as user curves

Filter Averages 2 to 64 input readings

#### **Control**

**Control loops** 2

Control type Closed loop digital PID with manual heater output or open loop

Tuning Autotune (one loop at a time), PID, PID zones

Control stability Sensor dependent—see Input Specification table

PID control settings

**Proportional (gain)** 0 to 1000 with 0.1 setting resolution

Integral (reset) 1 to 1000 (1000/s) with 0.1 setting resolution

**Derivative (rate)** 1 to 200% with 1% resolution

Manual output 0 to 100% with 0.01% setting resolution

Zone control 10 temperature zones with P, I, D, manual heater out, and heater range Setpoint ramping 0.1 K/min to 100 K/min

Safety limits Curve temperature, power up heater off, short circuit protection



<sup>&</sup>lt;sup>12</sup> Current source error has negligible effect on measurement accuracy

<sup>&</sup>lt;sup>13</sup> Diode input excitation current can be set to 1 mA — refer to the Model 325 user manual for details

<sup>&</sup>lt;sup>14</sup> Current source error is removed during calibration

<sup>&</sup>lt;sup>15</sup> Accuracy specification does not include errors from room temperature compensation



#### Loop 1 heater output

	25 Ω setting	50 Ω setting		
Туре	Variable DC current source			
D/A resolution	16-	-bit		
Max power	25	W		
Max current	1 A	0.71 A		
Voltage compliance	25 V	35.4 V		
Heater load range	20 $\Omega$ to 25 $\Omega$	40 $\Omega$ to 50 $\Omega$		
Heater load for max power	25 Ω	50 Ω		
Ranges	2 (2.5 V	V/25 W)		
Heater noise (<1 kHz)	1 μA + 0.01	% of output		
Grounding	Output referenced to chassis ground			
Heater connector	Dual banana			

#### Loop 2 heater output

Model 325 Temperature Controller

	25 Ω setting	50 Ω setting		
Туре	Variable DC v	oltage source		
D/A resolution	16-	-bit		
Max power	1 W	2 W		
Max voltage	5 V	10 V		
Current compliance	0.2 A			
Heater load range	≥ 25 Ω	≥ 50 Ω		
Heater load for max power	25 Ω	50 Ω		
Ranges	1			
Heater noise (<1 kHz)	50 μV + 0.01% of output			
Grounding	Output referenced to chassis ground			
Heater connector	Detachable terminal block			

#### **Front panel**

**Display** 2-line × 20-character, liquid crystal display with 5.5 mm character height

Number of reading displays 1 to 4

Display units K,  $^{\circ}$ C, V, mV,  $\Omega$ 

Reading source Temperature, sensor units

Display update rate 2 rdg/s

Temp display resolution  $0.001^\circ$  from  $0^\circ$  to  $99.999^\circ$ ,  $0.01^\circ$  from  $100^\circ$  to  $999.99^\circ$ ,  $0.1^\circ$ above 1000°

Sensor units display resolution Sensor dependent; to 5 digits

Other displays Setpoint, Heater Range, and Heater Output (user selected)

Setpoint setting resolution Same as display resolution (actual resolution is sensor dependent)

Heater output display Numeric display in percent of full scale for power or current **Heater output resolution 1%** 

Display annunciators Control Input, Remote, Autotune

Keypad 20-key membrane, numeric and specific functions

Front panel features Front panel curve entry, keypad lock-out

#### Interface

**IEEE-488** interface

**Features** SH1, AH1, T5, L4, SR1, RL1, PP0, DC1, DT0, C0, E1

Reading rate To 10 rdg/s on each input

Software support LabVIEW<sup>™</sup> driver (see www.lakeshore.com)

Serial interface

**Electrical format** RS-232C

**Baud rates** 9600, 19200, 38400, 57600 Connector 9-pin D-style, DTE configuration Reading rate To 10 rdg/s on each input

#### General

Ambient temperature 15 °C to 35 °C at rated accuracy, 5 °C to 40 °C at reduced accuracy

**Power requirement** 100, 120, 220, 240 VAC, +6%, -10%, 50 or 60 Hz, 85 VA

Size 216 mm W  $\times$  89 mm H  $\times$  368 mm D (8.5 in  $\times$  3.5 in  $\times$  14.5 in), half rack

Weight 4.00 kg (8.82 lb) Approval CE mark, RoHS

## Ordering information

Part number **Description** 

325 2 diode/resistor inputs temperature controller, includes one

dual banana jack heater output connector (106-009), two 6-pin DIN plug sensor input mating connectors (G-106-233), one 2-pin terminal block (106-735), a calibration certificate

and a user's manual

325-T1 Model 325 with one diode/RTD and one thermocouple input

Model 325 with two thermocouple inputs 325-T2

#### Please indicate your power/cord configuration:

1 100 V-U.S. cord (NEMA 5-15)

2 120 V-U.S. cord (NEMA 5-15)

3 220 V—Euro cord (CEE 7/7)

4 240 V—Euro cord (CEE 7/7) 5 240 V—U.K. cord (BS 1363)

6 240 V-Swiss cord (SEV 1011)

**7** 220 V—China cord (GB 1002)

#### **Accessories**

112-177 Temperature controller cable, 3 m (10 ft)—IN STOCK

112-178 Temperature controller cable, 6 m (20 ft) 112-180 Temperature controller cable, 10 m (33 ft)

117-017 1 m (3.3 ft long) IEEE-488 (GPIB) computer interface cable

assembly

RM-1/2 Kit for mounting one ½ rack temperature controller in a

482.6 mm (19 in) rack, 90 mm (3.5 in) high Kit for mounting two ½ rack temperature controllers in a

RM-2 482.6 mm (19 in) rack, 135 mm (5.25 in) high

G-106-735 Terminal block, 2-pin

G-106-233 Sensor input mating connector (6-pin DIN plug)

106-009 Banana plug, dual

CAL-325-CERT Instrument calibration with certificate

CAL-325-DATA Instrument recalibration with certificate and data Model 325 temperature controller manual 119-041













## Model 224 Temperature Monitor



#### **Model 224 features**

- Lake Shore's most capable cryogenic temperature monitor
- Equipped with 12 sensor channels for maximum monitoring capabilities
- Precisely measures in both higher temperature and cryogenic applications—down to 300 mK
- Ideal for multi-sensor lab uses, particularly for monitoring Cernox® sensors
- Ethernet, USB and IEEE-488 computer interfaces
- Proven, intuitive user interface
- Customizable display enables you to label individual input channels
- CE certification
- Full 3 year standard warranty





#### Introduction

The Lake Shore Model 224 temperature monitor offers precision measurement in a wide range of cryogenic and higher-temperature applications with the ability to easily monitor up to 12 sensor channels. It provides better measurement performance in applications where researchers need to ensure accuracy and precision in their low cryogenic temperature monitoring. Used with Lake Shore's Cernox® sensors, the Model 224 enables reliable and repeatable temperature measurement over a broad range and as low as 300 mK.

Cernox thin-film RTD sensors offer high sensitivity and low magnetic field-induced errors at cryogenic temperatures. The Model 224 has been optimized for use with these well-respected temperature sensors, and features many of the same advanced capabilities of Lake Shore Model 336 temperature controller, including its proven high-precision input circuitry.

In addition to Cernox, the Model 224 supports other NTC RTDs, PTC RTDs such as platinum sensors, and diodes such as the Lake Shore DT-670 Series. In cryogenic applications, the monitor is an ideal addition to any university or commercial low-temperature research lab requiring measurement flexibility using multiple sensors and sensor types. Used with silicon diodes, it provides accurate measurements in cryo-cooler and cryo-gas production applications from 1.4 K to above room temperature. Connected to PTC RTDs (platinum and rhodium-iron sensors), the Model 224 works well in cryogenic applications at liquid nitrogen temperatures.



You can set up different sensor types and responses on each input to support simultaneous measurement of various critical points in a system. Examples include monitoring multiple cryogenic refrigeration systems (e.g., liquid nitrogen Dewars, He-4 cryostats, and closed-cycle refrigerators), multiple stages within systems operating at different temperature levels, thermal gradient profiling, redundant measurements of critical values, leak detection, and other cryogenic applications where you need accurate readings at multiple points. Alarm thresholds can be configured independently for each input, and alarm events can activate the unit's relay outputs for hard-wired triggering of other systems or audible annunciators. Relays can be activated on high, low, or both alarms for any input.

#### **Configure each input independently**

Because the Model 224 features 12 independently configurable 6-pin DIN inputs, you can set it up for a different sensor on each input and run a number of different measurements simultaneously for various critical points in a system. Two inputs (A and B) are dedicated and non-scanned, updated at 10 rdg/s. The remaining 10 are scanned channels—inputs C and D can have up to five input devices each. These scanned channels are read anywhere from 1 to 10 rdg/s, depending on how many are being used at once.



Press any of the 4 input buttons (A, B, C or D) to view or change the parameters for each channel in the input display mode

The Model 224 features four high-resolution, 24-bit analog-to-digital converters for fast measurements. Optical isolation of input circuitry reduces line noise—interference that can skew low-level measurements—while providing repeatable sensor measurements.

Current reversal eliminates thermal electromotive force (EMF) errors when using resistance sensors. Also, nine excitation currents enable temperature measurements down to 300 mK when you use the appropriate NTC RTDs. When autoranging is enabled, the range will be automatically selected so that the excitation voltage is below 10 mV. This keeps the power dissipated in the sensor at a minimum, yet still at enough of a level to provide accurate measurements.

#### Monitor locally or remotely—from anywhere

For local monitoring, the front panel of the Model 224 features a bright liquid crystal display with an LED backlight that shows up to 12 readings simultaneously, or, you can even display a single sensor input to see greater detail at a glance.

Plus, monitoring can be done over a network. Using the Ethernet port on the Model 224, you can keep an eye on temperatures and log measurement data remotely via a networked local PC or even remotely over a TCP/IP Internet connection from anywhere. A chart recorder utility embedded in the Ethernet module enables real-time charting of temperatures using a convenient graphical interface. You can also interface with the temperature monitor or link it to a data acquisition system via its serial USB or parallel IEEE-488 ports.



#### Intuitive, configurable display

The Model 224 front panel features a 23-key keypad and intuitive user interface for easy navigation of the temperature monitor's functions.

For added convenience, you can also custom label each sensor input, eliminating the guesswork in remembering or determining the location to which a sensor input is associated.

A:	Sample Space	B:	Rad Shield
	5.2014 <sub>K</sub>		27.8645k
C	Second Rad Shield	D:	4 K Shield
	26.2431 <sub>K</sub>		25.8957k

Custom display modes can show multiple configurations of channels. The display above shows the 4 main inputs with their custom labels, while the one below shows all 12 channels plus 4 additional settings

A:	0.5410 K	B: 1.4072	K
C1:	4.6103 K	D1: 40.0471	
C2:	5.2014 K	D2: 77.2116	K
C3:	27.8645 K	D3: 80.2037	K
C4:	35.6387 K	D4: 88.7783	AIS .
C5:	36.4864 K	D5: 89.6519	
C5:	36.0622 K Mn	D5: 89.2277	K Mn
C5:	40.7284 K Mx	D5: 93.8939	K Mx

#### **Stores response curves**

Like the Lake Shore Model 336, the Model 224 includes standard temperature sensor calibration curves for silicon diodes, platinum RTDs, and  $Rox^{TM}$  (ruthenium oxide) RTDs.

The monitor's non-volatile memory enables users to store up to 39 200-point CalCurves for Lake Shore calibrated sensors or user curves. Lake Shore also offers curve handler software, which allows you to upload and manipulate temperature sensor calibration data.

#### Sensor Selection


#### Sensor temperature range (sensors sold separately)

		Model	Useful range	Magnetic field use
Negative	Cernox®	CX-1010	0.3 K to 325 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
temperature	Cernox®	CX-1030-HT	0.3 K to 420 K <sup>1, 2</sup>	$T > 2 K \& B \le 19 T$
coefficient	Cernox®	CX-1050-HT	1.4 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
RTDs	Cernox®	CX-1070-HT	4 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1080-HT	20 K to 420 K <sup>1</sup>	$T > 2 K \& B \le 19 T$
	Germanium	GR-300-AA	0.35 K to 100 K <sup>2</sup>	Not recommended
	Germanium	GR-1400-AA	1.8 K to 100 K <sup>2</sup>	Not recommended
	Rox™	RX-102	0.3 K to 40 K <sup>2</sup>	$T > 2 K \& B \le 10 T$
	Rox™	RX-103	1.4 K to 40 K	$T > 2 K \& B \le 10 T$
	Rox™	RX-202	0.3 K to 40 K <sup>2</sup>	$T > 2 K \& B \le 10 T$
Diodes	Silicon diode	DT-670-SD	1.4 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-670E-BR	30 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-414	1.4 K to 375 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-421	1.4 K to 325 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-470-SD	1.4 K to 500 K	$T \geq 60~K~\&~B \leq 3~T$
	Silicon diode	DT-471-SD	10 K to 500 K	$T \geq 60~K~\&~B \leq 3~T$
	GaAlAs diode	TG-120-P	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-PL	1.4 K to 325 K	$T > 4.2 \text{ K & B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-SD	1.4 K to 500 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
Positive	100 Ω platinum	PT-102/3	14 K to 873 K	$T > 40$ K & B $\leq 2.5T$
temperature	100 Ω platinum	PT-111	14 K to 673 K	$T > 40$ K & B $\leq 2.5T$
coefficient	Rhodium-iron	RF-800-4	1.4 K to 500 K	$T > 77 \text{ K & B} \le 8 \text{ T}$
RTDs	Rhodium-iron	RF-100T/U	1.4 K to 325 K	$T > 77 \text{ K & B} \le 8 \text{ T}$

- <sup>1</sup> Non-HT version maximum temperature: 325 K
- $^{2}\;$  Low temperature specified with self-heating error:  ${\leq}5\;\text{mK}$

#### Model 224 rear panel

- 1 Sensor input connectors
- 2 Line input assembly
- 3 Terminal block
- 4 Ethernet interface
- 5 USB interface
- 6 IEEE-488 interface

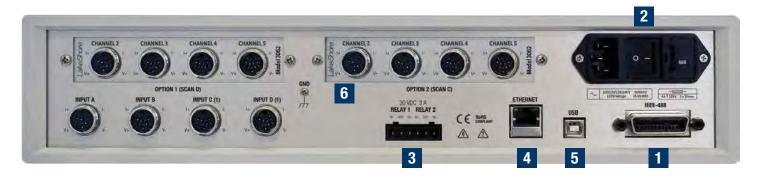



#### Ideal applications

- Labs with multiple temperature sensors
- Applications where both cryogenic and higher temperature readings are required
- Monitoring of simple Dewars and LN cryostats (>4.2 K)
- Closed-cycle refrigerators (CCRs) at 3 K to 4 K
- Pumped He-4 (1.4 K) and He-3 (300 mK) systems
- Temperature monitoring where superconducting magnets are used, such as in mass spectrometer and particle accelerator equipment



- Low magnetic field-induced errors
- A temperature range of 100 mK to 420 K (model dependent)

Model 224 Temperature Monitor

- High sensitivity at low temperatures and good sensitivity over a broad range
- Excellent resistance to ionizing radiation
- Bare die cryogenic temperature sensor with fast characteristic thermal response times: 1.5 ms at 4.2 K, 50 ms at 77 K
- Broad selection of models to meet your thermometry needs
- **Excellent stability**
- A variety of packaging options

These thin-film resistance cryogenic sensors offer significant advantages over diodes and conventional RTD sensors. The smaller package size makes them useful in a wide range of experimental mounting schemes, and they are also available in a chip form.






## Model 224 Specifications

#### **Input specifications**

	Sensor temperature coefficient	Input Range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy (at 25 °C)	Measurement temperature coefficient
NTC RTD		0 $\Omega$ to 10 $\Omega$	1 mA <sup>3</sup>	0.1 mΩ	0.15 mΩ	$\pm 0.002~\Omega~\pm 0.06\%$ of rdg	$(0.01 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
10 mV		0 $\Omega$ to 30 $\Omega$	300 μA <sup>3</sup>	0.1 mΩ	$0.45~\text{m}\Omega$	$\pm 0.002~\Omega~\pm 0.06\%$ of rdg	$(0.03 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 100 $\Omega$	100 μA <sup>3</sup>	1 mΩ	$1.5~\text{m}\Omega$	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	$(0.1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 300 $\Omega$	30 μA³	1 mΩ	$4.5~\text{m}\Omega$	$\pm 0.01~\Omega~\pm 0.04\%$ of rdg	$(0.3 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$
	Negative	0 $\Omega$ to 1 k $\Omega$	10 μA³	10 mΩ	15 m $\Omega$ +0.002% of rdg	$\pm 0.1~\Omega~\pm 0.04\%$ of rdg	$(1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 3 k $\Omega$	3 μA³	10 mΩ	45 m $\Omega$ +0.002% of rdg	$\pm 0.1~\Omega~\pm 0.04\%$ of rdg	$(3 \text{ m}\Omega + 0.0015\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 10 k $\Omega$	1 μA³	100 mΩ	150 m $\Omega$ +0.002% of rdg	$\pm 1.0~\Omega~\pm 0.04\%$ of rdg	$(10 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 30 $k\Omega$	300 nA <sup>3</sup>	100 mΩ	450 m $\Omega$ +0.002% of rdg	$\pm 2.0~\Omega~\pm 0.04\%$ of rdg	$(30 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 100 $k\Omega$	100 nA <sup>3</sup>	1 Ω	$1.5 \Omega + 0.005\%$ of rdg	$\pm 10.0~\Omega~\pm 0.04\%$ of rdg	$(100 \text{ m}\Omega + 0.002\% \text{ of rdg})/^{\circ}\text{C}$
Diode	Negative	0 V to 2.5 V	$10  \mu A \pm 0.05\%^{4,5}$	10 μV	10 μV	$\pm 80~\mu V~\pm 0.005\%$ of rdg	$(10 \mu V + 0.0005\% \text{ of rdg})/^{\circ}C$
	iveyalive	0 V to 10 V	$10  \mu A \pm 0.05\%^{4,5}$	100 μV	20 μV	$\pm 320~\mu V~\pm 0.01\%$ of rdg	$(20 \mu V + 0.0005\% \text{ of rdg})/^{\circ}C$
PTC RTD		0 $\Omega$ to 10 $\Omega$	1 mA <sup>3</sup>	0.1 mΩ	$0.2~\text{m}\Omega$	$\pm 0.002~\Omega~\pm 0.01\%$ of rdg	$(0.01 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 30 $\Omega$	1 mA <sup>3</sup>	$0.1~\text{m}\Omega$	$0.2~\text{m}\Omega$	$\pm 0.002~\Omega~\pm 0.01\%$ of rdg	$(0.03 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 100 $\Omega$	1 mA <sup>3</sup>	1 mΩ	2 mΩ	$\pm 0.004~\Omega~\pm 0.01\%$ of rdg	$(0.1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
	Positive	0 $\Omega$ to 300 $\Omega$	1 mA <sup>3</sup>	1 mΩ	2 mΩ	$\pm 0.004~\Omega~\pm 0.01\%$ of rdg	$(0.3 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 1 $k\Omega$	1 mA <sup>3</sup>	10 mΩ	20 mΩ	$\pm 0.04~\Omega~\pm 0.02\%$ of rdg	$(1 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
		0 $\Omega$ to 3 $k\Omega$	1 mA <sup>3</sup>	10 mΩ	20 mΩ	$\pm 0.04~\Omega~\pm 0.02\%$ of rdg	$(3 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$
		0 Ω to 10 kΩ	1 mA <sup>3</sup>	100 mΩ	200 mΩ	$\pm 0.4~\Omega~\pm 0.02\%$ of rdg	$(10 \text{ m}\Omega + 0.001\% \text{ of rdg})/^{\circ}\text{C}$

- <sup>3</sup> Current source error is removed during calibration
- <sup>4</sup> Current source error has negligible effect on measurement accuracy
- <sup>5</sup> Diode input excitation can be set to 1 mA

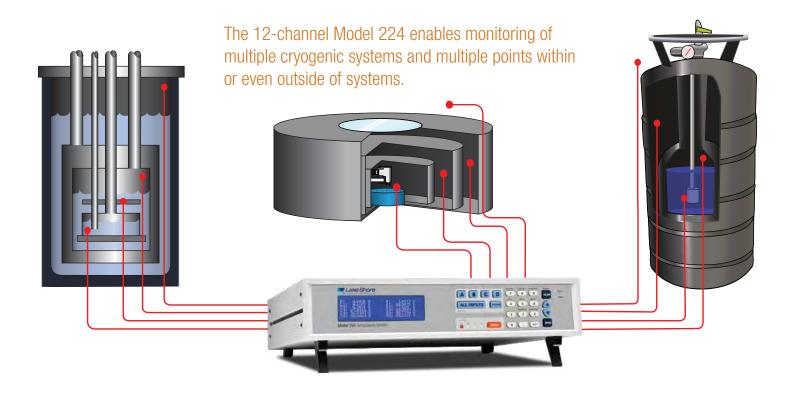



#### **Sensor input configuration**

	Diode/RTD
Measurement type	4-lead differential
Excitation	Constant current with current reversal for RTDs
Supported sensors	RTDs: Cernox®,100 Ω platinum,1000 Ω platinum, germanium, carbon-glass, and Rox™ Diodes: silicon, GaAlAs
Standard curves	DT-470, DT-670, DT-500-D, DT-500-E1, PT-100, PT-1000, RX-102A, RX-202A
Input connector	6-pin DIN

#### **Thermometry**

Number of inputs 12 (2 dedicated; 10 scanned)

Input configuration Inputs can be configured independently from the front panel to accept any of the supported input types

Isolation Sensor inputs optically isolated from other circuits but not from each other A/D resolution 24-bit

Input accuracy Sensor dependent, refer to Input Specifications table

Measurement resolution Sensor dependent, refer to Input Specifications table

Maximum update rate 10 rdg/s on each non-scanned input; 5 rdg/s when configured as 100 k $\Omega$  NTC RTD with reversal on; 2 rdg/s on each scanned input; update rate is dependent on the number of channels enabled (typically from 10 rdg/s for 1 channel to 2 rdg/s for all 10 scanned channels)

Autorange Automatically selects appropriate NTC RTD or PTC RTD range

**User curves** Room for 39 200-point CalCurves<sup>™</sup> or user curves

SoftCal<sup>™</sup> Improves accuracy of DT-470 diode to ±0.25 K from 30 K to 375 K; improves accuracy of platinum RTDs to ±0.25 K from 70 K to 325 K; stored as user curves

Math Maximum and minimum

Filter Averages 2 to 64 input readings

#### Front panel

Display 8-line by 40-character (240  $\times$  64 pixel) LCD display module with LED backlight

Number of reading displays 1 to 16

Display units K, °C, V, mV, Ω

Reading source Temperature, sensor units, max, and min

Display update rate 2 rdg/s

Temperature display resolution 0.0001° from 0° to 99.9999°, 0.001° from 100° to 999.999°, 0.01° above 1000°

Sensor units display resolution Sensor dependent, to 6 digits

Other displays Input name

**Display annunciators** Alarm

LED annunciators Remote, Ethernet status, alarm

Keypad 23-key silicone elastomer keypad

Front panel features Front-panel curve entry, display contrast control, and keypad lockout

#### **Interface**

Capabilities SH1, AH1, T5, L4, SR1, RL1, PP0, DC1, DT0, C0, E1

Reading rate To 10 rdg/s on each input

LabVIEW™ driver (see www.lakeshore.com) Software support

**Function** Emulates a standard RS-232 serial port

**Baud rate** 57.600 Connector B-type USB

Reading rate To 10 rdg/s on each input Software support LabVIEW™ driver (see www.lakeshore.com)

**Function** TCP/IP, web interface with built-in utilities

Connector **RJ-45** 

Reading rate To 10 rdg/s on each input

Software support LabVIEW™ driver (see www.lakeshore.com)

Number 12, high and low for each input Data source Temperature or sensor units

Settings Source, high setpoint, low setpoint, deadband, latching or

non-latching, audible on/off, and visible on/off

**Actuators** Display annunciator, beeper, and relays

Number

**Contacts** Normally open (NO), normally closed (NC), and common (C)

**Contact rating** 30 VDC at 3 A

Operation Activate relays on high, low, or both alarms for any input,

or manual mode

Connector Detachable terminal block

#### General

Model 224 Temperature Monitor

Ambient temperature 15 °C to 35 °C at rated accuracy: 5 °C to 40 °C at reduced accuracy

Power requirement 100, 120, 220, 240 VAC, ±10%, 50 or 60 Hz, 35 VA Size 435 mm W  $\times$  89 mm H  $\times$  368 mm D (17 in  $\times$  3.5 in  $\times$  14.5 in), full rack

Weight 7.6 kg (16.8 lb) Approval CE mark, RoHS

## Ordering information

#### Part number **Description**

224-12 Temperature monitor with 12 diode/RTD inputs—includes

twelve 6-pin DIN plug sensor input mating connectors (G-106-233), one 6-pin terminal block (106-737), a calibration

certificate and a user's manual

#### Please indicate your power/cord configuration:

1 100 V-U.S. cord (NEMA 5-15)

2 120 V-U.S. cord (NEMA 5-15)

3 220 V—Euro cord (CEE 7/7)

4 240 V—Euro cord (CEE 7/7)

**5** 240 V—U.K. cord (BS 1363)

6 240 V-Swiss cord (SEV 1011)

**7** 220 V—China cord (GB 1002)

#### **Accessories**

RM-1 Kit for mounting one full rack instrument G-106-233 Sensor input mating connector (6-pin DIN plug)

106-737 Terminal block (6-pin)

CAL-224-CERT Instrument calibration with certificate

CAL-224-DATA Instrument recalibration with certificate and data 119-062 Model 224 temperature monitor user manual











- CE certification
- Full 3 year standard warranty





#### Introduction

The Model 218 is our most versatile temperature monitor. With eight sensor inputs, it can be used with nearly any diode or resistive temperature sensor. It displays all eight channels continuously in K, °C, V or  $\Omega$ . The measurement input was designed for the demands of cryogenic temperature measurement, however, the monitor's low noise, high resolution, and wide operating range make it ideal for noncryogenic applications as well.

#### Sensor input reading capability

The Model 218 has eight constant current sources (one for each input) that can be configured for a variety of sensors. The inputs can be configured from the front panel or via a computer interface, and are grouped in two sets of four. Each set of four inputs is configured for the same sensor type (i.e., all 100  $\Omega$  platinum or all silicon diodes).

Two high-resolution A/D converters increase the update rate of the Model 218. It can read sensor inputs more quickly than other scanning monitors because it does not have to wait for current source switching. The result is 16 new readings per second, allowing all inputs to be read twice each second. Inputs can be turned off to obtain a higher reading rate on fewer sensors.

#### **Temperature response curves**

The Model 218 has standard temperature sensor response curves for silicon diodes and platinum RTDs. It can support a wide variety of temperature sensors because a unique 200-point user curve can be stored for each of the eight inputs. CalCurves™ for Lake Shore calibrated sensors can be stored as user curves.

Model 218 Temperature Monitor

The built in SoftCal<sup>™1</sup> algorithm can also be used to generate improved curves for DT-670 diodes and platinum RTDs that are stored as user curves.

#### Interface

The Model 218 is available with both parallel (IEEE-488, 218S only) and serial (RS-232C) computer interfaces. Each input has a high and low alarm which offer latching and non-latching operation. The eight relays on the Model 218S can be used with the alarms to alert the operator of a fault condition or perform simple on-off control. The Model 218S includes two analog voltage outputs. The user may select the scale and data sent to the output, including temperature, sensor units, or linear equation results. Under manual control, the analog voltage output can also serve as a voltage source for other applications.

#### Interface features of the Model 218S and Model 218E


	218S	218E
Numeric keypad		
Front panel curve entry	•	•
Alarms		
RS-232C interface		
IEEE-488 interface		
Two analog voltage outputs		
Eight relays		

#### **Display**

The eight display locations on the Model 218 are user configurable. Sources for readout data are temperature units, sensor units, and results of the math function. Input number and data source are always displayed for convenience. The display is updated twice each second.

<sup>1</sup> The Lake Shore SoftCal<sup>™</sup> algorithm for silicon diode and platinum RTD sensors is a good solution for applications requiring more accuracy than a standard sensor curve but not in need of traditional calibration. SoftCal uses the predictability of a standard curve to improve the accuracy of an individual sensor around a few known temperature reference points.

#### Model 218 rear panel

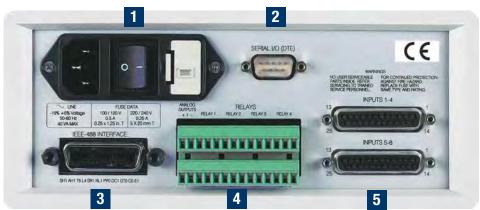

- Line input assembly
- 2 RS-232C or printer interface
- 3 IEEE-488 interface (218S)
- 4 Terminal block with relays and analog voltage outputs (218S only)



#### **Sensor Selection**

#### Sensor temperature range (sensors sold separately)

		Model	Useful range	Magnetic field use
Diodes	Silicon diode	DT-670-SD	1.4 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-670E-BR	30 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-414	1.4 K to 375 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-421	1.4 K to 325 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-470-SD	1.4 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-471-SD	10 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	GaAlAs diode	TG-120-P	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-PL	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-SD	1.4 K to 500 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
Positive temperature	100 Ω platinum	PT-102/3	14 K to 873 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
coefficient RTDs	100 Ω platinum	PT-111	14 K to 673 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
	Rhodium-iron	RF-800-4	1.4 K to 500 K	$T > 77 \text{ K \& B} \le 8 \text{ T}$
	Rhodium-iron	RF-100T/U	1.4 K to 325 K	$T > 77 \text{ K & B} \le 8 \text{ T}$
Negative	Cernox®	CX-1010	2 K to 325 K <sup>5</sup>	$T > 2 K \& B \le 19 T$
temperature	Cernox®	CX-1030-HT	3.5 K to 420 K <sup>3, 6</sup>	$T > 2 K \& B \le 19 T$
coefficient RTDs <sup>2</sup>	Cernox®	CX-1050-HT	4 K to 420 K <sup>3, 6</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1070-HT	15 K to 420 K <sup>3</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1080-HT	50 K to 420 K <sup>3</sup>	$T > 2 K \& B \le 19 T$
	Germanium	GR-300-AA	1.2 K to 100 K <sup>4</sup>	Not recommended
	Germanium	GR-1400-AA	4 K to 100 K <sup>4</sup>	Not recommended
	Rox™	RX-102A	1.4 K to 40 K <sup>5</sup>	$T > 2 K \& B \le 10 T$

<sup>&</sup>lt;sup>2</sup> Single excitation current may limit the low temperature range of NTC resistors

Silicon diodes are the best choice for general cryogenic use from 1.4 K to above room temperature. Diodes are economical to use because they follow a standard curve and are interchangeable in many applications. They are not suitable for use in ionizing radiation or magnetic fields.

**Cernox®** thin-film RTDs offer high sensitivity and low magnetic field-induced errors over the 2 K to 420 K temperature range. Cernox sensors require calibration.

**Platinum RTDs** offer high uniform sensitivity from 30 K to over 800 K. With excellent reproducibility, they are useful as thermometry standards. They follow a standard curve above 70 K and are interchangeable in many applications.

<sup>&</sup>lt;sup>3</sup> Non-HT version maximum temperature: 325 K

<sup>&</sup>lt;sup>4</sup> Low temperature limited by input resistance range

<sup>&</sup>lt;sup>5</sup> Low temperature specified with self-heating error: ≤ 5 mK

 $<sup>^{\</sup>rm 6}~$  Low temperature specified with self-heating error:  $\leq$  12 mK



#### Typical sensor performance—see Appendix F for sample calculations of typical sensor performance

	Example Lake Shore sensor	Temperature	Nominal resistance/ voltage	Typical sensor sensitivity <sup>7</sup>	Measurement resolution: temperature equivalents	Electronic accuracy: temperature equivalents	Temperature accuracy including electronic accuracy, CalCurve™, and calibrated sensor
Silicon diode	DT-670-SD	1.4 K	1.644 V	-12.49 mV/K	1.6 mK	±26 mK	±38 mK
	with 1.4H	77 K	1.028 V	-1.73 mV/K	11.6 mK	±152 mK	±174 mK
	calibration	300 K	0.5597 V	-2.3 mV/K	8.7 mK	±94 mK	±126 mK
		500 K	0.0907 V	-2.12 mV/K	9.4 mK	±80 mK	±130 mK
Silicon diode	DT-470-SD-13	1.4 K	1.6981 V	-13.1 mV/K	1.5 mK	±26 mK	±38 mK
	with 1.4H	77 K	1.0203 V	-1.92 mV/K	10.5 mK	±137 mK	±159 mK
	calibration	300 K	0.5189 V	-2.4 mV/K	8.4 mK	±88 mK	±120 mK
		475 K	0.0906 V	-2.22 mV/K	9.1 mK	±77 mK	±127 mK
GaAlAs diode	TG-120-SD	1.4 K	5.391 V	-97.5 mV/K	0.2 mK	±13 mK	±25 mK
	with 1.4H	77 K	1.422 V	-1.24 mV/K	16.2 mK	±359 mK	±381 mK
	calibration	300 K	0.8978 V	-2.85 mV/K	7 mK	±120 mK	±152 mK
		475 K	0.3778 V	-3.15 mV/K	6.4 mK	±75 mK	±125 mK
100 Ω platinum RTD	PT-103	30 K	3.66 Ω	0.19 Ω/Κ	10.5 mK	±25 mK	±35 mK
500 Ω full scale	with 1.4J	77 K	20.38 Ω	0.42 Ω/Κ	4.8 mK	±20 mK	±32 mK
000 11 1411 00410	calibration	300 K	110.35 Ω	0.39 Ω/Κ	5.2 mK	±68 mK	±91 mK
		500 K	185.668 Ω	0.378 Ω/Κ	5.3 mK	±109 mK	±155 mK
Cernox®	CX-1050-SD-HT8	4.2 K	3507.2 Ω	-1120.8 Ω/K	45 μK	±1.4 mK	±6.4 mK
0011101	with 4M	77 K	205.67 Ω	-2.4116 Ω/K	20.8 mK	±75.6 mK	±91.6 mK
	calibration	300 K	59.467 Ω	-0.1727 Ω/K	290 mK	±717 mK	±757 mK
		420 K	45.03 Ω	-0.0829 Ω/K	604 mK	±1.43 K	±1.5 K
Germanium	GR-300-AA	1.2 K	600 Ω	-987 Ω/K	51 μK	±0.3 mK	±4.5 mK
	with 0.3D	1.4 K	449 Ω	-581 Ω/K	86 µK	±0.5 mK	±4.7 mK
	calibration	4.2 K	94 Ω	-27 Ω/K	1.9 mK	±5.2 mK	±10.2 mK
		100 K	3 Ω	-0.024 Ω/K	2.10 K	±4.25 K	±4.27 K
Germanium	GR-1400-AA	4 K	1873 Ω	-1008 Ω/K	50 μK	±0.8 mK	±5.0 mK
dormaniani	with 1.4D	4.2 K	1689 Ω	-862 Ω/K	58 μK	±0.9 mK	±5.1 mK
	calibration	10 K	253 Ω	-62 Ω/K	807 μK	±3.2 mK	±8.2 mK
		100 K	3 Ω	-0.021 Ω/K	2.40 K	±4.86 K	±4.88 K
Carbon-glass	CGR-1-2000	4.2 K	2260 Ω	-2060 Ω/K	25 μK	±0.5 mK	±4.5 mK
(no longer available)	with 4L	77 K	21.65 Ω	-0.157 Ω/K	319 mK	±692 mK	±717 mK
( Singor available)	calibration	300 K	11.99 Ω	-0.015 Ω/K	3.33 K	±7 K	±7.1 K

Model 218 Temperature Monitor

## Specifications

#### **Thermometry**

Number of inputs 8

Input configuration Inputs separated into two groups of four (each group must be the same sensor type) – inputs can be configured from the front panel to accept any of the supported input types

Input accuracy Sensor dependent—refer to Input Specifications table

 $\textbf{\textit{Measurement resolution}} \ \textbf{Sensor dependent} \textbf{--} \textbf{refer to Input Specifications table}$ 

Maximum update rate 16 readings per s total

**User curves** Room for 8 (1 per input) 200-point CalCurves<sup>™</sup> or user curves

SoftCal™ Improves accuracy of DT-470 diode to ±0.25 K from 30 K to 375 K; improves accuracy of platinum RTDs to ±0.25 K from 70 K to 325 K; stored as user curves

 $\textbf{Math} \ \text{Maximum, minimum, and linear equation } (Mx+B) \ \text{or} \ M(x+B)$ 

Filter Averages 2 to 64 input readings

#### **Front panel**

Display 4 line by 20 character backlit LCD display

Number of reading displays 1 to 8  $\,$ 

Display units K, °C, V, and  $\Omega$ 

Reading source Temperature, sensor units, max, min, and linear equation

Display update rate All displayed inputs twice in 1 s

Temp display resolution 0.001° from 0° to 99.999°, 0.01° from 100° to 999.99°, 0.1° above 1000°

Sensor units display resolution Sensor dependent to 5 digits

**Display annunciators** Remote operation, alarm, data logging, max, min, and linear

Keypad Membrane keypad, 20-key, numeric and specific functions

Front panel features Front panel curve entry and keypad lock-out

<sup>&</sup>lt;sup>7</sup> Typical sensor sensitivities were taken from representative calibrations for the sensor listed

<sup>&</sup>lt;sup>8</sup> Non-HT version maximum temperature: 325 K





#### Input specifications

	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy
Diode	negative -	0 V to 2.5 V	10 μA ±0.05% <sup>9</sup>	100 μV	20 μV	$\pm 200 \mu V \pm 0.01\%$ of rdg
	negauve	0 V to 7.5 V	10 μA ±0.05% <sup>9</sup>	100 μV	20 μV	$\pm 350~\mu V~\pm 0.02\%$ of rdg
PTC RTD		0 $\Omega$ to 250 $\Omega$	1 mA ±0.3% <sup>10</sup>	10 mΩ	2 mΩ	$\pm 0.06~\Omega~\pm 0.02\%$ of rdg
	positive	0 $\Omega$ to 500 $\Omega$	1 mA ±0.3% <sup>10</sup>	10 mΩ	2 mΩ	$\pm 0.06~\Omega~\pm 0.02\%$ of rdg
		0 $\Omega$ to 5000 $\Omega$	1 mA ±0.3% <sup>10</sup>	100 mΩ	20 mΩ	$\pm 0.4~\Omega~\pm 0.04\%$ of rdg
NTC RTD	negative	0 $\Omega$ to 7500 $\Omega$	10 μA ±0.05% <sup>9</sup>	100 mΩ	50 mΩ	$\pm 0.8~\Omega~\pm 0.04\%$ of rdg

<sup>&</sup>lt;sup>9</sup> Current source error has negligible effect on measurement accuracy

**Cryogenic Instruments** 

#### **Sensor input configuration**

	Diode/RTD
Measurement type	4-lead differential
Excitation	8 constant current sources
Supported sensors	Diodes: Silicon, GaAlAs RTDs: 100 Ω Platinum, 1000 Ω Platinum, Germanium, Carbon-Glass, Cernox®, and Rox™
Standard curves	DT-470, DT-500D, DT-670, CTI-C, PT-100, and PT-1000
Input connector	25-pin D-sub

#### Interface

IEEE-488.2 interface (218S)

**Features** SH1, AH1, T5, L4, SR1, RL1, PP0, DC1, DT0, C0, E1

Reading rate To 16 rda/s LabVIEW<sup>™</sup> driver Software support

Serial interface

RS-232C **Electrical format** Max baud rate 9600 baud Connector 9-pin D-sub

Reading rate To 16 readings per s (at 9600 baud)

**Printer capability** Support for serial printer through serial interface port used

with data log parameters

**Alarms** 

Number 16: high and low for each input

Data source Temperature, sensor units, and linear equation Settings Source, high setpoint, low setpoint, deadband, latching or non-latching, and audible on/off

**Actuators** Display annunciator, beeper, and relays (218S)

Relays (218S)

Number

Contacts Normally open (NO), normally closed (NC), and common (C)

**Contact rating** 30 VDC at 5 A

Operation Each input may be configured to activate any or all of the eight relays-relays may be activated on high, low, or both alarms

for any input, or manually

Detachable terminal block Connector

Analog voltage output (218S)

Number

Scale User selected Update rate To 16 rda/s

Temperature, sensor units, and linear equation Data source

Range  $\pm 10 V$ Resolution 1.25 mV **Accuracy** ±2.5 mV

Min load resistance 1 kΩ (short-circuit protected)

**Data logging** 

Channels

Operation Data log records can be stored in memory or sent to the printer:

stored data may be displayed, printed, or retrieved by

computer interface

**Data memory** Maximum of 1500 single reading records, non-volatile

#### General

Ambient temperature 15 °C to 35 °C at rated accuracy, 10 °C to 40 °C at reduced accuracy Power requirement 100, 120, 220, 240 VAC, (+6%, -10%), 50 or 60 Hz, 18 VA Size 216 mm W  $\times$  89 mm H  $\times$  318 mm D (8.5 in  $\times$  3.5 in  $\times$  12.5 in), half rack Weight 3 kg (6.6 lb) Approval CE mark, RoHS

#### Ordering information

#### Part number **Description**

**218S** Standard temperature monitor (8 inputs, IEEE-488 and serial

interface, alarms, relays, corrected analog output, data logging)—includes two 25-pin D-sub sensor input plugs (G-106-253), two 25-pin D-sub sensor input shells (G-106-264), two 14-pin relay/analog output conectors (106-772), a

calibration certificate and a user's manual

218E Economy temperature monitor (8 inputs, serial interface,

alarms, data logging)—includes same accessories as the

**218S** 

#### Please indicate your power/cord configuration:

1 100 V-U.S. cord (NEMA 5-15)

2 120 V-U.S. cord (NEMA 5-15)

3 220 V—Euro cord (CEE 7/7)

240 V—Euro cord (CEE 7/7)

5 240 V—U.K. cord (BS 1363)

6 240 V-Swiss cord (SEV 1011)

**7** 220 V—China cord (GB 1002)

#### **Accessories**

4005 1 m IEEE-488 (GPIB) computer interface cable assembly—

includes extender which allows connection of IEEE cable and

relay terminal block simultaneously

RM-1/2 Kit for mounting one half rack instrument RM-2 Kit for mounting two half rack instruments

DB-25 plug, qty 1 G-106-253 G-106-264 DB-25 hood; qty 1

106-772 Terminal block mating connector, 14-pin connector, 218S only

8000 The CalCurve<sup>™</sup> breakpoint table from a calibrated sensor

loaded on a CD-ROM for customer uploading

8002-05-218 The breakpoint table from a calibrated sensor stored in a

NOVRAM for installation at the customer location

CAL-218-CERT Instrument calibration with certificate

CAL-218-DATA Instrument recalibration with certificate and data 119-007 Model 218 temperature monitor manual



<sup>&</sup>lt;sup>10</sup> Current source error is removed during calibration



# Model 211 Temperature Controller



Model 211 Temperature Monitor

#### **Model 211 features**

- Operates down to 1.2 K with appropriate sensor
- One sensor input
- Supports diode and RTD sensors
- 0 V to 10 V or 4 mA to 20 mA output
- Large 5-digit LED display
- RS-232C serial interface and alarm relays
- CE certification
- Full 3 year standard warranty





#### Introduction

The Lake Shore single-channel Model 211 temperature monitor provides the accuracy, resolution, and interface features of a benchtop temperature monitor in an easy to use, easily integrated, compact instrument. With appropriate sensors, it measures from 1.2 K to 873 K, including temperatures in high vacuum and magnetic fields. Alarms, relays, user-configurable analog voltage or current output, and a serial interface are standard features on the Model 211. It is a good choice for liquefied gas storage and monitoring, cryopump control, cryo-cooler, and materials science applications, and when you need greater accuracy than thermocouples allow.

#### Sensor input reading capability

The Model 211 temperature monitor supports diode temperature sensors and resistance temperature detectors (RTDs). It can be configured for the type of sensor in use from the instrument front panel. Ensuring high accuracy and 5-digit measurement resolution are 4-lead differential measurement and 24-bit analog-to-digital conversion.

The Model 211 converts voltage or resistance to temperature units based on temperature response curve data for the sensor in use. Standard temperature response curves for silicon diodes and platinum RTDs are included in instrument firmware. It also provides non-volatile memory for one 200-point temperature response curve, which can be entered via the serial interface.

#### **Interface**

With an RS-232C serial interface and other interface features, the Model 211 is valuable as a stand-alone monitor and is easily integrated into other systems. Setup and every instrument function can be performed via serial interface or the front panel. Temperature data can be read up to seven times per second over computer interface; the display is updated twice each second. High and low alarms can be used in latching mode for error limit detection and in non-latching mode in conjunction with relays to perform simple on-off control functions. The analog output can be configured for either 0 to 10 V or 4 to 20 mA output.

#### Sensor Selection

#### Sensor temperature range (sensors sold separately)

		Model	Useful range	Magnetic field use
Diodes	Silicon diode	DT-670-SD	1.4 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-670E-BR	30 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-414	1.4 K to 375 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-421	1.4 K to 325 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-470-SD	1.4 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	Silicon diode	DT-471-SD	10 K to 500 K	$T \ge 60 \text{ K \& B} \le 3 \text{ T}$
	GaAlAs diode	TG-120-P	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-PL	1.4 K to 325 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
	GaAlAs diode	TG-120-SD	1.4 K to 500 K	$T > 4.2 \text{ K \& B} \le 5 \text{ T}$
Positive temperature	100 Ω platinum	PT-102/3	14 K to 873 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
coefficient RTDs	100 Ω platinum	PT-111	14 K to 673 K	$T > 40 \text{ K \& B} \le 2.5 \text{ T}$
	Rhodium-iron	RF-800-4	1.4 K to 500 K	$T > 77 \text{ K & B} \le 8 \text{ T}$
	Rhodium-iron	RF-100T/U	1.4 K to 325 K	$T > 77 \text{ K & B} \le 8 \text{ T}$
Negative	Cernox <sup>®</sup>	CX-1010	2 K to 325 K <sup>4</sup>	$T > 2 K \& B \le 19 T$
temperature	Cernox®	CX-1030-HT	3.5 K to 420 K <sup>2, 5</sup>	$T > 2 K \& B \le 19 T$
coefficient RTDs1	Cernox®	CX-1050-HT	4 K to 420 K <sup>2, 5</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1070-HT	15 K to 420 K <sup>2</sup>	$T > 2 K \& B \le 19 T$
	Cernox®	CX-1080-HT	50 K to 420 K <sup>2</sup>	$T > 2 K \& B \le 19 T$
	Germanium	GR-300-AA	1.2 K to 100 K <sup>3</sup>	Not recommended
	Germanium	GR-1400-AA	4 K to 100 K <sup>3</sup>	Not recommended
	Rox™	RX-102A	1.4 K to 40 K <sup>4</sup>	$T > 2 K \& B \le 10 T$


- Single excitation current may limit the low temperature range of NTC resistors
- 3 Low temperature limited by input resistance range
- $^4$  Low temperature specified with self-heating error:  $\leq 5$  mK
- Non-HT version maximum temperature: 325 K <sup>5</sup> Low temperature specified with self-heating error: ≤ 12 mK

**Silicon diodes** are the best choice for general cryogenic use from 1.4 K to above room temperature. Diodes are economical to use because they follow a standard curve and are interchangeable in many applications. They are not suitable for use in ionizing radiation or magnetic fields.

**Cernox®** thin-film RTDs offer high sensitivity and low magnetic field-induced errors over the 2 K to 420 K temperature range. Cernox sensors require calibration.

Platinum RTDs offer high uniform sensitivity from 30 K to over 800 K. With excellent reproducibility, they are useful as thermometry standards. They follow a standard curve above 70 K and are interchangeable in many applications.

#### Model 211 rear panel

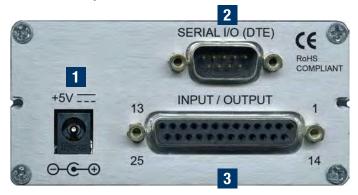

- 1 Power input connector
- 2 Serial (RS-232C) I/O (DTE)
- 3 Analog output



#### **Display**

The Model 211 has a 6-digit LED display with measurements available in temperature units K, °C, °F, or sensor units V or  $\Omega$ .

#### **Specifications**

#### **Sensor input configuration**

	Diode/RTD
Measurement type	4-lead differential
Excitation	8 constant current sources
Supported sensors	Diodes: silicon, GaAlAs
	RTDs: 100 $\Omega$ platinum, 1000 $\Omega$ platinum, germanium,
	carbon-glass, Cernox®, and Rox™
Standard curves	DT-470, DT-670, CTI-C, PT-100, and PT-1000
Input connector	Shared 25-pin D-sub

#### **Thermometry**

Model 211 Temperature Monitor

Number of inputs 1

Input configuration Input can be configured from the front panel to accept any of the supported input types

Isolation Measurement is not isolated from chassis ground

A/D resolution 24-bit

Input accuracy Sensor dependent—refer to Input Specifications table

Measurement resolution Sensor dependent—refer to Input Specifications table Maximum update rate 7 rdg/s

User curve One 200-point CalCurve™ or user curve in non-volatile memory

#### Front panel

Display 5-digit LED

Number of reading displays 1 **Display units** K,  $^{\circ}$ C,  $^{\circ}$ F, V, and  $\Omega$ 

Reading source Temperature and sensor units

Display update rate 2 rdg/s

Temp display resolution 0.001° from 0° to 99.999°, 0.01° from 100° to 999.99°, 0.1° above 1000°

Sensor units display resolution Sensor dependent to 5 digits

Display annunciators K, °C, °F, and V/ $\Omega$ 

Keypad 4 full travel keys, numeric and specific functions

Front panel features Display brightness control, keypad lock-out

#### Typical sensor performance—see Appendix F for sample calculations of typical sensor performance

	Example Lake Shore sensor	Temperature	Nominal resistance/ voltage	Typical sensor sensitivity <sup>6</sup>	Measurement resolution: temperature equivalents	Electronic accuracy: temperature equivalents	Temperature accuracy including electronic accuracy, CalCurve™, and calibrated sensor
Silicon diode	DT-670-SD	1.4 K	1.644 V	-12.49 mV/K	1.6 mK	±29 mK	±41 mK
	with 1.4H	77 K	1.028 V	-1.73 mV/K	11.6 mK	±175 mK	±197 mK
	calibration	300 K	0.5597 V	-2.3 mV/K	8.7 mK	±111 mK	±143 mK
		500 K	0.0907 V	-2.12 mV/K	9.4 mK	±99 mK	±149 mK
Silicon diode	DT-470-SD-13	1.4 K	1.6981 V	-13.1 mV/K	1.5 mK	±28 mK	±40 mK
	with 1.4H	77 K	1.0203 V	-1.92 mV/K	10.5 mK	±157 mK	±179 mK
	calibration	300 K	0.5189 V	-2.4 mV/K	8.4 mK	±105 mK	±137 mK
		475 K	0.0906 V	-2.22 mV/K	9.1 mK	±94 mK	±144 mK
GaAlAs diode	TG-120-SD	1.4 K	5.391 V	-97.5 mV/K	0.2 mK	±15 mK	±27 mK
	with 1.4H	77 K	1.422 V	-1.24 mV/K	16.2 mK	±512 mK	±534 mK
	calibration	300 K	0.8978 V	-2.85 mV/K	7 mK	±186 mK	±218 mK
		475 K	0.3778 V	-3.15 mV/K	6.4 mK	±135 mK	±185 mK
100 Ω platinum RTD	PT-103	30 K	3.66 Ω	0.19 Ω/K	10.5 mK	±320 mK	±330 mK
500 Ω full scale	with 1.4J	77 K	20.38 Ω	0.42 Ω/K	4.8 mK	±153 mK	±165 mK
	calibration	300 K	110.35 Ω	0.39 Ω/Κ	5.2 mK	±210 mK	±232 mK
		500 K	185.668 Ω	0.378 Ω/K	5.3 mK	±257 mK	±303 mK
Cernox <sup>®</sup>	CX-1050-SD-HT <sup>7</sup>	4.2 K	3507.2 Ω	-1120.8 Ω/K	45 μK	±2.0 mK	±7.0 mK
	with 4M	77 K	205.67 Ω	-2.4116 Ω/K	20.8 mK	±366 mK	±382 mK
	calibration	300 K	59.467 Ω	-0.1727 Ω/K	290 mK	±4.8 K	±4.8 K
		420 K	45.03 Ω	-0.0829 Ω/K	604 mK	±9.9 K	±9.9 K
Germanium	GR-300-AA	1.2 K	600 Ω	-987 Ω/K	51 μK	±0.6 mK	±5.3 mK
	with 0.3D	1.4 K	449 Ω	-581 Ω/K	86 µK	±1 mK	±5 mK
	calibration	4.2 K	94 Ω	-27 Ω/K	1.9 mK	±16 mK	±20 mK
		100 K	3 Ω	-0.024 Ω/K	2.10 K	±2.5 K	±2.5 K
Germanium	GR-1400-AA	4 K	1873 Ω	-1008 Ω/K	50 μK	±1.1 mK	±5.1 mK
	with 1.4D	4.2 K	1689 Ω	-862 Ω/K	58 μK	±1.2 mK	±5.2 mK
	calibration	10 K	253 Ω	-62 Ω/K	807 μK	±1.8 mK	±6.3 mK
		100 K	3 Ω	-0.021 Ω/K	2.40 K	±2.9 K	±2.9 K
Carbon-glass	CGR-1-2000	4.2 K	2260 Ω	-2060 Ω/K	25 μK	±0.6 mK	±4.6 mK
(no longer available)	with 4L	77 K	21.65 Ω	-0.157 Ω/K	319 mK	±410 mK	±435 mK
	calibration	300 K	11.99 Ω	-0.015 Ω/K	3.33 K	±4.2 K	±4.2 K

<sup>&</sup>lt;sup>6</sup> Typical sensor sensitivities were taken from representative calibrations for the sensor listed

Non-HT version maximum temperature: 325 K



#### **Input specifications**

Sensor type	Sensor temperature coefficient	Input range	Excitation current	Display resolution	Measurement resolution	Electronic accuracy	Instrument temperature coefficient
Silicon diode	negative	0 V to 2.5 V	10 μA ±0.05% <sup>8</sup>	100 μV	20 μV	$\pm 200~\mu V~\pm 0.01\%$ of rdg	$\pm 10~\mu V~\pm 5$ PPM of rdg/°C
GaAlAs diode	negative	0 V to 7.5 V	10 μA ±0.05% <sup>8</sup>	100 μV	20 μV	$\pm 350~\mu V~\pm 0.02\%$ of rdg	$\pm 20~\mu V \pm 5~PPM~of~rdg/^{\circ}C$
100 $\Omega$ platinum RTD, 250 $\Omega$ full scale	positive	0 $\Omega$ to 250 $\Omega$	1 mA ±0.3%9	$10 \text{ m}\Omega$	$2\ m\Omega$	$\pm 0.004~\Omega~\pm 0.02\%$ of rdg	$\pm 0.2~\text{m}\Omega~\pm 5$ PPM of rdg/°C
100 $\Omega$ platinum RTD, 500 $\Omega$ full scale	positive	0 $\Omega$ to 500 $\Omega$	1 mA ±0.3%9	10 mΩ	2 mΩ	$\pm 0.004~\Omega~\pm 0.02\%$ of rdg	$\pm 0.2 \text{ m}\Omega \pm 5 \text{ PPM of rdg/°C}$
1000 Ω platinum RTD	positive	0 $\Omega$ to 5000 $\Omega$	1 mA ±0.3%9	100 mΩ	20 mΩ	$\pm 0.06~\Omega~\pm 0.04\%$ of rdg	$\pm 2.0~\text{m}\Omega~\pm 5~\text{PPM}$ of rdg/°C
Cernox® RTD	negative	0 $\Omega$ to 7500 $\Omega$	10 μA ±0.05%	$100~\text{m}\Omega$	$50~\text{m}\Omega$	$\pm 0.8~\Omega~\pm 0.04\%$ of rdg	$\pm 20$ m $\Omega$ $\pm 15$ PPM of rdg/°C

8 Current source error has negligible effect on measurement accuracy

<sup>9</sup> Current source error is removed during calibration

#### Interface

Serial interface

RS-232C
Max baud rate
Connector
Reading rate

PS-232C
9600 baud
9-pin D-sub
Up to 7 rdg/s

**Alarms** 

**Number** 2, high and low **Data source** Temperature

Settings High setpoint, Low setpoint, Dead band, Latching or Non-latching

**Actuators** Display message, relays

Relays

Number 2

Contacts Normally Open (NO), Normally Closed (NC), and Common (C)

Contact rating 30 VDC at 1 A

**Operation** Activate relays on high or low input alarm or manual

**Connector** Shared 25-pin D-sub

Analog output

Isolation Output is not isolated from chassis ground

**Update rate** 7 readings per s **Data source** Temperature

	Voltage	Current
Range	0 V to 10 V	4 mA to 20 mA
Accuracy	±1.25 mV	±5.0 μA
Resolution	0.3 mV	0.6 μΑ
Min load resistance	500 Ω	NA
Compliance voltage	NA	10 V
Load regulation	NA	$\pm 0.02\%$ of reading 0 to 500 $\Omega$

	Temperature	Sensor units (fixed by type)
Scales:	0 K to 20 K	Diodes: $1 V = 1 V$
	0 K to 100 K	100 $\Omega$ platinum: 1 V = 100 $\Omega$
	0 K to 200 K	1000 Ω platinum: 1 V = 1000 Ω
	0 K to 325 K	NTC resistor: 1 V = 1000 $\Omega$
	0 K to 475 K	
	0 K to 1000 K	

**Settings** Voltage or current, scale **Connector** Shared 25-pin D-sub

#### General

Ambient temperature 15 °C to 35 °C at rated accuracy, 10 °C to 40 °C at reduced accuracy

Power requirements Regulated +5 VDC at 400 mA

Size  $96 \text{ mm W} \times 48 \text{ mm H} \times 166 \text{ mm D}$  ( $3.8 \text{ in} \times 1.9 \text{ in} \times 6.5 \text{ in}$ )

Mounting Panel mount into  $91 \text{ mm W} \times 44 \text{ mm H}$  ( $3.6 \text{ in} \times 1.7 \text{ in}$ ) cutout

Weight 0.45 kg (1 lb)
Approvals CE mark, RoHS



**2111** Single 1/4 DIN panel-mount adapter,  $105 \text{ mm W} \times 132 \text{ mm H} (4.1 \text{ in} \times 5.2 \text{ in})$ 



**2112** Dual 1/4 DIN panel-mount adapter, 105 mm W  $\times$  132 mm H (4.1 in  $\times$  5.2 in)



Comes standard with interchangeable input plugs

**Power requirements** 100 to 240 VAC, 50 or 60 Hz, 0.3 A max

Output +5 V at 1.2 A

Size 40.5 mm W × 30 mm H × 64 mm D (1.6 in × 1.2 in × 2.5 in)

**Weight** 0.15 kg (0.33 lb)



#### Part number Description

211S Model 211 single channel temperature monitor—includes 100

to 240 V, 6 W universal power supply with interchangeable input plugs (109-132), one DB-25 sensor input mating connector (G-106-253), one sensor input mating connector shell (G-106-264), a calibration certificate and a user's manual

shell (G-106-264), a calibration certificate and a user's manuments.

Model 211S with all accessories except the power supply

#### **Accessories**

**109-132** 100-240 VAC power supply with interchangeable plugs for US,

UK, Europe, Australia, and China application
Single 1/4 DIN panel-mount adapter
Dual 1/4 DIN panel-mount adapter

**8000** CalCurve<sup>™</sup>, CD-ROM (included with calibrated sensor)

**G-106-253** DB-25 plug, qty 1 DB-25 hood; qty 1

CAL-211-CERT Instrument recalibration with certificate Instrument recalibration with certificate and data 119-043 Model 211 temperature monitor manual






## 230 Series Temperature Transmitters/Monitors

230 Series Temperature Transmitters/Monitors

The Model 231, Model 231P, and Model 234 each support a different sensor type

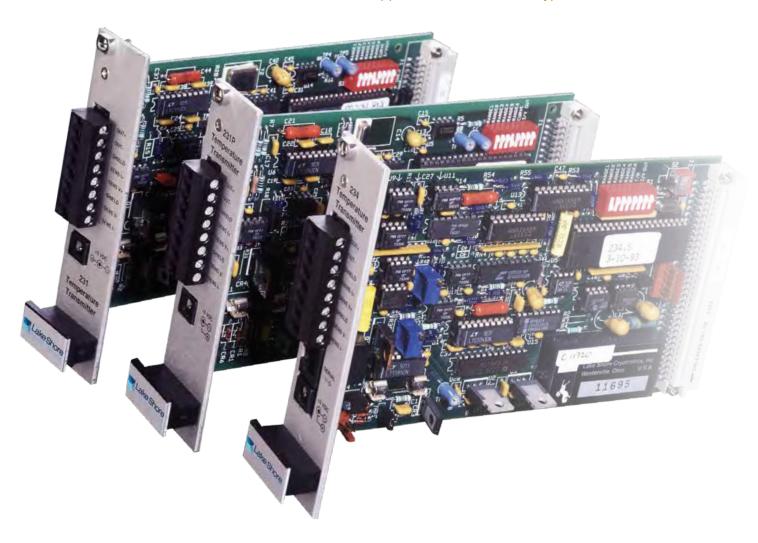

#### 230 Series features

- Sensor input fully isolated from power supply potential
- Different models support various sensor types
- 4-lead differential measurement
- Output range of 4 mA to 20 mA or 0 mA to 20 mA (0 V to 10 V)
- Available rack-mount case holds up to 12 units
- Full 3 year standard warranty



#### **Model 231 features**

Operates from 1.4 K to 500 K with appropriate diode

#### **Model 231P features**

Operates from 1.4 K to 800 K with appropriate PTC RTD

#### **Model 234 features**

- Operates from 100 mK to 420 K with appropriate NTC RTD
- Includes serial interface



#### Model 231

The Model 231 operates with either silicon diode or gallium-aluminum-arsenide (GaAlAs) diode sensors.

Excited with a 10 µA current source from the Model 231, the sensors produce a voltage that depends on temperature. A microcontroller reads the voltage through an A/D converter and translates it into temperature using a temperature response curve. The Model 231 includes two standard curves for DT-470 and DT-670 diode sensors. It also supports a single CalCurve™ option for calibrated sensors (TG-120 diodes require a CalCurve™).

#### Model 231P

The Model 231P uses a PT-100 Series platinum sensor. The Model 231P excites the sensor with a 500  $\mu$ A current to produce a measurable signal. Either the standard platinum curve (IEC 751) or a CalCurve<sup>TM</sup> is used for temperature conversion.

#### Model 234

The Model 234 operates with Cernox®, carbon-glass, germanium, or other negative temperature coefficient (NTC) resistance temperature sensors. The Model 234 excites the sensor with a constant voltage of 10 mV or less to minimize the effects of sensor self-heating at low temperatures.

The Model 234 employs an analog control circuit to maintain a constant voltage signal across the sensor. A series of reference resistors convert the resulting sensor current to a voltage. A microcontroller reads the voltage with an A/D converter, calculates sensor resistance, and converts the resistance to temperature by table interpolation (requires a CalCurve™ for temperature conversion). The sensor excitation voltage is reversed each reading to compensate for thermal voltages and offsets.

Once one of the 230 Series obtains temperature data, it transmits a current of 4 mA to 20 mA. The current output changes linearly with sensor temperature. Output scale depends on the selected temperature range. Several switch-selected ranges are available. Highest accuracy and sensitivity are achieved when the output is set for a narrow temperature band. A 0 mA to 20 mA output is also available to convert output to a voltage scaled from zero. A 500  $\Omega,\,\pm0.02\%$  output load resistor produces the maximum full-scale output of 10 V.

Circuitry for the Model 230 Series is powered by a single +5 VDC supply applied either from the front panel connector or the power pins on the VME bus connector. The outputs are isolated so several transmitters can be run off the same supply without interference.

Mechanical mounting is easy because the 230 Series is built on a standard size VME card. It fits directly into a single height (3U) VME card holder. The transmitter does not use the electrical bus format, only its physical shape and power supply.

The Model 234 includes a RJ11 serial interface.

Model 234 measurement scales, excitation, resolution, and accuracy

	Sensor resistance	Sensor excitation voltage	Resolution	Accuracy ±(% rdg + Ω)
0	1 $\Omega$ to 6 $\Omega$	5 mV	$0.0003~\Omega$	0.5 + 0.0006
1	4.5 $\Omega$ to 12.5 $\Omega$	5 mV	0.0001 Ω	0.1 + 0.0013
2	9 $\Omega$ to 60 $\Omega$	10 mV	0.001 Ω	0.1 + 0.006
3	45 Ω to 125 Ω	5 mV	0.001 Ω	0.1 + 0.013
4	90 $\Omega$ to 360 $\Omega$	10 mV	0.003 Ω	0.1 + 0.036
5	290 k $\Omega$ to 1.25 k $\Omega$	10 mV	0.01 Ω	0.1 + 0.13
6	900 k $\Omega$ to 3.6 k $\Omega$	10 mV	0.03 Ω	0.1 + 0.36
7	$2.9~\text{k}\Omega$ to $12.5~\text{k}\Omega$	10 mV	0.1 Ω	0.1 + 1.3
8	9 kΩ to 36 kΩ	10 mV	0.3 Ω	0.1 + 3.6
9	$29~\text{k}\Omega$ to $300~\text{k}\Omega$	10 mV	6.8 Ω	0.5 + 30


Thermometry	231	231P	234
lumber of inputs	1	1	1
Measurement type	4-lead differential	4-lead differential	4-lead differential
Sensor type	Silicon diode, GaAlAs diode	Platinum	Carbon-glass, germanium, Cernox®
Sensor temperature coefficient	Negative	Positive	Negative
Sensor units	Volts (V)	Ohms (Ω)	Ohms (Ω)
nput range	0 V to 5 V	0 Ω to 312 Ω	1 Ω to 300 kΩ
Sensor excitation	10 μA ±0.1% DC current	500 μA ±0.02% DC current	Constant voltage pinned at 5 mV or 10 mV dependent on resistance range
Jpdate rate	5 rdg/s	5 rdg/s	4 rdg/s (2 rdg/s on Scale 0 only)
recision curve storage	One curve loaded at Lake Shore	One curve loaded at Lake Shore	One curve, loaded at Lake Shore or in the field via serial interface
Example Lake Shore sensor	DT-470-C0	PT-103	CGR-1-1000 with 1.4L calibration
Sensor temperature range	1.4 K to 475 K	30 K to 800 K	1.4 K to 325 K
tandard curve	Lake Shore Curve 10	IEC 751	Requires calibrated sensor and CalCurve™
Typical sensor sensitivity <sup>1</sup>	-30 mV/K at 4.2 K	0.19 Ω/K at 30 K	-700 Ω/K at 4.2 K
,	40 mK at 77 K	11 mK at 77 K	±0.12 mK at 30 K
	32 mK at 300 K	13 mK at 300 K	±6.6 mK at 77 K
		14 mK up to 800 K	±67 mK at 300 K
Measurement resolution <sup>1</sup>			
Sensor units	76.3 μV	4.8 mΩ	Range dependent
emperature equivalence	2.5 mK at 4.2 K	22 mK at 30 K	±0.04 mK at 4.2 K
	40 mK at 77 K	11 mK at 77 K	±0.12 mK at 30 K
	32 mK at 300 K	13 mK at 300 K 14 mK up to 800 K	±6.6 mK at 77 K ±67 mK at 300 K
Electronic measurement accuracy <sup>1</sup>		14 IIIK up to 600 K	±07 IIIK at 300 K
ensor units	$\pm 75 \mu\text{V} \pm 0.01\%$ of reading	$\pm 0.05 \Omega \pm 0.05\%$ of reading	Range dependent (see table on 128)
	±7.5 μν ±0.01% of reading ±0.07 K at 4.2 K		±2 mK at 4.2 K
emperature accuracy	±0.07 K at 4.2 K ±0.16 K at 77 K	±0.2 K at 30 K ±0.15 K at 77 K	±2 mK at 4.2 K ±8 mK at 10 K
	±0.10 K at 77 K ±0.12 K at 300 K	±0.13 K at 77 K ±0.3 K at 300 K	±18 mK at 77 K
	_0.12 N dt 000 N	±0.7 K up to 800 K	±1.2 K at 300 K
Measurement temperature coefficie	nt		
Sensor units (% of reading/°C umbient)	0.0006% of resistance rdg /°C	0.002% of resistance rdg /°C	0.0125% of resistance rdg /°C
Temperature equivalence	3 mK/°C at 4.2 K	0.4 mK/°C at 30 K	±0.18 mK/°C at 4.2 K
	3 mK/°C at 77 K	1 mK/°C at 77 K	±0.8 mK/°C at 10 K
	1.2 mK/°C at 300 K	6 mK/°C at 300 K	±18 mK/°C at 77 K
		18 mK/°C at 800 K	±100 mK/°C at 300 K
Magnetic field use			
Silicon diode	Recommended for $T \ge 60 \text{ K}$ and $B \le 3 \text{ T}$	NA	NA
GaAlAs diode	Recommended for $T \ge 4.2 \text{ K}$ and $B \le 5 \text{ T}$	NA	NA
Platinum	NA	Recommended for $T > 40 \text{ K}$ and $B \le 2.5 \text{ T}$	NA
Carbon-glass	NA	NA	Recommended for $T \ge 2 \text{ K}$ and $B \le 19 \text{ T}$
	231	231P	234
General			
Ambient temperature range	15 °C to 35 °C	15 °C to 35 °C	15 °C to 35 °C
Power requirements	±5 (±0.25) VDC, 500 mA, 2.5 W	$\pm 5 \ (\pm 0.25) \ VDC, 500 \ mA, 2.5 \ W$	234: ±5 (±0.25) VDC, 500 mA, 2.5 W
Enclosure type	See diagrams	See diagrams	See diagrams
Mounting	VME end panel and back plane:	VME end panel and back plane:	VME end panel and back plane:
	transmitters do not use	transmitters do not use	transmitters do not use
	electrical bus format, only its	electrical bus format, only its	electrical bus format, only its
Nime.	physical shape and power	physical shape and power	physical shape and power
Size	100 mm H × 160 mm D × 30.5 mm W	100 mm H × 160 mm D × 30.5 mm W	234: 100 mm H × 160 mm D × 30.5 mm
Weight	0.25 kg (0.5 lb)	0.25 kg (0.5 lb)	0.25 kg (0.5 lb)


<sup>&</sup>lt;sup>1</sup> See Appendix F for sample calculations of typical sensor performance



Output	231	231P	234	
Number of outputs	1	1	1	
Output type			·	
1 21	Current source, isolated from power source—output or sensor can be grounded, but not both (all models) 4 mA to 20 mA or 0 mA to 20 mA (for 0 V to 10 V with provided 500 Ω 0.02%, 25 ppm resistor) (all models)			
Output range				
Output compliance	10 V (500 Ω max load)	10 V (500 Ω max load)	10 V (500 Ω max load)	
Output temperature ranges				
Range 1	0 K to 20 K	0 K to 20 K	0 K to 10 K	
Range 2	0 K to 100 K	0 K to 100 K	0 K to 20 K	
Range 3	0 K to 200 K	0 K to 200 K	0 K to 100 K	
Range 4 Range 5	0 K to 325 K 0 K to 475 K	0 K to 325 K 0 K to 475 K	0 K to 200 K 0 K to 300 K	
Range 6	0 K to 1000 K	0 K to 1000 K	75 K to 325 K	
4 mA to 20 mA output	0 K to 1000 K	0 K to 1000 K	70 K to 020 K	
Output resolution				
Current	1.22 µA (0.006% of full scale)	1.22 µA (0.006% of full scale)	1.22 µA (0.006% of full scale)	
Temperature equivalence				
Range 1	1.5 mK	Not used	0.8 mK	
Range 2	7.6 mK	7.6 mK	1.5 mK	
Range 3	15.3 mK	15.3 mK	7.6 mK	
Range 4	24.8 mK	24.8 mK	15.3 mK	
Range 5	36.2 mK	36.2 mK	22.9 mK	
Range 6	76.3 mK	76.3 mK	19.1 mK	
Output accuracy				
Current	±2 μA (±0.01% of full scale)	$\pm 2 \mu\text{A} (\pm 0.01\% \text{of full scale})$	$\pm 5 \mu\text{A} (\pm 0.025\% \text{of full scale})$	
Temperature equivalence				
Range 1	2.5 mK	Not used	3.1 mK	
Range 2	12.5 mK	12.5 mK	6.2 mK	
Range 3	25 mK	25 mK	31.2 mK	
Range 4	41 mK 59 mK	41 mK 59 mK	62.5 mK 93.7 mK	
Range 5 Range 6	125 mK	125 mK	78.1 mK	
Output temperature coefficient	120 11110	120 11110	70.1 11110	
Current (%/°C ambient)	±0.0055% of output current per °C	±0.0055% of output current per °C	±2 μΑ/°C (±0.01%/°C)	
Temperature equivalence	_0.0000 % of output outfolk por 0	_0.000076 01 04444 04110111 por 0	μι σ (_σ.σ., νω σ)	
Range 1	1 mK/°C	Not used	±1 mK/°C	
Range 2	6 mK/°C	6 mK/°C	±2 mK/°C	
Range 3	12 mK/°C	12 mK/°C	±10 mK/°C	
Range 4	18 mK/°C	18 mK/°C	±20 mK/°C	
Range 5	26 mK/°C	26 mK/°C	±30 mK/°C	
Range 6	55 mK/°C	55 mK/°C	±25 mK/°C	
0 mA to 20 mA output (0 V to 10 V w Output resolution	with 500 $\Omega$ , 0.02% load resistor)			
Voltage	0.6 mV	0.6 mV	0.61 mV	
Temperature equivalence				
Range 1	1.2 mK	Not used	0.6 mK	
Range 2	6.1 mK	6.1 mK	1.2 mK	
Range 3	12.2 mK 19.8 mK	12.2 mK 19.8 mK	6.1 mK 12.2 mK	
Range 5	29 mK	29 mK	18.3 mK	
Range 6	61 mK	61 mK	15.2 mK	
Output accuracy				
Voltage	3 mV (0.03% of full scale)	3 mV (0.03% of full scale)	$\pm 4.5$ mV ( $\pm 0.025\%$ of full scale $\pm 0.02\%$ resistor accuracy)	
Temperature equivalence	G mal/	Notd	4.5 mg/	
Range 1 Range 2	6 mK 30 mK	Not used 30 mK	4.5 mK 9.0 mK	
Range 3	60 mK	60 mK	45.0 mK	
Range 4	98 mK	98 mK	90.0 mK	
Range 5	143 mK	143 mK	135.0 mK	
Range 6	300 mK	300 mK	112.5 mK	
Output temperature				
coefficient Voltage (% output/°C ambient)	±0.008%/°C	±0.008%/°C	±1.25 mV/° C (±0.01%/°C ±0.0025%/°C of load resistor)	
Temperature equivalence			ioda resistot)	
Range 1	2 mK/°C	Not used	±1.2 mK/°C	
Range 2	8 mK/°C	8 mK/°C	±2.5 mK/°C	
Range 3	16 mK/°C	16 mK/°C	±12 mK/°C	
Range 4	26 mK/°C	26 mK/°C	±25 mK/°C	
Range 5	38 mK/°C	38 mK/°C	±36 mK/°C	
Range 6	80 mK/°C	80 mK/°C	±30 mK/°C	



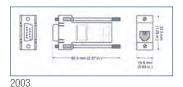

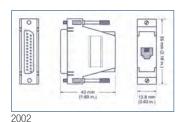



#### Multiple card enclosure (2308-12)

The 2308-12 VME card case holds up to 12 temperature transmitters. A +5 VDC power supply with universal input is provided with the case. Wall mount power supplies are not necessary with a 2308-12.

Card slots 12

Output voltage +5 VDC, 100 mV peak to peak ripple

Output current 6 A (max)

Input power Universal 85 to 265 VAC, 47 to 440 Hz, 60 W

**Ambient temp range** 15 °C to 35 °C (59 °F to 95 °F) **Enclosure mounting** Bench or full (19 in) rack

**Size**  $450 \text{ mm W} \times 178 \text{ mm H} \times 260 \text{ mm D} (17.7 \text{ in} \times 7 \text{ in} \times 10.25 \text{ in})$ 

Weight 5.5 kg (12 lb)

#### **Power connections**

Lake Shore temperature transmitters are powered by a +5 VDC supply if the transmitter card is ordered without a rack or plug-in supply. The voltage must be regulated to within  $\pm 0.25$  VDC. Each transmitter draws up to 500 mA from the supply.

## Ordering information

#### Part number Description

231 Transmitter card for use with silicon diode—includes sensor/

output mating connector (106-739), 500 ohm, 0.02% PPM output resistor (103-626), calibration certificate and user's

manual

231-115 231 transmitter with a 115 VAC (50/60 Hz) wall plug-in power

supply

231-230 231 transmitter with 230 VAC wall plug-in power supply
231P Transmitter for use with platinum resistor—includes sensor/

output mating connector (106-739), 500 ohm, 0.02% PPM output resistor (103-626), calibration certificate and user's

manual

231P-115 231P transmitter with 115 VAC (50/60 Hz) wall plug-in power

supply

231P-230 231P transmitter with 230 VAC wall plug-in power supply 234 Transmitter for use with carbon-glass, germanium, and

Cernox®—includes sensor/output mating connector (106-739), 500 ohm, 0.02% PPM output resistor (103-626),

calibration certificate and user's manual 234 transmitter with 115 VAC (50/60 Hz) wall plug-in power

gunnly

234-230 234 transmitter with 230 VAC wall plug-in power supply

#### **Accessories**

234-115

2001 RJ11 4 m (14 ft) modular serial cable
2002 RJ11 to DB25 adapter—connects RJ11 cable
to a 25-pin RS-232C serial port on rear of computer
2003 RJ11 to DB9 adapter—connects RJ11 cable to a 9-pin

RS-232C serial port on rear of computer

**2308-12** VME rack and power supply (holds up to 12 transmitters)

2308-BP VME rack blank panel

 8002-231
 231 CalCurve™ data, field installed

 8002-231P
 231P CalCurve™ data, field installed

 8002-234
 234/234D CalCurve™ data, field installed

 CAL-231-CERT
 Instrument recalibration with certificate

 CAL-231P-CERT
 Instrument recalibration with certificate

 CAL-234-CERT
 Instrument recalibration with certificate





# Model 121 Programmable Current Source



#### **Model 121 features**

- 6 decades of output current, selectable in 13 ranges
- Programmable current output, 100 nA to 100 mA
- Low-noise output
- Large 3 digit LED display
- Simple user interface
- Current reversal feature
- USB interface enables integration with automated test systems
- DIN panel mountable package
- Detachable output terminal block
- CE certification
- Full 3 year standard warranty



An accurate, stable source of current is key

devices, where the voltage drop across the device can be dependent upon temperature,

magnetic field, and other parameters. The

instrument's wide output range is of great

whose resistance can vary with temperature

current reversal feature enables compensation

Basic device QC ("good/bad" verification)

LED brightness testing (constant device

by as much as 6 orders of magnitude. The

for thermal EMF, important for accurately

measuring resistors at very low excitation

levels. Example applications include:

value when used with RTD-type sensors

to ensuring consistent operation of these



#### Introduction

The Model 121 programmable current source is a precision instrument suitable for bench-top use or panel-mounted operation in labs, test facilities, and manufacturing environments. It provides a low noise, highly-stable source of current up to 100 mA, with convenient manual selection through 13 pre-set output levels, each representing a ten-fold change in power when attached to a resistive load. A "user" setting allows the current output to be defined anywhere within the operating range of the unit, from 100 nA to 100 mA.

Fully automated operation is also possible via the instrument's USB computer interface, through which the Model 121 can be commanded to output any desired current at any time. Thus, application-specific test patterns can be created.

The instrument operates at 5 VDC, and power is supplied by the external AC wall-mount supply provided with the standard Model 121. The supply will automatically conform to any AC line voltage ranging from 100 VAC to 240 VAC, 50 to 60 Hz.

#### **Applications**

The Model 121 current source is ideally suited for testing, measuring, and operating resistive and semiconductor devices, such as:

- Lake Shore Cernox® temperature sensors
- Other resistance temperature detectors (RTDs) such as platinum sensors
- Diode temperature sensors, including Lake Shore DT-670s
- LED devices
- Hall sensors used for magnetic field measurement

rs, including

Temperature sensor calibration (determine

current)


Temperature measurement (using a voltmeter readout)

resistance at fixed calibration points)

- Magnetic sensor calibration and measurement
- Semiconductor device measurements (IV curves for diodes, transistors, etc.)
- Circuit prototyping (fixed current source)
- Small scale electro-chemical applications

Whether operating over a wide range of environmental conditions, establishing precise sensor calibrations or simply testing devices for conformance, the Model 121 provides a convenient and reliable alternative to simple voltage-based circuits, and a very affordable alternative to more expensive multi-function current sources. It can be readily integrated into automated test systems using its built-in USB computer interface and offers a highly readable, simple-to-use operator display.


#### Model 121 rear panel

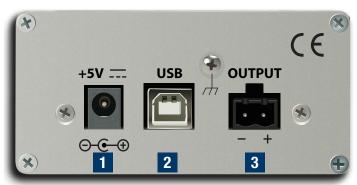


- 1 Power input connector
- 2 USB interface
- 3 Output current




#### **Model 121 possible applications**

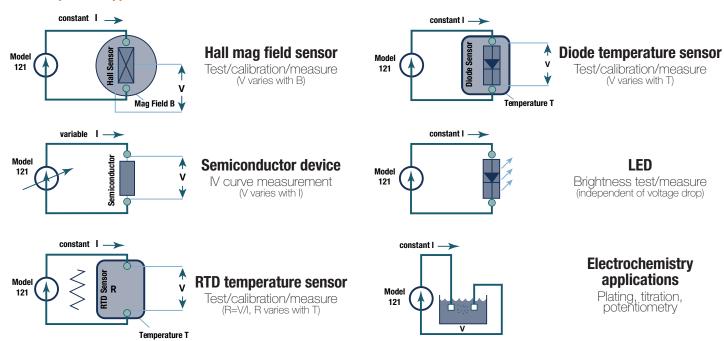

#### **Application range of Model 121**

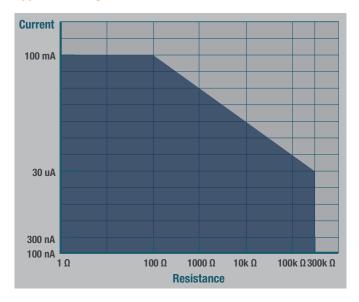

#### Using the Model 121 with a resistive device/sensor

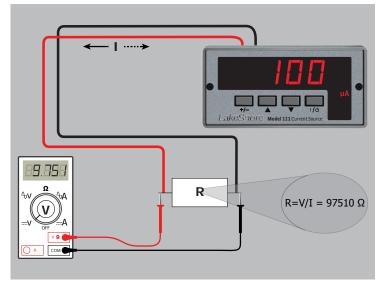



#### **Output**

Type: Bipolar, DC current source

Current ranges: 13 fixed ranges of 100 nA, 300 nA, 1  $\mu$ A, 3  $\mu$ A, 10  $\mu$ A, 30  $\mu$ A, 100  $\mu$ A, 300  $\mu$ A, 1 mA, 3 mA, 10 mA, 30 mA, 100 mA, plus a user programmable range

Accuracy: 0.05% on 10  $\mu A$  range, 0.5% on 100 nA and 300 nA ranges, 0.1% on all other ranges

Compliance voltage: ±11 V up to 30 mA, 10 V up to 100 mA

AC current ripple: Less than 0.1% on 100 nA and 300 nA ranges, Less than 0.01% on all other ranges in properly shielded system

Current ripple frequency: Dominated by the switching power supply line frequency/ harmonics

Temperature coefficient: 0.03% of range/°C for the 100 nA range, 0.01% of range/°C for all other ranges

**Line regulation:** Less than 0.01% change in output for 5% change in the DC input voltage **Load regulation:** Less than 0.01% change in output current over the full range scale **Stability (24 h):**  $\pm 0.05\%$  on 100 nA range,  $\pm 0.01\%$  per day on all other fixed ranges

**Connections:** Detachable terminal block

Maximum load:  $300~k\Omega$ Maximum lead length: 50~ft

#### **User setting**

#### **Programming**

Operation: Output current settable via computer interface

Resolution: 3 significant digits

Accuracy:  $\pm 0.5\%$  of 100 nA and 300 nA ranges,

±0.25% of all other ranges

Maximum current: 100 mA Minimum current: 100 nA

Front panel

Display: LED – 3 digits plus sign Display units: mA, μA, and nA
Display update rate: 2 rdg/s

Display anunciators: mA,  $\mu$ A, nA, and compliance

Keypad: 4 full travel keys

Keypad functions: Range Up, Range Down, Current Polarity, Output Inhibit

Interface

USB

Function: Emulates a RS-232 serial port

Baud rate: 57,600

Connector: B-type USB connector

Reading rate: To 10 rdg/s

Software support: LabVIEW<sup>™</sup> driver (see www.lakeshore.com)

#### General

Ambient temperature: 15 °C to 35 °C at rated accuracy; 5 °C to 40 °C at reduced accuracy

Power requirement: +5 VDC  $\pm5\%$  at 400 mA, barrel plug 5.5 mm 0D  $\times$  2.1 mm ID  $\times$  9.9 mm L, center pin positive

Size:  $96 \text{ mm W} \times 48 \text{ mm H} \times 166 \text{ mm D} (3.8 \text{ in} \times 1.9 \text{ in} \times 6.5 \text{ in})$ Mounting: Panel mount into  $91 \text{ mm W} \times 44 \text{ mm H} (3.6 \text{ in} \times 1.7 \text{ in})$  cutout

Weight: 0.45 kg (1 lb) Approval: CE mark

# Power supply (109-132) Comes standard with interchangeable input plugs Power requirements 100 to 240 VAC, 50 or 60 Hz, 0.3 A max Output +5 V at 1.2 A Size 40.5 mm W × 30 mm H >. mm D (1.6 in × 1.2 in × 2.5 ) Weight 0.15 kg (0.33 lb)

#### Ordering information

#### Part number Description

121 Programmable current source—includes one 100 V to 240 V, 10 W power supply with universal input interchangeable input

10 W power supply with universal input interchangeable input plugs (109-132), calibration certificate, and user manual Programmable current source—no power supply. Includes

calibration certificate and user manual

#### **Accessories**

121N

CAL-102-CERT Model 121 recalibration with certificate





# Model 625 Superconducting Magnet Power Supply



- Linear regulation minimizes noise
- Ripple  $\leq$  0.007% of maximum current (into a 1 mΩ load)
- 1 mA per hour stability
- Parallel operation to ±120 A
- Full 3 year standard warranty





#### Introduction

The Model 625 superconducting magnet power supply is the ideal supply for small to medium sized superconducting magnets used in high sensitivity materials research applications. The Model 625 is a practical alternative to both the larger, one size fits all, superconducting magnet supplies and the endless adaptations of generic power supplies. By limiting output power, Lake Shore was able to concentrate on the performance requirements of the most demanding magnet users. The resulting Model 625 provides high precision, low noise, safety, and convenience.

Precision in magnetic measurements is typically defined as smooth continuous operation with high setting resolution and low drift. Achieving these goals while driving a challenging load, such as a superconducting magnet, requires a unique solution. The Model 625 delivers up to 60 A at a nominal compliance voltage of 5 V, with the supply acting as either a source or a sink in true 4-quadrant operation. Its current source output architecture with analog control enables both smooth operation and low drift. A careful blending of analog and digital circuits provides high setting resolution of 0.1 mA and flexible output programming.

Lake Shore chose linear input and output power stages for the moderate 300 W output of the Model 625. Linear operation eliminates the radiated radio frequency (RF) noise associated with switching power supplies, allowing the Model 625 to reduce the overall noise in its output and the noise radiated into surrounding electronics.

Safety should never be an afterthought when combining stored energy and liquid cryogens in a superconducting magnet system. The Model 625 incorporates a variety of hardware and firmware protection features to ensure the safety of the magnet and supply.

Instrument users have come to rely on Lake Shore for convenience and ease of use. The Model 625 includes the features necessary to conveniently manage a superconducting magnet, such as a persistent switch heater output, calculated field reading, current ramping, and quench detection. Computer interfaces are also integrated for automation of the magnet system. The Model 625 is truly an excellent one-box solution for controlling a superconducting magnet.

Model 625 Superconducting Magnet Power Supply

#### **Output architecture**

True 4-quadrant output capability of the Model 625 is ideal for the charge and discharge cycling of superconducting magnets for both positive and negative fields. Tightly integrated analog control of the 4-quadrant output provides smooth current change with very low overshoot on output change. The Model 625 has the ability to charge and discharge magnets up to a 5 V rate.

True 4-quadrant operation eliminates the need for external switching or operator intervention to reverse the current polarity, significantly simplifying system design. The transition through zero current is smooth and continuous, allowing the user to readily control the magnetic field as polarity changes. At static fields, output current drift is also kept low by careful attention in the analog control circuits and layout. The high stability and low noise of the Model 625 make it possible in many situations to run experiments without going into persistent mode. This can help to reduce the time necessary to gather data.

The Model 625 output architecture relies on low noise, linear input and output stages. The linear circuitry of the Model 625 permits operation with less electrical noise than switch-mode superconducting magnet power supplies. One key benefit of this architecture is CE compliance to the electromagnetic compatibility (EMC) directive, including the radiated emissions requirement.

#### **Output programming**

The Model 625 output current is programmed internally via the keypad or the computer interface, externally by the analog programming input, or by the sum of the external and internal settings. For the more popular internal programming, the Model 625 incorporates a proprietary digital-to-analog converter (DAC) that is monotonic over the entire output range and provides a resolution of 0.1 mA.

#### Model 625 rear panel

- Positive and negative outputs
- Analog I/O
- Line input assembly
- IEEE-488 interface
- Serial (RS-232C) I/O (DTE)
- PSH output
- Digital I/O





The Model 625 generates extremely smooth and continuous ramps with virtually no overshoot. The digitally generated constant current ramp rate is variable between 0.1 mA/s and 99.999 A/s. To assure a smooth ramp rate, the power supply updates the high-resolution DAC 27 times per second. A low-pass filter on the output DAC smooths the transitions at step changes during ramping. Ramping can also be initiated by the trigger input.

The output compliance voltage of the Model 625 is settable to a value between 0.1 V and 5 V, with a 100  $\mu$ V resolution. The voltage setting is an absolute setting, so a 2 V setting will limit the output to greater than -2.0 V and less than +2.0 V.

#### **Output readings**

The Model 625 provides high-resolution output readings. The output current reading reflects the actual current in the magnet, and has a resolution of 0.1 mA. The output voltage reading reports the voltage at the output terminals with a resolution of 100  $\mu V$ . A remote voltage reading is also available to more accurately represent the magnet voltage by bypassing voltage drops in the leads connecting the power supply to the magnet. All output readings can be prominently displayed on the front panel and read over the computer interface.

#### **Protection**

Managing the stored energy in superconducting magnets necessitates several different types of protection. The Model 625 continuously monitors the load, line voltage, and internal circuits for signs of trouble. Any change outside of the expected operating limits triggers the supply to bring the output to zero in a fail-safe mode. When line power is lost, the output crowbar (SCR) will activate and maintain control of the magnet, discharging at a rate of 1 V until it reaches zero.

Quench detection is necessary to alert the user and to protect the magnet system. The Model 625 uses a basic and reliable method for detecting a quench. If the current changes at a rate greater than the current step limit set by the operator, a quench is detected and the output current is safely set to zero.

The remote inhibit input allows an external device to immediately set the output current to zero in case of a failure. This input is normally tied to an external quench detection circuit, the fault output of a second power supply, or an emergency shutdown button. The fault output is a relay contact that closes when a fault condition occurs. The contact closure alerts other system components of the fault.

#### **Parallel operation**

If an application requires more output current than a single Model 625 can provide, two supplies can be connected in parallel for 120 A/5 V operation. Each unit is programmed for half of the total output current, operates independently, and retains 0.1 mA resolution at 60 A operation. When the units are properly configured, either unit can detect a fault, protect itself, and issue a fault output signaling the other unit to automatically enter the proper protection mode.

#### Persistent switch heater output


The integrated persistent switch heater output is a controlled DC current source capable of driving most switch heaters. It sources from 10 mA to 125 mA with a setting resolution of 1 mA and selectable compliance voltage of 12 V or 21 V. The minimum load that the persistent switch heater can drive is 10  $\Omega$ . Persistent mode operation is integrated into the instrument firmware to prevent misoperation of the magnet.

#### **Interfaces**

The Model 625 includes IEEE-488 and RS-232C computer interfaces that provide access to operating data, stored parameters, and remote control of all front panel operating functions. In addition, the Model 625 includes a trigger function that is used to start an output current ramp. When the trigger is activated, either by an external trigger or by computer interface command, the power supply will begin ramping to the new setpoint.

The Model 625 provides two analog outputs to monitor the output current and voltage. Each output is a buffered, differential, analog voltage representation of the signal being monitored. The current monitor has a sensitivity of 1 V = 10 A, while the voltage monitor has a sensitivity of 1 V = 1 V.

#### **Current change using internal programming**

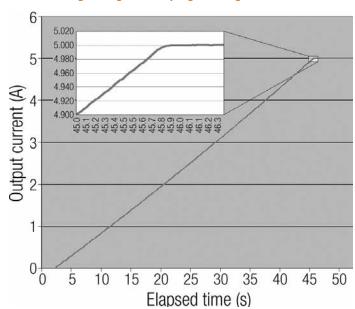

This plot illustrates an actual 5 A current change into an  $8.6\,H$  superconducting magnet. A smooth, 95 mA/s ramp is shown with minimal overshoot highlighted in the detail area. (Output current monitor measured at  $58.88\,Hz$  rate with a HP 34401-data multiplied by  $10\times$  to obtain output current results.)



#### Display and keypad

The Model 625 incorporates a large 8-line by 40-character vacuum fluorescent display. Output current, calculated field in tesla or gauss, output voltage, and remote voltage sense readings can be displayed simultaneously. Five LEDs on the front panel provide quick verification of instrument status, including ramping, compliance, fault, PSH status, and computer interface mode. Error conditions are indicated on the main display along with an audible beeper. Extended error descriptions are available under the Status key.



Current and voltage settings, current and voltage readings, ramp rate, voltage sense, and persistent switch heater status and instrument status displayed simultaneously



The instrument can be set up to show calculated field along with output field setting, current ramp rate, the output current reading, the output current setting, the output voltage setting, the voltage compliance setting, and the remote voltage sense reading

The keypad is arranged logically to separate the different functions of the instrument. The most common functions of the power supply are accessed using a single button press. The keypad can be locked to either lock out all changes or to lock out just the instrument setup parameters allowing the output of the power supply to be changed.

#### **Specifications**

#### **Output**

Type Bipolar, 4-quadrant, DC current source

Current generation Linear regulation with digital setting and analog control

Current range ±60 A

Compliance voltage ±5 V maximum (nominal, both source and sink)

Maximum power 300 W Load reactance 0 H to 100 H

Current ripple (max) 4 mA RMS at 60 A, (0.007%) into 1 m $\Omega$  load (significantly reduced

into a reactive load or at lower current)

Current ripple frequency Dominated by line frequency and its harmonics

Temperature coefficient ±15 ppm of full scale/°C

Line regulation 15 ppm/6% line change

Source impedance  $25 \Omega$ 

Stability (1 h) 1 mA/h (after warm-up)

Stability (24 h) 10 mA/24 h (typical, dominated by temperature coefficient and line regulation)

Isolation Output optically isolated from chassis to prevent ground loops Parallel operation 2 units can be paralleled for ±120 A, ±5 V operation

Protection Quench, line loss, low line voltage, high line voltage, output over voltage, output over current, over temperature, and remote inhibit (on critical error conditions, magnet discharges at 1 V nominal)

#### **Output programming**

#### Internal current setting

Resolution 0.1 mA (20-bit)

Settling time 600 ms for 1% step to within 0.1 mA into a resistive load

**Accuracy**  $\pm 10$  mA  $\pm 0.05\%$  of setting Operation Keypad, computer interface Protection **Current setting limit** 

#### Internal current ramp

0.1 mA/s to 99.999 A/s (compliance limited) Ramp rate

Update rate 27.7 increments/s

Ramp segments

Operation Keypad, computer interface, and trigger input

**Protection** Ramp rate limit

#### **External current programming**

6 V = 60 ASensitivity Resolution Analog

**Accuracy** ±10 mA ±1% of setting

Bandwidth (3 dB) 40 Hz, 2-pole, low-pass filter (10 Hz pass band,

compliance limited)

Input resistance >50 k $\Omega$ 

Operation Voltage program through rear panel

Connector Shared 15-pin D-sub Internally clamped at 6.1 V Limits

#### Compliance voltage setting

Range 0.1 V to 5.0 V Resolution 100 uV

Accuracy ±10 mV ±1% of reading

#### Readings

#### **Output current**

Resolution 0.1 mA

 $\pm 1$  mA  $\pm 0.05\%$  of reading Accuracy

Update rate 2.5 readings/s display, 10 readings/s interface

Compensation Compensated for lead resistance and 25  $\Omega$  source resistance

**Output voltage (at supply terminals)** Resolution 100 μV

> Accuracy  $\pm 1$  mV  $\pm 0.05\%$  of reading

2.5 readings/s display, 5 readings/s interface Update rate



## **Cryogenic Instruments**

#### Remote voltage (at magnet leads)

Resolution 100 μV

**Accuracy** ±1 mV ±0.05% of reading

**Update rate** 1.25 readings/s Input resistance  $>50 k\Omega$ 

Connector Shared 15-pin D-sub

#### Persistent switch heater output (PSHO)

Current range 10 mA to 125 mA

Compliance voltage (minimum) 12 V or 21 V selectable

Heater resistance (minimum)  $10~\Omega$ 

Setting resolution 1 mA

Accuracy ±1 mA

Operation On/Off with lockout delay of 5 s to 100 s

Protection Open or shorted heater detection, error message if off and on output

currents differ **Connector BNC** 

#### Front panel

Display type 8-line by 40-character, graphic vacuum fluorescent display module Display readings Output current, calculated field (T or G), output voltage, and remote

Display settings Output current, calculated field, compliance voltage, and ramp rate

**Display annunciators** Status and errors

LED annunciators PSHO on, remote, compliance limit, fault, and ramping

Keypad type 26 full travel keys

Keypad functions Direct access to common operations, menu driven setup

#### Interface

IEEE-488.2 interface

**Features** SH1, AH1, T5, L4, SR1, RL1, PP0, DC1, DT1, C0, E1

Reading rate To 10 readings/s

National Instruments LabVIEW™ driver Software support

(see www.lakeshore.com)

Serial interface

Electrical format RS-232C

9600, 19200, 38400, 57600 **Baud rates** 

Reading rate To 10 readings/s

Connector 9-pin D-sub

**Output current monitor** 

Sensitivity 60 A = 6 V**Accuracy** ±1% of full scale

1 mV Noise Source impedance 20 Ω

Connector Shared 15-pin D-sub

**Output voltage monitor** 

Sensitivity 1 V = 1 VAccuracy ±1% of full scale Noise 1 mV

Source impedance  $20 \Omega$ 

Connector Shared 15-pin D-sub

**Fault output** 

Relay (closed on fault) Type Relay contact 30 VDC at 1 A Connector Shared 25-pin D-sub

Remote inhibit input

Type TTL or contact closure Connector Shared 25-pin D-sub

**Trigger input** 

TTL or contact closure Type Connector Shared 25-pin D-sub

Ambient temperature 15 °C to 35 °C Cooling Air cooled with internal 2-speed fan Warm-up 30 minutes at output current setting

Line power 100, 120, 220, 240 VAC +6% -10%, single phase, 50 or 60 Hz, 850 VA Size 483 mm W  $\times$  178 mm H  $\times$  520 mm D (19 in  $\times$  7 in  $\times$  20.5 in), rack mount

(integrated rack mount ears) Weight 27.2 kg (60 lb)

Calibration schedule 1 year

#### Ordering information

#### Part number Description

625 Superconducting magnet power supply—includes two front

handles (6241), two rear handles/protectors (6242), one output shorting bar and terminal fasteners (6243), one 25-pin D-sub digital I/O mating connector (6251), one 15-pin D-sub analog I/O mating connector (6252), a calibration certificate

and a user's manual

625-DUAL Two Model 625s and one 6263 dual supply interconnect cable

#### Please indicate your power/cord configuration:

1 100 V-U.S. heavy duty cord (NEMA 5-15)

2 120 V—U.S. heavy duty cord (NEMA 5-15)

3 220 V—Euro cord (CEE 7/7) 4 240 V—Euro cord (CEE 7/7)

5 240 V-U.K. cord (BS 1363) 6 240 V-Swiss cord (SEV 1011) 7 220 V—China cord (GB 1002)

#### **Accessories**

1 m (3.3 ft) long IEEE-488 (GPIB) computer interface cable 117-017

assembly

6261 10 ft magnet cable kit, AWG 4 20 ft magnet cable kit, AWG 4 6262

6263 Dual supply interconnect cable kit including magnet cables

and safety interlock cable

CAL-625-CERT Instrument recalibration with certificate CAL-625-DATA Instrument recalibration with certificate and data

119-037 Model 625 user manual











# Accessories

144 Cryogenic Accessories

145 Wire

149 Cable


152 Solder

**155 Epoxy** 

156 Grease

157 Varnish

158 Miscellaneous Accessories

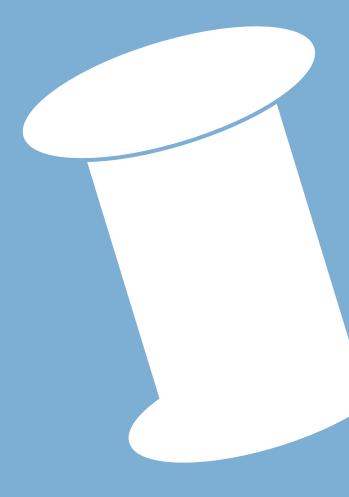



# Cryogenic Accessories

Lake Shore offers a complete line of accessories for sensor installation and general-purpose cryogenic use.

#### **Cryogenic wire**

Used to minimize heat leak into the sensor and cryogenic system, cryogenic wire has a much lower thermal conductivity (and higher electrical resistivity) than copper wire.

The most common type of cryogenic wire is phosphor bronze. This wire is available in one-, two-, and four-lead configurations. Four-lead configurations are available as Quad-twist™ (two twisted pairs) or Quadlead™ (ribbon). Wire gauge is 32 or 36 AWG, with polyimide or polyvinyl formal (Formvar®) used to insulate the wires.

Other common cryogenic wires and coaxial cables include manganin, nichrome heater wire, and HD-30 heavy-duty copper wire. For high-frequency signals, Lake Shore provides various coaxial cables: ultra miniature coaxial cables and semi-rigid coaxial with a stainless steel center conductor.

#### **Solders**

The most common electrical connections are solder joints. Solder can also be used to install various sensors to improve thermal heat sinking. Common solders are indium solder and 90/10 Pb/Sn. Indium solder is used for various applications including sensor installation to provide excellent thermal contact with the sample. 90/10 Pb/Sn solder is used for applications requiring a higher temperature (liquidus point of 575 K and solidus point 458 K). Ostalloy® 158 solder is used as a seal for demountable vacuum cans and electric feedthroughs in cryogenic systems.

#### Varnish, thermal grease, and epoxy

Thermal greases and epoxies are used to install and fasten sensors, while providing thermal contact and/or electrical insulation, with the sample. Epoxy can be used for mechanical attachment and joints.

The most common varnish for cryogenic installations is VGE-7031 varnish. It has good chemical resistance, bonds to a variety of materials, and has a fast tack time. Stycast® 2850FT is composed of a black epoxy resin, and has a thermal expansion coefficient that is matched to copper. A silver-filled, low-temperature conducting epoxy provides excellent strength, along with electrical and thermal conductivity.

Thermal grease, Apiezon N and Apiezon H, is suitable for enhancing thermal contact, especially for sensors inserted into cavities. Apiezon N is for low temperature applications, while H is for high temperature.

#### **Miscellaneous**

Lake Shore also supplies heat sink bobbins, a beryllium oxide heat sink chip, and a four-lead resistance sample holder. Cartridge heaters and vacuum feed through products are also available.



# Wire

Abbreviations used in	this section
American wire gauge	AWG
Single lead wire	SL
Duo-Twist <sup>™</sup> wire	DT
Quad-Twist <sup>™</sup> wire	QT
Quad-Lead <sup>™</sup> wire	QL



#### **Specifications**

		Phosphor bronze	Copper	Nichrome	Manganin
Melting range		1223 K to 1323 K	1356 K	1673 K	1293 K
Coefficient of the	rmal expansion	1.78 × 10 <sup>-5</sup>	20 × 10 <sup>-6</sup>	_	19 × 10 <sup>-6</sup>
Chemical composition (nominal)		94.8% copper, 5% tin, 0.2% phosphorus	_	80% nickel, 20% chromium	83% copper, 13% manganese, 4% nickel
Electrical resistiv (at 293 K)	ity	11 μΩ·cm	1.7 μΩ·cm	120 μΩ·cm	48 μΩ·cm
Thermal	0.1 K	NA	9	NA	0.006
conductivity (W/(m·K))	0.4 K	NA	30	NA	0.02
(**/(111 13))	1 K	0.22	70	NA	0.06
	4 K	1.6	300	0.25	0.5
	10 K	4.6	700	0.7	2
	20 K	10	1100	2.6	3.3
	80 K	25	600	8	13
	150 K	34	410	9.5	16
	300 K	48	400	12	22

	AWG	Resi	Resistance (Ω/m)		Diameter			Number of leads	Name	Insulated	Insulation type	Insulation	Insulation								
		4.2 K	77 K	305 K	(mm)	current air (A)				diameter (mm)		thermal rating (K)	breakdown voltage (VDC)								
Phosphor							,	1	SL-32	0.241	Polyimide										
bronze	32	3.34	3.45	4.02	0.203	4.2	3.1	2	DT-32	0.241	Polyimide	493	400								
	32	3.34	3.43	4.02	0.203	4.2	3.1	4	QT-32	0.241	Polyimide	493	400								
												4	QL-32	0.241	Polyimide						
		36 8.56	8.83 10				1	SL-36	0.152	Formvar®	368	250									
	36			10.3	0.127	107 0.6	1.4	2	DT-36	0.152	Polyimide	493	400								
	30	0.00	0.03		10.5	10.5	10.5	0.127	3 U.121	2.0	121 2.0	2.6 1.4	1.4	1.4	1.4	4	QT-36	0.152	Formvar®	368	250
																				4	QL-36
Nichrome	32	33.2	33.4	34	0.203	2.5	1.8	1	NC-32	0.241	Polyimide	493	400								
Copper	30	0.003	0.04	0.32	0.254	10.2	8.8	1	HD-30	0.635	Teflon®	473	250								
	34	0.0076	0.101	0.81	0.160	5.1	4.4	2	CT-34	0.254	Teflon®	473	100								
Manganin	30	8.64	9.13	9.69	0.254	4.6	4.3	1	MW-30	0.295	Heavy Formvar®		400								
	32	13.5	14.3	15.1	0.203	3.8	3.5	1	MW-32	0.241	Heavy Formvar®	378	400								
	36	34.6	36.5	38.8	0.127	2.6	2.5	1	MW-36	0.152	Heavy Formvar®	1	250								



# Phosphor bronze wire

Phosphor bronze wires (QL, QT, DT and NM), are suitable for almost all cryogenic applications. The low magnetoresistance of these wires make them the ideal choice for magnetic field use.

#### **Physical properties**

Melting range: 1223 K to 1323 K (950 °C to 1050 °C) Coefficient of thermal expansion:  $1.78 \times 10^{-5}$ Thermal conductivity: 48 W/(m · K) at 293 K

Electrical resistivity (annealed):  $1.15 \times 10^{-7} \,\Omega \cdot m$  at 293 K

Specific heat: 376.4 J/(kg · K)

Stress relief temperature (1 h): 423 K to 498 K (150 °C to 225 °C) **Chemical composition:** Nominal 94.8% copper, 5% tin, 0.2% phosphorus

#### Single strand cryogenic wire—SL-32, SL-36

- Phosphor bronze wire
- Non-ferromagnetic
- Single strand
- 32 and 36 AWG
- Polyimide insulation (SL-32)
- Formvar® insulation, clear (SL-36)

Lake Shore non-magnetic (NM) single lead (SL) wire is a phosphor bronze (CuSnP alloy) wire. This wire has a relatively low temperature dependence of its resistance from room temperature to helium temperatures, SL-32 can be used for sensor installations requiring stronger and more 'rugged' leads.SL-36 wire is recommended for general sensor installation.

#### Ordering information

Part number WSL-32-100 WSL-32-250 WSL-36-500

Description 32 AWG, 30 m (100 ft) 32 AWG, 76 m (250 ft) 36 AWG, 152 m (500 ft)





#### Insulation

#### Polyvinyl formal (Formvar®)

Magnet wire is insulated with vinyl acetal resin, as a smooth uniform film. Formvar® has excellent mechanical properties such as abrasion resistance and flexibility. The film will stand excessive elongation without rupture. When stressed during winding, Formvar® has a tendency to craze upon contact with solvents such as toluol, naphtha, and xylol, therefore, it should be given an annealing preheat prior to varnish application. Formvar® can be removed mechanically during terminal preparation. Formvar® is rated to 3525 VAC for 32 AWG, 2525 VAC for 36 AWG.

#### Polyimide (ML)

ML is a film coated insulation made with polyimide resin. It is a Class 220 thermal life insulation with exceptional resistance to chemical solvents and burnout. It will operate at temperatures in excess of 493 K (220 °C) for intermittent duty. ML is unaffected by prolonged exposure to varnish solvents and is compatible with virtually all systems. Polyimide insulation is rated to 3525 VAC for 32 AWG, 2525 VAC for 36 AWG.

Note: At Lake Shore, we strip both Formvar® and polyimide mechanically using an Eraser Rush Model RT-2 mechanical stripper.

147



# Duo-Twist<sup>™</sup> cryogenic wire—DT-32, DT-36

- Phosphor bronze wire
- Non-ferromagnetic
- Single twisted pair (2 wires)
- Color coded leads
- Minimizes pickup noise
- 32 and 36 AWG
- Polyimide insulation

Duo-Twist™ is a single twisted pair (2 leads) of 32 or 36 AWG phosphor bronze wire twisted at 3.15 twists per centimeter (8 twists per inch). This wire is a good choice when any possibility of pickup noise to a diode sensor or sample by induced currents through the leads needs to be minimized.

#### Ordering information

Part number WDT-32-25 WDT-32-100 WDT-32-500 WDT-36-25 WDT-36-100 WDT-36-500

Description 32 AWG, 7.6 m (25 ft) 32 AWG, 30 m (100 ft) 32 AWG, 152 m (500 ft) 36 AWG, 7.6 m (25 ft) 36 AWG, 30 m (100 ft) 36 AWG, 152 m (500 ft)







#### Quad-Twist<sup>™</sup> cryogenic wire—QT-32, QT-36

- Phosphor bronze wire
- Non-ferromagnetic
- 2 twisted pairs (4 wires), color coded
- Minimizes pickup noise
- Polyimide insulation (QT-32)
- Formvar® insulation (QT-36)

Quad-Twist<sup>™</sup> is 2 twisted pairs (4 leads) of 32 or 36 AWG phosphor bronze wire. Each pair incorporates 3.15 twists per centimeter (8) twists per inch), and the 2 pairs are entwined at 1.57 twists per centimeter (4 twists per inch). This wire is a good choice when pickup noise to a diode sensor or sample by induced currents through the leads needs to be minimized. Use one twisted pair for sensor excitation and the other twisted pair for sensor output voltage to minimize pickup of electromagnetic noise.

#### Ordering information

Part number WQT-32-25 WQT-32-100 WQT-32-500 WQT-36-25 WQT-36-100 WQT-36-500

Description 32 AWG, 7.6 m (25 ft) 32 AWG, 30 m (100 ft) 32 AWG, 152 m (500 ft) 36 AWG, 7.6 m (25 ft) 36 AWG, 30 m (100 ft) 36 AWG, 152 m (500 ft)







# Quad-Lead<sup>™</sup> cryogenic wire—QL-32, QL-36

- Phosphor bronze wire
- Non-ferromagnetic
- Four color coded leads
- 32 and 36 AWG
- Polyimide insulation

The Quad-Lead™ wire is a 4-wire "ribbon cable", which makes heat sinking and dressing leads much easier than working with individual wires. Noninductive (bifilar) windings are simple to make for heat sinks and heaters using the Quad-Lead™ wire. In addition, the wire is color coded for easy lead identification, and can be split to yield 2 wire pairs. Quad-Lead™ wire is also useful in standard 4-lead measurements in magnetic field applications due to its low magnetoresistance.

Note: The Quad-Lead™ wires are formed into a "ribbon cable" using Bond Coat 999 bonding film. Wire separation can be accomplished mechanically through the use of a razor blade or other tool equipped with a sharp, flat blade.

# Ordering information

Part number WQL-32-25 WQL-32-100 WQL-32-500 WQL-36-25 WQL-36-100 WQL-36-500

Description 32 AWG, 7.6 m (25 ft) 32 AWG, 30 m (100 ft) 32 AWG, 152 m (500 ft) 36 AWG, 7.6 m (25 ft) 36 AWG, 30 m (100 ft) 36 AWG, 152 m (500 ft)











#### Nichrome heater wire—NC-32

- Nominal 80% nickel, 20% chromium
- Non-ferromagnetic
- 32 AWG
- Polyimide insulation

This high resistance wire is typically used for heater requirements. The relatively large wire size provides sufficient surface area to dissipate the heat generated within the wire with only a moderate rise in wire temperature

Note: We have had poor experience with heaters made using wire smaller than 32 AWG and supplying 25 W or more power. A possible alternative is one of the Lake Shore cartridge heaters, see page 160.

#### Ordering information

Part number WNC-32-100 WNC-32-250

Description

32 AWG, 30 m (100 ft) 32 AWG, 76 m (250 ft)







#### Twisted lead wire—CT-34

- Silver-plated copper, 34 AWG
- Teflon® insulation

These low resistance twisted pair wires are ideal for extending the lead length of Lake Shore cryogenic Hall generators.

# Ordering information

Part number WCT-YB-34-25 WCT-YB-34-50 WCT-YB-34-100 WCT-RB-34-25 WCT-RB-34-50 WCT-RB-34-100 Description Yellow/blue, 7.6 m (25 ft) Yellow/blue, 15 m (50 ft) Yellow/blue, 30 m (100 ft)

Red/black, 7.6 m (25 ft) Red/black, 15 m (50 ft) Red/black, 30 m (100 ft)









#### Heavy duty lead wire—HD-30

- 30 AWG
- Seven 38 AWG silver-plated twisted copper strands
- Black etched Teflon® for adhesion to ероху

This more rugged wire is useful as a lead wire to resistance heaters in cryogenic environments where low resistance to the heater is required or desired.

# Ordering information

Part number WHD-30-100

Description 30 AWG, 30 m (100 ft)









# Manganin wire—MW-30, MW-32, MW-36

- Nominal 83% copper, 13% manganese, and 4% nickel
- Non-ferromagnetic
- 30, 32, and 36 AWG
- Heavy Formvar® insulation

Lake Shore manganin wire is often used for cryostat wiring or heater requirements in nonmagnetic applications.

# Ordering information

Part number WMW-30-100 WMW-30-500 WMW-32-100 WMW-32-500 WMW-36-100

WMW-36-500

Description

30 AWG, 30 m (100 ft) 30 AWG, 152 m (500 ft) 32 AWG, 30 m (100 ft) 32 AWG, 152 m (500 ft) 36 AWG, 30 m (100 ft) 36 AWG, 152 m (500 ft)







149

# **Coaxial Cable**

#### **Specifications**

	Туре С	Type SC	Type SS	Type SR
Dimensions				
Center conductor—AWG (diameter)	32 (0.2032 mm [0.008 in])	32 (0.2032 mm [0.008 in])	32 (0.2032 mm [0.008 in])	37 (0.1143 mm [0.004 in])
Dielectric/insulating material (diameter)	0.56 mm (0.022 in)	0.406 mm (0.016 in)	0.406 mm (0.016 in)	0.38 mm (0.015 in)
Shield (diameter)	0.025 mm (0.001 in) thickness	0.711 mm (0.028 in)	0.711 mm (0.028 in)	0.51 mm (0.02 in)
Drain wire (parallel to conductor)	32 AWG (0.203 mm [0.008 in])	NA	NA	NA
Jacket outer dimension	0.7874 mm × 1.016 mm (0.031 in × 0.039 in)	1.0 mm (0.04 in)	1.0 mm (0.04 in)	0.51 mm (0.02 in)
Material				
Center conductor	Silver-plated copper	Stranded copper <sup>1</sup>	304 stainless steel <sup>2</sup>	Carbon steel <sup>3</sup>
Dielectric/insulating material	Gore-Tex® expanded PTFE	Teflon® FEP	Teflon® FEP	Teflon® PTFE
Shield	Aluminized polyester <sup>4</sup>	Braided gold-plated copper <sup>5</sup>	304 braided stainless <sup>6</sup>	304 stainless steel <sup>7</sup>
Drain wire	Silver-plated copper	NA	NA	NA
Jacket material	FEP	Teflon® FEP	Teflon® FEP	NA
Jacket color	Blue	Gold	Gray	NA
Electrical properties				
Resistance Ω/m (Ω/ft)				
Center conductor at 293 K (20 °C)	0.541 (0.165)	0.282 (0.086)	23.62 (7.2)	4.30 (1.31)
Shield at 296 K (23 °C)	NA	0.085 (0.026)	3.61 (1.1)	8.63 (2.63)
Drain wire at 296 K (23 °C)	0.541 (0.165)	NA	NA	NA
Center conductor maximum DC voltage	150 V	600 V	600 V	700 V
Center conductor maximum DC current	150 mA	200 mA	200 mA	200 mA
Temperature range	10 mK to 400 K	<1 K to 400 K	10 mK to 473 K	10 mK to 400 K
Characteristic impedance	50 Ω (±5 Ω)	35 Ω at 10 MHz	40 $\Omega$ at 10 MHz 50 $\Omega$ (±2 $\Omega$ )	
Nominal capacitance at 5 kHz	79 pF/m (24 pF/ft)	154.2 pF/m (47 pF/ft)	173.9 pF/m (53 pF/ft)	95.14 pF/m (29 pF/ft)

<sup>&</sup>lt;sup>1</sup> 65 strands of 50 AWG

<sup>&</sup>lt;sup>2</sup> 64 strands of 50 AWG 304 SS wire

<sup>&</sup>lt;sup>3</sup> Silver-plated copper-clad carbon steel (0.103 mm outer diameter carbon steel covered by 0.0057 mm thick copper cladding covered by 0.001 mm thick silver plating

<sup>&</sup>lt;sup>4</sup> Aluminized polyester laminated tape, spirally applied at a 40–50% overlap, aluminum side in

 $<sup>^{\</sup>text{5}}$  12  $\times$  3 matrix of 42 AWG wire

 $<sup>^{\</sup>rm 6}$  12 imes 4 matrix of 44 AWG wire

<sup>&</sup>lt;sup>7</sup> A seamless tubular metal jacket serves as the outer conductor/shield



# Ultra miniature coaxial cable – Type C, SC, SS

- Very flexible
- Long flex life
- Available in three configurations:

**C** – solid copper center conductor, drain wire, and aluminized/polyester shield

**SC** – stranded copper conductors SS - stranded 304 stainless steel conductors

Ultra miniature coaxial cable is for use when a strong and flexible cable is needed. Type C and SC are recommended when low conductor resistance is a prime consideration. Type SC and type SS are mechanically the most flexible, due to their braided construction. Type SS is recommended for use when both shielding and low thermal losses are important.

For technical specifications on types SS, C, SC and SR, see page 149.

#### Thermal conductivity of copper—units are W/(m · K)

	4 K	20 K	30 K	77 K	300 K
$RRR^8 = 20$	122	719	870	502	397
RRR = 100	460	2460	2070	533	407

<sup>8</sup> RRR = residual resistance ratio


$$\frac{R_{273\,K}}{R_{4.2\,K}} = RRR$$

	Normal attenuation (dB/m)				
	<b>C</b> <sup>(1)</sup>	SC	SS		
1 MHz	0.092	0.108	0.569		
5 MHz	0.167	0.240	1.272		
10 MHz	0.224	0.344	1.799		
15 MHz	0.257	0.421	2.850		
20 MHz	0.294	0.486	2.545		
50 MHz	0.427	0.769	4.031		
100 MHz	0.623	1.090	5.694		
500 MHz	1.312	2.453	12.749		
1 GHz	1.886	3.488	18.048		
2 GHz	2.625	_	_		
5 GHz	_	7.968	40.526		

Type C has a bandwidth to at least 3 GHz—above that, the aluminum/polyester becomes a less effective

# Ordering information

Part number	Description
CC-C-25	Solid copper, 7.6 m (25 ft)
CC-C-50	Solid copper, 15 m (50 ft)
CC-C-100	Solid copper, 30 m (100 ft)
CC-C-500	Solid copper, 152 m (500 ft)
CC-SC-25	Stranded copper, 7.6 m (25 ft)
CC-SC-50	Stranded copper, 15 m (50 ft)
CC-SC-100	Stranded copper, 30 m (100 ft)
CC-SC-500	Stranded copper, 152 m (500 ft)
CC-SS-25	Stranded stainless, 7.6 m (25 ft)
CC-SS-50	Stranded stainless, 15 m (50 ft)
CC-SS-100	Stranded stainless, 30 m (100 ft)
CC-SS-500	Stranded stainless, 152 m (500 ft)
USD VISA	PayPal

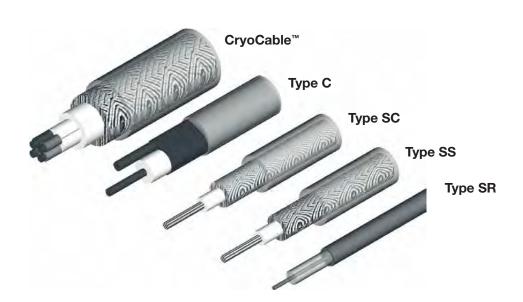



# Semi-rigid coaxial cable—type SR

- Easily bent, coiled, stripped, machined, soldered, or connected without impairing performance
- Solid center conductor provides the optimum geometrical surface for transmission
- Low standing wave ratio (SWR) with a dielectric controlled to exacting tolerances
- Low thermal conductivity ( $\approx 0.4 \text{ W/(m\cdot K)}$  at 4.2 K)9
- Matching minimizes reflective power loss
- Provides shielding isolation for virtually no extraneous signal pickup
- Tubular outer conductor offers minimum size and maximum conductor integrity; stainless steel jacket can be soldered directly to circuit boards
- 37 AWG, silver-plated copper-weld steel center conductor
- Thermal conductivity at low temperatures is dominated by the copper cladding around the center conductor

This cable transmits and receives high-speed, high-frequency microwave signals. Typically used for transmission lines in cryogenic-vacuum test systems.

#### To remove the outer conductor:

- Score jacket
- Bend at score until shield kinks, fatigues, and breaks
- Slide off outer conductor

Extreme caution must be used in this process to avoid damage to the cable

SR coaxial cable frequency response specifications				
	Insertion loss dB/m (dB/ft)	Power CW (20 °C, sea level, W)		
0.5 GHz	4.43 (1.35)	7.6		
1.0 GHz	6.27 (1.91)	5.3		
5.0 GHz	14.09 (4.30)	2.4		
10.0 GHz	20.01 (6.10)	1.7		
20.0 GHz	28.45 (8.67)	1.2		

# Ordering information

Part number CC-SR-10

Description Semi-rigid, 3 m (10 ft)








# CryoCable<sup>™</sup>—type CYRC

- Robust: the NbTi wire cores are strong and fatigue resistant, and the cable overbraid of 304 stainless steel adds significant strength and crush resistance
- Low heat leak due to all metal alloy and Teflon® construction
- **Solderable:** the CuNi wire surface is easy to solder with conventional rosin fluxes
- **Cryo-compatible:** all Teflon® (PFA) insulation is heat strippable for ease of preparation

A robust, 4-wire cable for use in cryogenic environments to room temperature is now available. The cable is designed around 32 AWG (203 µm) diameter superconductive wires consisting of a NbTi core (128 µm diameter) and a Cu-10% Ni jacket.

Minimum bend radius: 15 mm (0.6 in) Critical temperature: 9.8 K Critical field: 10 T

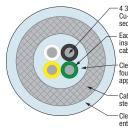

4 32 AWG wires: Nb-48wt%Ti core with Cu-10wt%Ni jacket, CuNi to NbTi cross sectional area ratio = 1.5:1

Each wire overcoated with Teflon® (PFA) insulation 0.003 in (75 µm) thick; wires cabled with approx. 25 mm twist pitch

-Clear Teflon® (PFA) extruded over the four-wire cable to an overall diameter of approx. 1.2 mm (0.048 in)

Cable overbraided with 304 stainless

-Clear Teflon® (PFA) extruded over the entire cable; finished cable has an overall diameter of 2.4 mm  $\pm 0.2$  mm (0.095 in  $\pm 0.008$  in)

	Temperature (K)			
	295	77	4.2	
Wire resistance — per wire ( $\Omega/m$ )	9.2	8.4	010	
Overbraid resistance (Ω/m)	0.90	0.64	0.62	
Thermal conductivity — entire cable assembly $(\Omega/(m \cdot K))$	7.6	2.8	0.17	

<sup>10</sup> Superconducting

Field	Critical current (per wire)
3 T	35 A
5 T	25 A
7 T	15 A
9 T	6 A

# Ordering information

Part number CRYC-32-25 CRYC-32-50 CRYC-32-100 Description

CryoCable<sup>™</sup>, 7.6 m (25 ft) CryoCable<sup>™</sup>, 15 m (50 ft) CryoCable<sup>™</sup>, 30 m (100 ft)









# Solder

	Indium foil	High temperature solder	Ostalloy <sup>®</sup>
Melting point	430 K	Solidus 548 K Liquidus 575 K	343.16 K
Electrical thermal conductivity	84 W/(m·K) at 293 K	35 W/(m·K) at 293 K	18.6 W/(m·K) at 293 K
Resistivity	$9 \times 10^{-4}  \Omega \cdot m$ at 293 K	$204\times10^{\text{-9}}\Omega\text{-m}$ at 293 K	_
Tensile strength	2.61 MPa to 3.55 MPa	30 MPa	<del>-</del>
Density	7.3 g/cm <sup>3</sup>	10.75 g/cm <sup>3</sup>	9.67 g/cm <sup>3</sup>
Composition	99.99% pure Indium	90% Pb 10% Sn	49.5% Bi, 27.3% Pb, 13.1% Sn, 10.1% Cd


#### Indium foil/solder

- Foil form
- Exceptional pressure seal
- Extremely malleable
- 99.99% pure
- Acts as a metallic seal against corrosion
- Flexible sensor mounting material for low stress at cryogenic temperatures

Indium can be used to create solder "bumps" for microelectronic chip attachments and also as gaskets for pressure and vacuum sealing purposes. When used as a washer between a silicon diode or other temperature sensors and refrigerator cold stages, indium foil increases the thermal contact area and prevents the sensor from detaching due to vibration. It also may be used as a sealing gasket for covers, flanges and windows in cryogenic applications.

Indium, a semiprecious, nonferrous metal, is softer than lead, and extremely malleable and ductile. It stays soft and workable down to cryogenic temperatures. It is an excellent choice for cryogenic pumps, high vacuum systems and other unique joining and sealing applications. Indium lends itself to this application due to its characteristic "stickiness" or "tackiness" and ability to conform to many irregular surfaces.

Note: Indium foil becomes a superconductor at 3.38 K (-270 °C), below which the thermal conductivity decreases.

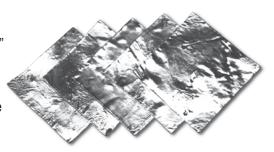

#### **Specifications**

Melting point: 430 K (157 °C)

Thermal conductivity at 293 K (20 °C):  $84 \text{ W/(m} \cdot \text{K)}$ Superconducting transition: 3.38 K (-270 °C)

Volume resistivity ( $\Omega$ ·m):  $8.27 \times 10^{-4}$  at 273 K (0 °C);  $9.00 \times 10^{-4}$  at 293 K (20 °C);  $30.11 \times 10^{-4}$  at 455 K (182 °C)

Thermal expansion coefficient:  $24.8 \times 10^{-6}$  at 300 K(27 °C)

Magnetism: Diamagnetic

**Dimensions:**  $0.127 \text{ mm} \times 50.8 \text{ mm} \times 50.8 \text{ mm}$  (0.005 in  $\times$  2 in  $\times$  2 in)

Tensile strength: 2.61 MPa to 3.55 MPa (380 psi to 515 psi)

Specific heat: 290 J/(kg  $\cdot$  K) at 293 K

# Ordering information

Part number	Description
IF-5	5 indium foil sheets,
	$0.127 \text{ mm} \times 50.8 \text{ mm} \times 50.8$
	mm (0.005 in $\times$ 2 in $\times$ 2 in)
ID-10-31	10 indium disks,
	7.925 mm diameter $\times$ 0.127 mm)
	$(0.312 \text{ in diameter} \times 0.005 \text{ in})$
ID-10-56	10 indium disks,
	14.27 mm diameter $\times$ 0.127 mm
	$(0.562 \text{ in diameter} \times 0.005 \text{ in})$
USD VISA	PayPal

#### High temperature solder

- 90% Pb, 10% Sn
- Good for connecting hardware
- Solidifies quickly

This solder has a higher lead content than normal electronics solder, and can be used for connecting hardware for use at cryogenic temperatures. Its higher melting point also makes it perfect for soldering leads to silicon diode, platinum, or rhodium-iron temperature sensors for operation up to 500 K (227 °C).

#### **Specifications**

**Solidus:** 548 K (275 °C) Liquidus: 575 K (302 °C) **Density:** 10.75 g · cm<sup>-3</sup> Diameter: 0.787 mm (0.031 in)

#### Ordering information

Part number **SLT-10** 

Description 90% Pb, 10% Sn solder,

3 m (10 ft)







# Ostalloy® 158 solder

- Does not shrink, but exhibits expansion upon solidification
- Low melting temperature 343 K (70 °C), requiring only a simple melting pot and a gas or electric heat source
- Reusable many times
- Oxide separated easily in hot water
- Solidifies quickly
- Creates almost no dross because of its low melting temperature

This is a low melting point solder, nearly identical to what is commonly called Wood's Metal. An alloy of bismuth, tin, lead, and cadmium, it is an eutectic alloy with a sharply defined melting point of 343.16 K (70 °C). Ostalloy® 158 has proven itself in production processes—there is no equal to be found to its special advantages.

Mainly used as sealing for demountable vacuum cans and electric feedthroughs in cryogenic testing facilities. Good for soldering any items which cannot be subjected to high temperatures. Ostalloy® 158 solder is used for tool fixturing, holding small parts to be machined, tube shaping and bending, nesting fixturing dies, and internal and external support of thin walled tools and parts. This solder is not recommended for general temperature sensor lead attachment due to its low joint strength.



#### **Specifications**

Composition of Ostalloy® 158 solder: 49.5% Bi, 27.3% Pb, 13.1% Sn, 10.1% Cd

# Ordering information

Part number S0SY-16

Description Ostalloy® 158 solder, 454 g (16 oz)









# Epoxy, Grease, and Varnish

#### **Specifications**

		Conductive epoxy	Stycast® epoxy	Apiezor	ı grease	Varnish	
				Type N	Type H		
Maximum temperat	ture	573 K	403 K	303 K	513 K	423 K	
Glass transition ten	perature	≥353 K	359 K	_	_	_	
Thermal conductivi	ty						
	1 K	_	0.0065 W/(m · K)	0.001 W/(m · K)	_	0.034 W/(m · K)	
	4.2 K	_	0.064 W/(m · K)	0.005 W/(m · K)	_	0.062 W/(m · K)	
	77 K	_	_	_	_	0.22 W/(m · K)	
	100 K	_	_	0.11 W/(m · K)	_	0.24 W/(m · K)	
	300 K	2.5 W/(m · K)	1.3 W/(m · K)	0.19 W/(m · K)	0.22 W/(m · K)	0.44 W/(m · K)	
Thermal expansion	(1/K)	$>$ 360 K: 150 $\times$ 10 <sup>-6</sup> $<$ 360 K: 43 $\times$ 10 <sup>-6</sup>	29 × 10 <sup>-6</sup>	0.00072	0.00072	_	
Volume resistivity		At 298 K: 0.0001 to 0.0004 (Ω·cm)	298 K: $5 \times 10^{14}$ (Ω·m) 394 K: $1 \times 10^{10}$ (Ω·m)	$2 \times 10^{16}$ ( $\Omega$ ·m)	$4.6 \times 10^{13}$ ( $\Omega$ ·m)	>1 × 10 <sup>15</sup> (Ω·m)	
Shelf-life (298 K ma	ix)	12 months from date of manufacture	12 months from date of manufacture	_	_	12 months from date on the can when stored at room temperature	
Pot life		4 days, ~1 day working time	45 min, ~20 min working time	_	_	_	
Cure schedule		323 K: 12 h 353 K: 90 min 393 K: 15 min 423 K: 5 min 448 K: 45 s	298 K: 16 to 24 h 318 K: 4 to 6 h 338 K: 1 to 2 h	NA	NA	5 min to 10 min drying time	
Dielectric strength		NA (conductive)	14.4 kV/mm	_	_	(Dry) 118 kV/mm	
Dielectric constant		NA (conductive)	(1 mHz): 5.01	_	_	_	
Vapor pressure		_	<13.3 Pa (0.1 Torr) at 298 K	$2.27 \times 10^{-7} \text{ Pa}$ (1.7 × 10 <sup>-9</sup> Torr) at 293 K	$2.27 \times 10^{-7}  \text{Pa}$ (1.7 × 10 <sup>-9</sup> Torr) at 293 K	Partial	
Outgassing		_	TML: 0.25% CVCM: 0.1%	TML: <1% CVCM: <0.1%	TML: <1% CVCM: <0.1%	_	



# **Epoxy**

#### Low temperature conductive epoxy

- Excellent low temperature thermal and electrical conductivity
- Low viscosity
- Thixotropic
- No resin bleed during curing
- Low weight loss
- Low volatility

This epoxy is used to permanently attach test samples or temperature sensors to sample holders. It is a 100% solid, two component, low temperature curing, silver-filled epoxy which features very high electrical and thermal conductivity combined with excellent strength and adhesive properties.

Note: Epoxy must be cured at a minimum of 50 °C for 12 h to achieve proper electrical and physical properties. Curing at 175 °C for 45 s will achieve optimum properties.

ESF-2-5 and ESF-2-10 can be used to 300 mK and below. Results may vary based on application and materials used.

#### **Specifications**

Maximum operating temperature: 573 K (300 °C) Thermal conductivity: 300 K (27 °C)—2.5 W/(m · K) Thermal expansion coefficient (K-1): Above 360 K (85 °C)—150  $\times$  10<sup>-6</sup>; below 360 K (85 °C)—43  $\times$  10<sup>-6</sup> Volume resistivity (Ω-cm) at 298 K (25 °C): 0.0001 to 0.0004

Shelf life (25 °C [298 K] max): 12 months from date of manufacture

Pot life: 4 days, about 1 day working time Cure schedule: 323 K (50 °C)-12 h; 353 K (80 °C)-90 min; 393 K (120 °C)—15 min; 423 K (150 °C)—5 min; 448 K (175 °C)-45 s

#### Ordering information

#### Part number Description

ESF-2-5 Low temperature conductive epoxy, 5 packets, 2 g each ESF-2-10 Low temperature conductive

epoxy, 10 packets, 2 g each



# Stycast® epoxy 2850-FT, catalyst 9

- Mixed and applied from two-part flexible packets
- Excellent low temperature properties
- Permanent mounting
- Exceptional electrical grade insulation properties
- Low cure shrinkage
- Low thermal expansion
- Resistance to chemicals and solvents

Stycast® is the most commonly used, highly versatile, nonconductive epoxy resin system for cryogenic use. The primary use for Stycast® is for vacuum feedthroughs or permanent thermal anchors. Lake Shore uses this product in vacuum tight lead-throughs with excellent thermal cycle reliability. Stycast® is an alternative to Apiezon® N grease when permanent sensor mounting is desired. (Can place stress on sensor—see Appendix C.)

Note: Can be chemically removed with methylenechloride (several hour soak). A commercial stripper is available from Miller-Stephenson Co. (http://www.millerstephenson.com/) part number MS-111.

Shipped as a Dangerous Good.

#### **Specifications**

Maximum operating temperature: 403 K (130 °C) Glass transition temperature: 359 K (86 °C)

Thermal conductivity:

1 K (272 °C)-0.0065 W/(m · K) 4.2 K (269 °C)-0.064 W/(m · K) 300 K (27 °C)-1.3 W/(m · K)

Thermal expansion coefficient (1/K):  $29 \times 10^{-6}$ Volume resistivity [Ω·m]

298 K (25 °C)-5 × 1014 394 K (121 °C)-1 × 1010

Shelf life (25 °C [298 K] max): 12 months from date of manufacture

Pot life: 45 minutes, about 20 minutes working time Cure schedule: 298 K (25 °C)—16 h to 24 h

318 K (45 °C)—4 h to 6 h 338 K (65 °C)—1 h to 2 h

Dielectric strength: 14.4 kV/mm Dielectric constant (1 MHz): 5.01

Vapor pressure at 298 K (25 °C): <13.3 Pa (0.1 Torr) **CVCM:** 0.01% **Outgassing TML:** 0.25%

# Ordering information

Part number

Description

ES-2-20 Stycast® epoxy, 20 packets,

2 g each









# Grease

# Apiezon® grease—Types N and H

**Accessories** 

- Stable
- Nonpermanent sensor mounting
- Chemically inert
- Nontoxic
- Easily applied and removed
- Excellent lubrication properties

Apiezon® grease is well-suited for cryogenic use because of its low vapor pressure and high thermal conductivity. It is often used for nonpermanent mounting and thermal anchoring of cryogenic temperature sensors as well as for lubricating joints and o-rings.

Apiezon® N: this general purpose grease enhances thermal contact and provides a temporary mounting method for temperature sensors. It is pliable at room temperatures and solidifies at cryogenic temperatures, which makes it easy to apply and remove the sensor (without damage) at room temperature. The grease is not an adhesive and will not necessarily hold a sensor or wires in place without some mechanical aid, such as a spring clip or tape. It is very good for sensors inserted into holes. Contains a high molecular weight polymeric hydrocarbon additive which gives it a tenacious, rubbery consistency allowing the grease to form a cushion between mating surfaces.

Apiezon® H: this grease will withstand temperatures up to 523 K (250 °C) without melting. It is designed for general purposes where operating temperatures necessitate the use of a relatively high melting point grease.

Note: Can be removed using Xylene with an isopropyl alcohol rinse.



#### **Specifications**

	Type N	Type H		
Approximate melting point:	316 K (43 °C)	523 K (250 °C)		
Thermal conductivity:				
293 K (20 °C)	0.19 W/(m·K)	0.22 W/(m·K)		
1 K (-272 °C)	0.001 W/(m·K)	_		
4.2 K (-269 °C)	0.005 W/(m·K)	_		
100 K (-173 °C)	0.15 W/(m·K)	_		
300 K (27 °C)	0.44 W/(m·K)	_		
Volume resistivity:	$2 \times 10^{16} \Omega \cdot m$	$4.6 \times 10^{13} \Omega$ ·m		
Thermal expansion coefficient (K <sup>-1</sup> ):	0.00072			
Vapor pressure at	2.67 × 10 <sup>-7</sup> Pa	3.60 × 10 <sup>-7</sup> Pa		
293 K (20 °C):	$(2 \times 10^{-9} \text{ Torr})$	$(2.7 \times 10^{-9}  \text{Torr})$		
Solvent system:	n: Hydrocarbons or chlorinated solvents			

# Ordering information

Part number
GAN-25
GAH-25
GAH-25
Apiezon® N grease, 25 g tube
Apiezon® H grease, 25 g tube
USD VISA
PayPar

# Varnish

#### VGE-7031 varnish

- Clear modified phenolic
- Can be air-dried or baked
- Use up to 470 K for 1 to 2 hour maximum
- Varnish will not outgas after baking
- Can be used in vacuum (1.33  $\times$  $10^{-6} \text{ Pa} [9.98 \times 10^{-9} \text{ Torr}]$
- Superior electrical properties
- Excellent chemical resistance
- May be applied by dipping, roller coating, brushing, or spraying
- Moderately good, low stress adhesive
- Enhances thermal contact

VGE-7031 insulating varnish and adhesive possesses electrical and bonding properties which, when combined with its chemical resistance and good saturating properties. make it an excellent material for cryogenic temperatures. As an adhesive, VGE-7031 bonds a variety of materials, has fast tack time, and may be air-dried or baked. It is excellent for laminating many types of materials, and may be applied to parts to be bonded and either baked shortly after applying or allowed to air dry and baked after the parts are stored and assembled hours, days, or even weeks later. It is also an electrically insulating adhesive at cryogenic temperatures, and is often used as a calorimeter cement. VGE-7031 is compatible when dry with a wide variety of materials, including cotton, Dacron® polyester fiber, nylon glass tapes, laminates, Mylar® polyester film, mica products, polyester products, vinyl products, wire enamels, paints, rayon, plastics, and metals. When soaked into cigarette paper, it makes a good, high thermal conductivity, low electrical conductivity heat sinking layer.

Note: May be thinned to the desired application viscosity with a 50:50 mix of denatured alcohol and toluene.

The solvents in the varnish have a tendency to craze Formvar® wire insulation. The wire cannot be disturbed during curing of the varnish (typically 12 to 24 hours at room temperature).

Classified as hazardous cargo by the U.S. Government. UPS Ground shipment only.

Available in continental U.S. only.

#### **Specifications**

Maximum operating temperature: 423 K (150 °C)

Thermal conductivity:

1 K (-272 °C)-0.034 W/(m·K) 4.2 K (-269 °C)---0.062 W/(m·K) 77 K (-196 °C)-0.22 W/(m·K)

100 K (-173 °C)-0.24 W/(m·K) 300 K (27 °C)-0.44 W/(m·K)

Percent solids by weight: 18 to 20%

Viscosity at 298 K (25 °C): 1.3 kg/(m·s) (1300 cP) Specific gravity at 298 K (25 °C): 0.88

Flash point, closed cup: 269 K (-4 °C)

Shelf life: 12 months from date on the can when stored at room temperature

Drying time (25 µm film, tack free): 5 min to 10 min at 298 K (25 °C); 2 min to 5 min at 398 K (125 °C) Solvent system: Xylene, alcohol, acetone

#### Ordering information

Part number VGF-7031

Description

Insulating varnish and adhesive, 0.47 L (1 pt) can










# Miscellaneous Accessories

#### Heat sink bobbins

Heat sink bobbins for cryostat lead wires are gold-plated OFHC or ETP copper for removing heat flowing down sensor leads. The small bobbin holds 4 to 8 phosphor bronze or manganin wires, and the large bobbin holds up to 40, depending on wire gauge and number of wraps. 4 or 5 wraps are usually sufficient, using VGE-7031 varnish or Stycast® epoxy for potting the wires. Do not use copper or other high conductivity wires.

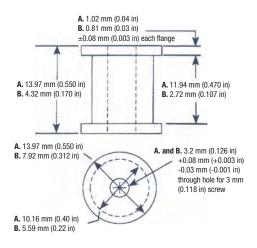


# Ordering information Part number HSB-40 Large heat sink bobbin (use "A" dimensions) Small heat sink bobbin (use "B" dimensions) USD VISA (SP) PayPail

#### Beryllium oxide heat sink chip

Beryllium oxide heat sink chips can be used to heat sink electrical leads or samples at low temperature with good electrical isolation. They can also be used as a buffer layer to take up expansion mismatch between an object with large expansion coefficient (e.g., copper, epoxy) and an object with a low expansion coefficient (e.g., a DT-670-SD diode sensor). One side is fully metallized with molybdenum/ manganese, followed by nickel and gold. It is easily soldered with In/Ag solders. Sn/Pb solders can pull up metallization under some circumstances. The other side has two 1.27 mm (0.05 in) by 4.06 mm (0.16 in) electrically isolated solder pads. The thermal conductivity is several times that of copper in the liquid nitrogen region but about 1000 times lower at liquid helium temperature. The magnetic susceptibility is about that of nonmagnetic stainless steel.

**Note:** Due to metallization irregularities and surface dirt, it is not recommended that these chips isolate more than 100 V.

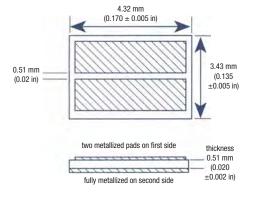


# Ordering information Part number Description HSC-4 Heat sink chip (package of 10)

#### 19-pin vacuum feedthrough


This hermetically sealed glass-to-metal electronic connector is designed to meet the dimensional requirements of MIL-C-26482 and is furnished with a silicone o-ring to seal against the mating connector plug shell. It is commonly used to pass electrical signals into a vacuum chamber from the outside.

Note: The VFT19-FMC threads should be sealed with Teflon® tape or epoxy if a vacuum seal is important.


#### VFT19-FMC

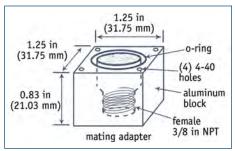


#### VFT19

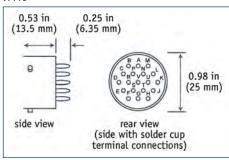

#### VFT19-MC

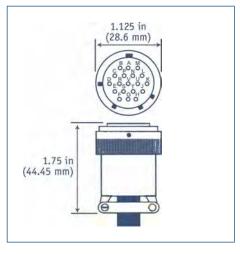

#### VFT19-F

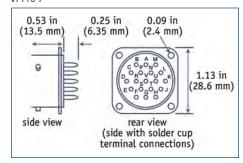

#### **Specifications**

Shell: Mild steel


Contacts: High nickel iron alloy Finish: Fused tin over cadmium

# Ordering information

Ordering information							
Part number	Description						
VFT19	19-pin vacuum feedthrough						
VFT19-F	19-pin vacuum feedthrough in						
	flange						
VFT19-FMC	Mating adapter for mounting						
	VFT19-F to 3/8 NPT hole pipe						
	feedthrough						
VFT19-MC	Mating connector plug to VFT19						
USD VISA	PayPal 9-F						

# 4-lead resistance sample holder

- 4 pre-tinned and drilled solder pads
- Plug-in convenience (4-pin plug)
- Mating socket included

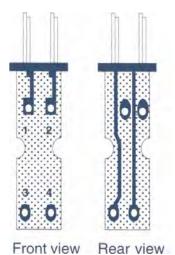

#### **Specifications**

Temperature range: 4.2 K to 373 K (-269 °C to 100 °C)

Current: 1 A at 100 VDC

Insertion force: 227 g (8 oz) per pin

**Dimensions:** 5.1 mm wide  $\times$  27.9 mm long (0.2 in wide

 $\times$  1.1 in long)

Hole diameter: 0.8 mm (0.03125 in)

Hole spacing: 2.5 mm (0.1 in) between holes 1 & 2 and 3 & 4; 15.2 mm (0.6 in) between holes 1 & 4 and 2 & 3

Mating connector: Black thermoplastic Sockets: Phosphor bronze with gold over nickel Socket diameter: 0.41 mm to 0.51 mm (0.016 in to 0.020 in square)

Socket depth: 2.03 mm to 6.35 mm (0.080 in to 0.25 in)

# Ordering information

Part number 700RSH

Description

4-lead resistance sample holder and mating connector; 200 cycle minimum when used below room temperature










#### Cartridge heaters

Precision-wound nickel-chromium resistance wire

**Accessories** 

- Efficient magnesium oxide insulation
- CSA component recognition
- 2 solid pins
- High-temperature rated INCOLOY® case\*
- \*The nickel lead wires and INCOLOY (iron/nickel/chromium alloy) case make these heaters unsuitable for magnetically sensitive locations

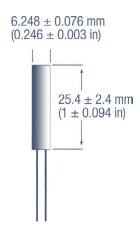

Lake Shore cartridge heaters can be used with all of our temperature controllers. Heaters have wattage ratings in dead air. In cryogenic applications, these cartridge heaters can handle many times the rated value if properly heat sunk or in liquid.

#### **Specifications**

**Diameter:** 6.248 mm  $\pm 0.076$  mm (0.246 in  $\pm 0.003$  in) recommended to fit hole of 6.35 mm (0.25 in) Insulation between leads and case: Magnesium oxide\*\* **Leads:** Nickel, 0.635 mm (0.025 in) diameter  $\times$  50.8 mm

Dielectric strength of insulation is reduced when hot, forming leakage current



#### Electrical tape for use at cryogenic temperatures

- Excellent tape for use at cryogenic temperatures—does not degrade with time like masking tape
- CHR Industries electrical tape
- Yellow polyester film

#### Specifications at 25 °C

Backing: Polyester film

Temperature class (upper limit): 403 K (130 °C)

**Total thickness:** 0.064 mm (0.0025 in)

Dielectric breakdown: 5 kV Insulation resistance: >1 M $\Omega$ **Breaking strength:** 55 N (12.5 lb) **Elongation:** 100% at break

# Ordering information

Part number T3M-72

Description 1 roll cryogenic tape  $12.7 \text{ mm} \times 65.8 \text{ m}$ 

 $(0.5 \text{ in} \times 72 \text{ yd})$ 


USD VISA





# Ferrite bead for high frequency filtering

RF pickup can affect an experiment by upsetting the instrument reading, by being rectified by a diode thermometer to appear as an offset, or by transmitting through the system wiring to pollute the experimental environment. A ferrite bead will reduce the effect of RF pickup on instrument leads by acting like a high impedance (resistance) to high frequency noise. DC and slow moving signals are not affected. The bead can be clamped around existing wiring for ease of installation.

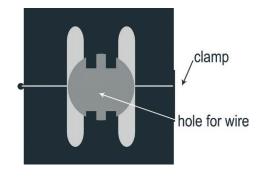

#### **Specifications**

Material: Fair-Rite® 43

Impedance with wire passed once through bead: 110  $\Omega$  at 25 MHz. 225  $\Omega$  at 100 MHz

Impedance with wire passed twice through bead: 440  $\Omega$  at 25 MHz, 900  $\Omega$  at 100 MHz

Construction: 2 halves of a ferrite bead held by a plastic

**Overall dimensions:** 22.1 mm  $\times$  23.4 mm  $\times$  32.3 mm  $(0.87 \text{ in} \times 0.92 \text{ in} \times 1.27 \text{ in})$ 

Cable opening diameter: 10.2 mm (0.4 in) Temperature range: 288 K to 308 K (15 °C to 35 °C)

Weight: 0.046 kg (0.1 lb)

# Ordering information

Part number 2071

Description Ferrite bead





# **Appendices**

164 Appendix A: Overview of Thermometry

166 Appendix B: Sensor Characteristics

176 Appendix C: Sensor Packaging and Installation

189 Appendix D: Sensor Calibration Accuracies

198 Appendix E: Temperature Measurement System

207 Appendix F: PID Temperature Control

210 Appendix G: Sensor Temperature Response

**Data Tables** 

217 Appendix H: Common Units and Conversions

219 Appendix I: Cryogenic Reference Tables





# Appendix A: Overview of Thermometry

#### **General thermometry and temperature scales**

Thermodynamically speaking, temperature is the quantity in two systems which takes the same value in both systems when they are brought into thermal contact and allowed to come to thermal equilibrium. For example, if two different sized containers filled with different gasses at different pressures and temperatures are brought into thermal contact, after a period of time, the final volumes, pressures, entropies, enthalpies, and other thermodynamic properties of each gas can be different, but the temperature will be the same.

Thermodynamically, the ratio of temperature of two systems can always be determined. This allows a thermodynamic temperature scale to be developed, since there is an implied unique zero temperature. Additionally, it allows the freedom to assign a value to a unique state. Therefore, the size of a temperature unit is arbitrary.

The SI temperature scale is the Kelvin scale. It defines the triple point of water as the numerical value of 273.16, i.e., 273.16 K. The unit of temperature in this scale is the kelvin (K).

Another scale is the Rankine scale, where the triple point of water is defined as the value 491.688 °R (degrees Rankine). On the Rankine scale, temperature is 9/5 the Kelvin temperature.

The Kelvin and Rankine scales are both thermodynamic, however, other non-thermodynamic scales can be derived from them. The Celsius scale has units of °C (degrees Celsius) with the size of the unit equal to one Kelvin.

$$T(^{\circ}C) = T(K) - 273.15$$
 Eqn. 1

While the Fahrenheit scale is defined as

$$T(^{\circ}F) = T(^{\circ}R) - 459.67$$
 Eqn. 2

Additionally,

$$T(^{\circ}C) = [T(^{\circ}F) - 32]*(5/9)$$
 Eqn. 3

Both Celsius and Fahrenheit are non-thermodynamic temperature scales, i.e., the ratio of temperature is not related to thermodynamic properties (a 50 °F day is not two times "hotter" than a 25 °F day!) These scales are used for their pragmatic representation of the range of temperature that is experienced daily.

At the most basic level, a thermometer is a device with a measurable output that changes with temperature in a reproducible manner. If we can explicitly write an equation of state for a thermometer without introducing any unknown, temperature-dependent quantities, then we call that thermometer a **primary** thermometer. These include the gas thermometer, acoustic thermometer, noise thermometer, and total radiation thermometer. A secondary thermometer has an output that must be calibrated against defined fixed temperature points. For example, a platinum resistance temperature detector (RTD) is based on the change in resistance of a platinum wire with temperature.

Since primary thermometers are impractical (due to size, speed, and expense), secondary thermometers are used for most applications. The common practice is to use secondary thermometers and calibrate them to an internationally recognized temperature scale based on primary thermometers and fixed points. The most recent efforts in defining a temperature scale have resulted in the International Temperature Scale of 1990 (ITS-90) and the Provisional Low Temperature Scale of 2000 (PLTS-2000).

The ITS-90 is defined by 17 fixed points and 4 defining instruments. It spans a temperature range from 0.65 K to 10,000 K. For cryogenic purposes the three defining instruments are helium vapor pressure thermometry, gas thermometry, and platinum resistance thermometry.

For temperature below 1 K there is the Provisional Low Temperature Scale of 2000 (PLTS 2000). The PLTS-2000 is defined by a polynomial relating the melting pressure of He3 to temperature from the range 0.9 mK to 1 K. The pressure to temperature relationship is based on primary thermometers such as Johnson noise and nuclear orientation. Realization of the PLTS-2000 requires a helium-3 melting pressure thermometer (MPT). For the best realization of PLTS-2000, an MPT with an absolute pressure standard is used. This is a costly and time consuming method. Another method is to use the MPT as an interpolating instrument in conjunction with superconducting fixed points.

Few, if any, individuals or laboratories can afford the expense of maintaining the equipment necessary for achieving the ITS-90 and PLTS-2000. It is more customary to purchase thermometers calibrated by a standards laboratory. Even then, this thermometer is typically two or three times removed from primary thermometers.

Normally the temperature scale, once defined, is transferred from the primary thermometers to secondary thermometers maintained by government agencies, such as the National Institute of Standards and Technology (NIST), the National Physical Laboratory (NPL), or the Physikalisch-Technische Bundesanstalt (PBT). The most common of these secondary thermometers is the resistance thermometer, which is normally a high purity platinum or a high purity rhodium-iron alloy. Standards grade platinum resistance thermometers are referred to as standard platinum resistance thermometers (SPRT) while rhodiumiron resistance thermometers are referred to as RIRTs. Both materials are highly stable when wire-wound in a strain-free configuration. These standards grade resistance thermometers are maintained for calibrating customers' thermometers in a convenient manner. A standards laboratory would maintain a temperature scale on a set of resistance thermometers calibrated by that government agency. This is extremely expensive and time consuming. Thus, primary standards would not be used in day-to-day operation. Instead, the standards laboratory would calibrate a set of working standards for that purpose. These are the standards used to calibrate thermometers sold to customers. Each step in the calibration transfer process introduces a small additive error in the overall accuracy of the end calibration.

In addition to the sensor calibration process, there is also a class of sensors where the manufacturing process is highly reproducible. All of these sensors have a similar output to temperature response curve to within a specified tolerance. Industrial grade platinum thermometers and silicon diodes are examples of sensors that are interchangeable, i.e., their output as a function of temperature (R vs. T or V vs. T) is so uniform that any sensor can be interchanged with another—without calibration—and the temperature reading will still be accurate. The level of accuracy is specified by tolerance bands. With silicon diodes it is possible for a sensor to be interchangeable to within 0.25 K.

#### **References:**

Schooley, James F. **Thermometry**. Boca Raton, Florida: CRC Press Inc., 1986.

Quinn, T.J. Temperature. Academic Press, 1983.

Callan, H.B. **Thermodynamics and an introduction to Thermostatistics, Second Edition**, New York: Wiley, 1985.

Mangum, B. W. and G. T. Furukawa. **Guidelines for Realizing the International Temperature Scale of 1990 (ITS-90)**. NIST Technical Note 1265, 1990.

#### **Fixed points**

Repeatable temperature points are referred to as fixed points. These are simply points that occur reproducibly at the same temperature. There are numerous examples of fixed points. These include boiling points, freezing points, triple points, superconducting transition points, and superfluid transition points.

Figure 1 shows a typical pressure-temperature phase diagram. Matter can exist in three states: solid, liquid, and gas. The pressure-temperature diagram intuitively makes sense. If we heat matter to a high enough temperature, it becomes gaseous. If we subject matter to a high enough pressure, it becomes a solid. At combinations of pressure and temperature in between these limits, matter can exist as a liquid. The boundaries that separate these states of matter are called the melting (or freezing) curve. the vaporization (or condensation) curve and the sublimation curve. The intersection of all three curves is called the triple point. All three states of matter can coexist at that pressure and temperature. When we say the freezing point or boiling point of a substance is reproducible, it is implied that we are measuring that point at the same nominal pressure as in previous measurements. As is shown in the diagram, there is not a single freezing point or a single boiling point. There are an infinite number of freezing points and boiling points which form the boundaries between the solid and liquid states of matter. There is, however, a single triple point, which makes it inherently reproducible. There is only one combination of pressure and temperature for a substance that allows the triple point to be obtained.

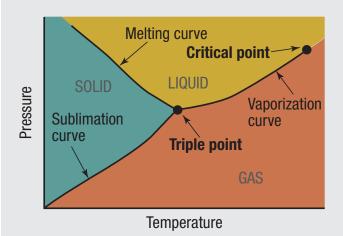


Figure 1 – generic pressure vs. temperature curve



# Appendix B: Sensor Characteristics

#### Types of temperature sensors

Any temperature dependent parameter can be used as a sensor if it fits the requirements of the given application. These parameters include resistance, forward voltage (diodes), thermal EMFs, capacitance, expansion/contraction of various materials, magnetic properties, noise properties, nuclear orientation properties, etc. The two most commonly used parameters in cryogenic thermometers are voltage (diodes) and resistance. There are distinct reasons for choosing diode thermometry or resistance thermometry.

#### **Diodes**

A diode temperature sensor is the general name for a class of semiconductor temperature sensors. They are based on the temperature dependence of the forward voltage drop across a p-n junction. The voltage change with temperature depends on the material. The most common is silicon, but gallium arsenide and gallium aluminum arsenide are also used.

Silicon diodes can be used from 1.4 K to 500 K. From 25 K to 500 K, a silicon diode has a nearly constant sensitivity of 2.3 mV/K. Below 25 K the sensitivity increases and is nonlinear. The temperature response curve is shown in Figure 1. Diode temperature sensors from Lake Shore (DT-670 Series) typically are mounted in a special semiconductor package (SD package). The semiconductor packaging is robust and allows for solder mounting for probes and circuits and easy installation and handling.

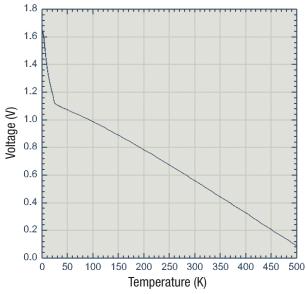


Figure 1 - Curve DT-670

Silicon diode sensors are typically excited with a constant 10 µA current. The output signal is fairly large: 0.5 V at room temperature and 1 V at 77 K. This can be compared to platinum where a 100  $\Omega$ PRT with a 1 mA excitation has only a 100 mV signal at 273 K. The straightforward diode thermometry instrumentation is shown in Figure 2.

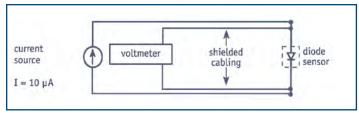


Figure 2 - Typical diode sensor instrumentation schematic

An important feature of silicon diodes is their interchange-ability. Silicon diodes from a particular manufacturer are interchangeable, or curve-matched over their whole range. This is typically defined in terms of tolerance bands about a standard voltage-temperature response curve. They are classified into different tolerance bands with the best accuracy being about  $\pm 0.25$  K from 2 K to 100 K and  $\pm 0.3$  K from 100 K to 300 K.

The large temperature range, nearly linear sensitivity, large signal and simple instrumentation make the diode useful for applications that require a better accuracy than thermocouples. Also, because of the large signal, a diode can be used in a two-lead measurement with little lead resistance error. AC noise-induced temperature errors, to which resistors are immune (aside from heating effects), can be prevalent in diodes.

#### **Resistors**

Temperature sensors based on the changing resistance with temperature can be classified as positive temperature coefficient (PTC) or negative temperature coefficient (NTC). Platinum RTDs are the best example of PTC resistance sensors. Other PTC RTDs include rhodiumiron, nickel, and copper RTDs. Figure 3 shows a typical resistance sensor instrumentation schematic.

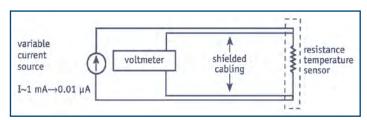


Figure 3 - Typical resistance sensor instrumentation schematic

A PTC RTD is typically metallic (platinum) and has a fairly linear temperature-resistance response. NTC RTDs are semiconductors or semi-metals (doped germanium, Cernox®). They have extremely nonlinear response curves, but are much more sensitive to temperature change.

**Positive temperature coefficient (PTC) RTDs**—The most common type of PTC RTD is platinum. Platinum RTDs are the industry standard due to their accuracy and reproducibility over a wide range of temperatures, as well as their interchangeability. Measurements in the range from -258 °C to 600 °C are routinely made with a high degree of accuracy using platinum RTDs. Industrial-grade platinum RTDs are wire-wound devices that are encapsulated in glass or ceramic, making them durable for general-purpose use.


Platinum RTDs follow a standard response curve to within defined tolerances (IEC 751). The industry standard for class B accuracy is specified as  $\pm 0.3$  K and  $\pm 0.75\%$  variation in the specified 0.00385 K $^{-1}$  temperature coefficient of resistance at 273 K. Below 70 K, a platinum RTD is still usable but requires an individual calibration.

Like all resistors, platinum RTDs can be measured by current excitation or voltage measurement. Common configurations are two-, three-, and four-lead measurements. Two-lead measurements do not correct for lead resistance, so therefore can only be used in applications where the sensor is close to a temperature transmitter. Because their resistance change with temperature is linear over a wide range, a single current excitation (1 mA) can be used for the whole range.

**Negative temperature coefficient (NTC) RTDs**—NTC resistors are normally semiconductors with a very strong temperature dependence of resistance, which decreases with increasing temperature. It is not uncommon for the resistance to change five orders of magnitude over their useful temperature range. The three most common are germanium, Cernox®, and ruthenium oxide (Rox $^{\text{TM}}$ ) RTDs. Carbon-glass RTDs are still used, but they are generally being replaced by Cernox® for nearly all applications.

Cernox® is the trade name for zirconium oxy-nitride manufactured by Lake Shore Cryotronics, Inc. It is a sputter-deposited thin film resistor. Cernox shows good temperature sensitivity over a wider range (0.1 K to 420 K) and is highly resistant to magnetic field-induced errors and ionizing radiation. Cernox can be packaged in the same robust hermetically sealed SD package (Figure 4) that is used for diode temperature sensors. This makes Cernox more robust than other NTC RTDs.

Germanium and carbon-glass (Figures 5 and 6) have very large sensitivities, but more narrow operating ranges than Cernox. Germanium is very stable and is recognized as a secondary standard for T < 30 K. Both sensors are piezoresistive, so the sensing element must be mounted in a strain-free package, which

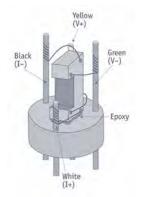



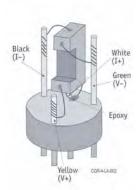


Figure 6 – Typical carbon-glass packaging

provides a very weak thermal link to their surroundings. Both sensors are sealed in a helium atmosphere, but at lower temperatures the pressure is very low and the gas eventually liquefies, reducing the thermal contact. The requirement of strain-free mounting also results in a very fragile sensor. Dropping a sensor from a height of a few centimeters can cause shifts in the calibration.

Ruthenium oxide is a generic name for a class of bismuth ruthenate thick-film resistors. They are epoxied to a BeO header, mounted, and sealed in gold-plated copper AA canisters. Unlike other NTC RTDs, Ruthenium oxide resistors are interchangeable and follow a standard curve. They can be used to below 50 mK and up to 40 K. Their sensitivity is negligible for T > 40 K.

For NTC RTD temperature sensors, up to 70% of the thermal connection to the sensor is through the leads. The large resistance change coupled with thermal considerations results in a requirement for a variable current source for measurement in which the current must be varied over several orders of magnitude (i.e., from about <0.01  $\mu$ A to 1 mA or above) as well as a voltmeter capable of measuring voltages near 1 mV.

#### **Capacitors**

Capacitors are also used for low temperatures, but usually not for temperature measurement. Capacitance temperature sensors have the advantage of being insensitive to magnetic fields, but they commonly experience calibration shifts after thermal cycling and the  ${\rm SrTiO}_3$  capacitors have been known to drift over time while at low temperatures. Phase shifts in the ferroelectric materials are probably the cause of the thermal cycling shifts. The time response of capacitance sensors is usually limited by the physical size and low thermal diffusivity of the dielectric material. The capacitance is measured by an AC technique.



**Appendices** 

#### **Thermocouples**

Thermocouples are only useful where low mass or differential temperature measurements are the main consideration. They must be calibrated in-situ because the entire length of the wire contributes to the output voltage if it traverses a temperature gradient. Errors of 5 K to 10 K can easily occur.

#### Sensor selection

The most important question to ask when selecting a temperature sensor and instrumentation system is "What needs to be measured?" A simple question, but it can be surprisingly easy to answer incorrectly. Some processes need extremely high resolution over a narrow temperature range. Other systems need only a gross estimate of the temperature but over a very wide range.

Design requirements dictate the choice of temperature sensor and instrumentation. Not all applications require the same choice. Even within an application, different temperature sensors can be required. Selecting the appropriate sensor requires prioritizing the most important design attributes. Some attributes are not exclusive to others: The most stable sensors also have a very slow response rate and can be expensive, while sensors with the highest sensitivity have the smallest range.

Design requirements can be classified into four categories:

**Quality of measurement**—This concerns measurement uncertainty, resolution, repeatability, and stability.

**Experimental design**—This is related to constraints due to the experiment (or cryogenic system). It concerns the physical size of the sensors, temperature range of operation, and power dissipation.

**Environmental constraints**—These are effects due to external conditions such as magnetic fields or ionizing radiation. Other external constraints would be vibration or ultra high vacuum.

**Utility requirements**—These are primarily requirements for cost, ease of use, installation, packaging, and long-term reliability.

#### **Quality measures**

#### **Accuracy versus uncertainty**

The term "accuracy" has been almost universally used in literature when presenting specifications, and is often used interchangeably with uncertainty. However, from a strict metrology viewpoint, a distinction does exist between accuracy and uncertainty. Accuracy refers to the closeness of agreement between the measurement and the true value of the measure quantity. Accuracy is a qualitative concept and should not have numbers associated with it. This can be understood since, in practice, one does not have a priori knowledge of the true value of the measured quantity. What one knows is the measured value and its uncertainty, i.e., the range of values which contain the true value of the measured quantity. The uncertainty is a quantitative result and the number typically presented in specifications.

In any proper measurement, an estimate of the measurement uncertainty should be given with the results of the measurement. There are often many sources that contribute uncertainties in a given measurement, and rigorous mathematical methods exist for combining the individual uncertainties into a total uncertainty for the measurement. Temperature sensors, installation, environment, instrumentation, thermal cycling, and thermal EMFs can all contribute to the measurement uncertainty.

A sensor calibration is a method to assign voltage or resistance measurements to a defined temperature scale (i.e., ITS-90 or PLTS-2000). The level of confidence at which this can be done (measuring voltage or resistance AND transferring those values to a defined temperature) is defined by the uncertainty of the calibration.

The uncertainty of the Lake Shore calibration is only one component in a customer measurement system.

It is possible to degrade this accuracy specification by as much as one or two orders of magnitude with improper installation and/or poor shielding and measurement techniques.

#### Repeatability (of the measurement)

The exact definition of repeatability is the closeness of the agreement between the results of successive measurements of the same measurand carried out under the same conditions of measurement (repeatable conditions).

Repeatability is a measure of how well a sensor repeats its measurement under the same conditions. This is often thought of as measurement performed over a period of time (seconds, minutes, hours) at the same temperature. This property is often called precision or stability of the measurement. This value is primarily an instrumentation specification. The sensors themselves are very stable under successive measurements. The stability of the instrument used to measure the sensor needs to be included.



#### Reproducibility

The definition of reproducibility is the closeness of agreement between the results of the measurements of the same measurand carried out under changed conditions of measurements. Often the changed conditions are thermal cycling or mounting (or unmounting) of the sensors. Temperature sensors are complex combinations of various materials bonded together. Aging, thermal cycling, mechanical shock from handling, etc. all affect the reproducibility.

Lake Shore quantifies the reproducibility under thermal cycling in two manners:

Short-term reproducibility: Changes in response values under repeated, successive cycles from ambient to liquid helium (4.2 K).

*Long-term stability:* Changes in response after 200 thermal shocks in  ${\rm LN_2}$  (77 K). Calibrations are performed prior to and after the thermal cycles.

Actual long-term stability for a specific sensor depends on the treatment of the sensor in terms of handling and thermal cycling. A single mechanical shock can cause an immediate calibration shift.

Users should include the short-term reproducibility value in their total uncertainty estimates.

#### Sensitivity and resolution

Sensitivity can be presented in a variety of ways. Typically, it is given in terms of the signal sensitivity, which is the change in a measured parameter per change in temperature ( $\Omega$ /K or V/K). These sensitivities can be a very strong function of temperature. Diodes have sensitivities that range from 2 mV/K to 180 mV/K. Resistor sensitivities can range from less than 0.001  $\Omega$ /K to 1,000,000  $\Omega$ /K, depending upon the device type and temperature.

For resistors, the above signal sensitivity (dR/dT) is geometry dependent (i.e., dR/dT scales directly with R), consequently, very often this sensitivity is normalized by dividing by the measured resistance to give a sensitivity,  $S_{\tau}$ , in change per kelvin

$$S_{T} = (1/R)(dR/dT),$$
 Eqn. 1

where T is the temperature in kelvin and R is the resistance in ohms. This is a common method to express the sensitivity of metal resistors like platinum RTDs.

When comparing different resistance sensors, another useful materials parameter to consider is the dimensionless sensitivity. The dimensionless sensitivity  $\mathbf{S}_{\mathrm{D}}$  for a resistor is a material-specific parameter given by

$$S_D = (T/R)^*(dR/dT) = d(InR)/d(InT)$$
 Eqn. 2

Equivalent definitions are made for diodes with resistance replaced by forward voltage and for capacitors with resistance replaced by capacitance.  $\mathbf{S}_{_{\!D}}$  is also the slope of the resistance versus temperature on a log-log plot, normally used to illustrate resistance versus temperature for negative temperature resistance sensors since their resistance varies by many orders of magnitude.  $\mathbf{S}_{_{\!D}}$  ranges from 0.2 to 6 for most common cryogenic temperature sensors, depending on temperature and sensor type.

Temperature resolution is the smallest temperature difference that can be determined by your measurement system and sensor choice. It is a combination of sensor sensitivity and instrument resolution ( $\Delta R$ ). It can be expressed as

$$\Delta T = \Delta R / (dR/dT)$$
 Eqn. 3  
(or as a ratio  $\Delta T/T = (\Delta R/R)/S_n$ )

Instrument manufacturers will either express the resolution of the measurement as fraction of full scale (i.e., 1 part per million) or as an absolute  $\Delta R$  (i.e., 1  $\Omega$  for 10,000  $\Omega$  scale). Do not confuse temperature resolution with display resolution; actual temperature resolution can be greater or less than the digital display resolution.

#### **Experimental design**

#### Range of use

Two factors limit the useful range of a sensor. First, the physical phenomena responsible for the temperature dependence of the property being measured must occur at a measurable level in both absolute signal and sensitivity to temperature change. Second, the materials used in construction of the temperature sensor must be appropriate to the temperature range of use. Materials such as epoxies, solders, and insulators that are very useful at low temperatures can break down at higher temperatures. Exposure to extreme temperatures (either high or low) can induce strains in the sensor due to changes in the packaging materials or in the leads; the resulting strain can cause a shift in the low temperature calibration for that sensor.



#### Physical size, construction, and thermal response times

As a general rule, larger sensors will be more stable, but they may have a longer thermal response time and may not fit into many experimental schemes. This can be somewhat deceptive, however, because the actual thermal response time depends integrally upon the physical construction of the sensor (i.e., the temperature sensing element) and its associated packaging. Strain-free mounting of sensor elements inside the package necessarily makes for poor thermal connection and longer thermal response times. The choice of package materials can also have a great effect on thermal response times at low temperatures.

Thermal response times are determined by physical construction material and mass of the temperature-sensing element. Strain-free mounted sensors tend to have longer thermal response times. Diode sensors that are mounted directly on a sapphire substrate will be in very good thermal contact with the surroundings and hence have short thermal response times. Thermal response times for various sensors are given in Table 1. The values listed are the 1/e response times.

**Table 1—Thermal response times** 

	401/	77 1/	070 1/
	4.2 K	77 K	273 K
DT-470-SD	<10 ms	100 ms	200 ms
DT-420	<10 ms	50 ms	NA
CX-XXXX-BC	1.5 ms	50 ms	135 ms
CX-XXXX-SD	15 ms	250 ms	0.8 s
CX-XXXX-AA	0.4 s	1 s	1 s
GR-200A-1000	200 ms	3 s	NA
CGR-1-1000	1 s	1.5 s	NA
PT-102	NA	1.75 s	12.5 s
PT-111	NA	2.5 s	20 s
TG-120-PL	100 ms	250 ms	3 s
RF-100-AA	0.8 s	3.6 s	14.5 s
RF-100-BC	2 ms	12 ms	35 ms

#### **Power dissipation**

Diode, resistance, and capacitance temperature sensors must all be energized electrically to generate a signal for measurement. The power dissipated within the temperature sensor must be appropriate for the temperature being measured; the joule heating within the temperature sensor causes an incremental temperature rise within the sensor element itself (self-heating). Consequently, this temperature rise must be kept negligible compared to the temperature of interest.

For diodes, a fixed excitation current of 10 µA is a compromise between power dissipation and noise immunity. The power dissipated is the product of voltage times current. Since the voltage increases with decreasing temperature, power also increases, resulting in a practical lower temperature limit for diode thermometers of slightly above 1 K.

Resistors, on the other hand, have a linear I-V relationship that allows (at a fixed temperature) the measurement of resistance at many different currents and voltages. Since positive temperature coefficient resistance temperature sensors vary relatively linearly with temperature, they can normally be measured by utilizing a fixed current chosen such that self-heating over the useful temperature range is minimized.

In the case of negative temperature coefficient resistance temperature sensors such as Cernox® or germanium RTDs, resistance can vary by as much as five orders of magnitude. To keep the joule heating low, their resistance must be measured either at a fixed voltage or with a variable current selected to keep the resulting measured voltage between 1 mV and 15 mV.

Table 2 gives some typical values of appropriate power levels to use with various temperature sensors in various ranges. These power dissipation levels should keep the temperature rise below 1 mK.

Table 2—Power (W)

	Cernox®, carbon-glass, germanium, Rox™	Platinum, rhodium-iron
0.02 K	10 <sup>-14</sup>	_
0.1 K	10 <sup>-10</sup>	_
1 K	10 <sup>-9</sup>	10-7*
2 K to 10 K	10 <sup>-8</sup>	10 <sup>-6</sup> *
10 K to 100 K	10 <sup>-7</sup> to 10 <sup>-6</sup>	10-5
273 K	10 <sup>-6</sup> (CGR, CX)	3 × 10 <sup>-5</sup>

<sup>\*</sup> Rhodium-iron only

#### **Environmental**

#### **Usefulness in magnetic fields**

Probably the most common harsh environment that temperature sensors are exposed to is a magnetic field. Magnetic fields cause reversible calibration shifts, which yield false temperature measurements. The shift is not permanent and sensors will return to their zero-field calibration when the field is removed.

The usefulness of resistance temperature sensors in magnetic fields depends entirely on the particular resistance temperature detector (RTD) chosen. The Lake Shore Cernox® thin-film resistance sensors are the recommended choice for use in magnetic fields. The Cernox® sensors are offered in a variety of packages and have a wider temperature range than carbon-glass. Ruthenium oxide RTDs are a good choice for temperature below 1 K and down to 50 mK or lower. Due to their strong magnetoresistance and associated orientation effect, germanium sensors are of little use in magnetic fields.



Depending on the desired accuracy, silicon diodes can be used effectively in certain temperature ranges (<0.5% error above 60 K in 1 T fields). However, special care must be taken in mounting the diode to ensure that the junction is perpendicular to field, i.e., current flow is parallel to the magnetic field. Diodes are strongly orientation dependent.

Capacitors are excellent for use in magnetic field environments as control sensors. They can be used in conjunction with another type of sensor (Cernox®, carbon-glass, germanium, etc.) to control temperature. The temperature is set using the other sensor before the field is turned on. Control is then accomplished with the capacitor. Table 3 (page 162) shows magnetic field dependence for some Lake Shore sensors.

#### **Usefulness in radiation**

lonizing radiation refers to a broad class of energetic particles and waves. The effects of radiation can produce temporary or permanent calibration shifts. The exposure can be measured using standard dosimetry techniques, but the actual absorbed dose will vary depending on the material. Due to extensive work performed on the effects of radiation on biological tissue and Si semiconductor devices, the dose is often expressed either in tissue equivalent dose or dose Si, i.e., grays (1 gray = 100 rad).

The data for neutron radiation is more difficult to interpret than gamma radiation data because effects occur due to both the neutrons and the associated background gamma radiation. In both cases it is difficult to calculate or measure the actual absorbed dose. The actual absorbed dose depends on dose rates, energy of the radiation, exposure dose, material being irradiated, etc. Figures 7a to 7e (pages 163 to 164) show data for various sensors.

#### Usefulness in ultra high vacuum systems

The bakeout procedure performed in most ultra high vacuum systems can be damaging to the materials used in the construction of a temperature sensor. Even if the sensor withstands the high bakeout temperature, the sensor's calibration may shift. Without the bakeout, (and possibly with it) some materials in the sensor (Stycast®, for example) may interfere with the high vacuum by acting as a virtual leak. There can be a considerable outgassing from various types of epoxies and ceramics, and some of these materials would not survive the high temperature bake. With proper packaging, diodes and Cernox® can be easily used in ultra high vacuum systems that require a high temperature bake out.

Specific factors to be aware of in an ultra high vacuum environment are:

- Check the compatibility of construction materials of the sensor with ultra high vacuum before using it in such an environment. This includes thermal grease, epoxies, and solders (e.g., Apiezon® N grease cannot be used in these systems due to vapor pressure).
- Solders may not be compatible. Welding may be required.
- Typical insulation used for cryogenic wire may be incompatible with high temperature bakeouts and ultra high vacuums due to thermal ratings and outgassing.

The Lake Shore SD package for diodes is considered UHV compatible. A special package exists for the Cernox® sensor that uses spot welded platinum leads.

A useful website with more information on outgassing properties of materials is found at <a href="http://outgassing.nasa.gov">http://outgassing.nasa.gov</a>.

#### Vibration (shock) environments

Subjecting a temperature sensor to vibrations can permanently shift the calibration, either slowly or catastrophically. Sensors such as germanium and carbon-glass are mounted in a strain-free manner, and mechanical shocks due to vibration will have the same effect on the sensor as dropping it. Other sensors including Cernox® and silicon diodes, due to their physical construction and packaging are less susceptible to vibration-induced errors.

#### Flight qualified

For special applications, Lake Shore will test and qualify sensors to flight standards. Silicon diode and Cernox® sensors, due to their characteristics, performance, construction, and packaging are ideally suited for many flight and large projects applications. Tests are performed to the required standards (for example MIL-STD-750 or MIL-STD-883). Some tests include burn-in lifetime tests, thermal shock, vibration, PIND, gross and fine leak (hermeticity), x-ray, and long and short-term stability.



Table 3—Typical magnetic field-dependent temperature errors ΔT/T (%) at B (magnetic induction)


		Magnetic flux density B					
	T(K)	2.5 T	8 T	14 T	19 T	Notes	
Cernox® 1050 (CX series)	2	1.3	3.1	3.9	5	Best sensor for use in magnetic field	
(	4.2	0.1	-0.15	-0.85	-0.8	(T > 1  K)	
	10	0.04	-0.4	-1.1	-1.5		
	20	0.04	0.02	-0.16	-0.2		
	30	0.01	0.04	0.06	0.11	-	
	77	0.002	0.022	0.062	0.11	-	
	300	0.003	0.004	0.004	0.006	-	
Carbon-glass resistors (CGR series)	4.2	-0.5	-2.3	-4.9	-6.6		
	10	-0.2	-1.1	-2.6	-3.8	_	
	25	0.02	0.22	0.54	0.79	_	
	45	0.07	0.48	1.32	2.2		
	88	0.05	0.45	1.32	2.3		
	306	< 0.01	0.22	0.62	1.1		
Rox™ 102A	2	-1.4	-7.9	-13	-17	Recommended for use over the 0.05 K to 40 K temperature range.	
	3	-1.5	-7	-14	-18	Consistent behavior between devices in magnetic fields.	
	4	-0.56	-6.7	-14	-18		
	8	-1.3	-6.1	-13	-21		
	16	-0.40	-3.4	-9.6	-16		
	23	-0.31	-2.2	-6.2	-11		
Rox™ 103A	2	0.58	1.5	2.2	2.6	Excellent for use in magnetic fields from 1.4 K to 40 K.	
	3	0.44	1.1	1.7	2.0	Predictable behavior.	
	4	0.27	0.95	1.4	1.7		
	8	0.11	0.49	0.71	0.80		
	16	0.018	0.076	0.089	0.040	-	
	23	0.0051	0.0058	-0.0060	0.095	-	
Rox™ 202A	2	-0.13	-2.2	-3.9	-5.2	Recommended for use over the 0.05 K to 40 K temperature range.	
	3	0.18	-0.68	-2.7	-3.7	Consistent behavior between devices in magnetic fields.	
	4	0.77	0.046	-1.8	-3.2		
	8	-0.023	0.16	-0.65	-3.0		
	16	0.03	0.16	-0.48	-1.5	_	
	23	-0.05	-0.08	-0.39	-0.92		
Platinum resistors	20	20	100	250	_	Recommended for use when $T \ge 40 \text{ K}$ .	
(PT Series)	40	0.5	3	6	8.8		
	87	0 .04	0.4	1	1.7		
	300	< 0.01	0.02	0.07	0.13		
Rhodium-iron	4.2	11	40	_		Not recommended for use below	
(RF Series)	40	1.5	12	30	47	77 K in magnetic fields.	
	87	0.2	1.5	4	6		
	300	< 0.01	0.1	0.4			
Capacitance CS-501 Series		ΔΤ/ ΔΤ/T(%)	Γ(%) <0.015 a <0.05 at 77 K	t 4.2 K and 18.7 and 305 K and	' T 18.7 T	Recommended for control purposes. Monotonic in C vs. T to nearly room temperature.	
Germanium resistors	2.0	-8	-60	_	_	Not recommended except at low B owing to large,	
(GR Series)	4.2	-5 to -20	-30 to -55	-60 to -75	_	orientation-dependent temperature effect.	
,	10	-4 to -15	-25 to -60	-60 to -75	_	· ·	
	20	-3 to -20	-15 to -35	-50 to -80	_	-	
Chromel-AuFe (0.07%)	10	3	20	30		Data taken with entire thermocouple in field,	
( , , , , , , , , , , , , , , , , , , ,	45	1	5	7	_	cold junction at 4.2 K; errors in hot junction.	
	100	0.1	0.8		_	• •	
		1	3	7		Useful when T $\geq$ 10 K.	
Type E thermocouples	10						
Type E thermocouples (chromel-constantan)	20	<1	2	4	_	Refer to notes for Chromel-AuFe (0.07%).	


	T(K)	1 T	2 T	3 T	4 T	5 T	Notes
Silicon diodes	4.2	-200	-300	-350	-400	-500	Strongly orientation dependent.
Junction parallel to field	20	-10	-20	-25	-30	-40	
(DT Series)	40	-4	-6	-8	-10	-12	
	60	-0.5	-1	-2	-3	-3.5	
	80	<0.1	-0.5	-0.8	-1.1	-1.5	
	300	<0.1	<-0.1	<-0.1	<-0.1	<-0.1	
Silicon diodes	4.2	-8	-9	-11	-15	-20	Strongly orientation dependent.
Junction perpendicular to field	20	-4	-5	-5	-5	-10	
(DT Series)	40	-1.5	-3	-4	-5	-5.5	
	60	-0.5	-1	-2	-3	-3.5	
	80	-0.1	-0.3	-0.5	-0.6	-0.7	
	300	< 0.1	0.2	0.5	0.6	0.6	
GaAlAs diodes	4.2	2.9	3.8	3.7	2.8	1	Shown with junction perpendicular (package base
(TG Series)	30	0.2	0.2	0.3	0.3	0.2	parallel) to applied field B. When junction is parallel
	78	< 0.1	<0.1	0.17	0.16	0.1	to B, induced errors are typically less than or on the
	300	-0.1	<0.1	<0.1	<0.1	<0.1	order of those shown.



#### Figure 7a—Gamma rays

Temperature shift as a function of temperature due to 10,000 Gy gamma radiation dose from a Cs-137 source. Dose rate was 0.5 Gy/min with irradiation performed at 298 K.

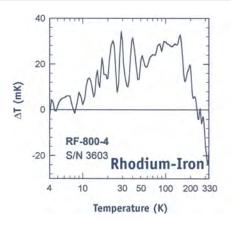

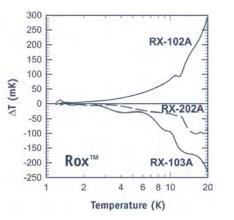

2000

1500


∆T (mK)


TG-120PL


S/N 8554

#### Figure 7b—Neutrons and gamma rays

Temperature shift as a function of temperature due to a  $2.5 \times 10^{12}$  neutron/cm² fluence from a nuclear pool reactor. The neutron flux was  $3.75 \times 10^7$  neutron/cm²/s with irradiation performed at 298 K (associated gamma ray dose of 29 Gy).

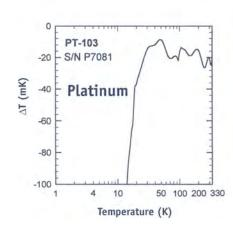


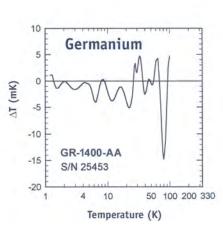

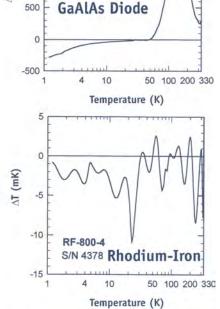


#### Figure 7d—Gamma rays


Temperature shift as a function of temperature due to 10,000 Gy gamma radiation dose from a Co-60 source. Dose rate was 40 Gy/min with irradiation performed at 4.2 K.

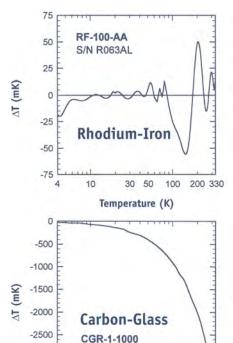



# Figure 7c—Neutrons and gamma rays

Temperature shift as a function of temperature due to a  $10^{14}$  neutron/cm<sup>2</sup> fluence from a nuclear pool reactor. The neutron flux was  $2 \times 10^{12}$  neutron/cm<sup>2</sup>/s with irradiation performed at 298 K (associated gamma ray dose of 116 Gy).

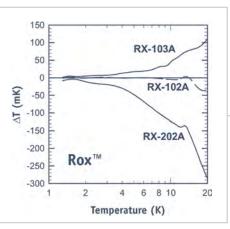

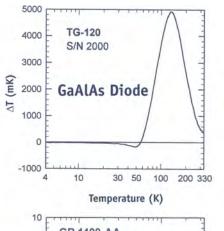

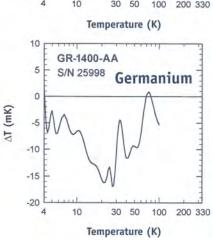

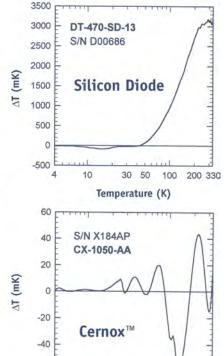

-60

#### Figure 7e—Neutrons and gamma rays

S/N C13672

30 50


Temperature (K)


100

200 330

-3000 -3500

Temperature shift as a function of temperature due to a  $2 \times 10^{12}$  neutron/cm² fluence from a nuclear pool reactor. The neutron flux was  $7.5 \times 10^7$  neutron/cm²/s with irradiation performed at 4.2 K (associated gamma ray dose of 23 Gy).

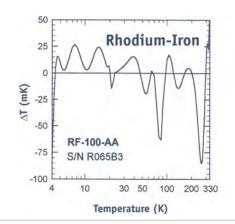

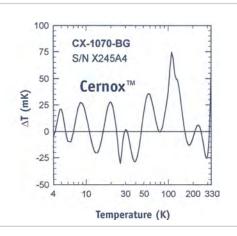

30 50

Temperature (K)

100

200 330



#### Utility

#### Interchangeability

It is very convenient and cost effective to have temperature sensors that match a standard curve, thus not requiring individual calibration. Such sensors are termed "interchangeable." In industry, interchangeable sensors make equipment design and manufacture simpler. Any monitoring equipment for those sensors can be identical. Time is saved in research settings since new calibrations do not have to be programmed into control and data acquisition equipment each time a new sensor is installed.

Some cryogenic temperature sensors exist at present which are interchangeable within a given tolerance band. Silicon diodes from Lake Shore are interchangeable. Series DT-670 diodes conform closely to a curve that Lake Shore calls Curve 670. The conformance is indicated by placing the diodes within tolerance bands. These sensors can be ordered by simply specifying a tolerance band. In this case, individual calibrations are not performed. If the greater accuracy is required, a calibration is necessary. Calibration can decrease the uncertainty by a factor of 10 or more. The DT-470 also follows a unique standard curve and is interchangeable with other DT-470s.

In addition to silicon diodes, platinum and ruthenium oxide RTDs both follow standard curves. Platinum RTDs match an industry standard curve (IEC 751) in terms of resistance versus temperature. Industrial platinum resistance temperature sensors are broken into Class B tolerances and Class A tolerances. Lake Shore offers only Class B sensors.

Ruthenium oxide RTD sensors also follow a standard curve. Like silicon diodes, this curve is unique to each manufacturer.

#### Signal size

For resistors, values lie between approximately 10  $\Omega$  and 100,000  $\Omega$ . Resistance measurements outside this range become more difficult to perform, especially at ultra-low temperatures. Keep in mind that for carbon-glass, Cernox®, and germanium sensors, there are several resistance ranges available to suit the appropriate temperature range(s). Because of their rapidly changing resistance and use at ultra-low temperature, it is necessary to use a small excitation current. The resulting voltage measurement can be in the nanovolt range in some cases. At these low voltages a variety of noise sources begin to affect the measurement.

Diode temperature sensors have a relatively large output (about 1 V) and a fixed current excitation of 10  $\mu$ A. This allows for simple instrumentation compared to NTC RTDs like Cernox®.

#### **Packaging**

Sensors come in various packages and configurations. Apart from the size considerations discussed previously, there are practical considerations as well. A cylindrical package is obviously better suited for insertion into a cylindrical cavity than a flat or square-shaped package. Lake Shore offers a variety of sensor packages and mounting adaptors as well as probe assemblies. The most common package is the SD package. It is a robust and reliable hermetically sealed flat package. With a metallized and insulated bottom, the SD package can be indium soldered to the experimental surface. It can also be mechanically clamped as well as varnished or epoxied. The SD package can also be mounted into adaptor packages like the CU bobbin.

Many RTDs like germanium and Cernox® are mounted in cylindrical AA canisters. This is a requirement for GRTs due to their strain-free mounting. Cernox® is also available in a SD package.

Many cryogenic sensors can be packaged into custom probes and thermowells. Lake Shore has many standard probe configurations and can manufacture special customer designed probes for various applications.



# Appendix C: Sensor Packaging and Installation

#### Installation

Once you have selected a sensor and it has been calibrated by Lake Shore, some potential difficulties in obtaining accurate temperature measurements are still ahead. The proper installation of a cryogenic temperature sensor can be a difficult task. The sensor must be mounted in such a way so as to measure the temperature of the object accurately without interfering with the experiment. If improperly installed, the temperature measured by the sensor may have little relation to the actual temperature of the object being measured.

Figure 1—Typical sensor installation on a mechanical refrigerator

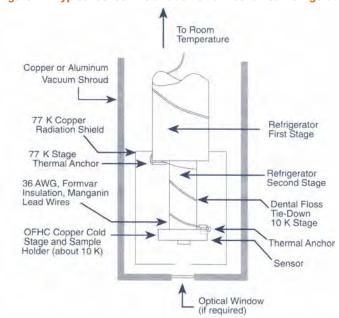


Figure 1 shows a typical sensor installation on a mechanical refrigerator. Note the additional length of lead wire wrapped around the refrigerator stages to minimize thermal conductance along the leads. If the optical radiation load through the window is large, the sample temperature will not necessarily be the same as that of the sensor in the block. A sensor placed in more intimate contact with the sample may be required.

#### General installation considerations

Even with a properly installed temperature sensor, poor thermal design of the overall apparatus can produce measurement errors.

#### **Temperature gradients**

Most temperature measurements are made on the assumption that the area of interest is isothermal. In many setups this may not be the case. The positions of all system elements—the sample, sensor(s), and the temperature sources—must be carefully examined to determine the expected heat flow patterns in the system. Any heat flow between the sample and sensor, for example, will create an unwanted temperature gradient. System elements should be positioned to avoid this problem.

#### **Optical source radiation**

An often overlooked source of heat flow is simple thermal or blackbody radiation. Neither the sensor nor the sample should be in the line of sight of any surface that is at a significantly different temperature. This error source is commonly eliminated by installing a radiation shield around the sample and sensor, either by wrapping super-insulation (aluminized Mylar®) around the area, or through the installation of a temperature-controlled aluminum or copper shield (see Figure 1).

#### 2-lead versus 4-lead measurement

4-lead measurements are recommended for all sensors. 2-lead measurements can be performed with diode sensors with a small increase in uncertainty. Refer to Appendix E: Temperature Measurement, for a detailed discussion.



#### **High temperature effects**

Below room temperature, the primary effect of using dissimilar materials bonded together in sensing elements or packages is stress-induced by different expansion coefficients. Above room temperature, additional problems can occur. Alloying, diffusion (Kirkendahl voids), chemical reactions, and corrosion (especially in the presence of moisture and chlorine) accelerate as the temperature increases. These factors can cause catastrophic failure with time, or a shift in the sensor calibration. Completely accurate de-rating data for all situations that could be encountered is impossible to compile.

#### **Conduction (lead attachment)**

Another source of heat flow that is often neglected is conduction through the electrical leads that run between the sensor and the ambient environment. 32- or 36-gauge, low thermal conductivity wire such as phosphor bronze or manganin is used to alleviate this problem. These leads must also be thermally anchored at several successive temperature points between ambient temperature and the sensor. Performing a 4-lead measurement will overcome the high lead resistance.

The physical mounting of the leads of a sensor is as important as the mounting of the sensor itself. Thermal contact to the active element in a cryogenic sensor occurs both through the sensor body and the electrical leads. In fact, for some sensors (e.g., germanium resistance thermometers) the primary thermal contact is through the leads. For accurate temperature readings, the sensor and its leads must be anchored so they are at the same temperature as the sample being measured. Table 1 shows typical heat sinking lengths.

Table 1—Wire heat-sinking length required to thermally anchor to a heat sink at temperature T to bring the temperature of the wire to within 1 mK of  $T_{lower}$ 

			Heat-sinking length (mm) for wire sizes					
	Tupper (K)	Tlower (K)	0.21 mm² (24 AWG)	0.032 mm <sup>2</sup> (32 AWG)	0.013 mm <sup>2</sup> (36 AWG)	0.005 mm <sup>2</sup> (40 AWG)		
Copper	300	80	160	57	33	19		
	300	4	688	233	138	80		
Phosphor bronze	300	80	32	11	6	4		
	300	4	38	13	7	4		
Manganin	300	80	21	4	4	2		
	300	4	20	7	4	2		
304 SS	300	80	17	6	3	2		
	300	4	14	5	3	2		

Note: values are calculated assuming wires are in a vacuum environment, and the thermal conductivity of the adhesive is given by the fit to the thermal conductivity of VGE-7031 varnish

There are a number of ways in which sensor leads can be properly anchored, with the choice usually determined by the needs and constraints of the particular application. Longer leads may be wound directly around a sensor adaptor or another anchor adjacent to the sample and varnished into place. The varnish serves two purposes: it physically holds the leads in place, and it increases the contact surface area between the wire and the sample, or sample holder. VGE-7031 varnish is widely used as a low-temperature adhesive and can be easily removed with methanol. As long as the leads are electrically insulated with an enamel-type coating, such as Formvar® (see caution note) or polyimide, the varnished-down leads provide a suitable thermal anchor (thermal short) to their surroundings. Leads with heavy insulation, such as Teflon®, minimize the potential for making a thermal short to the surroundings, resulting in more thermal conduction down the leads into the sensing element. Resulting temperature measurement errors can be significant.

#### TIP: maintain electrical isolation

To maintain good electrical isolation over many thermal cycles, a single layer of cigarette paper can be varnished to the thermal anchor first, and the wire then wound over the paper and varnished down. The actual sensor leads are then soldered to this thermally anchored lead wire after the sensor body is mounted. For a more permanent installation, replace the VGE-7031 varnish with a suitable epoxy such as Stycast® 2850-FT.

# Caution: varnish can cause crazing of Formvar<sup>®</sup> insulation.

One can make a separate thermal anchor to which the thermometer leads are attached. A typical technique for producing a physically compact anchor uses small gauge wire (32 AWG) insulated with Formvar®, polyimide, or a similar coating. The wire is wound around the sample in a bifilar manner or onto a separate bobbin and bonded with varnish. For most applications, a bonded length of 5 cm to 10 cm provides a sufficient thermal anchor unless poor practices elsewhere in the system permit excessive heat leaks down the leads. Copper wire may require several meters for heat sinking.



#### What you may need:

#### Wire

- ✓ Phosphor bronze
- ✓ Manganin
- ✓ Nichrome
- √ Conner
- Constantan
- ✓ Stainless steel coaxial cable

#### **Solders**

- 60/40 lead/tin
- √ 90/10 lead/tin
- ✓ Silver
- ✓ Ostalloy® 158 (Wood's metal)
- Indium-silver
- ✓ Indium

#### **Fluxes**

RMA Keep Clean flux Stay Clean flux Stay Silv® flux

#### **Insulating materials**

Ceramics Masking tape

- ✓ Polyester tape
  Kapton® films
  Teflon® tape
  Heat shrink tubing
  G-10
  Mylar® (polyester film)
  Fiberglass sleeving
- ✓ Epoxies
- ✓ VGE-7031 varnish
- ✓ Stycast® 2850 FT epoxy
- Cigarette paper

  Greases (Apiezon® N & H)

#### **Conducting materials**

- ✓ Silver filled epoxy Silver conductive paint
- ✓ Indium foil

#### **Fasteners**

- Dental floss
- Clamps Screws/bolts
- ✓ VGE-7031 varnish
- √ Stycast® 2850 FT epoxy

#### **Heat sinking**

- ✓ Copper bobbins
- ✓ Metallized ceramic chips

#### Other accessories

- ✓ Vacuum feedthroughs
- ✓ Cartridge heaters
- Lake Shore stocks these accessories as a convenience to our customers

#### Cryogenic accessories for installation

#### **Cryogenic wire**

Cryogenic wire is different from normal wire due to its low thermal conductivity and high electrical resistivity. The most common types of cryogenic wire are phosphor bronze and manganin. Phosphor bronze is a nonmagnetic copper alloy. Manganin wire has a lower thermal conductivity (a factor of about ½) and higher resistivity compared to phosphor bronze wire. Both are readily available in small gauges ranging from 32 to 42 AWG. Either polyimide or polyvinyl formal (Formvar®) is used to insulate the wires. The polyimide is a resin with a 220 °C thermal rating. It has exceptional resistance to chemical solvents and toxic heat. It also is unaffected by exposure to varnish solvent. Formvar® is a vinyl acetate resin rated at 105 °C. It has excellent mechanical properties such as abrasion resistance and flexibility. The film will withstand excessive elongation without rupture when stressed during winding. Formvar® has a tendency to craze upon contact with solvents such as toluol, naptha, and xylol. It should be given an annealing preheat prior to varnish application. The Formvar® insulation can be removed mechanically or chemically during terminal preparation.

Phosphor bronze wire is readily available in multifilar form with 2 or 4 wires. In bifilar form, the wires are twisted to minimize noise pickup. In quadfilar form, the wires are either straight or 2 twisted pairs twisted together. The latter form is most useful for standard 4-lead measurements. The wires are bonded together for ease in heat sinking while the twisting helps minimize noise pickup. Straight Quad-Lead™ wire can be bonded together with the help of VGE-7031 varnish. The bonding agent is soluble in alcohol.

Other types of common cryogenic wires include nichrome wire, which has a very high electrical resistivity making it excellent for heater windings. Ultra miniature flexible coaxial cables with 304 stainless steel or copper conductors are available for providing shielded leads when necessary. For low resistance, heavy duty lead wires and multifilar silver-plated twisted copper wire are available. Constantan wire is another copper alloy having just a little more copper content than manganin. As such, its resistivity is a little lower, while its thermal conductivity is a little higher. Evanohm® wire is a very high resistivity wire (about 5 times the resistivity of nichrome) with very small temperature dependence. This wire is also excellent for heater windings.

#### TIP: Making your own ribbon cable—ease of handling

Two nails should be hammered into a piece of wood at a distance of just over half the needed lead length. The wire is wrapped continuously from one nail to the other. With a rubber or plastic glove, apply a thin coating of VGE-7031 varnish along the entire length of the wires and allow to dry. Then the cable can be cut for full length. (Remember that the solvents in VGE-7031 varnish will attack Formvar® insulation.)



#### **Solders and fluxes**

The most common electrical connections are solder joints. There are a number of solder compounds available such as 60/40 tin/lead, silver, Wood's metal, cadmium/tin, and indium. They have varying melting points, and the melting points sometimes determine the upper temperature limit for a sensor. Care should be taken when using these solders, as the fumes are toxic. Also, many of these solders become superconducting at lower temperatures. The transition temperature should be checked if this could affect your experiment. (Read on to the fasteners section for more comments on solders.)

There are a number of fluxes that are used with these solders. Rosin Mildly Activated (RMA) soldering flux is an electronic grade rosin flux typically used for soldering wires to temperature sensors. Keep Clean flux is a mild acid flux used when RMA flux is not effective. It is strong enough to clean the oxidation off the surface and the solder to promote a good joint. It is very useful in situations where joints are repeatedly made and broken. Stay Clean flux is a corrosive acid flux used when neither of the above are useful. It is commonly used with stainless steel and platinum. Due to its highly corrosive nature, it must be cleaned off with methanol or water or it will continue to corrode the material. Stay Silv® flux is a high temperature flux for use with high temperature solders such as silver solder. It is not useful on aluminum, magnesium, or titanium.

It is often difficult to make electrical connection to many of the materials used for electrical leads in cryogenic applications. These lead materials include Kovar, copper, gold, phosphor bronze, manganin, constantan, platinum, stainless steel, and nichrome. Soldering these materials can be problematic. The small diameter wire complicates the problem by making it difficult to heat the wire uniformly, allowing the solder to flow. Choosing a proper flux and solder for the wire is crucial to making a reliable electrical connection with minimal effort.

Most of the sensors shipped by Lake Shore have undergone testing to ensure proper operation. Their electrical leads have been tinned. For these sensors, a standard electronic grade RMA flux is appropriate. This flux is also appropriate for Kovar, gold, and copper leads that have not been tinned. For other wire types, a more corrosive acid flux is needed. Stay Clean flux is recommended for untinned wire consisting of constantan, manganin, phosphor bronze, platinum, nichrome, or stainless steel.

Note: Care must be taken to thoroughly clean the residual Stay Clean flux off with water or methanol after use to prevent further corrosion.

Typically, standard 60/40 Sn/Pb solder can be used for applications ranging from 0.05 K to 350 K (liquidus point of 461 K and solidus point 456 K). This solder can be used with any of the above material types after tinning. If the application requires a higher temperature, then use 90/10 Pb/Sn solder (liquidus point of 575 K and solidus point 548 K). For very high temperatures up to 800 K, use Stay Silv® flux with cadmium-free silver solder (liquidus point of 922 K; solidus point of 891 K).

#### **Insulating materials**

When installing electrical leads at low temperatures, it is important to know what insulation materials can be used. Insulating materials that work well at cryogenic temperatures include ceramics, temporary masking tape, polyester film tape, Kapton® film, Teflon® tape and tubing, G-10, Mylar®, epoxies, varnishes, cigarette paper (used under VGE-7031 varnish), and greases.

The most common varnish for cryogenic work is VGE-7031 varnish. It has good chemical resistance, bonds to a variety of materials and has a fast tack time. It may be air-dried or baked. VGE-7031 varnish is compatible with cotton, Dacron® polyester fiber, nylon, glass tapes, laminates, Mylar® polyester film, mica products, polyester products, vinyl products, wire enamels, paints, rayon, plastics, and metals. The solvents in VGE-7031 varnish will attack Formvar® insulation, causing it to craze, but in most cases this will not be a problem after drying thoroughly.

Stycast® 2850FT and GT are composed of a black epoxy resin, filled with silica powder to give them a lower thermal expansion coefficient. The FT is roughly matched to copper, while the GT is roughly matched to brass. The result is a material that is very strong, adheres well to metals, and tolerates brief exposure up to 200 °C for soldering. The drawbacks are that it is essentially unmachinable, has a non-negligible magnetic susceptibility and a temperature-dependent dielectric constant at low temperatures, and is somewhat permeable to helium at room temperature.

Another useful insulator is Kapton® tape. It is a polyimide tape with a thin coating of Teflon® FEP on either or both sides of the film to provide adhesion. The principal advantages of this severed tape insulation is its uniform, pinhole free covering and thermal stability for continuous use up to 240 °C. It has exceptional cut-through resistance under extreme temperature and pressure conditions. This Kapton® insulation offers excellent moisture protection and, because it is smooth and thin, has a space advantage over glass, Dacron® glass, paper, and fiber-over-film constructions. It is compatible with all standard varnishes, and is highly resistant to solvent attack.



#### **Conducting materials**

Sometimes it is desired to make electrical contact between materials. The solders previously mentioned are electrically conducting, as are certain epoxies (silver-filled) and silver conductive paint.

#### **Fasteners**

A variety of materials are suitable for fastening sensors at low temperatures. These include dental floss (Dacron® fiber), screws, bolts, pins, springs, tape, pastes, solders, epoxies, and varnishes. You must consider coefficients of linear expansion when deciding upon a mounting scheme. If linear expansion coefficients are too mismatched, mountings will simply come loose, or in the worst case, damage the mounting surface or the sensor. Expansion coefficients should never differ by more than a factor of 3 between two materials being bonded together. Greases such as Apiezon™ N grease, H grease, and Cry-Con® grease can be used to increase the surface area of contact between a sensor and the mounting surface. VGE-7031 varnish accomplishes the same purpose, as does Stycast® 2850. Mounting the sensor with Stycast® is more permanent.

If the Stycast® is being used with diodes, the user should be aware that stress on the diode package can cause piezoresistive shifts in the calibration curve. In extreme cases, (e.g., by using hard solder between the SD package and copper), the package can crack. The best joint in almost all cases is made by pure indium, which remains malleable at all temperatures. The exceptions are for service temperatures over 125 °C or where strength is paramount. Indium can also corrode rapidly in the presence of moisture under thermal cycling conditions.

#### TIP: Where to buy flux & solder

RMA flux is available from most electronics supply stores as well as Kester Solder, 515 E. Touhy Avenue, Des Plaines, IL 60018

60/40 Sn/Pb solder is also available from most electronics supply stores both with and without RMA flux

Stay Clean soldering flux, Stay Silv® white brazing flux, and cadmium-free silver solder are available from J. W. Harris Company, Inc., 10930 Deerfield Road, Cincinnati, OH 45242

#### SD package installation

Three aspects of using a cryogenic temperature sensor are critical to its optimum performance. The first involves the proper mounting of the sensor package; the second relates the proper joining of sensor lead wires and connecting wires; the final concern is the thermal anchoring of the lead wires. Although the sequence in which these areas should be addressed is not fixed, all elements covered under each aspect should be adhered to for maximum operating capabilities of the sensor.

#### **Sensor mounting**

- The mounting area should be prepared and cleaned with a solvent such as acetone followed by an isopropyl alcohol rinse. Allow time for the solvents to evaporate before sensor mounting.
- The list below provides brief instructions on mounting a sensor using a number of different methods. The constraints of your application should dictate the most appropriate mounting method to follow.

Mechanical—The preferred method for mechanically mounting an SD sensor is using the Lake Shore spring loaded clamp. This clamp should be ordered at the time the sensor is ordered (-CO suffix on sensor part number). The clamp holds the SD sensor in contact with the surface and also allows the sensor to be changed or replaced easily. A thin layer of Apiezon® N Grease (0.055 mm) or a flat 100% indium preform should be used between the sensor and mounting surface to enhance thermal contact. The spring keeps the sensor from getting crushed.

Indium solder (100% In)—A low wattage heat source should be used, as the sensor temperature must never exceed 200 °C (147 °C for Cernox®). The mounting surface and sensor should be tinned with a rosin flux (RMA is recommended) prior to mounting the sensor. A thin, uniform layer of indium solder should be the goal. Clean both the sensor and mounting surface of residual flux using rosin residue remover. Once the surface area is dry, reheat the mounting surface to the melting point of the solder (156 °C). Press the sensor into position and allow it to warm to the melting point of the solder. Remove the heat source and allow sufficient time for the solder to solidify (typically 2 to 3 seconds) before removing it.

Apiezon® N grease—This is best used as a thermal conductor when the sensor is mounted in a hole or recess, and when the sensor is intended to be removed. The sensor should be surrounded with thermal grease and placed into the mounting position. When the temperature is lowered, the thermal grease will harden, giving good support and thermal contact.

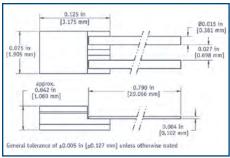


Figure 2 - SD Package

VGE-7031 varnish—Prepare varnish and apply a thin layer on the mounting surface. Press the sensor firmly against the varnish during curing to ensure a thin bond layer and good thermal contact. Varnish will air-dry in 5 to 10 minutes. Sufficient time must be allowed for the solvents in the varnish to evaporate. There is a small probability of ionic shunting across the sensor during the full cure period of the varnish (typically 12 to 24 hours).

**Stycast® 2850FT epoxy**—Prepare epoxy and apply a thin layer on the mounting surface. Press the sensor firmly into the epoxy during curing to assure a thin bond layer and good thermal contact. Epoxy will cure in 12 hours at 25 °C or in 2 hours at 66 °C.

Note: When using an electrically conductive adhesive or solder, it is important that the excess does not "creep-up" the edges of the sensor or come in contact with the sensor leads. There is a thin braze joint around the sides of the SD package that is electrically connected to the sensing element. Contact to the sides with any electrically conductive material will cause a short.

 Follow manufacturer's instructions for adhesive curing schedule. Never heat the sensor above 200 °C (147 °C for Cernox®).

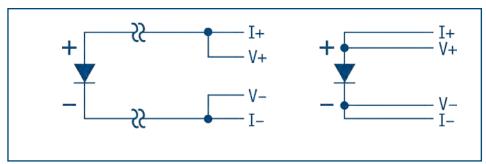


Figure 3 - 2-Lead versus 4-Lead Measurements

#### **Lead attachment**

- Although the SD sensor package (Figure 2) is a 2-lead device, measurements should preferably be made using a 4-wire configuration to avoid uncertainties associated with lead resistance.
  - **2-lead measurement scheme**—The leads used to measure the voltage are also the current carrying leads. The resulting voltage measured at the instrument is the sum of the temperature sensor voltage and the voltage drop across the 2 leads (see Figure 3).
  - **4-lead measurement scheme**—The current is confined to one pair of current leads with the sensor voltage measured across the voltage leads (see Figure 3).
- Lead polarity: for the silicon diode and for the GaAlAs diode, when viewed with the base down (the base is the largest flat surface) and the leads toward the observer, the positive lead (anode) is on the right and the negative (cathode) is on the left. For Cernox® there is no polarity.
- 3. Strip the insulation from the connecting wires by scraping delicately with a razor blade, fine sand paper, or steel wool. Phosphor bronze or manganin wire, in sizes 32 or 36 AWG, is commonly used as the connecting lead wire. These wires have low thermal conductivity and high resistivity, which help minimize the heat flow through the leads. Typical wire insulation is polyvinyl formal (Formvar®) or polyimide (ML). Formvar® insulation has better mechanical properties such as abrasion resistance and flexibility. Polyimide insulation has better resistance to chemical solvents, heat, and radiation.
- 4. Prepare the connecting wire ends with an RMA (rosin mildly active) soldering flux, and tin them with a minimal amount of 60% Sn 40% Pb solder. Use a low wattage soldering iron which does not exceed 200 °C.
- 5. Clean off residual flux with rosin residue remover. The sensor leads can be prepared in an identical manner.
- 6. Join one sensor lead with two of the connector wires. Apply the soldering iron to the connector wire above the joint area until the solders melt, then remove the iron. Repeat for the other set of connector wires and the other sensor lead. Heat sinking the SD sensor with a flat jaw alligator clip is good practice to eliminate heat buildup at the sensor element.
- 7. Avoid putting stress on the device leads, and leave enough slack to allow for the thermal contractions that occur during cooling which could fracture a solder joint or lead. Some epoxies and shrink-tubing can put enough stress on lead wires to break them.



#### **Heat sinking/thermal anchoring**

- Since the area being measured is read through the base of the sensor, heat flow through the connecting leads creates less of an offset between the sensor chip and the true sample temperature than with other types of packages. However, thermal anchoring of the connecting wires is necessary to ensure that the sensor and the leads are at the same temperature as the sample.
- Connecting wires should be thermally anchored at several temperatures between room temperature and cryogenic temperatures to guarantee that heat is not being conducted through the leads to the sensing element. Two different size copper bobbins are available from Lake Shore for heat sinking leads.
- 3. If connecting wires have a thin insulation such as Formvar® or polyimide, a simple thermal anchor can be made by winding the wires around a copper post, bobbin, or other thermal mass, A minimum of 5 wraps around the thermal mass should provide sufficient thermal anchoring, however, additional wraps are recommended for good measure if space permits. To maintain good electrical isolation over many thermal cycles, it is good practice to first varnish a single layer of cigarette paper to the anchored area then wrap the wire around the paper and bond in place with a thin layer of VGE-7031 varnish. Formvar® wiring insulation has a tendency to craze with the application of VGE varnish. If used, the wires cannot be disturbed until the varnish is fully cured and all solvents have evaporated (typically ≥24 hours).
- A final thermal anchor at the sample itself is good practice to ensure thermal equilibrium between the sample and temperature sensor.

#### CU, DI, CY, and CD package installation

Three aspects of using a cryogenic temperature sensor are critical to its optimum performance. The first involves the proper mounting of the sensor package; the second relates the proper joining of sensor lead wires and connecting wires; the final concern is the thermal anchoring of the lead wires. Although the sequence in which these areas should be addressed is not fixed, all elements covered under each aspect should be adhered to for maximum operating capabilities of the sensor.

#### **Sensor mounting**

The CU, DI, and CY packages (Figures 4 and 5) combine a standard SD sensor with a goldplated copper mounting bobbin. The mounting bobbin of these packages each has a hole designed for mounting with a #4-40 screw. The CD package is shown in Figure 6.

- A threaded hole in your mounting surface is necessary for mounting the sensor package.
   The hole in the sensor package will accommodate a #4-40 screw. A brass screw is recommended due to the thermal contractions/expansions of the final assembly.
- The threaded hole and surrounding surface should be cleaned with a solvent such as acetone followed by an isopropyl alcohol rinse. Allow time for the solvents to evaporate before sensor mounting.
- 3. Apply a small amount of Apiezon® N grease to the threads of the screw. To ensure good thermal contact between the sensor and mounting surface, use an indium washer/preform or a thin layer of Apiezon® N grease between the mounting surface and the sensor package. Note: An overabundance of grease will increase the thermal barrier. Keep the thickness to 0.05 mm or less.
- 4. Insert screw through sensor mounting bobbin and tighten screw firmly enough to hold sensor in place. Avoid overtightening (torque of 30 to 50 in·oz [0.2 to 0.35 N⋅m] should be sufficient).

#### **Lead attachment**

The SD sensor has been attached to the mounting bobbin and encapsulated in Stycast® epoxy. The 0.92 m (36 in) Polyimide (ML) insulated sensor leads are 36 AWG phosphor bronze wire which are thermally anchored to the bobbin. Teflon® tubing is used as a strain relief to reinforce the leads at the bobbin assembly. The difference between the CU package and the DI package is the connecting lead configuration. Standard lead configuration for the CU is a 4-lead device [red (I-), green (V-), black/dark blue (V+), clear (I+)] while standard lead configuration for the DI package is a 2-lead device [green = cathode (-), clear = anode (+)].

Figure 4—CU & DI package

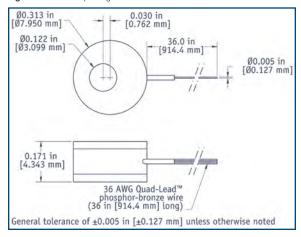


Figure 5—CY package


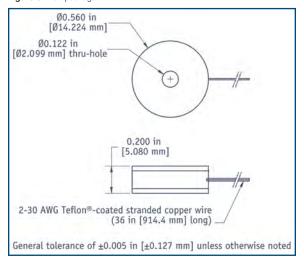


Figure 6—CD package

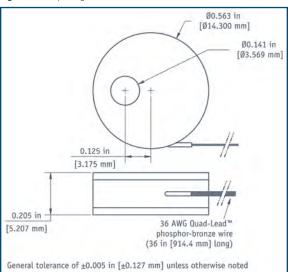

#### DI package—2-lead measurement scheme

The leads used to measure the voltage are also the current carrying leads. The resulting voltage measured at the instrument is the sum of the temperature sensor voltage and the voltage drop within the two current leads (see Figure 3).

#### CU package—4-lead measurement scheme

The current is confined to one pair of current leads with the sensor voltage measured across the voltage leads (see Figure 3).

Thirty-six inches of lead wire is attached during the production process. If additional connection wire is required, use the following instructions:

- Prepare the sensor leads with an RMA (rosin mildly active) soldering flux, and tin them with a minimal amount of 60% Sn 40% Pb solder. Use a low wattage soldering iron that does not exceed 200 °C. Clean off residual flux with rosin residue remover.
- 2. Strip the insulation from the connecting wires by scraping delicately with a razor blade, fine sand paper, or steel wool. (Phosphor bronze or manganin wire, in sizes 32 or 36 AWG, is commonly used as the connecting lead wire. These wires have low thermal conductivity, which help minimize the heat flow through the leads. Typical wire insulation is Formvar® or Polyimide (ML). Formvar® insulation has better mechanical properties such as abrasion resistance and flexibility. Polyimide insulation has better resistance to chemical solvents and burnout.) Follow the same procedure as Step 1 for preparing connecting wires.
- 3. DI package—join one sensor lead with two of the connector wires. Apply the soldering iron above the joint area until the solders melt, then remove the iron immediately. Repeat for the other connecting wires and the other sensor lead. Insulate the joints appropriately.

CU package—identify lead polarities and apply the soldering iron above the joint area until the solders melt, then remove the iron immediately. Leave enough slack to allow for the thermal contractions that occur during cooling, which could fracture a solder joint or lead. Insulating the soldering joint is recommended to prevent shorts. Use heat shrink tubing. Teflon® and Kynar® shrink tubings are more resistant to cracking at low temperatures than polydelefin.



#### **Heat sinking/thermal anchoring**

Depending on the application, sufficient heat sinking of the leads may already exist in the bobbin. Use the following procedure if additional heat sinking is recommended:

- Connecting wires should be thermally anchored at several temperatures between room temperature and cryogenic temperatures to guarantee that heat is not being conducted through the leads to the sensing element.
- 2. A simple thermal anchor can be made by winding the wires around a copper post, bobbin, or other thermal mass. A minimum of 5 wraps around the thermal mass should provide sufficient thermal anchoring, however, additional wraps are recommended for good measure if space permits. To maintain good electrical isolation over many thermal cycles, it is good practice to first varnish a single layer of cigarette paper to the anchored area then wrap the wire around the paper and bond in place with a thin layer of VGE-7031 varnish. Formvar® wiring insulation has a tendency to craze with the application of VGE varnish. If used, the wires cannot be disturbed until the varnish is cured and all solvents have evaporated (typically  $\geq$  24 hours).

#### Copper AA package

Three aspects of using a temperature sensor are critical to its optimum performance. The first involves the proper mounting of the sensor package; the second relates to the proper joining of sensor lead wires and connecting wires; the final concern is the thermal anchoring of the lead wires. Although the sequence in which these areas should be addressed is not fixed, all elements covered under each aspect should be adhered to for maximum operating capabilities of the sensor.

#### **Sensor mounting**

Shown in Figure 7, the copper AA package (or "can") is designed for mounting in a 3.2 mm ( $\frac{1}{2}$  in) hole.

- 1. A hole should be drilled 3.2 mm (1/8 in) diameter by 8.5 mm (0.335 in) deep minimum for the copper can.
- 2. Surface area should be cleaned with a solvent such as acetone followed by an isopropyl alcohol rinse. Allow time for the solvents to evaporate before sensor positioning.
- A small amount of Apiezon® N grease should be applied around the mounting surface and the sensor to enhance thermal contact.
- 4. Position the copper can so that it is fully submerged in the mounting hole.

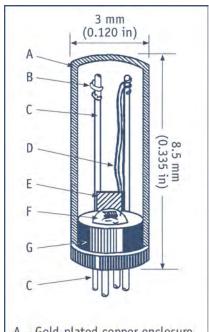


Figure 7—Copper AA package with Cernox® sensor shown. While internal connections are different for the other sensors, the overall package dimensions are the same.

- A Gold-plated copper enclosure B – Current contact zone
- C Phosphor bronze leads
- 0.20 mm (0.008 in) diameter D – Gold leads
- 0.05 mm (0.002 in) diameter
- E Sensing element F - Epoxy heat sink
- G Beryllium oxide base



#### **Lead configurations**

Four leads are attached with strain relief at the sensor. For Cernox®, germanium, and rhodium-iron sensors, each lead is 32 AWG (0.20 mm diameter) phosphor bronze wire, insulated with heavy build polyimide to an overall diameter of 0.24 mm (0.0095 in), 150 mm (6 in) long. For Rox™ sensors, each lead is 34 AWG (0.15 mm diameter) copper wire, insulated with heavy build polyurethane nylon to an overall diameter of 0.185 mm (0.0073 in), 15 cm (6 in) long. Thermal rating of the insulation is 220 °C. Leads are colorcoded at the base of each sensor.

Table 2—Key/color code

	Rox™*	Rox™* Cernox®		Rhodium- iron	
l+	_	White	White	White	
V+		White	Yellow	White	
<b>I</b> –	_	Black	Black	Black	
V—		Black	Green	Black	

\* The Rox™ ruthenium oxide RTD uses the copper AA package but is a 2-lead only device. The leads have no specific polarity. While the Rox™ is built as a 2-lead device, the sensor should be operated in a 4-lead measurement scheme to eliminate errors due to lead resistance, which can be significant.

#### **Extra lead attachment**

If extra long leads are to be attached, then it is recommended that a 4-lead measurement scheme be used with this sensor. Attaching four connecting wires to the sensor leads is recommended. Refer to Table 2 to determine sensor lead polarity.

- 1. Prepare the sensor leads and connecting lead wires with a RMA (rosin mildly active) soldering flux, and tin them with a minimal amount of 60% Sn/40% Pb solder. Use a low wattage soldering iron that will not exceed 200 °C. Clean off residual flux with rosin residue remover. The sensing element inside the package should be protected from excessive heat by putting a heat sink clip over the package.
- 2. Strip connecting wire insulation by delicately scraping with a razor blade, fine sand paper, or steel wool. Phosphor bronze or manganin wire, in sizes 32 or 36 AWG, is commonly used as the connecting lead wire. These wires have low thermal conductivity, which helps minimize the heat flow through the leads. Typical wire insulation is polyvinyl formal (Formvar®) or Polyimide (ML). Formvar® insulation has better mechanical properties such as abrasion resistance and flexibility. Polyimide insulation has better resistance to chemical solvents and burnout.
- 3. Prepare the connecting wire ends with a RMA (rosin mildly active) soldering flux, tin them with a minimal amount of 60% Sn 40% Pb solder. Use a low wattage soldering iron that will not exceed 200 °C.
- Clean off residual flux with rosin residue remover. The sensor lead can be prepared in an identical manner.
- Attach one sensor lead with the connector wire and apply the soldering iron above the joint area until the solders melt, then remove the iron immediately. Repeat for the other set of connector wires and the other sensor lead.
- Avoid putting stress on the device leads and leave enough slack to allow for the thermal contractions that occur during cooling that could fracture a solder joint or lead. This can be achieved with heat shrink tubing.



#### **Heat sinking/thermal anchoring**

- Since the heat flow through the connecting leads can create an offset between the sensor substrate and the true sample temperature, thermal anchoring of the connecting wires is necessary to assure that the sensor and the leads are at the same temperature as the sample.
- Connecting wires should be thermally anchored at several temperatures between room temperature and cryogenic temperatures to guarantee that heat is not being conducted through the leads to the sensing element.
- 3. If the connecting leads have a thin insulation such as Formvar® or polyimide, a simple thermal anchor can be made by winding the wires around a copper post, bobbin, or other thermal mass. A minimum of 5 wraps around the thermal mass should provide enough of an anchor, however, additional wraps are recommended for good measure if space permits. To maintain good electrical isolation over many thermal cycles, it is good practice to first varnish a single layer of cigarette paper to the anchored area, then wrap the wire around the paper and bond in place with a thin layer of VGE-7031 varnish. Formvar® wiring insulation has a tendency to craze with the application of VGE varnish. Once VGE varnish is applied, the wires cannot be disturbed until all solvents have evaporated and the varnish has fully cured (typically 12 to 24 hours).
- A final thermal anchor at the sample itself is a good practice to ensure thermal equilibrium between the sample and temperature sensor.

#### Bare chip installation

#### **General comments**

All of the possible permutations for mounting the chips have not been thoroughly tested. Also, in order to avoid possible adverse effects on stability and thermal mass, heat capacity thermal response times, etc., chips also are not protected by a coating over the active film. The customer must therefore assume some risk of damaging the chips during installation. The sensor and contact films on the Cernox® chips, however, are refractory materials and difficult to scratch. The material presented below includes the best techniques we know to help assure the successful application of unencapsulated chips.

- a. Use good fine-point tweezers. Grasp the chip by the edges at one end (at a contact pad end, if possible). This way, if the tweezers should scrape across the chip, the resistor will not be damaged. Alternately, the wires may be grasped with fingers or tweezers. In the latter case, the operator must develop a very light touch so the wires are not cut or damaged.
- b. If it is necessary to apply pressure to the chip, do so with a cotton swab over the contact area, or with harder objects only outside the patterned area. Do not rub the chip.
- c. Some dirt particles will not hurt the sensor reading, but conducting particles and moisture may, especially if halogen (e.g., chlorine, etc.) contaminants are present. If it is deemed necessary to clean the chips, place a few into a watch glass and rinse with appropriate solvents. (A watch glass is used because it has a curved surface and the sensor will touch only at its corners. It also has a shallow sloped surface, and the rinse liquids can be easily decanted.) Finish with a rinse of pure isopropyl alcohol. Decant the liquid and dry under a light bulb ( $\approx 50~{}^{\circ}\text{C}$ ). For chips with leads, hold the sensor by the leads and immerse it in isopropyl alcohol for a few seconds.  $\text{CO}_2$  snow cleaning can also be very effective, as can ultraviolet/ozone treatments.

#### **Attaching leads**

There are several ways to apply electrical leads to the contact pads, which are gold over contact metal (not wetted easily with solder). In all cases, clamp the sensor chip by the edges and, if possible, do not rely on hand control to position and attach the wires. A clamp can be made from a small, smooth-jawed alligator clip (Figure 8) by cutting off the jaw on the side to which the wire is normally soldered and then fastening that side of the clip to a plate. Another method uses tape to hold the sensors (Figure 9). Kapton® tape and its adhesive will withstand epoxy cure temperatures (165 °C) and the adhesive will not come off on the chip. Do not use Scotch® tape.

The best way by far to connect the chip is to use a thermosonic gold ball bonder. The bonding is clean, uses no flux, and can be done at or near room temperature. The ball attachment at the pad also provides a robust way of making a flying lead that can be attached at the other end later (50 µm diameter gold wire).

Another way is to use silver-loaded conducting epoxy. Make sure the wire and the pads are clean. Use a flexible wire, 40 AWG or smaller, so undue stress will not be applied to the pads. Use a needle to apply small amounts of epoxy to the pads and to the ceramic substrate as well. If the epoxy must be heated in order to cure, a temperature of up to 200 °C could be tolerated by the chip (not Cernox®). This should be done before calibrating, however, since the calibration may shift slightly (shift may amount to 1% of reading at temperatures above 50 K and 0.05% at 4.2 K and below).

#### **Mounting sensor chips**

There are several means of attaching a chip to a substrate. It is possible for strain-induced shifts in calibration to occur. Therefore, keep in mind that the greater the expansion difference between the sensor substrate, the bonding substance and the mating piece, the more likely a strain-induced shift in the calibration may occur. If the joint is stable, this shift probably will be reproducible, and an in-situ calibration may remove the uncertainty. The only substance we have found capable of relieving stress during use is pure indium. This will only work with metallized substrates and in systems that can be heated if the joint is to be soldered.

If it is deemed advisable to use an indium solder joint for reasons of strain, and the mating piece cannot be soldered, a "buffer" layer of metallized BeO or sapphire can be used. Solder the chip to the buffer with indium, and use Stycast® 2850FT/catalyst 9 or equivalent epoxy to attach the buffer to the mating piece.

Stycast® 2850FT or another low expansion, nonconducting epoxy can be used for direct mounting as well. If epoxy is used to completely encapsulate the chip, stress-induced calibration shifts of up to 0.5 K can occur at lower temperatures.

If a greased mounting is desired (Apiezon® N or equivalent), the sensor could be inserted into a hole lined with cigarette paper or tied to a greased surface with thread or dental floss, with paper over it to avoid abrasion. The leads must be insulated with plastic sleeving, fiberglass sleeving, epoxy, or other technique.

VGE-7031 varnish is also a good mounting adhesive and is more easily removed than epoxy. It can be soaked into cigarette paper for a more reliable insulating layer for the leads. The substrate of the sensor is already insulating.


#### Attaching cable wires to sensor leads


The lead wires on a chip sensor are necessarily small in diameter. 50 µm diameter gold wire has a break strength of about 25 g, and 62 µm (42 AWG) copper wire has a rated tensile strength of about 150 g, but the actual break strength is lower because the weak point is usually at the point of attachment or damage from handling (e.g., tweezer marks). The copper wire will only withstand 2 or 3 sharp 90° bends with a 10 g weight attached. The wire will also peel out of silver-loaded epoxy at a smaller force than the rated break strength. However, with reasonable care, loss from damaged leads is negligible.

#### **Soldering**

Both gold and copper wires will dissolve in In and Pb/Sn solders, but gold dissolves much faster. Gold can be successfully soldered by using a temperature controlled iron set just above the solder's melting point. The wire or other attachment point is tinned, and the gold wire stuck into the solder as the iron is removed. If the gold alloy is any length beyond the solder bead, the joint will be greatly weakened, but it is not difficult to repeatedly make successful joints.

Copper wire does not require the precautions above, but repeated soldering will gradually shorten the wire. Keep in mind that heat sinking may be necessary in some situations, but the joints on the chip, if any, will usually be well heat sunk through the chip.

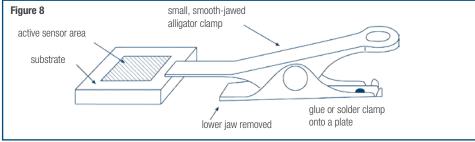

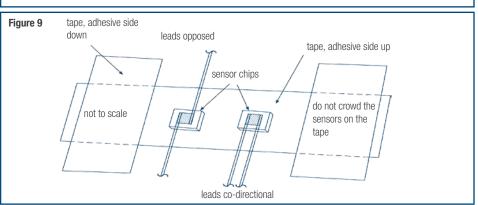





#### **Attachment**

The two most important requirements are that the attachment points of the fine sensor wires should be immobile under all operating conditions, and the sensor leads should have some slack to take up contraction upon cooling. If the leads are connected to a cable, the cable should be attached so it cannot twist at the end. 4-wire (kelvin) cabling schemes down to the sensor leads are preferred for resistance sensors. The lower the resistance of the sensor, the more necessary this becomes.

The following sequence is usually the best:

- 1. Fix the end of the wire or cable in place, with the ends pretinned.
- Apply an insulating layer on the mounting surface if it is a conductor. The uninsulated sensor leads can be kept separate using small Teflon® sleeving or by making channels out of the cigarette paper, Kapton® film, etc. used for the insulator. (See Figure 10.)
- Mount the sensor as desired.
- 4. Adjust the sensor leads into contact with the proper cable wire and solder the joint. It is best to do this by pushing or training the leads into place. (See Figure 11.) Grasping the wire while trying to solder it is inviting wire damage. It is unnecessary to twist the sensor leads around the cable wires. Slack can be built into the leads by using two pairs of tweezers to put an "s-curve" into the wire before soldering.

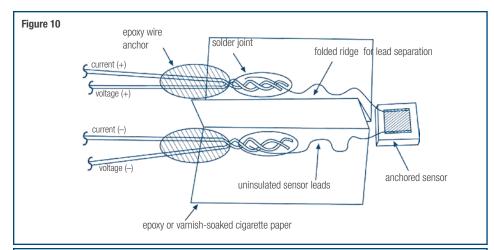

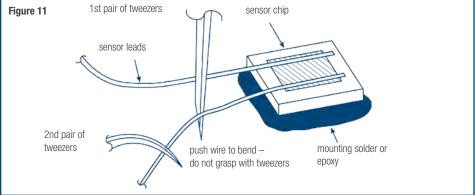

#### Cryogenic accessories

Recommended for proper installation and use of Lake Shore sensors—see Accessories section for more information

#### Stycast® epoxy 2850FT

Permanent attachment, excellent low temperature properties, electrical insulator, low cure shrinkage

#### Apiezon® N grease

Low viscosity, easy to use, solidifies at cryogenic temperatures, excellent lubricant

#### VGE-7031 varnish

Nonpermanent attachment, excellent thermal conductor, easy to apply and remove

#### Indium solder

99.99% pure, excellent electroplating material, foil form

#### 90% Pb 10% Sn solder

Greater lead content, for higher temperature applications greater than 200  $^{\circ}\text{C}$ 

#### Soldering flux

Variety of types

#### Phosphor bronze wire

Available in single, dual, and quad strands, no magnetic attraction, low thermal conduction

#### Manganin wire

Low thermal conductivity, high resistivity, no magnetic attraction

#### Heat sink bobbin

Gold-plated oxygen-free high-conductivity (OFHC) copper bobbins



## Appendix D: Sensor Calibration Accuracies

# Understanding what's available: Uncalibrated—Good SoftCal™—Better Calibrated—Best

The accuracy¹ of a sensor relates to how closely the measurement of resistance (or voltage) can be converted to temperature relative to the recognized international temperature scales (ITS-90 and PLTS-2000). Understanding how the accuracy of temperature sensors is specified begins with the definition of the response curve (e.g., voltage vs. temperature, resistance vs. temperature) for a particular sensor.

Temperature sensors either follow a known standard response within a given tolerance, or they must be calibrated against known standards. Details on calibration procedure are defined in this section. More information on the measurement system and uncertainty analysis is found in Appendix E: Temperature Measurement System.

It is convenient to have temperature sensors that match a standard curve and do not need an individual calibration. Such sensors are interchangeable. Interchangeable sensors follow the same response curve to within a given accuracy and can be interchanged routinely with one another.

Some cryogenic temperature sensors exist currently which are interchangeable within several tolerance bands. The Lake Shore DT-670 series silicon diodes are one example. These conform to five defined accuracy bands about a single curve (Curve 670) and can be ordered by simply specifying the tolerance band required for the experimental accuracy required. In this case, individual calibrations are not performed. However, if increased accuracy is required, full calibration over a specified range may be selected. This provides a fully characterized sensor that can be relied upon for much more accurate measurements, at the cost of being considered "interchangeable" with other sensors from that product line.

In addition to diodes, both platinum and ruthenium oxide sensors also follow a standard curve of resistance versus temperature. Platinum sensors follow an industry standard curve (IEC 751). Lake Shore offers platinum available in Class B tolerance band. If greater temperature accuracy is required, these sensors can be individually calibrated or a SoftCal<sup>™</sup> can be utilized to increase the accuracy of the temperature measurement.

Ruthenium oxide RTDs are also interchangeable. Like silicon diodes, they are interchangeable within a manufacturer lot. Two tolerance bands for ruthenium oxide are defined by Lake Shore.

Table 1, Table 2, and Table 5 summarize Lake Shore temperature sensor accuracies. They are categorized into Good, Better, and Best for each sensor type. The following pages explain the advantages of investing in SoftCal™ or a full calibration from Lake Shore to obtain improved accuracy.

Good	Uncalibrated	<ul> <li>Silicon diodes follow standard curve</li> <li>Platinum resistors follow standard curve</li> <li>Ruthenium oxide (Rox™) resistors follow standard curve</li> <li>GaAlAs diode, carbon-glass, Cernox®, germanium, and rhodium-iron sensors can be purchased uncalibrated but must be calibrated by the customer</li> </ul>
Better	SoftCal™	An abbreviated calibration (2-point: 77 K and 305 K; or 3-point: 77 K, 305 K, and 480 K) which is available for platinum sensors
Best	Calibration	■ All sensors can be calibrated in the various pre-defined temperature ranges. Lake Shore has defined calibration ranges available for each sensor type. The digits represent the lower range in kelvin, and the letter corresponds to high temperature limit, where:  A = 6 K B = 40 K D = 100 K L = 325 K M = 420 K H = 500 K J = 800 K
		For example: The calibration range "1.4L" would result in a sensor characterized from 1.4 K to 325 K

<sup>&</sup>lt;sup>1</sup> The use of the terms accuracy and uncertainty throughout this catalog are used in the more general and conventional sense as opposed to following the strict metrological definitions. For more information, see Appendix B: Accuracy versus Uncertainty.



#### Uncalibrated—Good

With the purchase of an uncalibrated sensor you will receive:

#### Silicon diodes

- Curve 670 data (DT-670)
- Installation instructions

#### **Platinum**

- Standard IEC-751 data
- Installation instructions

#### **Ruthenium oxide**

- Curve data (102, 103, or 202)
- Installation instructions

#### **Thermocouple**

Reference data

## Cernox®, germanium, GaAlAs, carbon-glass, capacitance

- Thermal cycling data—resistance, voltage, or capacitance readings at helium, nitrogen, and room temperature
- Installation instructions

#### Table 1—Uncalibrated sensors: typical accuracy (interchangeability)

							T	emperatur	e						
	0.05 K	0.5 K	1.4 K	2 K	4.2 K	10 K	20 K	25 K	40 K	70 K	100 K	305 K	400 K	500 K	670 K
Silicon diode															
DT-470-SD, Band 11		_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±0.5 K	±1.0 K	±1.0 K					
DT-470-SD, Band 11A	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±1% of temp	±1% of temp	±1% of temp	_				
DT-470-SD, Band 12		_	_	±0.5 K	±0.5 K	±0.5 K	±0.5 K	±1.0 K	±2.0 K	±2.0 K					
DT-470-SD, Band 12A	_	_	_	±0.5 K	±0.5 K	±0.5 K	±0.5 K	±1% of temp	±1% of temp	±1% of temp	_				
DT-470-SD, Band 13	_	_	_	±1.0 K	±1.0 K	±1.0 K	±1.0 K	±1% of temp	±1% of temp	±1% of temp	_				
DT-471-SD	_	_	_	_	_	±1.5 K	±1.5 K	±1.5 K	±1.5 K	±1.5 K	±1.5 K	±1.5% of temp	±1.5% of temp	±1.5% of temp	_
DT-414	_	_	_	±1.5 K	±1.5 K	±1.5 K	±1.5 K	±1.5% of temp	_	_	_				
DT-421	_	_	_	_	_	_	±2.5 K	±2.5 K	±2.5 K	±2.5 K	±2.5 K	±1.5% of temp	_	_	_
DT-670-SD, Band A	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±0.5 K	±0.5 K	±0.5 K	_				
DT-670-SD, Band B	_	_	_	±0.5 K	±0.5 K	±0.5 K	±0.5 K	±0.5 K	±0.33% of temp	±0.33% of temp	_				
DT-670-SD, Band C	_	_	_	±1.0 K	±1.0 K	±1.0 K	±1.0 K	±1.0 K	±0.5% of temp	±0.5% of temp	_				
DT-670-SD, Band D	_	_	_	_	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±0.50 K	±0.2% of temp	±0.2% of temp	_
DT-670-SD, Band E	_	_	_	_	_	_	_	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±0.25% of temp	±0.25% of temp	±0.25% of temp	_
Platinum															
PT-102	_	_	_	_	_	_	_	_		±1.3 K	±1.2 K	±0.5 K	±0.9 K	±1.4 K	±2.3 K
PT-103	_	_	_	_	_	_	_	_	_	±1.3 K	±1.2 K	±0.5 K	±0.9 K	±1.4 K	±2.3 K
PT-111	_			_	_		_	_	_	±1.3 K	±1.2 K	±0.5 K	±0.9 K	±1.4 K	±2.3 K
Rox™															
RX-102A-AA	±10 mK	±25 mK	±50 mK	±75 mK	±125 mK	±300 mK	±1.25 K	±1.5 K	±4.0 K	_	_	_	_	_	_
RX-102A-AA-M	±5 mK	±20 mK	±25 mK	±40 mK	±75 mK	±200 mK	±500 mK	±750 mK	±1.5 K	_	_	_	_	_	_
RX-202A-AA	±15 mK	±30 mK	±100 mK	±125 mK	±250 mK	±1 K	±2.5 K	±3 K	±5.0 K	_	_	_	_	_	_
RX-202A-AA-M	±10 mK	±25 mK	±50 mK	±75 mK	±150 mK	±500 mK	±1.0 K	±1.5 K	±2.0 K	_	_	_	_	_	_
RX-103A-AA		_	±150 mK	±180 mK	±400 mK	±1 K	±2.0 K	±2.5 K	±4.0 K	_	_	_	_	_	_
RX-103A-AA-M	_	_	±50 mK	±75 mK	±100 mK	±300 mK	±700 mK	±1 K	±1.5 K	_	_		_	_	_



#### SoftCal<sup>™</sup>—Better

SoftCal<sup>™</sup> is only available with platinum resistors. With the purchase of SoftCal<sup>™</sup> you will receive:

- Interpolation table and breakpoint interpolation table
- 2-point calibration report (thermal cycling data at LN<sub>2</sub> and room temperature K) OR
- 3-point calibration report (thermal cycling data at LHe, LN<sub>2</sub>, and either 305 K or 480 K)

The temperature characteristics of Lake Shore temperature sensors are extremely predictable, and exhibit excellent uniformity from device to device. The SoftCal™ feature (sensor specific interpolation/extrapolation techniques) allows an abbreviated calibration, based on two or three calibration points, to generate a resistance versus temperature or voltage versus temperature curve over the useful range of selected sensors with remarkable accuracy. In the case of

the Lake Shore platinum resistance sensors, the SoftCal™ procedure makes small adjustments to the IEC-751 curve so that the resulting curve matches the resistance versus temperature characteristic of the individual sensor more closely. SoftCal™ provides the means to generate accurate, inexpensive calibrations for selected Lake Shore sensors to use with either Lake Shore temperature controllers and monitors or the customer's own readout electronics.

Table 2—SoftCal™ (2- and 3-point soft calibration sensors): typical accuracy

	70 K	305 K	400 K	475 K	500 K	670 K
Platinum						
PT-102-2S <sup>2</sup>	±0.25 K	±0.25 K	±0.9 K	±1.3 K	±1.4 K	±2.3 K
PT-103-2S <sup>2</sup>	±0.25 K	±0.25 K	±0.9 K	±1.3 K	±1.4 K	±2.3 K
PT-111-2S <sup>2</sup>	±0.25 K	±0.25 K	±0.9 K	±1.3 K	±1.4 K	±2.3 K
PT-102-3S <sup>3</sup>	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±1.4 K	±2.3 K
PT-103-3S <sup>3</sup>	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±1.4 K	±2.3 K
PT-111-3S <sup>3</sup>	±0.25 K	±0.25 K	±0.25 K	±0.25 K	±1.4 K	±2.3 K

<sup>&</sup>lt;sup>2</sup> 2S (2-point at 77 K and 305 K)

 $<sup>^{\</sup>rm 3}$  3S (3-point at 77 K, 305 K, and 480 K)



#### Calibrated—Best

Lake Shore calibrations include the following:

- Certificate of calibration
- Calibration data plot
- Calibration test data
- Polynominal fit equation and fit comparisons (temperature as a function of resistance or voltage)
- Interpolation table (resistance or voltage as a function of temperature)
- Breakpoint interpolation table
- Instrument breakpoint table

Lake Shore provides precision temperature calibrations for all sensor types, and Lake Shore calibrations are traceable to internationally recognized temperature scales.

Above 0.65 K, calibrations are based on the International Temperature Scale of 1990 (ITS-90). The ITS-90 scale became the official international temperature scale on January 1, 1990; it supersedes the International Practical Temperature Scale of 1968 (IPTS-68) and the 1976 Provisional Temperature Scale (EPT-76). Internally, this scale is maintained on a set of germanium, rhodium-iron, and platinum standards grade secondary thermometers calibrated at the U.S. National Institute of Standards and Technology (NIST) or Great Britain's National Physical Laboratory (NPL), or another recognized national metrology laboratory. Working standard thermometers are calibrated against, and routinely intercompared with, these secondary standards.

For temperatures below 0.65 K, Lake Shore calibrations are based on the Provisional Low Temperature Scale of 2000 (PLTS 2000) adopted by the Comité International des Poids et Mesures in October 2000. Internally, this scale is maintained on a set of germanium and rhodiumiron resistance thermometers calibrated at the U.S. National Institute of Standards and Technology, Great Britain's National Physical Laboratory, or Germany's Physikalisch-Technische Bundesanstalt (PTB). Working standard thermometers are calibrated against, and routinely intercompared with, these secondary standards along with a nuclear orientation thermometer and superconducting fixed points sets.

#### **Calibration method**

Lake Shore performs comparison calibrations measuring the resistance or forward voltage of both the sensor under test and the working standard thermometer. All measurements are performed in a four-lead fashion to eliminate lead resistance.

The sensors to be calibrated are mounted, along with appropriate known standards, in a copper block designed to accommodate a variety of sensor styles. This block is enclosed within a quasi-adiabatic copper radiation shield, which, in turn, is thermally isolated within an outer vacuum jacket.

Constant temperature of the block is achieved by an appropriately mounted heater and precision temperature controller. The electrical, mechanical, and thermal designs of the calibration probe provide extremely stable and uniform temperatures within the copper block.

The calibration process above 4.2 K is computer controlled and the calibration data collected automatically. Data points are usually not at integer temperatures since the primary concern is temperature stability near a data point rather than the specific value. The precise temperature for each data point is subsequently determined. The typical number of data points collected is listed in Table 4 (page 194).

Calibration data is provided for each calibration, together with a calibration data plot and polynomial fits to that raw data, along with a computer generated smoothed interpolation table which is listed as a function of temperature. For resistance sensors, the raw data is given as temperature (T) and resistance (R); the interpolation table shows T, R, dR/dT and dimensionless sensitivity d(log R)/d(log T). For diode sensors, the raw data is given as forward voltage (V) and temperature (T), and the interpolation table presents T, V, and dV/dT.

The specific techniques for generating and controlling calibration temperatures vary, depending on the temperature involved.

Calibrations performed over a wide temperature span frequently entail the consecutive use of a variety of procedures and equipment. In these cases, data points are routinely overlapped to assure integrity of the calibration. The sections that follow describe the specific techniques used for the various temperature ranges.

#### Calibration method—1.2 K to 330 K

Temperatures from 1.2 K to 4.2 K are achieved by filling a He-4 subpot attached to the copper sensor block and pumping on the subpot through a vacuum regulator valve. Temperatures above 4.2 K are achieved by applying controlled power to a heater while the entire probe assembly remains immersed in liquid helium. In either case, the sensors themselves are maintained in a vacuum.



Extreme care is taken to ensure that the sensor block is thermally stable before calibration data is collected. The computer examines successive and interposed measurements of both the known standards and the sensors being calibrated at each data point to verify temperature stability.

Once temperature has stabilized, an appropriate DC excitation current is applied to the thermometer, and the resulting voltage is measured. In the case of resistance sensors, currents from 0.01 mA to 5 mA are selected as required. Sensor voltage is maintained between 1 mV and 3 mV for Cernox®, carbon-glass, germanium, and Rox™ elements up to 300 k $\Omega$ . Higher resistances are measured using a fixed current of 0.01 mA. Sensor power is held between 1 mW and 10 mW for platinum and rhodium-iron resistors.

For resistors, successive voltage readings taken with the current applied in opposite polarities are averaged together to eliminate thermal EMFs from the data. The resistance of the sensing element is determined and reported to five significant figures at each temperature.

Diode thermometers are normally excited with a 10  $\mu$ A current ( $\pm 0.1\%$ ) and the resulting forward voltage reported to five significant figures.

#### Calibration method—below 1.2 K

Calibration temperatures below 1.2 K are produced in a dilution refrigerator. Techniques similar to those for higher temperatures are followed to ensure reliable calibration data. The need for increased care at these lower temperatures, however, requires greater involvement on the part of a skilled system technician and less reliance on automation.

Sensors are measured with a Lake Shore Model 370 AC resistance bridge operated at 13.7 Hz. Germanium and Rox sensors are maintained at a nominal excitation voltage of 20  $\mu$ V RMS (0.05 K to 0.1 K) or 63  $\mu$ V RMS (0.1 K to 1.2 K). Cernox sensors are maintained at a nominal excitation voltage of 20  $\mu$ V RMS from 0.1 K to 0.5 K and 63  $\mu$ V RMS from (0.5 K to 1.2 K).

#### **Accuracy considerations**

The uncertainty associated with a sensor calibration is the net result of each step in the calibration process. A temperature scale disseminated by national standards laboratories is transferred to secondary thermometers maintained by Lake Shore. Those thermometers are used to calibrate in-house working standard thermometers which are then used to calibrate commercial thermometers. Each step introduces an uncertainty that depends on the instrumentation used in the calibration and the specific temperature dependent characteristics of the sensor type calibrated. Other considerations such as calibration block uniformity and stability must also be accounted for. As a result, the calibration accuracy varies with both temperature range and sensor type. Table 3 summarizes the uncertainties associated with the raw data for Lake Shore calibrations.

## Note: The values are the expanded uncertainty based upon a 95% (2 $\sigma$ ) confidence limit with respect to ITS-90.

In practice, however, the uncertainty of subsequent measurements performed with a calibrated sensor should include an additional uncertainty related to the short-term reproducibility of the sensor.

A summary of total calibration uncertainty for selected Lake Shore sensors at specific temperatures is given in Table 5. Errors in each case are expressed in millikelvin deviation from ITS-90 or PLTS-2000. The values in this table reflect the combination of all calibration uncertainties, and the short-term reproducibility upon temperature cycling. It should be noted that at a given temperature, uncertainties are highest for sensors with lowest normalized sensitivity [(1/R)(dR/dT) or (T/R)(dR/dT)] due to the low signal-to-noise ratio.

Lake Shore's calibration facility and procedures for diode and resistance sensor calibrations are traceable to recognized national metrology laboratories and are in compliance with ISO 9001. See page 197 regarding recalibration information.

Table 3—Calibration uncertainty for Lake Shore calibration for selected sensors4

	Germanium		Cernox®				Platinum			Rox™			Diode
		1010	1030	1050	1070	1080	100 Ω	25 Ω	102A	103A	202A	27 Ω	
1.4 K	±4 mK	±4 mK	±4 mK	±4 mK	_	_	_	_	±4 mK	±4 mK	±4 mK	±5 mK	±7 mK
4.2 K	±4 mK	±4 mK	±4 mK	±4 mK	±4 mK	_	_	_	±4 mK	±6 mK	±5 mK	±5 mK	±5 mK
10 K	±4 mK	±5 mK	±5 mK	±4 mK	±4 mK	_	_	_	±10 mK	±15 mK	±12 mK	±7 mK	±6 mK
20 K	±8 mK	±10 mK	±9 mK	±8 mK	±8 mK	±8 mK	±9 mK	±10 mK	±35 mK	±35 mK	±28 mK	±13 mK	±9 mK
30 K	±9 mK	±13 mK	±11 mK	±9 mK	±9 mK	±9 mK	±9 mK	±9 mK	±76 mK	±61 mK	±46 mK	±14 mK	±31 mK
50 K	±11 mK	±18 mK	±14 mK	±12 mK	±12 mK	±11 mK	±10 mK	±10 mK	_	_	_	±13 mK	±37 mK
100 K	±20 mK	±29 mK	±22 mK	±17 mK	±16 mK	±14 mK	±11 mK	±12 mK	_	_	_	±12 mK	±32 mK
300 K	_	±78 mK	±60 mK	±46 mK	±45 mK	±36 mK	±24 mK	±24 mK	_	_	_	±25 mK	±35 mK
400 K	_	±124 mK	±94 mK	±74 mK	±72 mK	±60 mK	±45 mK	±45 mK	_	_	_	±45 mK	±49 mK
500 K	_	_		_	_	_	±51 mK	±51 mK	_	_	_	_	±54 mK

<sup>&</sup>lt;sup>4</sup> All uncertainties are with respect to ITS-90 and represent an approximate 95% configience interval using a coverage factor k=2.



#### Lake Shore calibrations include:

- **1. Certificate of calibration**—This states the traceability of the calibrations performed by Lake Shore to international temperature scales and standards.
- **2. Calibration data**—The measured test data (resistance or forward voltage) is plotted as a function of the temperature. A straight-line interpolation is shown between the data points as a visual aid to the behavior of the sensor.
- 3. Calibration data plot—This table contains the actual calibration data recorded during the calibration of the temperature sensor. The indicated temperatures are those measured using the standard thermometers maintained by Lake Shore, while the voltage or resistance values are the measurements recorded on the device being calibrated.

Table 4—Number of calibration data points

Range (K)	Typical number of data points	Interpolation calibration printout interval
0.050-0.100	6	0.005
0.100-0.300	9	0.010
0.300-0.500	5	0.020
0.500-1.00	7	0.050
1.00-2.00	18	0.10
2.00-5.00		0.20
5.00-10.0	40	0.50
10.0–30.0		1.0
30.0-40.0		2.0
40–100		5.0
100–300	28	5.0
300–380		5.0
340-480 (silicon diodes)	10	5.0
340–480 platinum and rhodium- iron resistors (400 K upper limit)	15	5.0
480-800 platinum sensors only	2	5

**4. Curve fit**—A curve fit is given for each sensor, allowing temperature to be calculated from the measurement of the forward voltage (diodes) or the resistance. One of two curve fit types are used: the first curve fit type is a polynomial equation based on the Chebychev polynomials; the second curve fit type is based on a cubic spline routine. Cubic spline routines are preferred when fitting a rapidly varying function or when smoothing is not desired. In general, the differences between the spline technique and the polynomial fits will be considerably less than the measurement uncertainties.

#### Chebychev polynomial fits

A polynomial equation based on the Chebychev polynomials has the form

$$T(X) = \sum a_n t_n(X)$$
 Eqn. 1

where T(X) represents the temperature in kelvin,  $t_n(X)$  is a Chebychev polynomial,  $a_n$  represents the Chebychev coefficient, and the summation is performed from 0 to the order of the fit. The parameter X is a normalized variable given by

$$X = ((Z-Z_1)-(Z_1-Z))/(Z_1-Z_1).$$
 Eqn. 2

For diodes, Z is simply the voltage V. For resistors, Z is either the resistance R or  $Z = log_{10}(R)$  depending on the behavior of the resistance with temperature.  $Z_L$  and  $Z_U$  designate the lower and upper limit of the variable Z over the fit range.

The Chebychev polynomials can be generated from the recursion relation

$$t_{_{n+1}}(X)=2Xt_{_n}(X)-t_{_{n-1}}(X) \label{eq:tn}$$
 Eqn. 3 where  $t_n(X)=1,\,t_*(X)=X$ 

Alternately, these polynomials are given by

$$t_n(X) = \cos [n \cdot \arccos(X)].$$
 Eqn. 4

All the necessary parameters for using equations 1 through 4 to calculate temperatures from either resistance or voltage are given in the calibration report. This includes the Chebychev coefficients,  $\boldsymbol{Z}_L$  and  $\boldsymbol{Z}_U$ , and also the definition of Z. Depending on the sensor being calibrated and the calibration range, several different fit ranges may be required to span the full temperature range adequately.

The use of Chebychev polynomials is no more complicated than the use of the regular power series, and they offer significant advantages in the actual fitting process. The first step is to transform the measured variable, either R or V, into the normalized variable using equation 2. Equation 1 is then used in combination with equation 3 or 4 to calculate the temperature.

An interesting and useful property of the Chebychev fits is evident in the form of the Chebychev polynomial given in equation 4. The cosine function requires that  $[t_n(X)] \leq 1$ , so no term in equation 1 will be greater than the absolute value of the coefficient. This property makes it easy to determine the contribution of each term to the temperature calculation and where to truncate the series if the full accuracy of the fit is not required.



Table 5—Calibrated sensors: typical accuracy<sup>5</sup>

						Ţ	emperatu	re					
	0.05 K	0.1 K	0.3 K	0.5 K	1 K	1.4 K	4.2 K	10 K	20 K	77 K	300 K	400 K	500 K
Silicon diode													
DT-670-SD/C0	_	_	_	_	_	±12 mK	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	±45 mK	±50 mł
DT-670-CU/CO/LR/CY/ET/BO	_	_	_	_	_	±12 mK	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	_	_
DT-414	_	_	_	_	_	_	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	_	_
DT-421	_	_	_	_	_	±12 mK	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	_	_
DT-470-SD/C0	_	_	_	_	_	±12 mK	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	±45 mK	±50 mk
DT-470-B0/BR/CU/CY/ET/LR/MT	_	_	_	_	_	±12 mK	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	_	_
DT-471-SD/CO	_	_	_	_	_	_	_	±12 mK	±14 mK	±22 mK	±32 mK	±45 mK	±50 mk
DT-471-B0/BR/CU/CY/ET/LR/MT	_	_	_	_	_	_	_	±12 mK	±14 mK	±22 mK	±32 mK	_	_
GaAIAs diode													
TG-120-P	_	_	_	_	_	±12 mK	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	_	_
TG-120-PL	_	_	_	_	_	±12 mK	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	_	_
TG-120-SD/C0	_	_	_	_	_	±12 mK	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	±45 mK	±50 mK
TG-120-CU	_	_	_	_	_	±12 mK	±12 mK	±12 mK	±14 mK	±22 mK	±32 mK	_	_
Cernox®													
CX-1010-AA/CD/C0/CU/LR/ET/MT/SD	_	±3 mK	±3.5 mK	±4.5 mK	±5 mK	±5 mK	±5 mK	±6 mK	±9 mK	±25 mK	±75 mK	_	_
CX-1010-BC	_	_	_	_	_	±5 mK	±5 mK	±6 mK	±9 mK	±25 mK	±75 mK	_	_
CX-1030-AA/CD/C0/CU/LR/ET/MT/SD	_	_	±3 mK	±4 mK	±5 mK	±5 mK	±5 mK	±6 mK	±9 mK	±25 mK	±75 mK	_	_
CX-1030-BC		_	_	_	_	±5 mK	±5 mK	±6 mK	±9 mK	±25 mK	±75 mK		_
CX-1050-AA/BC/CD/C0/CU/LR/ET/MT/SD	_	_	_	_	_	±5 mK	±5 mK	±6 mK	±9 mK	±16 mK	±40 mK		_
CX-1070-AA/BC/CD/C0/CU/LR/ET/MT/SD		_	_	_	_		±5 mK	±6 mK	±9 mK	±16 mK	±40 mK		_
CX-1080-AA/BC/CD/C0/CU/LR/ET/MT/SD		_	_	_	_	_			±9 mK	±16 mK	±40 mK		
CX-1030-C0/SD-HT		_	±3 mK	±4 mK	±5 mK	±5 mK	±5 mK	±6 mK	±9 mK	±16 mK	±40 mK	±65 mK	
CX-1050-CO/SD-HT		_				±5 mK	±5 mK	±6 mK	±9 mK	±16 mK	±40 mK	±65 mK	_
CX-1030-00/SD-HT					_	±3111K	±5 mK	±6 mK	±9 mK	±16 mK	±40 mK	±65 mK	
CX-1080-C0/SD-HT			_	_		_	2011110		±9 mK	±16 mK	±40 mK	±65 mK	_
Carbon-glass	_				_	_	_	_	±9 IIIK	±10 IIIK	±40 IIIK	±05 IIIK	
CGR-1-500, CGR-1-500-CD	_		_	_	_	±4 mK	±4 mK	±5 mK	±8 mK	±25 mK	±105 mK	_	_
CGR-1-1000, CGR-1-1000-CD						±4 mK		±5 mK		±25 mK	±105 mK		
		_	_	_	_		±4 mK		±8 mK				
CGR-1-2000, CGR-1-2000-CD		_	_	_	_	±4 mK	±4 mK	±5 mK	±8 mK	±25 mK	±105 mK	_	_
Rox™	014	0.5(	417	4.5	5.5 ···· /	F17	4016	401/	071/				
RX-102A-AA/CD	±3 mK	±3.5 mK	±4 mK	±4.5 mK	±5.5 mK	±5 mK	±16 mK	±18 mK	±37 mK	<del>-</del>	_	_	_
RX-103A-AA/CD		-	4	45	-	±5 mK	±17 mK	±22 mK	±38 mK	<del>-</del>	_		_
RX-202A-AA/CD	±3 mK	±3.5 mK	±4 mK	±4.5 mK	±5.5 mK	±5 mK	±16 mK	±18 mK	±37 mK	_	_	_	_
Rhodium-iron													
RF-100T-AA/CD/BC/MC	_	_	_	_	_	±11 mK	±11 mK	±12 mK	±14 mK	±15 mK	±25 mK	_	_
RF-100U-AA/CD/BC		_	_	_	_	±11 mK	±11 mK	±12 mK	±14 mK	±15 mK	±25 mK		_
RF-800-4	_	_	_	_	_	±7 mK	±7 mK	±8 mK	±10 mK	±13 mK	±23 mK	±41 mK	±46 mK
Platinum													
PT-102		_	_	_	_	_	_	_	±10 mK	±12 mK	±23 mK	±40 mK	±46 mk
PT-103				_				_	±10 mK	±12 mK	±23 mK	±40 mK	±46 mk
PT-111	_	_	_	_	_	_	_	_	±10 mK	±12 mK	±23 mK	±40 mK	±46 mk
Germanium													
GR-50-AA/CD	±5 mK	±5 mK	±5 mK	±5 mK	±6 mK	±6 mK	±6 mK	_	_	_	_		_
GR-300-AA/CD		_	±4 mK	±4 mK	±4 mK	±4 mK	±4 mK	±4 mK	±8 mK	±25 mK	_		_
GR-1400-AA/CD	_	_	_	_	_	±4 mK	±4 mK	±4 mK	±7 mK	±15 mK	_	_	_

 $<sup>^5</sup>$ All accuracies are: 2  $\sigma$  figures; [(calibration uncertainty) $^2$  + (reproducibility) $^2$ ] $^{0.5}$ ; for additional information, please see Appendix D.

The Chebychev polynomial fit is a smoothing fit and often yields a better representation of the calibration, as it can eliminate some random errors. Along with each set of Chebychev coefficients, a deviation table is given to show how well the polynomial fits the measured test data. This table gives the measured resistance or voltage, the measured temperature, and the temperature calculated from the fit equation. The last column gives the difference in millikelvin (0.001 K) between the measured value and the calculated value. A root mean square (RMS) deviation is given as an indication of the overall quality of the fit and as an indication of the accuracy with which the equation represents the calibration data. Chebychev polynomial fits are provided for all resistance temperature sensor calibrations.

#### Cubic spline fit

Some device types (e.g., GaAlAs diode thermometers) have either a fine structure that is undesirably smoothed by a Chebychev polynomial fit or else a rapidly varying response with temperature. For these devices, a cubic spline fit is provided. A cubic spline fit creates a cubic equation for each interval between calibration points. At each calibration point, the method requires that the cubic equations on either side of the calibration point match in value, first derivative (slope), and second derivative (curvature) at the calibration point. For this fit method, a table is provided listing temperature (T), forward voltage (V), and curvature (C) for each calibration point. In use, the voltage V is measured at the unknown temperature T. Using the provided table, the bracketing calibration points V(k) and V(k+1) are determined and the following quantities are defined:

dV=V(k+1)-V(k), dT=T(k+1)-T(k), dx=V-V(k),	Eqn. 5
from which $S(0)=T(k)$ ,	Eqn. 6
$S(1)=(dT/dV)-dV \cdot (2 \cdot C(k)+C(k+1))/6,$	Eqn. 7
S(2)=C(k)/2, and	Eqn. 8
$S(3)=(C(k+1)-C(k))/(6 \cdot dV)$ are derived.	Eqn. 9

Finally, the temperature is calculated as

$$T=S(0)+S(1) \cdot dx+S(2) \cdot dx^2+S(3) \cdot dx^3$$
. Eqn. 10

A major difference between the Chebychev polynomial fit and the cubic spline fit is that the cubic spline fit provides no smoothing. The curve fit produced by this method passes through each calibration point exactly, so there are no error terms to report.

**5. Interpolation table**—A complete interpolation table is provided over the calibration range of the sensor. This table lists the temperature, the resistance (resistance sensors) or voltage (diode sensors), the sensitivity (dR/dT or dV/dT), and, in the case of resistors, a normalized dimensionless sensitivity [d(log R)/d(logT) = (T/R)·(dR/dT)]. The interpolation table lists resistance or voltage as a function

of temperature, which is the reverse of the curve fit, which gives temperature as a function of sensor units. A cubic spline routine is used to calculate the resistance or voltage at a predetermined set of temperatures. For resistors, the interpolation table is calculated from the smoothed data produced by the Chebychev curve fit. For diodes, however, the interpolation table is calculated from the raw data in order to maintain the fine structure of the sensors' temperature response. Consequently, slight differences between the polynomial equations and the interpolation table are expected. These differences may be on the order of the RMS deviations for the polynomial fits. For resistors, these differences are typically about one tenth the calibration uncertainty. For diodes, the differences may be on the order of the calibration uncertainty in the regions of high curvature and one tenth the calibration uncertainty in the linear regions.

6. Breakpoint table—Lake Shore temperature instruments provide a seamless solution for measuring temperature sensors and converting the measurement into temperature units. The conversion from sensor units to temperature units requires the entry of the temperature response curve into the instrument. For calibrated sensors, this is accomplished through the use of a breakpoint table. With each calibration, Lake Shore provides breakpoint table formats to optimize the performance of the sensor when used with a Lake Shore instrument. The formats provided are compatible with any Lake Shore instrument produced over the last twenty years that accepts user curves. Software is also provided to install the breakpoint table file into most instruments using USB, Ethernet, IEEE-488, or RS-232 interfaces (instrument dependent).

In addition to the breakpoint table and software mentioned above, the CalCurve™ service provides the user with additional alternatives for installing a temperature response curve into a Lake Shore instrument. When the sensor and instrument are ordered together, a factory installed CalCurve service can be provided. A CalCurve can be done in the field when additional or replacement sensors are installed. In this case, curve data is loaded into a non-volatile memory that can be installed into the instrument by the user.

If the sensor is used with customer provided equipment (e.g., voltmeter, current source, and computer) then the curve fit (Chebychev or cubic spline) described in number 4 above should be used. The breakpoint tables are not necessary in this case.

Caution: Proper calculation of a breakpoint table is based upon the interpolation method utilized by the specific instrument for which it is intended. The use of the breakpoint table in an instrument that uses a different interpolation method can cause significant conversion errors.

#### Lake Shore calibration services

- Recalibration
- Calibration report on CD-ROM
- CalCurve<sup>™</sup>
- Certificate of conformance

- Expanded interpolation table
- Second copy of calibration report

#### Recalibration

The stability of a temperature sensor over time is dependent on both its operating environment and history of use. These environmental effects contribute to the degradation of calibration over time:

- Ionizing radiation
- Thermal shock
- Thermal stress from continuous exposure to high temperatures (relative to the sensor materials)
- Mechanical shock
- Improper use
- Corrosion (a serious problem for systems of dissimilar metallurgies in the presence of moisture and chemical agents such as salts—this includes integrated circuits and other electronics)
- Electrical stress/electromagnetic interference (EMI)/electrostatic discharge (ESD)

There are no specific published regulations or guidelines that establish requirements for the frequency of recalibration of cryogenic temperature sensors. There are certainly military standards for the recalibration of measuring devices. However, these standards only require that a recalibration program be established and then adhered to in order to fulfill the requirements.

Temperature sensors are complex assemblies of wires, welds, electrical connections, dissimilar metallurgies, electronic packages, seals, etc., and hence have the potential for drift in calibration. Like a voltmeter, where components degrade or vary with time and use, all of the "components" of a temperature sensor may also vary, especially where they are joined together at material interfaces. Degradation in a sensor materials system is less apparent than deterioration in performance of a voltmeter.

Lake Shore sensor calibrations are certified for one year. Depending upon the sensor type and how it is used, it is recommended that sensors be recalibrated in the Lake Shore Calibration Service Department periodically. Certainly, recalibration before important experiments would be advisable.

Sensors stored at a consistent temperature have been shown to remain stable over a long period of time.

#### 8000 CalCurve™

The 8000 CalCurve<sup>™</sup> on CD-ROM is provided free of charge at the time of order to any customer who orders a calibrated sensor. The 8000 is the calibration breakpoint interpolation data. Also on the CD is a PC executable program to load the data into a Lake Shore instrument. Once the data is loaded into the instrument, the user can calculate and display temperature with the data. The following information is included with the 8000 CalCurve<sup>™</sup>:

- Raw data
- Coefficients
- Interpolation table
- Instrument breakpoints
- A program for installing curves into instrument
- Instructions describing all file formats and contents

There is a charge to load previously stored calibration curves.



## Appendix E: Temperature Measurement System

The goal is to measure the temperature of some system. The ability to do so accurately and with the required resolution depends on a variety of factors. The calibration report from Lake Shore (or any calibration facility) is only the first step in determining the accuracy of the temperature measurement in the end-user's system.

A more quantifiable term than accuracy is total uncertainty of the measurement. This is simply the measurement itself and an estimate of all the errors of the measurement. Smaller errors are considered more accurate. The first step in estimating the errors in a customer system is the calibration itself. Essentially, a calibration is a series of resistance or voltage measurements of an unknown sensor and a corresponding measurement of an established temperature. By accounting for all the uncertainties of the measurement (installation, instrumentation, etc.) a total uncertainty is estimated. The actual accuracy a customer can expect will depend on this and other factors:

## **1. Design errors:** Can the system be measured by the sensor?

These are errors of design and happen prior to sensor installation. For example, whether or not the sensor can be mounted on or near the sample to be measured could be a design error. If it is too far away, there can be thermal lags and offsets due to thermal conductance of the sample. Another example would be using too large a sensor to measure small samples. The thermal mass of the sensor could bias the temperature of the sample.

Design errors also apply to the physical construction of the sensor. This affects the reproducibility of the sensor over thermal cycling. Some sensors are more fragile than others and more prone to physical damage (for example carbon-glass RTDs).

## **2. Installation and environment errors:** Does the interaction of the sensor and system disturb the measurement?

This would include installation errors and environmental effects. If leads are not properly heat sunk, they will introduce a heat load into the sensor. This affects the sensor's measurement and can also affect the sample. It can bias the reading of temperature as well as directly affect the temperature if the heat leak is great enough. Other interactions include thermal radiation, magnetic fields, and radiation.

## **3. Operation and instrumentation errors:** Does the instrumentation introduce errors to the measurement?

Instrumentation is a crucial component to the total quality of the measurement. The choice of 2-lead or 4-lead measurements, excitation currents, instrument resolution, and accuracy all affect the measurement. Additionally, grounding errors, RF noise coupling, and thermal EMFs can introduce noise to the measurement.

Error terms can be classified into two classes:

**Type A**, (or random): Errors that can be evaluated by statistical methods.

**Type B**, (or systematic): Errors that can be evaluated by other means.

Most random errors are the result of instrumentation: uncertainty in the current source and voltage measurements. Other random errors are the actual assignment of a temperature (transferring ITS-90 or PLTS-2000), and interpolation errors. Design, installation, and environmental errors are systematic. For example, sensors in magnetic fields will create an offset to the measurement. This offset can be estimated from prior information or directly measured by other means (isothermal measurements with and without field). RF noise can also cause both random errors (adds to current noise) and systematic errors since at ultra-low temperatures the added noise can self-heat the sensors causing a systematic offset.

#### Installation

2-lead vs. 4-lead installations can lead to significant measurement errors. Even with a properly installed temperature sensor, poor thermal design of the overall apparatus can produce measurement errors. Installation issues are addressed in Appendix C: Sensor Packaging and Installation, along with detailed installation instructions for specific Lake Shore sensors.

#### **Environmental concerns**

Temperature sensors can be affected by changes in the environment. Examples include magnetic fields, ionizing radiation, or changes in the pressure, humidity, or chemistry of the environment. The most common are magnetic field and radiation-induced errors. These effects have been discussed previously. These environmental effects will create a systematic bias in the temperature measurement.

#### Instrumentation

#### 2-lead versus 4-lead

The measurement of resistance and diode temperature sensors requires passing a current through the temperature sensor to produce a sensor voltage that can be measured. The simplest resistance or voltage measurement configuration is a current source connected to the temperature sensor with a voltmeter connected to the current leads as shown in Figure 1. The current source can be represented as an ideal current source ( $I_s$ ) in parallel with a shunt resistance,  $R_s$ . The voltmeter, normally a digital multimeter (DMM) can be modeled as an ideal voltmeter ( $V_{io}$ ) in parallel with an input impedance,  $R_{io}$ .

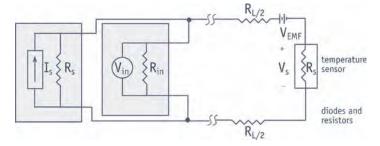


Figure 1—2-lead resistance measurement

The dominant source of error in a 2-lead resistance measurement is usually the resistance of the lead wires connecting the current source to the temperature sensor. In a cryogenic environment, the flow of heat down the leads of the cryostat is of critical concern due to the potential for sensor element heating. Normally, wire of small diameter and significant resistance per unit length is preferred to minimize this heat flow. Consequently, the resulting lead resistance can become a significant percentage of the resistance measured. The wire also has its own temperature sensitivity of resistance.

The equivalent error the lead resistance represents depends on the sensor type and sensor sensitivity. The 100  $\Omega$  platinum RTD has a nominal resistance of 100  $\Omega$  at 273.15 K (0 °C). The IEC 751 standard for the temperature sensitivity for platinum RTDs is 0.385  $\Omega/K$  between 273.15 K and 373.15 K (0 °C to 100 °C). Both the magnitude of the resistance and the temperature sensitivity are relatively small numbers, especially when the lead resistance may be several ohms. A 10  $\Omega$  lead resistance would result in a positive 26 K error in this temperature range (10  $\Omega/0.385$   $\Omega/K=26$  K). The effect of lead resistance becomes even greater as the temperature decreases, since the temperature sensitivity (dR/dT) of platinum sensors decreases with decreasing temperature.

Additionally, it is not uncommon for the internal lead resistance of the current leads (parasitic resistance) of a germanium or carbon-glass sensor to be as much as 10% to 20% or more of the sensor 4-lead resistance. Consequently, the 4-lead calibrated resistance-temperature data is of little use for a 2-lead measurement and the temperature error associated with 2-lead resistance measurements for germanium

and carbon-glass is almost always extremely large. The parasitic resistance for Cernox® temperature sensors, due to having common current and voltage contact, is extremely small. Even still, the low temperature error due to lead resistance can be at least 3 mK for  $100~\Omega$  of lead resistance. Since lead wire has its own temperature dependence, the error could be much larger. Table 1 shows typical error with 2-lead measurement.

In order to eliminate the effects of lead resistance, a 4-lead measurement (Figure 2) is normally used. Two of the leads, I+ and I-, supply current to the sensor, while the other two leads, V+ and V-, are used to eliminate the effect of lead resistance by measuring the voltage at the sensor voltage leads (4-lead sensor) or directly at the device leads (2-lead sensor). The reason this measurement scheme works is that the IR drop in the current leads is not measured, and the voltage drop in the voltage leads is extremely small due to the very small current required by the voltmeter (picoamperes or less) to make the voltage measurement.

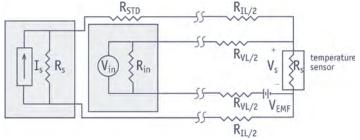


Figure 2—4-lead resistance measurement

A diode temperature sensor measurement requires a fixed 10  $\mu$ A current source and a voltmeter. As with resistance measurements, the dominant source of error in a 2-lead diode measurement is often the lead resistance. A 100  $\Omega$  lead resistance will result in a 1 mV voltage error at a current of 10  $\mu$ A. The Lake Shore DT-400 Series silicon diode temperature sensors have an average sensitivity of approximately –26 mV/K below 30 K, resulting in a temperature error of 40 mK (1 mV/26 mV/K = 0.038 K); above 30 K the sensitivity is approximately –2.3 mV/K, resulting in error exceeding 400 mK (1 mV/2.3 mV/K = 0.435 K). Consequently, unless the lead resistance can be reduced in magnitude or the resultant error can be tolerated, a 4-lead measurement is recommended.

Table 1—Typical errors for Cernox® 1070 resistor with lead resistance at 100  $\Omega$  (50  $\Omega$  each lead)

	R (Ω)	dR/dT (Ω/K)	ΔT (mK)
4.2 K	52444	-33321	-3.0
20 K	2851.6	-200.70	-498
100 K	317.44	-4.3908	large
300 K	68.949	-0.3052	very large



#### **Voltmeter input impedance**

The voltmeter input impedance is generally not a problem in 2- or 4-lead measurements. It is not uncommon for today's voltmeters to have a 109  $\Omega$  or 1010  $\Omega$  input impedance in the voltage ranges of interest, which is large when compared to the temperature sensor resistance. Consequently, virtually no current will be shunted from the temperature sensor into the voltage measurement circuitry at these input impedance levels. A voltmeter input impedance of 109  $\Omega$  would produce only a 0.0001% error in a 1000  $\Omega$  resistance measurement.

#### **Current source output impedance**

The output impedance of a good current source is also not ordinarily a problem in either 2- or 4-lead measurements, for the same reason. If the output impedance is not large compared to the sensor resistance, then a known series resistor should be placed in one of the current paths, and the current to the sensor should be measured by measuring the voltage across the known standard resistance.

#### Thermoelectric and zero offset voltages

Voltages develop in electrical conductors with temperature gradients when no current is allowed to flow (Seebeck effect). Thermoelectric voltages appear when dissimilar metals are joined and joints are held at different temperatures. Typical thermoelectric voltages in cryogenic measurement systems are on the order of microvolts.

This effect can be minimized by a few steps. The same material should be used for conductors whenever practical, and the number of connections, or joints, in the measurement circuit should be minimized. Low thermal EMF solder can also be used (cadmium-tin solder has a lower thermal EMF than tin-lead solder by a factor of ten).

In addition to thermal offset, the instrumentation can have a zero offset (the signal value measured with no input to the measuring instrument). The zero offset can drift with time or temperature and is usually included in the instrument specifications.

The total offset voltage can be measured by reversing the current. When reading the voltage with the current in the forward direction, the voltmeter will read:

$$V^+ = V_s + V_{EME}$$
 Eqn. 1

where  $\rm V_{\rm S}$  is the actual voltage reading of the sensor, and  $\rm V_{\rm EMF}$  is the lead thermal EMFs. When the current is reversed, the voltmeter will read

$$V^- = -V_S + V_{EMF}$$
 Eqn. 2

When the current is reversed, the voltage due to the sensor reverses sign while the thermal EMFs do not. The true voltage (V) across the sensor is

$$V = (V^+ - V^-)/2 = V_s$$
 Eqn. 3

By averaging the forward and reverse current voltage measurements, the error in the voltage measurement due to thermal EMFs is eliminated.

Diode measurements do not allow current reversal. The value of the offset voltage can be estimated by shorting the leads at the diode and measuring the offset voltage with zero excitation current at operating temperature.

Thermal EMFs in the sensor leads and connections do not have as big an effect on diode measurements as they do on resistance measurements, since the diode signal levels are much larger (typically a few tenths of a volt at room temperature to several volts at 4.2 K).

#### Grounding

Signal grounding is important to the stability and repeatability of measurements. A measurement system that includes sensors, instruments, cabling, and possibly computer interfacing requires careful grounding.

Improper grounding of instruments or grounding at multiple points can allow current flows which result in small voltage offsets. The current flow through ground loops is not necessarily constant, resulting in a fluctuating voltage. Current can flow in the ground loop as it acts as a large aperture for inductive pickup. Also, current can result if there is a potential difference due to multiple grounds.

As each instrument handles grounding differently, it is important to carefully read your instrument manual for grounding suggestions. The grounding and isolation is handled differently in the Model 370 than in other Lake Shore instruments, since it is used for ultra-low temperature measurements. Ideally, there should be one defined ground for the measurement, and the cryostat is the best choice. Realistically, however, there are many instruments, wiring, and pumps attached to the cryostat. Each instrument may have its own ground. Simply attaching ground straps may create more ground loops.

Books on grounding and shielding can help to identify and eliminate both ground loops and electromagnetic noise.

#### **Reducing AC signal interference (RF noise)**

Signal leads and cables are very susceptible to interference from unwanted AC signals in the RF frequency range. They act like antennas and pick up noise from computers, monitors, instrumentation, radio broadcasts, and other sources. Signals are either inductively coupled or capacitively coupled. The induced signals circulate as noise current in the measurement leads and distort measurements. There are other concerns when diodes are used as the sensing element, as discussed in the next section.

There are several ways to reduce the effect of AC signals. First, when possible, remove or shield the source of unwanted signals. Second, make each pair of signal leads as bad an antenna as possible. This can be accomplished by keeping them short and using twisted leads. Twisting reduces loop area to make leads that are prone to picking up noise smaller targets to electromagnetic signals. Twisting also helps to cancel unwanted signals in leads that are prone to transmit noise. In a typical 4-lead measurement, the current leads should be twisted together and the voltage leads should be twisted together. Third, put a conductive shield around all the leads to divert electric field signals and prevent capacitive coupling into the leads. Tie the shield to the ground closest in potential to the measurement. Many Lake Shore instruments provide a shield pin on the sensor connector for this purpose. The shield should be tied only at the instrument. Attaching at any other point can cause ground loops that were previously discussed.

In cases where shielding is not enough, filtering the unwanted signals can be considered. It is very difficult to add a filter to a measurement system without changing the measurement. One type of filter that has proven to work is a ferrite bead (see the Accessories section). The bead will act like a high impedance to unwanted high frequency signals and not affect the slow moving desired signals being measured. The Lake Shore 2071 ferrite bead can be clamped around existing wiring.

The greatest concern relates to leads external to the cryostat. Ideally, the cryostat itself acts as the shield for all wiring internal to it. However, it is still possible for cross-talk between different signal leads. In this application Lake Shore recommends Quad-Twist™ cryogenic wire, which has two twisted pairs of phosphor bronze wire that minimize noise pickup and allow proper heat sinking. In extreme cases coaxial cable may be needed, although it is much more difficult to heat sink.

#### Measurement errors in diode thermometers due to AC interference

Wiring techniques are especially important when using diode thermometers in a measurement system. Noise currents produce a shift in measurement. Because diodes have a nonlinear voltage response to the changing current, the shift is seen as a lower measured voltage corresponding to a higher measured temperature. The temperature error in noisy systems can be as high as several tenths of a kelvin.

The following equation can be used to estimate the temperature shift with DT-470 silicon diodes over the range 0 <  $V_{RMS}$  < 40 mV and 30 < T < 300 K. The temperature errors tend to decrease at temperatures below 30 K ( $\Delta T$  in K, T in K, and  $V_{RMS}$  in mV).

$$\Delta T = 2.7768 \cdot T^{-1.11953} \cdot V_{RMS}^{2.01803}$$
 Eqn. 4

There are two simple techniques that can be used to determine if this problem is present in the measuring system. The first is to connect a 10  $\mu$ F capacitor in parallel with the diode to act as a shunt for any induced AC currents. The capacitor must have low leakage current so

it does not alter the DC current through the diode. If the DC voltage reading across the diode increases with the addition of the capacitor, AC noise currents are present. The second method involves the measurement of the AC voltage across the diode. While an oscilloscope is the logical choice for looking at AC signals, many do not have the sensitivity required and often introduce unwanted grounds into the system and compound the problem. An AC voltmeter should be used.

Lake Shore instrumentation includes a 1  $\mu$ F capacitor across the current source in order to minimize the effects of noise related to power line frequency. A 0.1  $\mu$ F capacitor in parallel with a 30 pF to 50 pF capacitor at the voltage measurement input are used to minimize the effects of AC-coupled digital noise. The obvious disadvantage of the addition of AC filtering is that it slows down the response time of the measurement system.

#### **Effect of current source accuracy**

**Diode temperature sensors**—Measurement accuracy of diode sensors is not as strongly dependent upon the current source accuracy as is the case with resistance temperature sensors. Diode sensors possess a nonlinear forward current-voltage characteristic. Consequently, the forward voltage variation with changing current for diodes is smaller than for resistance temperature sensors, which have linear current-voltage characteristics.

Below 30 K, the sensitivity (dV/dT) of Lake Shore diode temperature sensors increases by an order of magnitude over sensitivities at higher temperatures. The slope (dV/dl) of the I-V curves (Table 2) stays relatively constant. Both characteristics further reduce the effect of any change in forward bias current on temperature measurement accuracy.

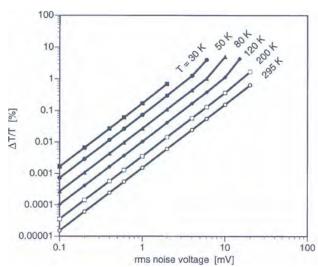


Figure 3—Calculated temperature reading shifts due to voltage noise across a Lake Shore Model DT-470 silicon diode temperature sensor



#### **Effect of current source accuracy**

Diode temperature sensors—Measurement accuracy of diode sensors is not as strongly dependent upon the current source accuracy as is the case with resistance temperature sensors. Diode sensors possess a nonlinear forward current-voltage characteristic. Consequently, the forward voltage variation with changing current for diodes is smaller than for resistance temperature sensors, which have linear current-voltage characteristics.

Below 30 K, the sensitivity (dV/dT) of Lake Shore diode temperature sensors increases by an order of magnitude over sensitivities at higher temperatures. The slope (dV/dl) of the I-V curves (Table 2) stays relatively constant. Both characteristics further reduce the effect of any change in forward bias current on temperature measurement accuracy.

Table 2 - Approximate dV/dl values for the DT-470 sensor

	Approximate dV/dI (Ω)
300 K	3000
100 K	1000
40 K	700
4.2 K	2800

Table 3—Equivalent temperature offsets for the DT-470 diode temperature sensors at selected current source uncertainties

	dV/dl (Ω)	dV/dT	Temperature offset (mK)					
		(mV/K)	dl(%)=0.05	dl(%)=0.1				
300 K	3000	-2.40	6.5	13				
100 K	1000	-2.04	2.5	5				
40 K	700	-1.74	2.0	4				
4.2 K	2800	-33.6	0.4	1				

Lake Shore diode current sources are typically set to 10  $\mu$ A  $\pm$ 0.1% or better and have a low-pass filter to minimize the effect of AC pickup in the current leads. Resultant errors due to current source inaccuracy are on the order of 10 mK or less for diode sensors.

If the output from a current source is not precisely 10  $\mu$ A, the resultant error in temperature can be calculated using this relationship between the dV/dT and dV/dI values:

$$\Delta T = ((dV/dI)/(dV/dT))\Delta I$$
 Eqn. 5

Note: dV/dl and dV/dT values are derived at the same temperature T.

In the above expression,  $R_d = dV/dI$  and  $R_s = V/I$  are the dynamic and static resistances of the temperature sensor. Note that the dynamic and static resistances of an ohmic sensor are equal. Results shown in Table 3.

**Resistance temperature sensors**—for resistance sensors, an error in current measurement is inversely related to the resultant measurement error of resistance:

$$R - \Delta R = V/(I + \Delta I)$$
 Eqn. 6 
$$\approx (V/I)(1 - \Delta I/I)$$
 
$$= R - R(\Delta I/I)$$

where I is the current setting,  $\Delta I$  is the variation from that setting, and  $\Delta R = R\Delta I/I$ .

The temperature error,  $\Delta T$ , due to current source uncertainty,  $\Delta I(\%)$ , is

$$\Delta T = \Delta R/(dR/dT)$$
Eqn. 7
$$= R(\Delta I/I)/(dR/dT)$$

$$\Delta T = \Delta I(\%)/[(100/R)(dR/dT)]$$

where  $\Delta I(\%) = 100 \Delta I/I$ 

All Lake Shore resistance current sources are typically set to 0.01%. For example (Table 4), temperature errors for a platinum resistance sensor near room temperature due to the current source can approach 36 mK and diminish to less than 10 mK below 100 K.

Table 4—Equivalent temperature offsets for selected resistance sensors at selected voltmeter and current source uncertainties

				Temperature	e offset (mK)
	T (K)	R (Ω)	dR/dT (Ω)/K	dV(%)=0.01 dl(%)=0.01	dV(%)=0.05 dl(%)=0.05
PT-100	300	110.452	0.388	28.5	142.5
	100	29.987	0.411	7.3	36.5
	40	5.938	0.291	2.0	10
CGR-1-1000	300	6.21021	-0.0047	132	660
	100	9.66389	-0.0465	20.8	104
	40	16.8227	-0.3211	5.2	26
	4.2	964.19	-842.21	0.1	1
CX-1050	300	50.1	-0.165	30.4	152
	100	154.62	-1.55	10.0	50
	40	376.1	-8.9	4.2	21
	4.2	4596	-1867	0.2	1
GR-200A-1000	100	4.95987	-0.0469	10.6	53
	40	18.7191	-0.844	2.2	11
	4.2	981.026	-451.3	0.2	1

#### **Effect of voltage measurement accuracy**

**Diode temperature sensors**—The effect of voltage measurement accuracy on resultant temperature measurement is not difficult to calculate, provided that diode sensitivity is known for the temperature of interest. The potential temperature error,  $\Delta T_{\nu}$  is

$$\Delta T_v = \Delta V/[dV/dT]$$
 Eqn. 8

Table 5 illustrates potential temperature error due to the voltage measurement.

Table 5—Equivalent temperature offsets for the DT-470 diode temperature sensor at selected voltmeter uncertainties

	Ţ		dV/dT	Temperature	offset (mK)
	(K)	(V)	(mV/K)	ΔV(%)=0.01	ΔV(%)=0.05
DT-470	300	0.51892	-2.40	21.6	108
	100	0.97550	-2.04	47.8	239
	40	1.08781	-1.74	62.5	313
	4.2	1.62602	-33.6	4.8	24

**Resistance temperature sensors**—for positive temperature coefficient resistors such as platinum or rhodium-iron, the potential temperature error,  $\Delta T_p$ , is

$$\Delta T_R = \Delta R / [dR/dT]$$
 Eqn. 9  
=  $[\Delta V/I] / [dR/dT]$ 

since from Ohm's law,  $\Delta V = I\Delta R$ .

But  $\Delta V$  (%) = 100 $\Delta V/V$ ; therefore

$$\Delta T_{R} = [V\Delta V (\%) / 100I] / [dR/dT]$$
 Eqn. 10  

$$= [\Delta V (\%) \cdot R/100] / [dR/dT]$$

$$\Delta T_{R} = \Delta V (\%) / [(100/R) (dR/dT)], \text{ and}$$

$$\Delta T_{R} = \Delta I(\%) / [(100/R) (dR/dT)]$$

The temperature offsets in Table 4 are calculated using both of the above equations.

This is not surprising, as we are dealing with Ohm's Law and a linear system.

#### **Self-heating**

Any difference between the temperature of the sensor and the environment the sensor is intended to measure produces a temperature measurement error or uncertainty. Dissipation of power in the temperature sensor will cause its temperature to rise above that of the surrounding environment. Power dissipation in the sensor is also necessary to make a temperature measurement. Minimization of the temperature measurement uncertainty thus requires balancing the uncertainties due to self-heating and output signal measurement.

Self-heating is really a combination of sensor design and instrumentation. The primary reason for self-heating offsets at low temperatures is the thermal boundary resistance between the active sensor element and its surroundings. The thermal boundary resistance has a very strong inverse cube relationship with temperature. This forces the instrumentation to be capable of sourcing a small excitation and measuring a small (voltage) signal. The optimum excitation power will be a function of sensor, resistance, and temperature.

Lake Shore temperature controllers each have different excitation currents for NTC RTDs which effectively defines the minimum temperature range of the instrument-sensor combination.

An estimate of the self-heating error including thermal resistance for select sensors and optimum excitation power is found in Table 6 (page 205).



#### Thermal (Johnson) noise

Thermal energy produces random motions of the charged particles within a body, giving rise to electrical noise. The minimum root mean square (RMS) noise power available is given by  $P_{\rm n}=4kT~\Delta f_{\rm n}$ , where k is the Boltzmann constant and  $\Delta f_{\rm n}$  is the noise bandwidth. Peak-to-peak noise is approximately five times greater than the RMS noise. Metallic resistors approach this fundamental minimum, but other materials produce somewhat greater thermal noise. The noise power is related to current or voltage noise by the relations:  $I=[P_{\rm n}/R_{\rm d}]^{0.5}$  and  $V=[P_{\rm n}~R_{\rm d}]^{0.5}$ . The noise bandwidth is not necessarily the same as the signal bandwidth, but is approximately equal to the smallest of the following:

- $\pi/2$  times the upper 3 dB frequency limit of the analog DC measuring circuitry, given as approximately 1/(4  $R_{eff}$   $C_{in}$ ) where  $R_{eff}$  is the effective resistance across the measuring instrument (including the instrument's input impedance in parallel with the sensor resistance and wiring) and  $C_{in}$  is the total capacitance shunting the input
- 0.55/t, where t, is the instrument's 10% to 90% rise time
- 1 Hz if an analog panel meter is used for readout
- One half the conversion rate (readings per second) of an integrating digital voltmeter

#### **Calibration uncertainty**

Commercially calibrated sensors should have calibrations traceable to international standards. About the best accuracy attainable is represented by the ability of national standards laboratories. Many laboratories provide calibrations for a fee. The calibration uncertainty typically increases by a factor of 3 to 10 between successive devices used to transfer a calibration.

#### **Calibration fit interpolation uncertainty**

Once a calibration is performed, an interpolation function is required for temperatures that lie between calibration points. The interpolation method must be chosen with care, since some fitting functions can be much worse than others. Common interpolation methods include linear interpolation, cubic splines, and Chebychev polynomials. Formulas based on the physics of the sensor material may give the best fits when few fit parameters are used.

Use of an interpolation function adds to the measurement uncertainty. The additional uncertainty due to an interpolation function can be gauged by the ability of the interpolation function to reproduce the calibration points. Each calibration can be broken up into several ranges to decrease the fitting uncertainties. Typical uncertainties introduced by the interpolation function are on the order of one tenth the calibration uncertainty.

#### Combining measurement uncertainties

Estimating the quality of a measurement involves the following steps:

1) identify the relevant sources of measurement uncertainty, 2) change the units of all uncertainties to temperature, and 3) combine all of the uncertainties using the root sum of squares method described later. Examples of source of measurement uncertainties affecting the accuracy, but not the precision of a measurement include offset voltages and calibration uncertainties.

The expected uncertainty of a measurement is expressed in statistical terms. As stated in the Guide to the Expression of Uncertainty in Measurement:

"The exact values of the contributions to the error of the measurement arising from the dispersion of the observations, the unavoidable imperfect nature of the corrections, and incomplete knowledge are unknown and unknowable, whereas the uncertainties associated with these random and systematic effects can be evaluated. ...the uncertainty of a result of a measurement is not necessarily an indication of the likelihood that the measurement result is near the value of the measurand; it is simply an estimate of the likelihood of nearness to the best value that is consistent with presently available knowledge."

The uncertainty is given the symbol u and has the same units as the quantity measured. The combined uncertainty  $u_{\rm c}$  arising from several independent uncertainty sources can be estimated by assuming a statistical distribution of uncertainties, in which case the uncertainties are summed in quadrature according to

$$u_{r} = \sqrt{u_{i}^{2} + u_{2}^{2} + ... + u_{i}^{2} + ... + u_{n}^{2}}$$
 Eqn. 11

Both random and systematic uncertainties are treated in the same way. Note that both sides of Equation 11 can be divided by the measurement quantity to express the measurement uncertainty in relative terms. Finding statistical data suitable for addition by quadrature can be a problem; instrument and sensor specifications sometimes give maximum or typical values for uncertainties. Two approaches may be taken when dealing with maximum uncertainty specifications. The conservative approach is to use the specification limit value in the combined uncertainty calculation. The less conservative approach is to assume a statistical distribution within the specification limits and assume the limit is roughly three standard deviations, in which case one third of the specification limit is used in uncertainty calculations. The manufacturer may be able to supply additional information to help improve uncertainty estimates. Practical recommendations and procedures for problems related to the estimation of measurement uncertainties are discussed in greater detail by Rabinovich.

Table 6 gives examples of uncertainty calculations for two types of temperature sensors, the DT-470-SD silicon diode sensor, and the CX-1050-AA Cernox® sensor.

When Lake Shore accounts for uncertainties in calibration measurements, all the above issues are taken into consideration, and their contributions are estimated.

#### **References:**

ISO/TAG 4/WG 3. Guide to the Expression of Uncertainty in Measurement, First Edition. Geneva, Switzerland: International Organization for Standardization, 1992.
 S. Rabinovich, Measurement Errors, College Park, Maryland: American Institute of Physics, 1993.

Table 6—Combined temperature measurement calculation examples

	DT-470-	-SD-11	CX-10	)50-AA
Temperature, T	80	K	4.	2 K
Mounting environment (N-greased to block)	vacu	ium	liquid	helium
Static Electrical Resistance, R <sub>s</sub>	101,525 Ω (s	tatic $R_s = V/I$ )	4920 Ω (sta	atic R <sub>s</sub> = V/I)
Dynamic Electrical Resistance, R <sub>d</sub>	1000 Ω (dynan	$\operatorname{R}_{d} = \operatorname{dV/dI}$	4920 Ω (dyna	$mic R_d = dV/dI$
Excitation current, I	10	μΑ	1	μΑ
Output voltage, V	1.015	525 V	4.92	2 mV
Dimensionless temperature sensitivity, S <sub>D</sub>	-0.1	521	-1	.71
Uncertainties due to:	Value used	Temperature uncertainty u <sub>T</sub> /T (PPM)	Value used	Temperature uncertainty u <sub>T</sub> /T (PPM)
Measurement instrumentation				
(Keithley Instruments 2000 DVM)				
Meter range full scale (FS)	10.00000 V		100.0000 mV	
Voltage accuracy specification (ppm)	±(30+5 FS/V)	521	±(50+35 FS/V)	445
Sensor self-heating				
Thermal resistance	$R_t = 1000 \text{ K/W}$	127	$R_{t} = 3500 \text{ K/W}$	4.1
Excitation uncertainty (Lake Shore Model 120-CS)				
Current accuracy specification	u/I = 0.05%	32	$u_{l}/I = 0.1\%$	585
Thermal noise		0.02		0.2
Thermal voltages and zero drift	10 μV	65	O <sup>1</sup>	0
Electromagnetic noise <sup>2</sup>	2 mV	1040	O <sup>1</sup>	0
Calibration uncertainty	0.250 K <sup>3</sup>	3130	4 mK	952
Interpolation uncertainty <sup>4</sup>	313		95.2	
Combined uncertainties (ppm)	3357		1206	

<sup>&</sup>lt;sup>1</sup> Eliminated by current reversal

 $<sup>^2 \</sup>quad \text{Assuming an AC voltage of 2 mV}_{\text{ms}} \text{ is read across the voltmeter terminals} \\ \text{—the voltage is converted to an approximate temperature shift}$ 

<sup>&</sup>lt;sup>3</sup> Calibration accuracy

<sup>&</sup>lt;sup>4</sup> Assumed to be one tenth the calibration uncertainty

#### **Estimating self-heating of temperature sensors**

Any difference between the temperature of the sensor and the environment the sensor is intended to measure produces a temperature measurement error or uncertainty. Dissipation of power in the temperature sensor will cause its temperature to rise above that of the surrounding environment. Power dissipation in the sensor is also necessary to make a measurement with most temperature sensors (exceptions include thermocouples and optical pyrometers). Minimization of the temperature measurement uncertainty thus requires balancing the uncertainties due to self-heating and output signal measurement. The possibility that other experimental considerations might impose more stringent limitations on the power that can be dissipated in the temperature sensor should also be considered.

Following are two approaches to dealing with the problem of self-heating:

 Choose an excitation that allows acceptable instrumentation measurement uncertainty and check to make sure self-heating is negligible at one or two points where it is likely to be most significant.

An easy way to check for self-heating is to increase the power dissipation and check for an indicated temperature rise. Unfortunately, this procedure will not work with non-linear devices such as semiconductor diodes. An indication of the self-heating error can be made by reading the diode temperature in both a liquid bath and in a vacuum at the same temperature, as measured by a second thermometer not dissipating enough power to self-heat significantly.

2. Measure the thermal resistance in the temperature range of interest and calculate the optimum operating point.

Examination leads to the conclusion that an increase in the sensor output voltage will result in a decreasing temperature uncertainty, so long as the voltage uncertainty remains constant. This is possible with an ohmic sensor by increasing the excitation current. Unfortunately, a larger excitation will dissipate more power in the temperature sensor, raising its temperature above the surroundings.


The self-heating depends on the excitation power according to the equation

$$\Delta T_{sh} = P_s R_t = I^2 R_e R_t = V^2 R_t / R_e$$
 Eqn. 12

where  $\Delta T_{sh}$  is the temperature rise due to self-heating,  $P_s$  is the power dissipated in the sensor, I is the excitation current,  $R_e$  is the electrical resistance, and  $R_t$  is the thermal resistance between the sensor and its environment. The thermal resistance is extremely difficult to calculate for all but the simplest cases and is best determined experimentally using the following procedure:

- Mount the sensor as it will be used on a temperature controlled block or directly in liquid
- 2. Record the output voltage as a function of excitation current (I-V curve) until significant self-heating is observed (when  $\rm R_e$  =V/I is no longer constant)
- Replot the data as sensor temperature reading versus power dissipated (T versus P),
- 4. Fit the data with a linear equation of the form  $T = T_0 + R_t P_s$  to find the thermal resistance,  $R_t$

Thermal resistance values determined from some commercial resistance temperature sensors in common mounting configurations are shown as a function of temperature in Figure 4. The thermal resistance varies with the environment in and around the sensor package (vacuum, gas, liquid), sensor mounting (solder, grease, clamp pressure, epoxy, etc.) and details of sensor construction. The thermal resistances shown in the figure should be used only as a guide with reference to the source papers and preferably measurement on the actual sensor in the temperature range and environment of use. See www.lakeshore.com for additional notes and papers.

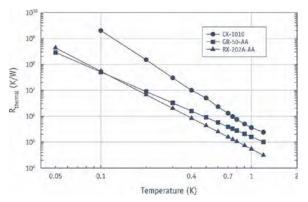

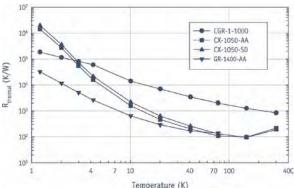


Figure 4—Thermal resistance data for various sensors as a function of T

## Appendix F: PID Temperature Control

#### Closed loop PID control

Closed loop PID control, often called feedback control, is the control mode most often associated with temperature controllers. In this mode, the controller attempts to keep the load at exactly the user entered setpoint, which can be entered in sensor units or temperature. To do this, it uses feedback from the control sensor to calculate and actively adjust the control (heater) output. The control algorithm used is called PID.

The PID control equation has three variable terms: proportional (P), integral (I), and derivative (D)—see Figure 1. The PID equation is:

HeaterOutput = 
$$P[e + I \int (e)dt + D\frac{de}{dt}]$$
 Eqn. 1

where the error (e) is defined as: e = Setpoint - Feedback Reading.

#### **Proportional (P)**

The proportional term, also called gain, must have a value greater than zero for the control loop to operate. The value of the proportional term is multiplied by the error (e) to generate the proportional contribution to the output: Output (P) = Pe. If proportional is acting alone, with no integral, there must always be an error or the output will go to zero. A great deal must be known about the load, sensor, and controller to compute a proportional setting (P). Most often, the proportional setting is determined by trial and error. The proportional setting is part of the overall control loop gain, as well as the heater range and cooling power. The proportional setting will need to change if either of these change.

#### Integral (I)

In the control loop, the integral term, also called reset, looks at error over time to build the integral contribution to the output:

Output(I) = 
$$PI \int (e)dt$$
. Eqn. 2

By adding integral to the proportional contribution, the error that is necessary in a proportional-only system can be eliminated. When the error is at zero, controlling at the setpoint, the output is held constant by the integral contribution. The integral setting (I) is more predictable than the proportional setting. It is related to the dominant time constant of the load. Measuring this time constant allows a reasonable calculation of the integral setting.


#### **Derivative (D)**

The derivative term, also called rate, acts on the change in error with time to make its contribution to the output:

Output(D) = PD 
$$\frac{de}{dt}$$
. Eqn. 3

By reacting to a fast changing error signal, the derivative can work to boost the output when the setpoint changes quickly, reducing the time it takes for temperature to reach the setpoint. It can also see the error decreasing rapidly when the temperature nears the setpoint and reduce the output for less overshoot. The derivative term can be useful in fast changing systems, but it is often turned off during steady state control because it reacts too strongly to small disturbances or noise. The derivative setting (D) is related to the dominant time constant of the load.

Figure 1—Examples of PID Control

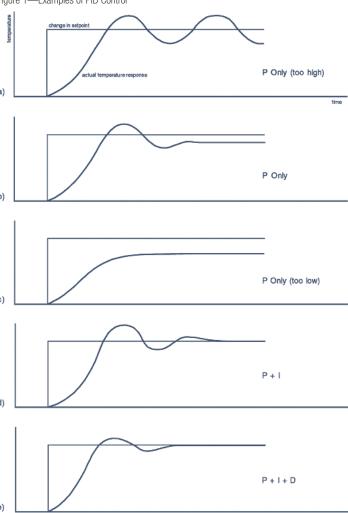



#### Tuning a closed loop PID controller

There has been a lot written about tuning closed loop control systems and specifically PID control loops. This section does not attempt to compete with control theory experts. It describes a few basics to help users get started. This technique will not solve every problem, but it has worked for many others in the field. It is also a good idea to begin at the center of the temperature range of the cooling system.

#### **Setting heater range**

Setting an appropriate heater output range is an important first part of the tuning process. The heater range should allow enough heater power to comfortably overcome the cooling power of the cooling system. If the heater range will not provide enough power, the load will not be able to reach the setpoint temperature. If the range is set too high, the load may have very large temperature changes that take a long time to settle out. Delicate loads can even be damaged by too much power.

Often there is little information on the cooling power of the cooling system at the desired setpoint. If this is the case, try the following: allow the load to cool completely with the heater off. Set manual heater output to 50% while in Open Loop control mode. Turn the heater to the lowest range and write down the temperature rise (if any). Select the next highest heater range and continue the process until the load warms up through its operating range. Do not leave the system unattended; the heater may have to be turned off manually to prevent overheating. If the load never reaches the top of its operating range, some adjustment may be needed in heater resistance or an external power supply may be necessary to boost the output power of the instrument.

The list of heater range versus load temperature is a good reference for selecting the proper heater range. It is common for systems to require two or more heater ranges for good control over their full temperature. Lower heater ranges are normally needed for lower temperature.

#### **Tuning proportional**

The proportional setting is so closely tied to heater range that they can be thought of as fine and coarse adjustments of the same setting. An appropriate heater range must be known before moving on to the proportional setting.

Begin this part of the tuning process by letting the cooling system cool and stabilize with the heater off. Place the instrument in closed loop PID control mode, then turn integral, derivative, and manual output settings off. Enter a setpoint above the cooling system's lowest temperature. Enter a low proportional setting of approximately 5 or 10 and then enter the appropriate heater range as described above. The heater display should show a value greater than zero and less than 100% when temperature stabilizes. The load temperature should stabilize at a temperature **below** the setpoint. If the load temperature and heater display swing rapidly, the heater range or proportional value may be set too high and should be reduced. Very slow changes in load temperature that could be described as drifting are an indication of a proportional setting that is too low (which is addressed in the next step).

Gradually increase the proportional setting by doubling it each time. At each new setting, allow time for the temperature of the load to stabilize. As the proportional setting is increased, there should be a setting in which the load temperature begins a sustained and predictable oscillation rising and falling in a consistent period of time. (Figure 1a). The goal is to find the proportional value in which the oscillation begins. Do not turn the setting so high that temperature and heater output changes become violent. In systems at very low temperature it is difficult to differentiate oscillation and noise. Operating the control sensor at higher than normal excitation power can help.

Record the proportional setting and the amount of time it takes for the load change from one temperature peak to the next. This time is called the oscillation period of the load. It helps describe the dominant time constant of the load, which is used in setting integral. If all has gone well, the appropriate proportional setting is **one half** of the value required for sustained oscillation. (Figure 1b).

If the load does not oscillate in a controlled manner, the heater range could be set too low. A constant heater reading of 100% on the display would be an indication of a low range setting. The heater range could also be too high, indicated by rapid changes in the load temperature or heater output less than 10% when temperature is stable. There are a few systems that will stabilize and not oscillate with a very high proportional setting and a proper heater range setting. For these systems, setting a proportional setting of one half of the highest setting is the best choice.

#### **Tuning integral**

When the proportional setting is chosen and the integral is set to zero (off), the instrument controls the load temperature below the setpoint. Setting the integral allows the control algorithm to gradually eliminate the difference in temperature by integrating the error over time. (Figure 1d). A time constant that is too high causes the load to take too long to reach the setpoint. A time constant that is too low can create instability and cause the load temperature to oscillate.

Note: The integral setting for each instrument is calculated from the time constant. The exact implementation of integral setting may vary for different instruments. For this example it is assumed that the integral setting is proportional to time constant. This is true for the Model 370, while the integral setting for the Model 340 and the Model 331 are the inverse of the time constant.



Begin this part of the tuning process with the system controlling in proportional only mode. Use the oscillation period of the load that was measured above in seconds as the integral setting. Enter the integral setting and watch the load temperature approach the setpoint. If the temperature does not stabilize and begins to oscillate around the setpoint, the integral setting is too low and should be doubled. If the temperature is stable but never reaches the setpoint, the integral setting is too high and should be decreased by half.

To verify the integral setting make a few small (2 to 5 degree) changes in setpoint and watch the load temperature react. Trial and error can help improve the integral setting by optimizing for experimental needs. Faster integrals, for example, get to the setpoint more quickly at the expense of greater overshoot. In most systems, setpoint changes that raise the temperature act differently than changes that lower the temperature.

If it was not possible to measure the oscillation period of the load during proportional setting, start with an integral setting of 50. If the load becomes unstable, double the setting. If the load is stable make a series of small setpoint changes and watch the load react. Continue to decrease the integral setting until the desired response is achieved.

#### **Tuning derivative**

If an experiment requires frequent changes in setpoint or data taking between changes in the setpoint, derivative should be considered. (Figure 1e). A derivative setting of zero (off) is recommended when the control system is seldom changed and data is taken when the load is at steady state.

A good starting point is one fourth the integral setting in seconds (i.e., ¼ the integral time constant). Again, do not be afraid to make some small setpoint changes: halving or doubling this setting to watch the effect. Expect positive setpoint changes to react differently from negative setpoint changes.

#### **Manual output**

Manual output can be used for open loop control, meaning feedback is ignored and the heater output stays at the user's manual setting. This is a good way to put constant heating power into a load when needed. The manual output term can also be added to the PID output. Some users prefer to set an output value near that necessary to control at a setpoint and let the closed loop make up the small difference.

NOTE: Manual output should be set to 0 when not in use.

Typical sensor performance sample calculation: Model 331S temperature controller operating on the 2.5 V input range used with a DT-670 silicon diode at 1.4 K

- Nominal voltage—typical value taken from Appendix G: Sensor Temperature Response Data Tables.
- Typical sensor sensitivity—typical value taken from Appendix G: Sensor Temperature Response Data Tables.
- Measurement resolution in temperature equivalents

Equation: Instrument measurement resolution/typical sensor sensitivity

 $10 \mu V / 12.49 mV/K = 0.8 mK$ 

The instrument measurement resolution specification is located in the Input Specifications table for each instrument.

Electronic accuracy in temperature equivalents

Equation: Electronic accuracy (nominal voltage)/typical sensor sensitivity

 $(80 \mu V + (0.005\% \cdot 1.644 V)) / 12.49 \text{ mV/K} = \pm 13 \text{ mK}$ 

The electronic accuracy specification is located in the Input Specifications table for each instrument.

Temperature accuracy including electronic accuracy, CalCurve<sup>™</sup>, and calibrated sensor

Equation: Electronic accuracy + typical sensor accuracy at temperature point of interest

 $13 \text{ mK} + 12 \text{ mK} = \pm 25 \text{ mK}$ 

The typical sensor accuracy specification is located in the Accuracy table for each instrument.

 Electronic control stability in temperature equivalents (applies to controllers only)

Equation: Up to 2 times the measurement resolution

 $0.8 \text{ mk} \cdot 2 = \pm 1.6 \text{ mK}$ 



### Appendices

## Appendix G: Sensor Temperature Response Data Tables

Silicon diode DT-670			
T (K)	V (volts)	dV/dT (mV/K)	
1.4	1.64429	-12.49	
4.2	1.57848	-31.59	
10	1.38373	-26.84	
20	1.19775	-15.63	
30	1.10624	-1.96	
50	1.07310	-1.61	
77.35	1.02759	-1.73	
100	0.98697	-1.85	
150	0.88911	-2.05	
200	0.78372	-2.16	
250	0.67346	-2.24	
300	0.55964	-2.30	
350	0.44337	-2.34	
400	0.32584	-2.36	
450	0.20676	-2.39	
500	0.09068	-2.12	

Silicon diode DT-470			
T (K)	V (volts)	dV/dT (mV/K)	
1.4	1.6981	-13.1	
4.2	1.6260	-33.6	
10	1.4201	-28.7	
20	1.2144	-17.6	
30	1.1070	-2.34	
50	1.0705	-1.75	
77.35	1.0203	-1.92	
100	0.9755	-2.04	
150	0.8687	-2.19	
200	0.7555	-2.31	
250	0.6384	-2.37	
300	0.5189	-2.4	
350	0.3978	-2.44	
400	0.2746	-2.49	
450	0.1499	-2.46	
475	0.0906	-2.22	

GaAlAs TG-120			
T (K)	V (volts)	dV/dT (mV/K)	
1.4	5.3909	-97.5	
4.2	4.7651	-214	
10	3.7521	-148	
20	2.5341	-97.5	
30	1.8056	-48.2	
50	1.4637	-2.82	
77.35	1.4222	-1.24	
100	1.3918	-1.48	
150	1.2985	-2.25	
200	1.1738	-2.64	
250	1.0383	-2.77	
300	0.8978	-2.85	
350	0.7531	-2.99	
400	0.6066	-2.97	
450	0.4556	-3.08	
475	0.3778	-3.15	

Cernox® CX-1010*			
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)
0.1	21389	-558110	-2.70
0.2	4401.6	-38756	-1.76
0.3	2322.4	-10788	-1.39
0.4	1604.7	-4765.9	-1.19
0.5	1248.2	-2665.2	-1.08
1	662.43	-514.88	-0.78
1.4	518.97	-251.77	-0.68
2	413.26	-124.05	-0.60
3	328.95	-58.036	-0.53
4.2	277.32	-32.209	-0.49
6	234.44	-17.816	-0.46
10	187.11	-8.063	-0.43
20	138.79	-3.057	-0.44
30	115.38	-1.819	-0.47
40	100.32	-1.252	-0.50
50	89.551	-0.929	-0.52
77.35	70.837	-0.510	-0.56
100	61.180	-0.358	-0.59
150	47.782	-0.202	-0.63
200	39.666	-0.130	-0.66
250	34.236	-0.090	-0.66
300	30.392	-0.065	-0.65

<sup>\*</sup>Cernox sensors do not follow a standard response curve — the listed values are typical, but can vary widely; consult Lake Shore to choose a specific range

Cernox® C>	(-1030*		
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)
0.3	31312	-357490	-3.43
0.4	13507	-89651	-2.65
0.5	7855.7	-34613	-2.20
1	2355.1	-3265.2	-1.39
1.4	1540.1	-1264.9	-1.15
2	1058.4	-509.26	-0.96
3	740.78	-199.11	-0.81
4.2	574.20	-97.344	-0.71
6	451.41	-48.174	-0.64
10	331.67	-19.042	-0.57
20	225.19	-6.258	-0.56
30	179.12	-3.453	-0.58
40	151.29	-2.249	-0.59
50	132.34	-1.601	-0.61
77.35	101.16	-0.820	-0.63
100	85.940	-0.552	-0.64
150	65.864	-0.295	-0.67
200	54.228	-0.184	-0.68
250	46.664	-0.124	-0.67
300	41.420	-0.088	-0.64
350	37.621	-0.065	-0.61
400	34.779	-0.050	-0.57
420	33.839	-0.045	-0.55

<sup>\*</sup>Cernox sensors do not follow a standard response curve — the listed values are typical, but can vary widely; consult Lake Shore to choose a specific range



Cernox® C	K-1050*		
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)
1.4	26566	-48449	-2.55
2	11844	-11916	-2.01
3	5733.4	-3042.4	-1.59
4.2	3507.2	-1120.8	-1.34
6	2252.9	-432.14	-1.15
10	1313.5	-128.58	-0.98
20	692.81	-30.871	-0.89
30	482.88	-14.373	-0.89
40	373.11	-8.392	-0.90
50	305.19	-5.507	-0.90
77.35	205.67	-2.412	-0.91
100	162.81	-1.488	-0.91
150	112.05	-0.693	-0.93
200	85.800	-0.397	-0.92
250	69.931	-0.253	-0.90
300	59.467	-0.173	-0.87
350	52.142	-0.124	-0.83
400	46.782	-0.093	-0.79
420	45.030	-0.089	-0.77

*Cernox sensors do not follow a standard response curve — the listed values are	
typical, but can vary widely; consult Lake Shore to choose a specific range	

Cernox® CX-1070*			
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)
4.2	5979.4	-2225.3	-1.56
6	3577.5	-794.30	-1.33
10	1927.2	-214.11	-1.11
20	938.93	-46.553	-0.99
30	629.90	-20.613	-0.98
40	474.89	-11.663	-0.98
50	381.42	-7.490	-0.98
77.35	248.66	-3.150	-0.98
100	193.29	-1.899	-0.98
150	129.60	-0.854	-0.99
200	97.626	-0.477	-0.98
250	78.723	-0.299	-0.95
300	66.441	-0.201	-0.91
350	57.955	-0.143	-0.86
400	51.815	-0.106	-0.81
420	49.819	-0.094	-0.80

<sup>\*</sup>Cernox sensors do not follow a standard response curve — the listed values are typical, but can vary widely; consult Lake Shore to choose a specific range

Cernox® CX-1080*				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	
20	6157.5	-480.08	-1.56	
30	3319.7	-165.61	-1.50	
40	2167.6	-79.551	-1.47	
50	1565.3	-45.401	-1.45	
77.35	836.52	-15.398	-1.42	
100	581.14	-8.213	-1.41	
150	328.75	-3.057	-1.40	
200	220.93	-1.506	-1.36	
250	163.73	-0.863	-1.32	
300	129.39	-0.545	-1.26	
350	106.98	-0.368	-1.20	
400	91.463	-0.261	-1.14	
420	86.550	-0.231	-1.12	

<sup>\*</sup>Cernox sensors do not follow a standard response curve — the listed values are typical, but can vary widely; consult Lake Shore to choose a specific range





Carbon-glass CGR-1-500				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)·(dR/dT)	
1.4	103900	-520000	-6.9	
4.2	584.6	-422.3	-3.1	
10	85.64	-13.39	-1.6	
20	36.21	-1.77	-0.98	
30	25.46	-0.653	-0.77	
50	18.05	-0.213	-0.59	
77.35	14.33	-0.098	-0.48	
100	12.75	-0.055	-0.43	
150	10.85	-0.027	-0.37	
200	9.79	-0.017	-0.34	
250	9.08	-0.012	-0.33	
300	8.55	-0.009	-0.33	

Carbon-glass CGR-1-2000			
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)
1.4	1401600	-8440000	-8.4
4.2	2260	-2060	-3.8
10	196.7	-39.1	-1.9
20	66.57	-4.05	-1.2
30	43.14	-1.35	-0.94
50	28.47	-0.401	-0.70
77.35	21.65	-0.157	-0.56
100	18.91	-0.094	-0.49
150	15.70	-0.044	-0.42
200	13.96	-0.027	-0.39
250	12.83	-0.019	-0.38
300	11.99	-0.015	-0.36

Carbon-glass CGR-1-1000				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)·(dR/dT)	
1.4	342900	-1900000	-7.8	
4.2	967.4	-802.77	-3.5	
10	104.9	-19.046	-1.8	
20	38.76	-2.183	-1.1	
30	25.88	-0.760	-0.88	
50	17.51	-0.233	-0.66	
77.35	13.51	-0.093	-0.53	
100	11.86	-0.057	-0.48	
150	9.92	-0.027	-0.41	
200	8.87	-0.017	-0.37	
250	8.18	-0.012	-0.36	
300	7.66	-0.009	-0.35	



Germanium GR-50-AA				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	
0.05	35000	-3642000	-5.2	
0.095	2725	-92170	-3.2	
0.2	364.6	-4043	-2.2	
0.3	164.0	-964.0	-1.8	
0.5	73.75	-202.9	-1.4	
1	33.55	-31.33	-0.93	
1.4	24.73	-13.15	-0.74	
2	19.32	-6.167	-0.64	
3	15.59	-2.334	-0.45	
4.2	13.66	-1.036	-0.32	
5	13.02	-0.624	-0.24	

Germanium GR-1400-AA				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	
1.4	35890	-94790	-3.7	
2	11040	-16670	-3.0	
3	3636	-2977	-2.5	
4.2	1689	-861.9	-2.1	
6	800.7	-278.9	-2.1	
10	252.8	-61.95	-2.5	
20	44.19	-5.41	-2.4	
30	17.18	-1.23	-2.1	
50	6.45	-0.211	-1.6	
77.35	3.55	-0.050	-1.1	
100	2.80	-0.021	-0.74	

Germanium GR-300-AA			
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)·(dR/dT)
0.3	35180	-512200	-4.4
0.5	5443	-34800	-3.2
1	875.7	-1901	-2.2
1.4	448.6	-581.3	-1.8
2	248.8	-187.4	-1.5
3	142.2	-60.14	-1.3
4.2	94.46	-26.56	-1.2
6	61.89	-12.42	-1.2
10	33.20	-3.97	-1.2
20	16.08	-0.750	-0.93
30	10.81	-0.383	-1.1
50	5.91	-0.148	-1.3
77.35	3.50	-0.050	-1.1
100	2.72	-0.024	-0.88





Rox™ RX-102A				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)·(dR/dT)	
0.05	63765	-2888654	-2.27	
0.1	19400	-266199	-1.37	
0.2	8277	-43000	-1.04	
0.3	5615	-16647	-0.89	
0.5	3701	-5478	-0.74	
1	2381	-1257	-0.53	
1.4	2005	-667	-0.47	
2	1726	-331	-0.38	
3	1502	-152	-0.30	
4.2	1370	-80.4	-0.25	
6	1267	-41.0	-0.19	
10	1167	-15.2	-0.13	
20	1089	-3.96	-0.07	
30	1063	-1.78	-0.05	
40	1049	-1.03	-0.04	

Rox™ RX-102B				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	
0.01	9856.38	-413888	-0.42	
0.02	7289.79	-170565	-0.47	
0.03	5975.92	-100138	-0.50	
0.04	5184.10	-62048	-0.48	
0.05	4676.87	-41480	-0.44	
0.1	3548.94	-12578	-0.35	
0.2	2813.75	-4116	-0.29	
0.3	2502.26	-2365	-0.28	
0.5	2187.50	-1056	-0.24	
1	1884.56	-350.8	-0.19	
1.4	1779.33	-197.7	-0.16	
2	1691.44	-114.5	-0.14	
3	1606.45	-63.53	-0.12	
4.2	1546.44	-40.04	-0.11	
6	1488.89	-26.05	-0.11	
10	1410.19	-15.43	-0.11	
20	1300.92	-7.82	-0.12	
30	1239.54	-4.83	-0.12	
40	1198.80	-3.41	-0.11	

Rox™ RX-103A			
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)
1.4	30745	-13571	-0.62
2	25089	-6551	-0.52
3	20712	-2943	-0.43
4.2	18149	-1559	-0.36
6	16131	-812	-0.30
10	14083	-337	-0.24
20	12289	-102	-0.17
30	11555	-53.0	-0.14
40	11137	-33.2	-0.12

Rox™ RX-202A				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	
0.05	69191	-3186379	-2.30	
0.1	21927	-256913	-1.17	
0.2	10962	-45681	-0.83	
0.3	8079	-18420	-0.68	
0.5	5904	-6434	-0.54	
1	4299	-1631	-0.38	
1.4	3820	-879	-0.32	
2	3443	-453	-0.26	
3	3127	-222	-0.21	
4.2	2929	-124	-0.18	
6	2765	-68.6	-0.15	
10	2582	-31.6	-0.12	
20	2389	-12.1	-0.10	
30	2299	-6.84	-0.09	
40	2243	-4.60	-0.08	



DI I' DT 400				
Platinum PT-100				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)·(dR/dT)	
20	2.2913	0.085	0.74	
30	3.6596	0.191	1.60	
50	9.3865	0.360	1.90	
77.35	20.380	0.423	1.60	
100	29.989	0.423	1.40	
150	50.788	0.409	1.20	
200	71.011	0.400	1.10	
250	90.845	0.393	1.10	
300	110.354	0.387	1.10	
400	148.640	0.383	1.00	
500	185.668	0.378	1.00	
600	221.535	0.372	1.00	
700	256.243	0.366	1.00	
800	289.789	0.360	1.00	

Rhodium-iron RF-800-4				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)-(dR/dT)	
1.4	1.5204	0.178	0.16	
4.2	1.9577	0.135	0.29	
10	2.5634	0.081	0.32	
20	3.1632	0.046	0.29	
30	3.5786	0.040	0.34	
50	4.5902	0.064	0.71	
77.4	6.8341	0.096	1.1	
100	9.1375	0.106	1.2	
150	14.463	0.105	1.1	
200	19.641	0.102	1.0	
250	24.686	0.101	1.0	
300	29.697	0.101	1.0	
350	34.731	0.101	1.0	
400	39.824	0.103	1.0	

Rhodium-iron RF-100				
T (K)	R (Ω)	dR/dT (Ω/K)	(T/R)·(dR/dT)	
1.4	6.892	0.489	0.10	
4.2	8.2053	0.418	0.21	
10	10.072	0.245	0.24	
20	11.858	0.137	0.23	
30	13.130	0.131	0.29	
50	16.724	0.242	0.72	
77.4	25.298	0.368	1.1	
100	34.123	0.403	1.2	
150	54.292	0.396	1.1	
200	73.692	0.381	1.0	
250	92.529	0.374	1.0	
300	111.19	0.373	1.0	

Thermocouple type E (T <sub>Ref</sub> = 273.15 K)			
T (K)	EMF (μV)	dV/dT (μV/K)	
3.2	-9834.9	1.59	
4.2	-9833.0	2.09	
10	-9813.3	4.66	
20	-9747.0	8.51	
30	-9643.8	12.1	
40	-9505.5	15.5	
50	-9334.2	18.7	
75	-8777.7	25.6	
100	-8063.4	31.4	
150	-6238.1	41.2	
200	-3967.4	49.3	
250	-1328.7	56.0	
300	1608.0	61.1	
350	4777.7	65.6	
400	8159.8	69.6	
500	15426	75.3	
600	23138	78.6	
670	28694	80.0	
700	31100	80.4	
800	39179	81.0	
900	47256	80.4	
1000	55247	79.3	
1100	63119	78.1	
1200	70842	76.3	
1270	76136	75.2	



Thermocouple type K ( $T_{Ref} = 273.15 \text{ K}$ )			
T (K)	EMF (μV)	dV/dT (μV/K)	
3.2	-6457.7	0.743	
4.2	-6456.9	0.916	
10	-6448.5	2.01	
10.5	-6447.4	2.12	
20	-6417.8	4.15	
30	-6365.1	6.39	
40	-6290.0	8.61	
50	-6193.3	10.7	
75	-5862.9	15.6	
100	-5417.6	19.9	
150	-4225.5	27.5	
200	-2692.8	33.5	
250	-897.60	38.0	
300	1075.3	40.6	
350	3135.8	41.5	
400	5200.0	40.8	
500	9215.6	40.3	
600	13325	41.7	
670	16264	42.2	
700	17533	42.4	
800	21789	42.6	
900	26045	42.4	
1000	30251	41.7	
1100	34373	40.7	
1200	38396	39.7	
1270	41153	39.0	
1300	42318	38.7	
1400	46131	37.5	
1500	49813	36.1	
1600	53343	34.5	
1640	54712	34.0	

Thermocouple type T ( $T_{Ref} = 273.15 \text{ K}$ )			
T (K)	EMF (μV)	dV/dT (μV/K)	
3.2	-6257.5	1.03	
4.2	-6256.2	1.40	
10	-6242.9	3.12	
20	-6199.2	5.58	
30	-6131.3	7.99	
40	-6040.0	10.2	
50	-5927.7	12.2	
75	-5573.6	16.0	
100	-5131.2	19.4	
150	-4004.3	25.6	
200	-2575.3	31.4	
250	-872.57	36.6	
300	1067.4	40.8	
350	3215.5	45.0	
400	5560.2	48.7	
500	10735	54.6	
600	16437	59.2	
670	20677	61.7	

(T <sub>Ref</sub> = 273.15 K)				
T (K)	EMF (μV)	dV/dT (μV/K)		
1.2	-5299.6	8.98		
2	-5292.0	10.1		
3.2	-5278.9	11.6		
4.2	-5266.8	12.6		
10	-5181.8	16.0		
20	-5014.0	17.0		
30	-4846.4	16.6		
40	-4681.5	16.5		
50	-4515.8	16.7		
75	-4084.6	17.8		
100	-3627.0	18.8		
150	-2645.2	20.4		
200	-1600.1	21.4		
250	-512.81	22.0		
300	597.44	22.4		
350	1696.3	21.8		
400	2805.7	22.7		
500	5135.3	23.4		
600	7470.7	23.4		

Thermocouple type Chromel-AuFe (0.07%)

Thermocouple type Chromel-AuFe (0.03%) $(T_{Ref} = 273.15 \text{ K})$			
T (K)	EMF (μV)	dV/dT (μV/K)	
3.5	-4671.4	16.1	
4.2	-4660.1	16.0	
10	-4570.7	14.9	
20	-4427.2	13.9	
30	-4290.7	13.5	
40	-4156.0	13.5	
50	-4020.0	13.7	
75	-3664.7	14.8	
100	-3281.4	15.9	
150	-2430.8	18.1	
200	-1480.7	19.8	
250	-471.53	20.4	
300	544.06	20.2	
350	1554.9	20.4	
400	2589.5	21.0	

$(T_{Ref} = 273.15 \text{ K})$			
T (K)	EMF (μV)	dV/dT (μV/K)	
4.2	-5075.3	15.3	
10	-4983.8	16.3	
20	-4811.6	18.1	
30	-4624.8	19.2	
40	-4431.5	19.4	
50	-4239.2	19.0	
75	-3785.8	17.4	
100	-3357.1	17.3	
150	-2436.2	19.4	
200	-1467.7	19.3	
250	-469.66	20.5	
300	503.22	17.8	
350	1493.9	23.5	

Thermocouple type Chromel-AuFe (0.15%)



## Appendix H: Common Units and Conversions

## **Temperature**

Fahrenheit to Celsius:  $^{\circ}$ C = ( $^{\circ}$ F-32)/1.8 Celsius to Fahrenheit:  $^{\circ}$ F = (1.8 ×  $^{\circ}$ C) + 32 Fahrenheit to Kelvin: convert  $^{\circ}$ F to  $^{\circ}$ C, then add 273.15 Celsius to Kelvin: add 273.15

#### Volume

1 liter (I) =  $1.000 \times 10^{-3}$  cubic meters (m<sup>3</sup>) = 61.02 cubic inches (in<sup>3</sup>)

#### Mass

1 kilogram (kg) = 1000 grams (g) = 2.205 pounds (lb)

#### **Force**

1 newton (N) = 0.2248 pounds (lb)

#### **Electric resistivity**

1 micro-ohm-centimeter ( $\mu\Omega$ ·cm)

- =  $1.000 \times 10^{-6}$  ohm-centimeter ( $\Omega$ -cm)
- =  $1.000 \times 10^{-8}$  ohm-meter ( $\Omega \cdot m$ )
- = 6.015 ohm-circular mil per foot ( $\Omega$ ·circ mil/ft)

#### **Heat flow rate**

- 1 watt (W) = 3.413 Btu/h
- 1 British thermal unit per hour (Btu/h) =  $0.2930 \, W$

#### A Note on SI

The values in this catalog are expressed in International System of Units, or SI (from the French Le Système International d'Unités). Whenever possible, the common CGS or British equivalent has been parenthetically included as well.

These common conversions and constants have been included as a reference.

Please refer to NIST Special Publication 811 "Guide for the Use of the International System of Units (SI)" for further standards and conversions.

## References:

Barry N. Taylor, NIST Special Publication 811, 1995 Edition, Guide for the Use of the International System of Units (SI), Washington, U.S. Government Printing Office, April 1995.

The NIST Physical Constants webpage (http://www.nist.gov/pml/data/physicalconst.cfm)

#### Length

	centimeter (cm)	meter (m)	inch (in)
centimeter (cm)	1	1.000 × 10 <sup>-2</sup>	0.3937
meter (m)	100	1	39.37
inch (in)	2.540	2.540 × 10 <sup>-2</sup>	1

1 micrometer (sometimes referred to as micron) =  $10^{-6}$  m

 $1 \text{ mil} = 10^{-3} \text{ in}$ 

#### Area

	cm²	m²	in²	circ mil
cm <sup>2</sup>	1	10-4	0.1550	$1.974 \times 10^{5}$
m²	10 <sup>4</sup>	1	1550	1.974 × 10 <sup>9</sup>
in²	6.452	$6.452 \times 10^{-4}$	1	$1.273 \times 10^{6}$
circ mil	5.067 × 10 <sup>-6</sup>	5.067 × 10 <sup>-10</sup>	$7.854 \times 10^{-7}$	1

#### **Pressure**

	pascal (Pa	) millibar (mbar)	torr (Torr)	atmosphere (atn	n) psi (lbf/in²)
pascal (Pa)	1	$1.000 \times 10^{-2}$	$7.501 \times 10^{-3}$	$9.868 \times 10^{-6}$	$1.450 \times 10^{-4}$
millibar (mbar)	1.000 × 10	<sup>2</sup> 1	7.502 × 10 <sup>-1</sup>	9.868 × 10 <sup>-4</sup>	1.450 × 10 <sup>-2</sup>
torr (Torr)	1.333 × 10	<sup>2</sup> 1.333 × 10 <sup>0</sup>	1	1.316 × 10 <sup>-3</sup>	1.934 × 10 <sup>-2</sup>
atmosphere (atm)	1.013 × 10	<sup>5</sup> 1.013 × 10 <sup>3</sup>	$7.600 \times 10^{2}$	1	1.470 × 10 <sup>1</sup>
psi (lbf/in²)	6.897 × 10	$6.895 \times 10^{1}$	$5.172 \times 10^{1}$	$6.850 \times 10^{-2}$	1
1 torr (Torr) =	133.332 1.33 0.001316 0.01934	pascal (Pa) millibar (mbar) atmosphere (atm) psi (lbf/in²)	1 pascal (Pa) =	$0.007501$ $9.87 \times 10^{-6}$	millibar (mbar) torr (Torr) atmosphere (atm) psi (lbf/in²)

## **Magnetic induction B**

	gauss (G)	kiloline/in²	Wb/m²	milligauss (mG)	gamma (y)
gauss (G)	1	$6.452 \times 10^{-3}$	10-4	1000	10 <sup>5</sup>
kiloline/in²	155.0	1	$1.550 \times 10^{-2}$	$1.550 \times 10^{5}$	$1.550 \times 10^{7}$
Wb/m²	104	64.52	1	10 <sup>7</sup>	10 <sup>9</sup>
milligauss (mG)	0.001	$6.452 \times 10^{-6}$	10-7	1	100
gamma (y)	10-5	6.452 × 10 <sup>-8</sup>	10-9	0.01	1

 $1 \text{ ESU} = 2.998 \times 10^6 \text{ Wb/m}^2$ 

## **Magnetomotive force**

	abampere-turn	ampere-turn	Gilbert (Gi)
abampere-turn	1	10	12.57
ampere-turn	0.1	1	1.257
Gilbert (Gi)	7.958 × 10 <sup>-2</sup>	0.7958	1

1 pragilbert =  $4\pi$  ampere-turn

1 ESU =  $2.655 \times 10^{-11}$  ampere·turn



## **Magnetic field strength H**

	abampere-turn/cm	ampere-turn/cm	ampere-turn/in	ampere-turn/m	oersted (Oe)
abampere-turn/cm	1	10	25.40	1000	12.57
ampere-turn/cm	0.1	1	2.540	100	1.257
ampere-turn/in	3.937 × 10 <sup>-2</sup>	0.3937	1	39.37	0.4947
ampere-turn/m	0.001	0.01	2.540 × 10 <sup>-2</sup>	1	1.257 × 10 <sup>-2</sup>
oersted (Oe)	7.958 × 10 <sup>-2</sup>	0.7958	2.021	79.58	1

<sup>1 0</sup>e = 1 Gi

## Energy, work, heat

	Btu	erg	J	cal	kW∙h
British thermal unit	1	$1.055 \times 10^{10}$	1055	252.0	$2.930 \times 10^{-4}$
erg	$9.481 \times 10^{-11}$	1	10 <sup>-7</sup>	$2.389 \times 10^{-8}$	$2.778 \times 10^{-14}$
joule (J)	$9.481 \times 10^{-4}$	10 <sup>7</sup>	1	0.2389	$2.778 \times 10^{-7}$
calorie (cal)	$3.968 \times 10^{-3}$	$4.186 \times 10^{7}$	4.186	1	$1.163 \times 10^{-6}$
kilowatt hour (kW·h)	3413	$3.6 \times 10^{13}$	$3.6 \times 10^{6}$	8.601 × 10 <sup>5</sup>	1

<sup>1</sup> electronvolt (eV) =  $1.602 \times 10^{-19}$  joules (J)

## **Fundamental physical constants**

Quantity	Symbol	Value*	Unit
speed of light in a vacuum	C, C <sub>0</sub>	299 792 458	m · s <sup>−1</sup>
magnetic constant	$\mu_{0}$	$4\pi\times10^{-7} = 12.566370614\times10^{-7}$	$N \cdot A^{-2}$
electric constant $1/\mu_0c^2$	ε <sub>0</sub>	$8.854187817 \times 10^{-12}$	F⋅m <sup>-1</sup>
characteristic impedance of vacuum $\sqrt{\mu_0}/\epsilon_0 = \mu_0 c$	Z <sub>o</sub>	376.730 313 461	Ω
Planck constant	h	6.6260693(11) × 10 <sup>-34</sup>	J·s
in eV · s		$4.13566743(35)\times 10^{-15}$	eV⋅s
h/2π	h	$1.05457168(18) \times 10^{-34}$	J·s
in eV $\cdot$ s		$6.58211915(56)\times 10^{-16}$	eV·s
elementary charge	е	1.602 176 53(14) × 10 <sup>-19</sup>	С
magnetic flux quantum h/2e	Φ <sub>0</sub>	2.067 833 72(18) × 10 <sup>-15</sup>	Wb
Avogadro constant	N <sub>A</sub> , L	6.0221415(10) × 10 <sup>23</sup>	mol <sup>-1</sup>
atomic mass constant $m_u = \frac{1}{12}m(^{12}C) = 1 u$	m <sub>u</sub>	1.660 538 86(28) × 10 <sup>-27</sup>	kg
Faraday constant N <sub>A</sub> e	F	96 485.3383(83)	C · mol⁻¹
molar gas constant	R	8.314472(15)	J · mol⁻¹ · K⁻¹
Boltzmann constant R/N <sub>A</sub>	k	1.380 650 5(24) × 10 <sup>-23</sup>	J · K <sup>-1</sup>
molar volume of ideal gas RT/p			
T = 273.15 K, p = 101.325 kPa	V <sub>m</sub>	$22.413996(39) \times 10^{-3}$	$\mathrm{m}^3\cdot\mathrm{mol}^{-1}$
T = 273.15 K, p = 100 kPa	V <sub>m</sub>	22.710 981(40) × 10 <sup>-3</sup>	$\mathrm{m}^3\cdot\mathrm{mol}^{-1}$
Stefan-Boltzmann constant (π²/60)k⁴/h³c²	σ	5.670 400(40) × 10 <sup>-8</sup>	W · m <sup>-2</sup> · K <sup>-4</sup>
electron volt: (e/C) J	eV	1.60217653(14) × 10 <sup>-19</sup>	J
Bohr magneton eh/2m <sub>e</sub>	$\mu_{\scriptscriptstyle B}$	$927.400949(80) \times 10^{-26}$	J·T⁻¹
in eV · T <sup>-1</sup> [ $\mu_{\text{B}}$ /(J · T <sup>-1</sup> )](e/C)		$5.788381804(39) \times 10^{-5}$	eV · T⁻¹

<sup>\*</sup> Values are shown in their concise form with uncertainty in parentheses. Numbers with uncertainty values are subject to revision. Refer to the NIST Reference on Constants, Units, and Uncertainty website for the latest values

<sup>1</sup> ESU =  $2.655 \times 10^{-9}$  ampere-turn/m

<sup>1</sup> pracersted =  $4\pi$  ampere-turn/m



## Appendix I: Cryogenic Reference Tables

## **Cryogenic heat flow calculations**

The heat flow  $\dot{\mathbf{Q}}$  conducted across small temperature differences can be calculated using the formula:

$$\dot{Q} = -KA \frac{dT}{dx} \cong -KA \frac{\Delta T}{I}$$
 Eqn. 1

where K is the thermal conductivity, A is the cross-sectional area,  $\Delta T$  is the temperature difference, and L is the length of the heat conduction path.

Thermal conduction across significant temperature differences should be calculated using thermal conductivity integrals.

Note that the thermal conductivity and the thermal conductivity integral of a material can depend strongly on composition and fabrication history. Without verification, the data in the accompanying figures should be used only for qualitative heat flow calculations.

Calculating the heat conduction through a body with its ends at greatly different temperatures is made difficult by the strong temperature dependence of the thermal conductivity between absolute zero and room temperature. The use of thermal conductivity integrals (called thermal boundary potentials by Garwin) allows the heat flow to be calculated as

$$\dot{Q} = -G(\theta_2 - \theta_3)$$
 Eqn. 2

where  $\theta$  is the integral of the temperature-dependent thermal conductivity, K, calculated as

$$\theta_1 = \int_0^{\tau_1} KdT$$
 Eqn. 3

and G is a geometry factor calculated as

$$\frac{1}{G} = \int_{x_{1}}^{x_{2}} \frac{dx}{A}$$
 Eqn. 4

where A(x) is the cross sectional area at position x along the path of heat flow. Note that G=A/L in the case of a body of length L and uniform cross-sectional area A.

Equation 1 is only applicable to bodies within which a common thermal conductivity integral function applies.

Reference: R. L. Garwin, Rev. Sci. Instrum. 27 (1956) 826.

Figure 1—Thermal conductivity of selected materials


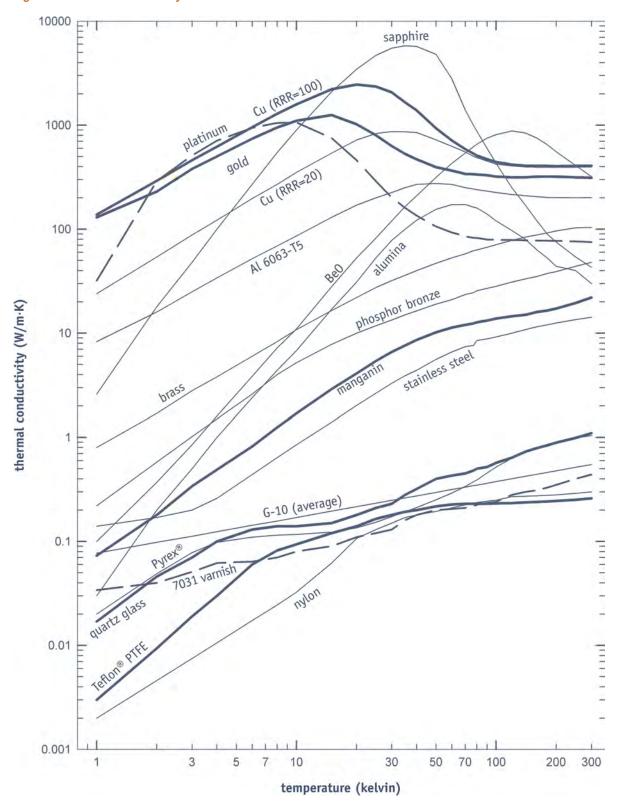




Figure 2—Thermal conductivity integral of selected materials



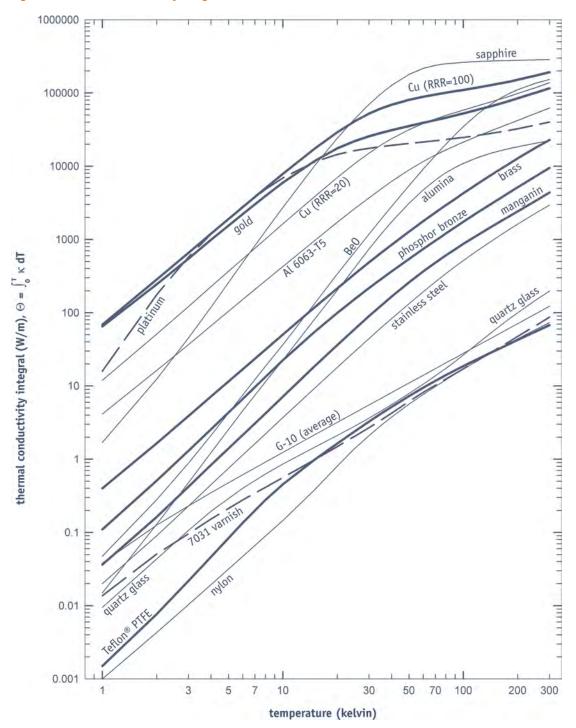

## Table 1—Thermodynamic properties for various cryogenic liquids

	Temperature (K) pressure			Latent heat of vaporization					
	Triple point	Normal boiling point	Critical point	Triple point (kPa)	Critical point (kPa)	Critical density (kg/m³)	L (J/g)	Density (g/ml)	
Helium	2.1768 <sup>a</sup>	4.222	5.1953	5.048	227.46	69.64	20.6	0.13	
Hydrogen	13.8	20.28	32.94	7.042	1283.8	31.36	441	0.07	
Neon	24.5561	27.09	44.44	43.35	2703	483.23	86	1.20	
Nitrogen	63.15	77.36	126.26	12.46	3399	313.11	199	0.81	
Oxygen	54.36	90.19	154.58	0.148	5043	436.14	213	1.14	
Argon	83.8	87.28	150.86	68.9	4906	535.70	162	1.40	
Krypton	115.76	119.77	209.39	73.2	5496	910.75	108	2.40	
Xenon	161.36	165.04	289.74	81.6	5821	1100	96	3.10	
CO <sub>2</sub>	216.58	_	304.21	518.16	7384	466.51	571	1.56	
Methane	90.69	111.63	190.55	11.7	4599	162.65	510	0.42	
Ethane	90.35	184.55	305.33	0.0011	4871	206.73	489	0.55	
Propane	85.47	231.07	369.85	0.1 × 10 <sup>-6</sup>	4248	220.49	425	0.58	
Ammonia	195.49	239.81	406.65	0.0662	11627	237.57	1371	0.68	

<sup>&</sup>lt;sup>a</sup> Triple point values for helium are those of the lambda point

Table 2—Gamma radiation-induced calibration offsets as a function of temperature for several types of cryogenic temperature sensors

	Radiation-induced offset (mK) at temperature					
	Model	4.2 K	20 K	77 K	200 K	300 K
Platinum <sup>b</sup>	PT-103	NA	-15	-10 <sup>d</sup>	10 <sup>d</sup>	10 <sup>d</sup>
Rhodium-iron <sup>b</sup>	RF-100-AA	2 <sup>d</sup>	15 <sup>d</sup>	15 <sup>d</sup>	5 <sup>d</sup>	5 <sup>d</sup>
Cernox®b	CX-1050-SD	-10	-10 <sup>d</sup>	-5 <sup>d</sup>	25 <sup>d</sup>	25 <sup>d</sup>
Carbon-glass <sup>b</sup>	CGR-1-1000	-30	-140	-700	-1300	-3400
Germanium <sup>b</sup>	GR-1400-AA	-5	-20	-25	NA	NA
Ruthenium oxide <sup>b</sup>	R0600	20	150	d	d	NA
GaAlAs diode <sup>b</sup>	TG-120P	-15	-25	2200	2500	400
Silicon diode <sup>b</sup>	DT-470-SD	25	1000	1300	1000	2700
Silicon diode <sup>b</sup>	DT-500P-GR-M	350	50	20	250	300
Silicon diode <sup>b</sup>	SI-410-NN	600	2000	300	450	1400
Platinum <sup>c</sup>	PT-103	NA	-50	5 <sup>d</sup>	50	75
Rhodium-iron <sup>c</sup>	RF-800-4	5 <sup>d</sup>	15 <sup>d</sup>	25	10 <sup>d</sup>	-15 <sup>d</sup>
Rhodium-iron <sup>c</sup>	RF-100-AA	-5 <sup>d</sup>	-5 <sup>d</sup>	5 <sup>d</sup>	-10 <sup>d</sup>	5 <sup>d</sup>
Carbon-glass <sup>c</sup>	CGR-1-1000	-25	-175	-1400	-4200	-6500
Germanium	GR-1400-AA	2 <sup>d</sup>	2 <sup>d</sup>	5 <sup>d</sup>	NA	NA
GaAlAs diodec	TG-120P	-50	-75	700	600	-250
Silicon diode <sup>c</sup>	DT-470-SD	+20	-200	1500	11000	18000
Silicon diode <sup>c</sup>	DT-500P-GR-M	10 <sup>d</sup>	10 <sup>d</sup>	-5 <sup>d</sup>	-5 <sup>d</sup>	-100

Sensors were irradiated in situ at 4.2 K with a cobalt-60 gamma source at a dose rate of 3,000 Gy/hr to a total dose of 10,000 Gy (1  $\times$  10 $^{6}$  rad)

Sensors were irradiated at room temperature with a cesium-137 gamma source at a dose of 30 Gy/hr to a total dose of 10,000 Gy (1 × 10<sup>e</sup> rad)

d Deviations smaller than calibration uncertainty

Table 3—Vapor pressure of some gases at selected temperatures in Pascal (Torr)

	4 K	20 K	77 K	150 K	Triple°point temperature
Water	f	f	f	$1.33 \times 10^{-4}  (10^{-7})$	273 K
Carbon dioxide	f	f	1.33 × 10 <sup>-5</sup> (10 <sup>-8</sup> )	1333 (10)	217 K
Argon	f	$1.33 \times 10^{-10}  (10^{-13})$	21332 (160)	h	84 K
Oxygen	f	$1.33 \times 10^{-10}  (10^{-13})$	19998 (150)	h	54 K
Nitrogen	f	1.33 × 10 <sup>-8</sup> (10 <sup>-11</sup> )	97325 (730)	g	63 K
Neon	f	4000 (30)	g	g	25 K
Hydrogen	1.33 × 10 <sup>-4</sup> (10 <sup>-7</sup> )	101,325 (760)	g	g	14 K

Note: estimates—useful for comparison purposes only (1 Torr = 133.3 Pa)

Table 4—Thermal contraction of selected materials between 293 K and 4 K

	Contraction (per 10⁴)
Teflon®	214
Nylon	139
Stycast® 1266	115
SP22 Vespel®	63.3
Stycast® 2850FT	50.8
Stycast® 2850GT	45
Al	41.4
Brass (65% Cu/35% Zn)	38.4
Cu	32.6
Stainless steel	30
Quartz a-axis	25
Quartz c-axis	10
Quartz mean, for typical transducer	15
Titanium	15.1
Ge	9.3
Pyrex®	5.6
Si	2.2

Table 5—Electrical resistivity of alloys (in μΩ·cm)

	Resistivity (295 K)	(4.2 K)
Brass	7.2	4.3
Constantan	52.5	44
CuNi (80% Cu/20% Ni)	26	23
Evanohm®	134	133
Manganin	48	43
Stainless steel	71 to 74	49 to 51

e Solid and vapor only at equilibrium below this temperature; no liquid

f Less than 10-13 Torr

g Greater than 1 atm

<sup>&</sup>lt;sup>h</sup> Above the critical temperature, liquid does not exist



## Table 6—Defining fixed points of the ITS-90

**Appendices** 

Temperature (T <sub>90</sub> /K)	Substance <sup>i</sup>	State <sup>j</sup>	Defining instrument			
0.65 to 3	ЗНе	Vapor pressure point	He vapor pressure			
3 to 5	Не	Vapor pressure point	thermometer			
13.8033	e-He <sub>2</sub>	Triple point				
~17	e-He <sub>2</sub> (or He)	Vapor pressure point or gas thermometer point		Constant volume		
~20.3	e-He <sub>2</sub> (or He)	Vapor pressure point or gas thermometer point		•		
24.5561	Ne	Triple point				
54.3584	02	Triple point			Platinum	
83.8058	Ar	Triple point			resistance	
234.3156	Hg	Triple point			thermometer	
273.16	H <sub>2</sub> 0	Triple point				
302.9146	Ga	Melting point				
429.7485	In	Freezing point				
505.078	Sn	Freezing point				
692.677	Zn	Freezing point				
933.473	Al	Freezing point				
1234.93	Ag	Freezing point				
1337.33	Au	Freezing point				Radiation
1357.77	Cu	Freezing point				

All substances except 3He are of natural isotopic composition; e-H, is hydrogen at the equilibrium concentration of the ortho- and para-molecular forms

## **Table 7—Saturated vapor pressure of helium**

Table 1 Oc	iturateu vapor	prossure or in	Ciiuiii		
T (K)	P (Pa)	T (K)	P (Pa)	T (K)	P (Pa)
5.1	211600	3.4	41590	1.7	1128
5	196000	3.3	36590	1.6	746.4
4.9	181000	3.2	32010	1.5	471.5
4.8	167000	3.1	27840	1.4	282.0
4.7	154300	3	24050	1.3	157.9
4.6	141900	2.9	20630	1.27	130.7
4.5	130300	2.8	17550	1.24	107.3
4.4	119300	2.7	14810	1.21	87.42
4.3	108900	2.6	12370	1.18	70.58
4.2	99230	2.5	10230	1.15	56.45
4.1	90140	2.4	8354	1.12	44.68
4	81620	2.3	6730	1.09	34.98
3.9	73660	2.2	5335	1.06	27.07
3.8	66250	2.1	4141	1.03	20.67
3.7	59350	2	3129	1	15.57
3.6	52960	1.9	2299		
3.5	47040	1.8	1638		

For complete definitions and advice on the realization of these various states, see "Supplementary Information for the ITS-90"



# **Customer Service**

228 Ordering Information

230 Lake Shore Limited Warranty

231 Sales Offices



## **Ordering Information**

## Shop online 24 hours a day

www.lakeshore.com is your one stop source for the most updated product information, pricing, and applications. Visit it to browse, compare, and purchase products. For your convenience, all the products found in this catalog can be purchased online using a credit card.

### To place an order

Lake Shore products are always being improved. We suggest contacting your local representative, Lake Shore headquarters, or our website for the most current product and pricing information. Please see page 231 to find your local sales representative or dealer.

For maximum efficiency, please be ready to provide us with the following items:

- Purchase order number
- Bill-to address and contact
- Ship-to address
- Name and phone number of purchasing agent
- Name and phone number of technical contact
- Tax status (if applicable)
- Name and model number of products ordered, including any options or accessories
- Line voltage (if applicable)

## **Pricing and quotations**

The price lists that may be included with this catalog are U.S. list prices at the time of approval for printing. Lake Shore reserves the right to change prices without notice. The price for an order is determined by the price prevailing at the time the order is received. Therefore, any prices included with this catalog are intended only for budgetary information. To obtain destination prices, formal quotations, pro forma invoices, or other information before ordering, contact Lake Shore or your local representative.

#### **Product changes**

Product information and illustrations in this catalog were current as of press time. Lake Shore, in a continuing effort to offer excellent products, reserves the right to change specifications, designs, and models without notice. Please visit www.lakeshore.com for the most updated information on products and services.

## Terms and conditions of sale

#### Scope

The following Terms and Conditions (Terms) apply to all sales made by Lake Shore Cryotronics, Inc. (Lake Shore), except where there is a conflict between the Terms and a contractual agreement contained in a purchaser's written Sales Agreement, in which case the Sales Agreement shall control.

#### **Payment**

Visa, MasterCard, and American Express—Lake Shore accepts Visa, MasterCard, and American Express from customers in Canada and the United States and from foreign companies when an American or Canadian financial institution issues the credit card. Credit card orders limited to \$50,000 USD or at the discretion of Lake Shore.

**Inside the U.S.**—Unless otherwise stated in writing, full payment must be received within 30 days of invoice date.

**Outside the U.S.**—Payment terms acceptable to Lake Shore Cryotronics, Inc., unless otherwise agreed upon in writing, are as follows:

- A) Payment with the order.
- B) Payment prior to shipment. Note: acceptable payment methods are bank draft (check) or telegraphic transfer (T/T). A sight draft is not acceptable.
- C) Irrevocable letter of credit. Note: There is an additional documentation processing fee if payment is made with a letter of credit. The purchaser shall be responsible for the additional fee at the time of Payment.

Product lines requiring operational training and verification (OT&V) will be billed in two installments:

- 1. When product is shipped; and,
- 2. When OT&V is performed.

Full payment for product lines requiring OT&V and for any product for which Lake Shore sends the purchaser an invoice under the above stated Payment terms must be received within 30 days of each invoice date, unless otherwise agreed upon in writing.

#### **Taxes**

Unless otherwise stated, all prices that are quoted to the Purchaser are exclusive of use, sales, and similar taxes. The Purchaser shall be responsible to pay Lake Shore the amount of all taxes for which Lake Shore is responsible to remit the tax on a product or service at the time of the Purchaser's Payment, as outlined above. Such taxes that do apply directly to the Purchaser's order shall be paid to Lake Shore by the purchaser, unless an exemption certificate is supplied to Lake Shore that is acceptable to pertinent taxing authorities.

Where Lake Shore is not responsible to remit the tax on sale to pertinent governmental entities, then the Purchaser shall be responsible to remit such tax. Value Added Tax (VAT) and other export taxes are the responsibility of the purchaser.

#### Minimum charge

There is a \$75 minimum charge for orders placed with a purchase order (P0) and a \$5,000 minimum charge for orders placed with a letter of credit (LC). There is no minimum order charge for other payment methods accepted by Lake Shore.



#### **Delivery**

Proposed delivery dates are based upon Lake Shore production schedules at the time of Lake Shore's receipt of the Purchaser's order. However, Lake Shore shall not be liable for any direct or indirect Purchaser damages for default in or delay in delivery. Transportation charges between the Lake Shore factory or warehouse and the Purchaser's delivery point are payable by the Purchaser as a separate invoice item unless otherwise set forth in the quotation and contract in writing.

#### **Shipment**

Hazardous materials: Lake Shore reserves the right to separate the hazardous materials from the rest of the purchaser's shipment and to send hazardous materials directly to the customer via the "Best Way" available, in accord with Lake Shore's discretion, without incurring any liability to Purchaser for default in or delay in delivery of the same to the Purchaser or for its decision to treat the goods as hazardous materials, regardless as to how the purchaser originally scheduled the transport of the goods purchased.

- Lake Shore cannot guarantee delivery of International hazardous materials shipments without conducting an investigation of each destination on a caseby-case basis because such materials are subject to a variety of controls and laws imposed by the various jurisdictions through which the products must pass.
- Purchasers who order hazardous materials and then refuse to accept the hazardous materials shipment will be charged for all shipping costs and a re-stocking fee.

**Risk of loss:** Unless otherwise indicated in writing, the title to Purchaser's goods and risk of loss for such goods pass to Purchaser when the products have been delivered to the Purchaser's transport agent at Lake Shore Cryotronics, Inc.'s plant dock or shipping department in Westerville, Ohio.

- Shipment is not insured by Lake Shore. If customer wishes insurance, this can added at customer's expense.
- The delivery receipt for the Purchaser's shipment will be prima facie evidence
  that the shipment was delivered in compliance with the terms of the sale
  agreement if there is no written notification of rejection of the goods sent by
  Purchaser to Lake Shore (whether such condition is immediately apparent or
  concealed). The Purchaser has the responsibility of verifying the condition of
  the goods on receipt;
- Purchaser must reject purchased goods (whether such condition is immediately apparent or concealed) by notifying Lake Shore in writing within five (5) days from the date of the Purchaser's receipt of goods. No goods may be returned for credit without prior written consent from Lake Shore.

#### Force majeure

Lake Shore shall not be liable for any damage or penalty for delay in or for failure in performance under its sale agreement with Purchaser, including any failure to give notice of delay in delivery due to the weather, acts of God, act of civil or military authority within any government, war, riot, concerted labor action, shortages of materials, or any other causes beyond the reasonable control of Lake Shore.

#### Operational training & verification (OT&V)

The Purchaser warrants that the site where the products are to be installed is in all respects suitable for the safe and lawful installation and operation of the products. The Purchaser shall obtain any certificates or other approvals required prior to

OT&V services and shall inform Lake Shore of all relevant safety, building, and electrical codes and other requirements relevant to the service. Purchaser shall indemnify Lake Shore against any liability or expense resulting from the Purchaser's failure to do so.

If Lake Shore has agreed to effect or supervise the Purchaser's OT&V of Lake Shore's products under the sale agreement, the Purchaser shall prepare the site prior to Lake Shore's arrival in accord with the sale and service agreement and the Purchaser shall provide all services (including labor) for efficient OT&V. If the Purchaser fails to perform its obligations under the OT&V sale and service agreement, Lake Shore shall be permitted to charge the Purchaser for Lake Shore's personnel costs due to the Purchaser's failure.

Upon completion of Purchaser's OT&V, Lake Shore may issue and Purchaser shall review and, if accurate, sign an "Acceptance of Product Certificate," (the Acceptance) within 10 days after the OT&V. The Acceptance shall verify OT&V services have been performed on the purchased equipment. The Acceptance shall be conclusive evidence of the Goods' and OT&V's conformity with the sale agreement contract.

If the Purchaser fails to return the Acceptance, or fails to respond to the issuance of the same to the Purchaser, then the Purchaser's failure to respond to the certificate within the time allotted shall constitute acceptance by the Purchaser.

#### Cancellation

The Purchaser may cancel orders for catalog items if the Purchaser pays Lake Shore's re-stocking charge. (See Returned Goods below.)

A Purchaser's order for custom fabricated or non-catalog products can not be cancelled under any condition, and the Purchaser must remit the purchase price of such product to Lake Shore; the Purchaser and Lake Shore may reach mutually agreed terms in writing regarding other costs or expenses that are attributable to the Purchaser's election to terminate shipment of a product. However, the Purchaser shall be responsible for all expense or costs incurred or obligated by Lake Shore in relation to such product prior to the date and time of Purchaser's notification to Lake Shore of Purchaser's intent. Purchaser shall further indemnify Lake Shore from all such expense or cost

## **Returned goods**

Goods may not be returned to Lake Shore for any reason, except with prior written authorization from Lake Shore. Unless otherwise agreed, authorized returned goods that are not properly rejected for compliance reasons are subject to a 15% re-stocking charge [there is a \$50.00 minimum charge for domestic US Purchasers and a \$60.00 minimum charge for International Purchasers on sensors and other temperature transducers], plus an assessment against the Purchaser of any additional expense required to return received material to first class salable condition.

#### **Export regulations**

Purchaser agrees to comply fully with all laws and regulations concerning the purchase and sale of products. In particular, Purchaser agrees to comply with the Export Administration Regulations of the United States in so far as they apply to the sale, re-sale and transport of products.

As part of Lake Shore's compliance with export regulations, Lake Shore collects and records end-user information to determine whether Lake Shore needs export licenses for its products. As such, all products must be delivered to the agreed upon contract ultimate destination, as shown on Lake Shore sale paperwork. Any in-transit diversion from the agreed upon ultimate destination is prohibited.



## Lake Shore Limited Warranty Statement

## WARRANTY PERIOD: THREE (3) YEARS

- 1. Lake Shore warrants that products manufactured by Lake Shore (the "Product") will be free from defects in materials and workmanship for three years from the date of Purchaser's physical receipt of the Product (the "Warranty Period"). If Lake Shore receives notice of any such defects during the Warranty Period and the defective Product is shipped freight prepaid back to Lake Shore, Lake Shore will, at its option, either repair or replace the Product (if it is so defective) without charge for parts, service labor or associated customary return shipping cost to the Purchaser. Replacement for the Product may be by either new or equivalent in performance to new. Replacement or repaired parts, or a replaced Product, will be warranted for only the unexpired portion of the original warranty or 90 days (whichever is greater).
- Lake Shore warrants the Product only if the Product has been sold by an authorized Lake Shore employee, sales representative, dealer or an authorized Lake Shore original equipment manufacturer (OEM).
- 3. The Product may contain remanufactured parts equivalent to new in performance or may have been subject to incidental use when it is originally sold to the Purchaser.
- 4. The Warranty Period begins on the date the Product ships from Lake Shore's plant.
- 5. This limited warranty does not apply to defects in the Product resulting from (a) improper or inadequate installation (unless OT&V services are performed by Lake Shore), maintenance, repair or calibration, (b) fuses, software, power surges, lightning and non-rechargeable batteries, (c) software, interfacing, parts or other supplies not furnished by Lake Shore, (d) unauthorized modification or misuse, (e) operation outside of the published specifications, (f) improper site preparation or site maintenance (g) natural disasters such as flood, fire, wind, or earthquake, or (h) damage during shipment other than original shipment to you if shipped through a Lake Shore carrier.
- 6. This limited warranty does not cover: (a) regularly scheduled or ordinary and expected recalibrations of the Product; (b) accessories to the Product (such as probe tips and cables, holders, wire, grease, varnish, feedthroughs, etc.); (c) consumables used in conjunction with the Product (such as probe tips and cables, probe holders, sample tails, rods and holders, ceramic putty for mounting samples, Hall sample cards, Hall sample enclosures, etc.); or, (d) non-Lake Shore branded Products that are integrated with the Product.
- 7. To the extent allowed by applicable law,, this limited warranty is the only warranty applicable to the Product and replaces all other warranties or conditions, express or implied, including, but not limited to, the implied warranties or conditions of merchantability and fitness for a particular

- purpose. Specifically, except as provided herein, Lake Shore undertakes no responsibility that the products will be fit for any particular purpose for which you may be buying the Products. Any implied warranty is limited in duration to the warranty period. No oral or written information, or advice given by the Company, its Agents or Employees, shall create a warranty or in any way increase the scope of this limited warranty. Some countries, states or provinces do not allow limitations on an implied warranty, so the above limitation or exclusion might not apply to you. This warranty gives you specific legal rights and you might also have other rights that vary from country to country, state to state or province to province.
- 8. Further, with regard to the United Nations Convention for International Sale of Goods (CISC,) if CISG is found to apply in relation to this agreement, which is specifically disclaimed by Lake Shore, then this limited warranty excludes warranties that: (a) the Product is fit for the purpose for which goods of the same description would ordinarily be used, (b) the Product is fit for any particular purpose expressly or impliedly made known to Lake Shore at the time of the conclusion of the contract. (c) the Product is contained or packaged in a manner usual for such goods or in a manner adequate to preserve and protect such goods where it is shipped by someone other than a carrier hired by Lake Shore.
- 9. Lake Shore disclaims any warranties of technological value or of non-infringement with respect to the Product and Lake Shore shall have no duty to defend, indemnify, or hold harmless you from and against any or all damages or costs incurred by you arising from the infringement of patents or trademarks or violation or copyrights by the Product.
- 10. This warranty is not transferrable.
- 11. Except to the extent prohibited by applicable law, neither Lake Shore nor any of its subsidiaries, affiliates or suppliers will be held liable for direct, special, incidental, consequential or other damages (including lost profit, lost data, or downtime costs) arising out of the use, inability to use or result of use of the product, whether based in warranty, contract, tort or other legal theory, regardless whether or not Lake Shore has been advised of the possibility of such damages. Purchaser's use of the Product is entirely at Purchaser's risk. Some countries, states and provinces do not allow the exclusion of liability for incidental or consequential damages, so the above limitation may not apply to you.
- 12. This limited warranty gives you specific legal rights, and you may also have other rights that vary within or between jurisdictions where the product is purchased and/or used. Some jurisdictions do not allow limitation in certain warranties, and so the above limitations or exclusions of some warranties stated above may not apply to you.
- 13. Except to the extent allowed by applicable law, the terms of this limited warranty statement do not exclude, restrict or modify the mandatory statutory rights applicable to the sale of the product to you.



## Sales Offices

## North America

#### **United States**

Lake Shore Cryotronics, Inc.

575 McCorkle Blvd. Westerville OH 43082 Tel: (614) 891-2244 Fax: (614) 818-1600

e-mail: sales@lakeshore.com

## **West Coast Sales**

CA, OR, and WA

Lake Shore Cryotronics, Inc.

Henderson, NV Contact: Vaden West Tel: (614) 891-2243 ext. 106 Direct line: (614) 212-1594 e-mail: vaden.west@lakeshore.com

## **Western Region Sales**

AK, AK, AZ, CO, HI, ID, MT, NM, KS, ND, NE, NV, OK, SD, TX, UT, WA and WY

Lake Shore Cryotronics, Inc.

Henderson, NV Contact: Vaden West Tel: (562) 436-8398

e-mail: vwest@lakeshore.com

## **Midwest Region Sales**

IA, IL, IN, MI, MN, MO, OH, PA, KY, WI, and WV

Lake Shore Cryotronics, Inc.

575 McCorkle Blvd.
Westerville OH 43082
Contact: Marshall Calhoun
Tel: (614) 891-2243 Ext: 164
e-mail: marshall.calhoun@lakeshore.com

#### **Southeast Region Sales**

AL, AR, GA, LA, SC, TN, MS, and FL

Lake Shore Cryotronics, Inc.

575 McCorkle Blvd. Westerville OH 43082 Contact: Marshall Calhoun Tel: (614) 891-2243 Ext: 164

e-mail: marshall.calhoun@lakeshore.com

## Northeast Region Sales

NY, NJ, PA, NC, VA, Washington DC, MD, and DE

Lake Shore Cryotronics, Inc.

575 McCorkle Blvd. Westerville OH 43082 Contact: Marshall Calhoun Tel: (614) 891-2243 Ext: 164

e-mail: marshall.calhoun@lakeshore.com

MA, ME, RI, VT, NH, and CT

Shain Associates, Inc. 958 Forest Street Marshfield MA 02050 Contact: Dave Shain Tel: (781) 771-2725 Fax: (614) 818-1600

e-mail: shaininc@comcast.net

## Canada

Phoenix Scientific, Inc.

24 Parfield Drive Toronto, Ontario M2J 1B9 Canada Contact: Mr. Albert Young

Tel: (647) 774-3649 Fax: (905) 829-0149

e-mail: sales@phoenix-scientific.ca

#### Mexico

**Valley Research Mexico** 

Canahutli 417 04369, Ciudad de Mexico Mexico

Contact: Rodolfo Carrera Tel: 512-453-0310 Fax: 512-453-0547

e-mail: rodolfo@valleyresearch.com

## South America

#### **Argentina/Brazil**

ANALOG-LAB Serviços e Tecnologia Ltda

Rua Nazira 74, Cotia-SP, CEP 06708-150, Brazil Tel: +55 11 4612 4387 e-mail: info@analoglab.com.br

## **Africa**

#### **Egypt**

**Advanced Biochemicals Company (ABCO)** 

65/8 Al-Mogawra Al-Tania, Al-Hay Al-Tamen, 6 October Cairo, Egypt

Contact: Dr. Hisham Abd El Hamid Tel: (202) 38384690 – 38384691

Fax: (202) 38384691 e-mail: abco@abco-eg.com

## **Rest of Africa**

Lake Shore Cryotronics, Inc.

575 McCorkle Blvd. Westerville OH 43082 Tel: (614) 891-2244 Fax: (614) 818-1600

e-mail: sales@lakeshore.com

## Asia-Pacific

#### Australia, New Zealand

#### **Coherent Scientific Pty. Ltd.**

116 Sir Donald Bradman Drive, Hilton, South Australia 5033 Australia

Tel: +61 8 8150 5200 Fax: +61 8 8352 2020

e-mail: sales@coherent.com.au

## **People's Republic of China**

### **LinkPhysics Corporation**

Room 1907, No.1, Lane 58, New Sports Plaza, Putuo District, Shanghai 200061 Contact: Tina Huang Tel: (86) 021-51088260 Fax: (86) 021-51088158

e-mail: xthuang@linkphysics.com

#### **Republic of China**

#### **Omega Scientific Taiwan**

3F., No. 127, Minshan Street, NeiHu District, Taipei City 11494, Taiwan (R.O.C.)

Contact: Steve Wang Tel: +886-2-2793-0033 Fax: +886-2-2791-5570

e-mail: omega001@ms3.hinet.net

#### India

## Specialise Instruments Marketing Company

305, Kails Industrial Complex, A-Wing, 3rd Floor, Building No.2, Parksite, Vikroli (West), Mumbai-400079 Contact: Sanjiv D.Deshpande Tel: +91 22 25171922/23 Mobile: +91 9820093150 e-mail: sanjiv@simc.co.in

#### Israel

#### El-Tan Technologies Ltd.

1 Coach Street Moshav-Salit, Israel 4588500 Contact: Tomer Shelef Tel: 972-9-740 5551

Fax: 972-9-740 6690 e-mail: tomer.s@el-tan.co.il

#### Japan

#### **Toyo Corporation**

1-6, Yaesu 1-chome, Chuo-ku, Tokyo 103-8284, Japan Contact: Tetsuro Nishida Tel: 81-3-3279-0771 Fax: 81-3-5205-2030 e-mail: lakeshore@toyo.co.jp

## Malaysia

#### Nanorian Technologies Sdn. Bhd.

No. 40 & 40-1 Jalan Kajang Perdana 3/2 Taman Kajang Perdana, 43000 Kajang,

Selangor, Malaysia

Contact: Mr. Mohd Yazid Ahmad

Tel: (+603) 8740 6787 Fax: (+603) 8741 6587

email: yazid@nanoriantech.com

#### S.Korea

#### **ASK Corporation**

Room #1101, Anyang Trade Center #161, Simindae-ro, Dongan-Ku, Anyang-city,

Kyunggi-do, Korea Contact: Henry Kim Tel: (82) 31-451-5600 Fax: (82) 31-451-5605 e-mail: ask@askcorp.co.kr

## **Thailand, Singapore**

## APP Systems Services Pte Ltd.,

11 Toh Guan Road East,
APP Enterprise Building, #03-01
Singapore 608603
Contact: Ronnie Tan
Tel: +65-64256611
Fax: +65-606616
email: ronnie.tan@appsystems.com.sg

## For all Asia-Pacific countries not listed

#### Lake Shore Cryotronics, Inc.

575 McCorkle Blvd. Westerville OH 43082 Contact: Nelson Chen Tel: (614) 891-2243 Ext: 107 Fax: (614) 818 1600

e-mail: nelson.chen@lakeshore.com

## Europe

### **Bulgaria & Republic of Macedonia**

#### **MultiM Systems E00D**

25, Nezabravka Street, District Izgrev, Park Hotel "Moskva", floor 2, room 204 BG - 1000, Sofia, Bulgaria Contact: Dragos Ofrim, PhD EE Tel/Fax: +359 2 870 3004 Mobile: +359 877 00 44 64 e-mail: office@ofrim.bg

## **Czech and Slovak Republic**

#### InterNET SRL

Piata Lahovari nr. 5A, et. 4, ap. 10, sector 1 R0 - 010464, Bucharest, Romania Contact: Dragos Ofrim, PhD EE Tel: +40.21.310.71.16/17 Fax: +40.21.312.16.62/63 e-mail: internet@inter-net.ro

#### **Denmark**

#### **AAGE Christensen A/S**

Skelmosevej 10, 2500 Valby, Denmark Contact: Martin Birch Jensen Tel: 45 3644 2444 Fax: 45 3644 2024

e-mail: mbj@aagechristensen.dk

#### **Estonia**

## Labochema Eesti OU

Aleksandri str.8 EE-51004 Tartu Contact: Mr. Juri Torma Tel: +372 7311131 e-mail: Jyri.torma@labochema.ee

#### **France and Belgium**

#### Cryoforum

52 rue Paul Doumer 78510 Triel sur Seine, France Contact: Phillipe Benoist Tel: 33 1 39 74 02 34 Fax: 33 1 39 27 75 10 e-mail: infos@cryoforum.com



#### Germany, Austria, Sweden & Finland

## **Cryophysics GmbH**

Dolivostrasse 9 D-64293 Darmstadt Germany Contact: Dr. Detlef Cieslikowski Tel: 49 6151 815 70

Fax: 49 6151 815 799 e-mail: info@cryophysics.de

## **Holland and Belgium**

#### Hositrad Holland BV,

Postbus 114
3870 CC Hoevelaken, Holland
Contact: Jurgen Tomassen
Tel: 31 33 253 7210
Fax: 31 33 253 5274
e-mail: info@hositrad.nl

#### **Hungary, Croatia, Slovenia**

#### **SciTEST Solutions Kft**

Szinyei Merse Pál street, no. 21, ap.5, district VI H-1063, Budapest, Hungary

Contact: Mrs. Renata Szabo Mobile: +40 756 642 128 Tel/Fax: +36 20 964 8257 e-mail: office@scitest.hu

## Latvia

## Labochema Latvija SIA

Dzerbenes str.27-260 LV-1006 Riga Contact: Mr. Ingus Rozeks Tel: +371 675 53686

e-mail: ingus.rozeks@labochema.lv

#### Lithuania

#### Labochema LT UAB

Vilkpedes str.22 LT-03151 Vilnius Contact: Mr. Simonas Rudys Tel: +370 5 2750815

e-mail: simonas.rudys@labochema.lt

#### **Poland**

### **Cryo-Tech International**

ul. Diamentowa 3/5 05-077 Warszawa - Wesola 4 Poland Contact: Zbigniew Joachimiak

Tel: 48 22 773 1847 Fax: 48 22 773 1602 e-mail: cryo@tech.x.pl

## Romania, Republic of Moldova, Serbia

#### InterNET SRL

Piata Lahovari nr. 5A, et. 4, ap. 10, sector 1 R0 - 010464, Bucharest, Romania Contact: Dragos Ofrim, PhD EE Tel: + 40.21.310.71.16/ 17 Fax: +40.21.312.16.62/ 63 e-mail: internet@inter-net.ro

#### Russia

#### **Cryotrade Engineering**

Office 36, Proezd Polesskiy, 16, building.1 Moscow, 125367, Russia Contact: Max Klenov Tel: +7 (495) 374-69-52 Fax: +7 (499) 374-69-52 e-mail: sales@cryotrade.ru

#### Italy

#### InterNET SRL

Piata Lahovari nr. 5A, et. 4, ap. 10, sector 1 R0 - 010464, Bucharest, Romania Contact: Dragos Ofrim, PhD EE Tel: + 40.21.310.71.16/ 17 Fax: +40.21.312.16.62/ 63 e-mail: internet@inter-net.ro

#### Spain, Portugal

## **TeraTorr Technologies SL**

C/Las Fábricas nº1, Local 12 28923, Alcorcón (Madrid) Tel: +34 910 298 111 e-mail: info@teratorr.com

#### **Switzerland**

#### **TECO René Koch**

Chemin du Paradis 10 CH-1807 Blonay, Switzerland Contact: René Koch

Tel: 41 21 943 77 90 Fax: 41 21 943 77 92

e-mail: info@teco-rene-koch.com

## **U.K.**, Ireland

#### **Elliot Scientific Ltd**

3 Allied Business Centre, Coldharbour Lane, Harpenden, Hertfordshire, AL5 4UT, UK Contact: lan Perry

Tel: +44 (0)1582 766300 Fax: +44 (0)1582 766340

e-mail: ian.perry@elliotscientific.com

#### For all European countries not listed

#### Lake Shore Cryotronics, Inc.

575 McCorkle Blvd. Westerville OH 43082 Contact: Shane Hritz Tel: (614) 891-2243 Ext: 112

Fax: (614) 818 1600

e-mail: shane.hritz@lakeshore.com



## **Cryogenic Sensors, Instruments, and Accessories**

**Temperature Sensors** 

AC Resistance Bridge

**Temperature Controllers** 

**Temperature Monitors** 


Temperature Transmitters

Programmable DC Current Source

Superconducting Magnet Power Supply

Cryogenic Accessories

**Reference Materials** 

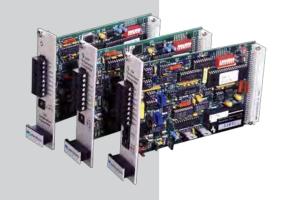





Lake Shore Cryotronics, Inc. 575 McCorkle Blvd. Westerville, OH 43082 Tel 614-891-2244 Fax 614-818-1600 info@lakeshore.com www.lakeshore.com

